




## INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

<b>(51) International Patent Classification <sup>4</sup> :</b>  <b>A61B 5/02</b>	<b>A1</b>	<b>(11) International Publication Number:</b> <b>WO 86/ 07248</b>  <b>(43) International Publication Date:</b> 18 December 1986 (18.12.86)
<b>(21) International Application Number:</b> PCT/US86/01193 <b>(22) International Filing Date:</b> 30 May 1986 (30.05.86)  <b>(31) Priority Application Number:</b> 742,088 <b>(32) Priority Date:</b> 5 June 1985 (05.06.85) <b>(33) Priority Country:</b> US  <b>(71) Applicant:</b> BOARD OF TRUSTEES OF UNIVERSITY OF ILLINOIS [US/US]; Administration Bldg., 506 S. Wright St., Urbana, IL 61801 (US).  <b>(72) Inventors:</b> ARAI, Makoto, R. ; 4950 North Marine Drive, Apartment 1209, Chicago, IL 60640 (US). McALPINE, Laura, E. ; 1543 West Fargo, Chicago, IL 60626 (US). GORDON, David ; 5141 South Kimbark, Chicago, IL 60615 (US).		<b>(74) Agent:</b> GRUBER, Lewis, S.; Marshall, O'Toole, Gerstein, Murray & Bicknell, Two First National Plaza, Suite 2100, Chicago, IL 60603 (US).  <b>(81) Designated States:</b> AT (European patent), AU, BB, BE (European patent), BG, BR, CF (OAPI patent), CG (OAPI patent), CH (European patent), CM (OAPI patent), DE (European patent), DK, FI, FR (European patent), GA (OAPI patent), GB (European patent), HU, IT (European patent), JP, KP, KR, LK, LU (European patent), MC, MG, ML (OAPI patent), MR (OAPI patent), MW, NL (European patent), NO, RO, SD, SE (European patent), SN (OAPI patent), SU, TD (OAPI patent), TG (OAPI patent).  <b>Published</b> <i>With international search report.</i>

**(54) Title:** METHODS AND APPARATUS FOR MONITORING CARDIOVASCULAR REGULATION USING HEART RATE POWER SPECTRAL ANALYSIS

**(57) Abstract**

A patient monitor (4) having an electrocardiographic signal source (2) and an electroplethysmographic respiratory signal source (3) provides inputs to an ECG trigger circuit (5) and an analog-to-digital interface respectively which in turn provide data and control signals to a personal computer (7) programmed to automatically correct the data for artifacts and analyze the spectral densities of the signals which are then shown on display (9).

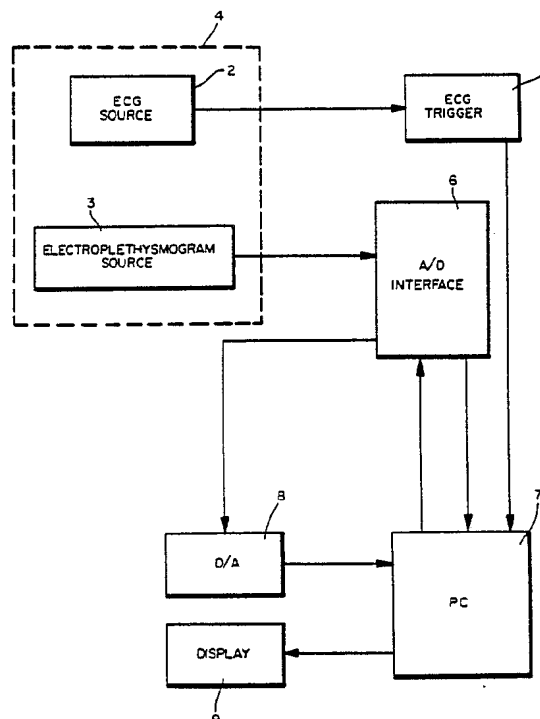

***FOR THE PURPOSES OF INFORMATION ONLY***

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	GA	Gabon	MR	Mauritania
AU	Australia	GB	United Kingdom	MW	Malawi
BB	Barbados	HU	Hungary	NL	Netherlands
BE	Belgium	IT	Italy	NO	Norway
BG	Bulgaria	JP	Japan	RO	Romania
BR	Brazil	KP	Democratic People's Republic of Korea	SD	Sudan
CF	Central African Republic	KR	Republic of Korea	SE	Sweden
CG	Congo	LI	Liechtenstein	SN	Senegal
CH	Switzerland	LK	Sri Lanka	SU	Soviet Union
CM	Cameroon	LU	Luxembourg	TD	Chad
DE	Germany, Federal Republic of	MC	Monaco	TG	Togo
DK	Denmark	MG	Madagascar	US	United States of America
FI	Finland	ML	Mali		
FR	France				

- 1 -

METHODS AND APPARATUS FOR MONITORING  
CARDIOVASCULAR REGULATION USING HEART  
RATE POWER SPECTRAL ANALYSIS

5

Background of the Invention

10           The present invention relates in general to  
methods and apparatus for monitoring cardiovascular  
regulation and in particular to methods and apparatus  
for heart rate spectral analysis.

          Changes in cardiovascular regulation  
15   associated with congestive heart failure include  
attenuation of activity in the parasympathetic division  
of the autonomic nervous system, enhancement of activity  
in the sympathetic division of the autonomic nervous  
system, cardiac catecholamine depletion, down regulation  
20   of the beta-receptor system, increased renin-angiotensin  
system activity, and alteration of baroreceptor  
function. All of these regulatory changes require  
either specific clinical manipulations, such as a stress  
test, a Valsalva maneuver, or the like, and/or invasive  
25   maneuvers, such as cardiac biopsy, plasma catecholamine  
measurement, or the like, in order to determine the  
extent of regulatory dysfunction and its impact upon the  
clinical state of the patient and upon prognoses for the  
patient. These procedures are time consuming, and  
30   generally do not permit the formation of a clinical  
judgment and subsequent action within the timeframe of  
the course of treatment for critically ill patients in  
an Intensive Care Unit.

          Fluctuations from heartbeat to heartbeat in  
35   measured properties of the circulatory system reflect  
both the presence of a variety of naturally occurring

- 2 -

physiological disturbances of the circulatory system homeostasis, and the dynamic response of cardiovascular control systems to these disturbances. For example, the cyclic variation in intrathoracic pressure which  
5 accompanies breathing mechanically affects the return of venous blood to the heart and also affects blood pressure in pulmonary vessels and in the aorta. The variation in intrathoracic pressure is also coupled to a cyclic variation in heart rate through a neural  
10 mechanism mediated by the central nervous system. Furthermore, the resulting cyclic variation in arterial blood pressure impinges on heart rate through a reflex, known as the baroreceptor reflex, which is mediated by the autonomic nervous system. Disturbances in  
15 cardiovascular homeostasis also occur with fluctuations in the resistance of peripheral blood vessels as vascular beds regulate local blood flow to match supply with demand. These fluctuations in peripheral resistance may perturb central blood pressure and,  
20 through the baroreceptor reflex, may also lead to a compensatory variation in heart rate.

Many types of medical instruments exist for studying heart rate variability. The instantaneous rate-meter is perhaps the earliest such instrument.  
25 This meter measures each RR interval through analog or digital circuitry and displays the instantaneous heart rate.

An improvement in the rate-meter is achieved by performing first order statistical evaluation on the  
30 RR-intervals. With mini- and micro-computer systems, histogram displays of RR-interval differences may be generated along with their mean and standard deviations.

Another technique for heart rate variability analysis involves the study of spectral content of the  
35 instantaneous heart rate time series. In one approach to spectral analysis in animals, the computations are



- 3 -

done on a computer. Akselrod, et al., Science, 213, 220-222 (1981) Hyndman, et al., Automedica, 1, 239-252 (1975). Such systems analyze data recorded on magnetic or punched tape. However, not only do these systems  
5 introduce additional errors during the recording process, they do not perform in real time. Furthermore, these systems are not multichannel in nature.

A Sparse Discrete Fourier Transform algorithm which may be implemented on a personal computer (CBM  
10 2016) and which may perform on-line monitoring of heart rate variability, based on a low pass filtered cardiac event series is disclosed in Rompelman, et al., IEEE Trans. Biomed. Engineering, BME-29, 503-510 (1982). A specialized hardware device also exists for low pass  
15 filtering the cardiac event series by a stepwise convolution to create the low pass filtered cardiac event series. Coenen, et al., Medical and Biological Engineering and Computing, 15, 423-430 (1977). Nevertheless, these instruments possess a limited band  
20 width and a limited frequency resolution capability.

There exists a need for an instrument which provides multi-channel spectral analysis of an instantaneous heart rate and of a respiratory activity time series. There also exists a need for an instrument  
25 wherein such calculations are performed in real time at the bedside.

#### Summary of the Invention

30 An apparatus according to the present invention corrects artifacts in a series of heartbeats. Means for collecting a series of heartbeat samples are coupled to means for determining a mean interval between heartbeats. Means for identifying a  
35 mean variance among the intervals between heartbeat samples are coupled to means for establishing an

- 4 -

acceptable of slewing rates as a function of the mean variance. Means for particularizing the absolute value of the slewing rate of a heartbeat sample relative to the mean interval are coupled to the means to  
5 determining and means for substituting the mean interval between heartbeats for all heartbeat interval samples having an absolute outside the range of acceptable slewing rates are coupled to the means for particularizing.

10 A method according to the present invention corrects artifacts in a series of heartbeats. A series of heartbeat interval samples is collected and an appropriate interval between heartbeats is determined. Variances in the intervals between heartbeats are  
15 identified and an acceptable range of slewing rates is established as a function of a mean variance. An absolute value of the slewing rate of a heartbeat sample relative to the mean interval is particularized. An appropriate interval is substituted for all heartbeat  
20 interval samples having an absolute value outside the range of acceptable slewing rates.

Apparatus according to the present invention calibrates a heart rate power spectrum monitor. Means for supplying a signal simulating a heart rate, means  
25 for generating a signal simulating a respiratory frequency fluctuation in heart rate and means for providing a signal simulating a low frequency fluctuation in heart rate are coupled to means for applying signals from these means to a heart rate power  
30 spectrum analyzer.

Apparatus according to the present invention performs heart rate fluctuation power spectral analysis. Means for providing an electrocardiogram signal and means for supplying electroplethysmogram  
35 signal are coupled to means for obtaining a heart rate fluctuation power spectrum from an electrocardiogram

- 5 -

signal and from an electroplethysmogram signal. Real time means for displaying a heart rate fluctuation power spectrum are coupled to the means for obtaining.

Apparatus according to the present invention  
5 trends heart rate fluctuation power spectral data.  
Means for providing an electrocardiogram signal and the means for supplying an electroplethysmogram signal are coupled to means for obtaining a heart rate fluctuation power spectrum from an electrocardiogram signal and from  
10 an electroplethysmogram signal. Means for storing heart rate fluctuation power spectral data are coupled to means for obtaining. Addressable means for transmitting stored heart rate fluctuation power spectral data are coupled to the means for storing and means for  
15 converting heart rate fluctuation power spectral data into graphic form are coupled to the addressable means for transmitting. Real time means for displaying heart rate fluctuation power spectra are coupled to the means for converting.

20 A method according to the present invention treats conditions related to malfunctions of the cardiovascular control system. A power spectrum of heart rate fluctuations in the patient are monitored. A level below about  $0.1 \text{ (beats/min.)}^2$  in the power  
25 spectrum of heart rate fluctuations is identified at a frequency between about 0.04 and about 0.10 Hz as indicative of cardiovascular instability. Procedures are applied to treat the condition and thereby to increase the level of heart rate fluctuations at a  
30 frequency between about 0.04 and about 0.10 Hz.

A method according to the present invention treats conditions related to malfunctions of the cardiovascular control system in a patient. A power spectrum of heart rate fluctuations is monitored in the  
35 patient. A marked increase to above about  $10 \text{ (beats/min.)}^2$  in heart rate fluctuations at a frequency

- 6 -

between about 0.04 to about 0.10 Hz is identified as indicative of cardiovascular stress. Procedures are applied to treat the condition and thereby to decrease the level of heart rate fluctuations between about 0.04  
5 and about 0.10 Hz.

Yet another method according to the present invention treats conditions related to malfunctions of the cardiovascular control system in a patient. A power spectrum of heart rate fluctuations in the patient is  
10 monitored. A ratio of the area under a heart rate power spectrum peak at a frequency between about 0.04 and 0.10 Hz to the area under a peak in the respiratory power spectrum centered at the mean respiratory rate about 0.1 Hz is identified as having an absolute value less than  
15 2.0 for longer than or equal to about one hour as indicating of cardiac instability. Procedures are applied to treat the condition and thereby to increase the ratio.

Still another method according to the present invention treats conditions related to malfunctions of the cardiovascular control system in a patient. A power spectrum of heart rate fluctuations in the patient is  
20 monitored. A ratio of the area under a heart rate power spectrum peak at a frequency between about 0.04 and 0.10 Hz to the area under a peak in the respiratory power spectrum centered at the mean respiratory rate about 0.1 Hz is identified as having an absolute value greater than or about 50 as indicating of cardiac instability. Procedures are applied to treat the condition and  
25 thereby to increase the ratio.  
30

#### Brief Description of the Drawings

Fig. 1 illustrates low frequency,  
35 mid-frequency and high frequency in the power spectrum of heart rate fluctuations in a dog according to the

- 7 -

prior art;

Fig. 2 illustrates aspects of the cardiovascular control system according to the prior art;

5 Fig. 3 is a block diagram of apparatus for heart rate fluctuation power spectral analysis according to the present invention;

Fig. 4 illustrates address buffers and address decoding in a data acquisition device according to the  
10 present invention;

Fig. 5 illustrates components according to the present invention for interfacing an ECG apparatus with a personal computer according to the present invention;

Fig. 6 illustrates a digital to analog  
15 converter according to the present invention;

Fig. 7 illustrates a ECG trigger according to the present invention;

Fig. 8 illustrates a portable calibrator according to the present invention;

20 Figs. 9A and B are halves of a flow chart for software applicable to an embodiment of the present invention on a IBM personal computer;

Fig. 10 illustrates a trend for a stable patient according to the present invention;

25 Fig. 11 illustrates a trend display for an unstable patient according to the present invention;

Fig. 12 is an illustration of an instantaneous heart rate according to the present invention;

Fig. 13 is an illustration of an instantaneous  
30 heart rate fluctuation spectrum of the sort obtainable from apparatus according to the present invention;

Fig. 14 is a stable patient's heart rate fluctuation power spectrum according to the present invention;

35 Fig. 15 is an unstable patient's heart rate fluctuation power spectrum according to the present

- 8 -

invention;

Fig. 16 depicts distributions in LFP data obtained according to the present invention for stable and for unstable patients;

5 Fig. 17 graphically depicts distributions of RFP data according to the present invention for stable and for unstable patients; and

Fig. 18 graphically depicts data for LFP/RFP ratios according to the present invention for stable and  
10 for unstable patients.

#### Detailed Description

Power spectral methods may be used to analyze  
15 the frequency content of fluctuations in heart rate and other hemodynamic parameters. Hyndman, et al., Nature, 233, 339-341 (1971); Sayers, Ergonomics, 16, 17-32 (1973). Short term (i.e., on a time scale of seconds to minutes) fluctuations in these parameters are  
20 concentrated in three principal spectral peaks as illustrated for a canine model in Fig. 1. Akselrod, et al., supra. One peak is centered at the respiratory frequency; this peak shifts with changes in the respiratory rate. The second identifiable spectral  
25 peak, the mid-frequency peak, occurs typically between 0.1 and 0.15 Hz. The oscillations associated with this second peak occur at 6-9 cycles per minute, a considerably lower frequency than the respiratory frequency, and are related to the frequency response of  
30 the baroreceptor reflex. The third peak of the spectrum typically occurs in the frequency band of 0.04 to 0.10 Hz. This low frequency peak is related to thermoregulatory fluctuations in vasomotor tone.

In one approach to the spectral analysis of  
35 heart rate, properties of the heart rate fluctuations in the conscious dog may be related to the activity of

- 9 -

three cardiovascular control systems - the parasympathetic nervous system, the sympathetic nervous system and the renin-angiotensin system. Akselrod, et al., Science, 213, 220-223 (1981). This model is  
5 further elaborated in Akselrod, et al., "Hemodynamic Regulation: Investigation by Spectral Analysis " (In Press). Heart rate fluctuations occurring at frequencies above roughly 0.1 Hz are mediated solely by the parasympathetic system. Blockade of the  
10 renin-angiotensin system leads to a dramatic increase in the amplitude of the low frequency peak. The effects of an autonomic blockade also exist in humans and changes in body posture alter sympathetic-parasympathetic balance as measured by the heart rate power spectrum.  
15 Pomeranz, et al., Am. J. Physiol., 248, H151-H153 (1985).

A simple model of the short term cardiovascular control system is illustrated in Fig. 2. Akselrod, et al., supra. In this model, heart rate  
20 is directly modulated by the sympathetic and parasympathetic nervous systems. Through a variety of receptors both these systems sense, fluctuations in cardiovascular parameters including arterial and venous pressures, vascular volumes, and correlates of blood  
25 flow and oxygenation. The parasympathetic system may respond over a wide frequency range while the sympathetic system may only respond at relatively low frequencies below roughly 0.1 Hz.

A hypothesis was proposed in Akselrod, et al.,  
30 Science, 213, 220-223 (1981), that fluctuations in vasomotor tone associated with the low frequency heart rate fluctuations are not solely related to thermoregulation but also reflect local adjustment to resistance in individual beds of blood vessels in order  
35 to match local blood flow to local metabolic demand. Such fluctuations in peripheral vasomotor tone lead to

- 10 -

fluctuations in central blood pressures which are in turn sensed by pressoreceptors. Stimulation of these pressoreceptors occasions an autonomically mediated baroreceptor reflex, which leads to compensatory fluctuations in heart rate at the corresponding frequency. In addition, the renin-angiotensin hormonal system senses blood pressure fluctuations and, through the elaboration of a substance called angiotensin II, plays the role of the guardian of the overall peripheral vascular resistance. Blockade of the renin-angiotensin system by a converting enzyme inhibitor, may remove this damping influence and may permit increased fluctuations in blood pressure and increased compensatory fluctuations in heart rate in the low frequency regime.

The critically ill infant or child prior to, during, and after cardiac surgery at times exhibits marked changes in heart rate, blood pressure, and peripheral perfusion. These changes may be of no clinical consequence or they may indicate the existence of a major unrecognized pathology whose first outward manifestation may be sudden cardiac arrest. To be able to quantify cardiovascular regulatory reserve permits objective assessment of a patient's cardiovascular stability as well as their response to medical and surgical interventions intended to improve cardiovascular function.

Spectral analysis of tape-recorded records of ECG and respiratory activity from patients with complex congenital heart diseases and myocarditis reveals peculiarities in low frequency heart rate fluctuations not seen in studies of healthy children and adults. In particular: (1) low levels of low frequency heart rate fluctuations are noted for critically ill patients in congestive heart failure, which levels revert to normal after surgical or medical treatment and (2) a marked increase in low frequency heart rate fluctuations is



- 11 -

observed in patients with otherwise undetected cardiac tamponade.

A transitional microprocessor-based monitoring instrument, which utilized a Z-80 microprocessor and a  
5 S-100 bus, was constructed along with a data acquisition system which interfaced the microprocessor with a Hewlett-Packard 78341 patient monitor.

A prototype system is described in Jerome C. Tu, "Microprocessor System for Real-Time Spectral  
10 Analysis Physiological Signals," Master of Department of Electrical Engineering and Computer Sciences, Science Thesis, Massachusetts Institute of Technology (1984). An electrocardiogram (ECG) was inputted into a the data acquisition system from a patient monitor for this  
15 prototype system.

In the data acquisition system, the analog voltage signal of the ECG was applied to the input of a variable frequency voltage-controlled oscillator in the data acquisition system. A counter coupled to the  
20 output of the VCO provided a digital representation of the voltage associated with the ECG peaks. The largest voltage peak, called the R voltage peak and associated in the ECG with ventricular contraction, was used to trigger a clock. Each R peak loaded the value of the  
25 clock into a holding register and restarted the clock. The value of the clock provided a measure of the heart rate as the inverse of the time between beats. (i.e., as the RR interval)

The regular respiratory signal of a patient on  
30 a ventilator was employed to obtain a respiratory spectrum and was similarly obtained through a VCO The respiratory frequency had to be manually entered in order to establish a fixed window for computing the power in the heart rate power spectrum in the  
35 respiratory peak.

Every 256 seconds the digitized ECG RR

- 12 -

intervals were inputted into the microprocessor from the data acquisition system. A smoothed heart rate "tachometer wave form" was created as follows: (1) the instantaneous heart rate time series was computed from the stored RR intervals; (2) a 1024 point time series of the instantaneous heart rate was computed from the stored instantaneous heart rate time series by sampling the latter at 4 Hz; (3) the mean heart rate computed from the 1024-point time series of instantaneous heart rate was subtracted from the smoothed series resulting in a "tachometer waveform". The heart rate power spectrum was computed from the heart rate "tachometer waveform" as follows: (1) a 1024-Point Fast Fourier Transform was computed using 1024 points of the tachometer cardiac tachometer waveform; and (2) the heart rate power spectrum was computed by squaring the absolute value of the previously calculated transform.

As new data was inputted into the computer's buffer, the results of the smoothed cardiac tachometer signal, power spectrum and integral of power spectrum were outputted onto a printer. Thus, for every 256-second time interval, a spectral representation of the preceding 256 seconds of instantaneous heart rate data was exhibited.

From the above data, the area under the low frequency peak (LFP) between 0.04 and 0.1 Hz and the area under the respiratory frequency peak (RFP) within a peak width window of 0.2 Hz were determined. Trend graphs of LFP, RFP, and LFP/RFP ratio were created. The 256 second data segments were rejected if, (1) the patient was not in sinus rhythm; (2) transients and/or artifact were present on the cardiac "tachometer wave form"; and (3) the LFP/RFP ratios were greater than 2 standard deviations from the mean for the study period.

The practical problems associated with this prototype monitoring instrument included the extremely

- 13 -

tedious calculations required for use of the prototype with free-breathing patients and the large amount of data (as much as 50%, in some instances) which had to be discarded due to the presence of motion artifacts.

5 These artifacts resulted from virtually any disturbance of the patient, even a disturbance so slight as holding the patient's hand. The prototype system had no capacity to identify or reject artifacts or to examine the data for dropped beats and premature triggers.

10           Upon reviewing clinical studies performed using the prototype, it was discovered that not only were attenuated low frequency heart rate fluctuations associated with a severely compromised regulatory reserve but also that the ratio of the power in the  
15 heart rate power spectrum at low frequency to the power at the respiratory frequency provided an even sharper discriminatory index between stable and critically ill patients. In addition it was noted that this ratio was markedly elevated in the setting of moderate to severe  
20 congestive heart failure, cardiac tamponade, and prior to the development of malignant ventricular arrhythmias.

A low value for LFP/RFP ( $<2$ ) which is sustained for greater than one hour or a value greater than or about 50 is associated with a clinical course  
25 characterized by cardiac arrest and/or profound hypotension. At times this ratio may be the only clinical indicator of cardiovascular instability. The LFP/RFP ratio provides a sensitive and specific index of cardiovascular instability and may provide a clinically  
30 important, continuous, non-invasive probe of cardiovascular stability.

In order to further examine the diagnostic value of the power spectrum of heart rate fluctuations and to overcome the difficulties with the prototype, a  
35 multipurpose microcomputer-based system, including data basing, instantaneous heart rate and respiratory

- 14 -

activity spectral monitor, was developed using a Hewlett Packard Series 200 Computer and Multiprogrammer as available from Hewlett-Packard. Advantages over the original design include: (1) error correcting routines which correct automatically for motion artifact and missed triggerings of the EKG, thus permitting a substantial increase (>30%) in available data; (2) automated trending of spectral densities along with the instantaneous heart rate and respiratory activity time series; and (3) a data basing program which permits accurate temporal correlation of spectral densities with virtually every clinical intervention, routine ventilatory changes, hemodynamic, fluid monitoring and laboratory results. Software incorporating these advantages is included herein as Appendix A.

In a further improvement, programs and a data acquisition system and programs were developed for use with an IBM PC or compatible personal computer. This improvement is illustrated in Figs. 3 through 12.

In Fig. 3, a block diagram of apparatus according to the present invention is illustrated. In Fig. 3, a source of an ECG signal 2 and a source of an electrolythsmogram signal 3 are contained within a patient monitor 4. A patient monitor for use with the present invention may be the System 2 Infant Monitor available from ARVEE, Incorporated, Battle Creek, Michigan. Source 2 is connected to an ECG trigger 5 which is in turn connected to a personal computer 7. Source 3 is connected to an analog to digital interface 6. Interface 6 is connected to analog converter 8 which is connected in turn to a personal computer 7. Personal computer 7 receives input from and provides output to interface 6. Personal computer 7 is connected to a display 9.

Source 2 receives input from pregelled electrodes adhered to the chest wall and thigh of the

- 15 -

patient. Source senses respiratory activity through a pair of electrodes by the impedance method. Personal computer 7 and display 9 are available as an IBM PC and a compatible display available from IBM, Incorporated, Armonk, New York. Elements 5, 6 and 8 are described below.

In a data acquisition device according to the present invention, address buffers and address decoding, as illustrated in Fig. 4, receive input from a PC bus 10. Nodes 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25 and 26 are respectively connected to address lines A0, A1, A2, A3, A4, A5, A6, A7, A8, A9, A10, A11, A12, A13, A14 and A15 in PC bus 10. A first address buffer 100 has address inputs A0, A1, A2, A3, A4, A5, A6 and A7 which are respectively connected to nodes 11-18. Buffer 100 also has two gate inputs, 1G and 2G, which are connected to ground along with a ground output GND of buffer 100. A power supply input  $V_{CC}$  of buffer 100 is connected to a node 102 at a potential of +5 volts.

A second address buffer 110 has address inputs A8, A9, A10, A11, A12, A13, A14 and A15 which are respectively connected to nodes 19-26. Buffer 110 also has two gate inputs, 1G and 2G, which are connected by way of a node 111 to ground. A ground GND output of buffer 110 is also connected to a common potential. Buffer 110 has a power supply input  $V_{CC}$  which is connected to a node 112 at a potential of +5 volts.

A status buffer 120 has address inputs A16, A17, A18 and A19 which are respectively connected to nodes 27, 28, 29 and 30. Nodes 27-30 are respectively connected to an address enable line AEN, a reset line RES, an input/output read line IOR and an input/output write line IOW in PC bus 10. Buffer 120 has two gate inputs, 1G and 2G, which are connected by way of a node 121 to ground. A ground output GND of buffer 120 is

- 16 -

also connected to ground by way of node 121. A power supply input  $V_{CC}$  of buffer 120 is connected to a node 122 at a potential of +5 volts.

According to the present invention, a data acquisition system board which is both reliable and compatible with a personal computer (PC) bus, preferably adheres to the timing requirements and the loading requirements supplied by the PC bus. This means that all connections to the PC bus should be buffered so that the load provided at any input or output of the bus is equivalent to 1 LS TTL load and high speed CMOS integrated circuits are provided for this purpose.

Because there are multiple devices attached to the address bus, address buffers are provided. This is done by buffers 100 and 110. Parts used for buffers 100, 110 and 120 are normally gated, but the gate enables, 1G and 2G, are tied to ground so that the gates are always enabled. Some of the status lines on the PC bus are buffered by a chip 120, in particular: the reset line RES; the read and write lines IOR and IOW respectively, for the input/output (IO) channels; and the address enable AEN.

An address decoder according to the present invention, as illustrated in Fig. 4, includes a three to eight line decoder 130. Decoder 130 has three line inputs A, B and C which are respectively connected to outputs B2, B3 and B4 of buffer 100. Decoder 130 has gate inputs G2A and G2B which are respectively connected to outputs B5 and B6 of buffer 100. A power supply  $V_{CC}$  input of decoder 130 is connected to a node 131 at a potential of +5 volts while a ground GND output of decoder 130 is connected to a common potential. Outputs Y0, Y1, Y2, Y3, Y4, Y5, Y6 and Y7 are connected to inputs of a NAND gate 140.

A NAND gate 151 has an input connected to each of outputs B8, B9 and B10 of buffer 110. An output B10

- 17 -

of buffer 110 is connected to an input of an inverter 152 which has an output connected to an input of NAND gate 151. Similarly, outputs B12, B13, B14 and B15 of buffer 110 are respectively connected to an input of each of inverters 153, 154, 155 and 156, each of which has an output connected to an input of NAND gate 151. NAND gate 151 has an output connected to an input of an inverter 157.

A NAND gate 158 has an input connected to an output of inverter 157 and has an output connected to an input of an inverter 159. An inverter 160 has an input connected to an output B7 of buffer 100 and has an output connected to an input of NAND gate 158. Likewise, an inverter 161 has an input connected to an output B16 of buffer 120 and has an output connected to an input of NAND 158. An output of inverter 159 is connected to a gate input G1 of decoder 130.

So that devices on the board are recognized at a particular IO channel address, address decoding is provided. In this particular case, a fixed address location, location hex 700 to 71F (a total of 32 channels), is used. The decoding of the fixed upper bytes in the address is provided by a combination of nine inverting gates, 152, 153, 154, 155, 156, 157, 159, 160 and 161, and NAND gates 151 and 158. These elements, in combination with decoder 140, provide chip enable signals which can be used to select one or another of the functional chips on our board. Each of the eight chip enable signals correspond to a block of four channels. For example, a chip select #0 from output to of decoder 130 corresponds to channels hex 700, 701, 702 and 703.

A logic network for driving a data buffer, as illustrated in Fig. 5, includes a NAND gate 171, an inverter 172 and a NAND gate 173. An output of inverter 172 is connected to a first input of NAND gate 173 while

- 18 -

an output of NAND gate 140 is connected by way of a node 174 to a second input of NAND 173 and to a first input of a NAND gate 175. A second input of NAND gate 175 is connected to an output of NAND gate 171.

5           In addition, a node 181 is connected to an output B0 of buffer 100. A node 182 is connected to an output B1 of buffer 100. Nodes 183 and 184 are respectively connected to output Y0 and output Y7 of decoder 130. Nodes 185 and 186 are respectively connected to an  
10           output of NAND gate 175 and an output of NAND gate 173. A node 187 is connected to an output B17 of buffer 120. A node 188 is connected an output B18 of buffer 120, to a first input of NAND gate 171 and to an input of inverter 172. A node 189 is connected to a second  
15           input of NAND 171 and to an output B19 of buffer 120.

          Additional chips are used to provide logic which drives a data buffer connected to a data bus. The data bus is bidirectional in order to both transmit data to and from devices on the board. In order that this be  
20           accomplished, one must determine at any time whether or not data is either being read from or written to the board. This logic is supplied by NAND gate 171, NAND gate 173, AND gate 175 and inverter 172 which translates the read and write signals for the input/output (IO)  
25           channel into an output enable and a transmit enable for a data buffer. The apparatus of Fig. 4 may be used to properly interface a device to the PC bus 10.

          As illustrated in Fig. 5, components according to the present invention for interfacing an ECG  
30           apparatus with a personal computer include a port expander 200. Port expander 200 has four sets of 8 nodes each, the four sets correspond to four ports A, B, C and D. The outputs for port A are A0, A1, A2, A3, A4, A5, A6 and A7. The inputs corresponding to port B are  
35           B0, B1, B2, B3, B4, B5, B6 and B7. Outputs corresponding to port C are C0, C1, C2, C3, C4, C5, C6



- 19 -

and C7. A set of outputs corresponding to port D includes D0, D1, D2, D3, D4, D5, D6 and D7. Expander 200 has a chip select input CS connected to node 184. Expander 200 also has a read input RD and a write input WR respectively connected to nodes 188 and 189. Expander 200 has two address inputs, AD0 and AD1 which are respectively connected to nodes 181 and 182. A reset RES input of expander 200 is connected to node 187. Inputs A0, A1, A2, A3, A4, A5, A6, A7 are respectively connected to nodes 291, 292, 293, 294, 295, 296, 297 and 298. Outputs D0-D7 are respectively connected to nodes 208, 207, 206, 205, 204, 203, 202 and 201 which define a data bus. A power supply input  $V_{CC}$  of expander 200 is connected to a node 209 at a potential of +5 volts. A ground GND output of expander 200 is connected to a common potential.

Port expander 200 is used to overcome the low speed of the data bus on both A/D converter 260 and a digital analog converter. This permits slowing down the read and write signals inasmuch as they may be provided artificially on port C of expander 200 or as chip select signals from address decoder 130. Port C of expander 200 is a bit addressable register which allows one to individually select or deselect bits without affecting any of the other bits. This is accomplished by sending a one byte command to expander 200. Because expander 200 is given the control function, the address of expander 200 is the highest address in the set of channels. In other words, expander 200 occupies IO channels hex 71C to hex 71F. The ports A, B and C on expander 200 are addresses 71C, 71D and 71E, respectively, and the control register internal for expander is at input/output I/O channel 71F.

A timer 220 according to the present invention has two address inputs, AD0 and AD1 respectively connected to nodes 181 and 182. Timer 220 also has a

- 20 -

read input RD connected to node 188, a write input WR connected to node 189 and a chip select input CS connected to node 184. A first gate input G0 is connected to the CO of expander 200 while a second gate  
5 input G1 and a third gate input G2 are both connected by way of a node 223 to output C1 of expander 200. Timer 220 has three clock inputs CLK0, CLK1 and CLK2, of which CLK1 is connected by way of node 222 to an output OUT0 of timer 220 and input CLK2 is connected to an output  
10 OUT1 of timer 221 by way of a node 31. An interrupt request line IRQ4 within PC bus 10 is also connected to node 31.

An output OUT2 is connected to a non-inverting input of an operational amplifier 224, an inverting  
15 input and a output of which are connected to a node 400.

A power supply input  $V_{CC}$  of timer 220 is connected to a node 221 which at a potential of +5 volts.

20 Timer 221 has seven outputs D0, D1, D2, D3, D4, D5, D6 and D7 which are respectively connected to nodes 208, 207, 206, 205, 204, 203, 202 and 201. A ground output of timer 220 is connected to a common potential.

25 Timer 220 includes three 16 bit timers which are addressed at hex locations 704, 705, 706, and 707. In other words, they are provided by chip select 1. The three clocks on timer 220 are connected in series which effectively converts it into a 48 bit counter. However,  
30 in the operation of the program, some of the bits in this counter are thrown away because the reset values are less than 65,536. The three clock registers are used in the following way. Counter 0, corresponding to input CLK 0, counts an onboard time base to be discussed  
35 later and provides an output which gives the minimum resolution of the heart rate counting. In other words,

- 21 -

it provides the counter time base for measuring the heart rate. Counter #1, corresponding to input CLK 1, counts the heart rate counter time base and provides as an output an interrupt at IRQ4. This signal drives the sampling of the respiratory signal at a constant frequency, and is also used to measure interbeat intervals. In the standard data collecting mode, where one is interested in measuring the respiratory signal at 4 hertz intervals, this means that the counter 0 is set to generate output pulses at 11 microsec. intervals and that these pulses are in turn counted by counter 1 to generate 4 hertz pulses which are used to drive data acquisition from the respiratory signal. The last counter register, counter #2, corresponding to input CLK2, is used to count the number of respiratory sampling pulses which have been supplied. This functions as an overflow counter and always has the reset value of 65,536. Thus the counter measuring interbeat intervals effectively overflows only every 65,536 respiratory sampling times, which is far in excess of what would be required to recover dropped beats which occur because the heart rate is not adequately detected.

A counter 240 has an input 1A connected to a clock line PC CLK in PC bus 10 by way of a node 32. Counter 240 has a first output 1QA connected to the CLK0 input of timer 220. Counter 240 has a second output 1QB and has a third output 1QC. A clear input CLR1 of counter 240 and a ground output GND of timer 240 are connected to a common potential by way of a node 242.

A data output buffer 280 has an output enable input OE connected to node 185 and has a transfer enable input TE connected to a node 186. Eight data inputs, A0, A1, A2, A3, A4, A5, A6 and A7, of buffer 280 are respectively connected to nodes 208, 207, 206, 205, 204, 203, 202 and 201. A power supply  $V_{CC}$  input of buffer

- 22 -

280 is connected to a source of potential at +5 volts. A ground GND output of buffer 280 is connected to a common potential. Outputs B0, B1, B2, B3, B4, B5, B6 and B7 of buffer 280 are respectively connected to data lines in PC bus 10 by way of nodes 33, 34, 35, 36, 37, 38, 39 and 40.

The time base for this clock system is provided by counter 240. Timer 220 counts only at a rate of 2.6 MHz megahertz which is exceeded by the IBM PC bus clock of 4.77 megahertz. The IBM PC bus clock is divided by 2 using counter 240 and the result used to provide a time base at 2.38 megahertz for timer 220. The 4.77 megahertz clock is also divided by 8 to provide a 596 kilohertz clock which is used to drive an analog to digital (A/D) converter. A/D converter 260 uses this clock signal in order to properly execute the successive approximation scheme to convert analog inputs into digital outputs.

A/D converter 260 has an output enable input OE connected to output C4 of expander 200. A/D converter 260 also has three inputs A, B and C which are respectively connected to outputs C5, C6 and C7 of expander 200. A clock input CLK of A/D converter 260 is connected to the LQC output of counter 240. An address latch enable ALE and a start input STR of A/D converter 260 are connected to a node 261. A power supply  $V_{CC}$  input and a reference voltage  $+V_{REF}$  input of A/D converter 260 are connected to a node 262 at a potential of +5 volts. A reference voltage  $-V_{REF}$  output and a ground GND output of A/D. converter 260 are connected to a common potential by way of a node 263. A/D converter 260 has seven outputs D0, D1, D2, D3, D4, D5, D6 and D7 which are respectively connected to inputs B0, B1, B2, B3, B4, B5, B6 and B7 of expander 200. In addition, A/D converter 260 has an end of count EOC output connected to a first input of the NAND gate 264, an output of

- 23 -

which is connected to an input of an inverter 265. A second input of NAND gate 264 is connected to an output of an inverter 266 which has an input connected to node 187. An output of inverter 265 is connected to node 261.

A/D converter 260 has a signal input IN connected to a node 267. An output of an operational amplifier 268 is connected to node 267 and to a first lead of a resistor 269. A second lead of resistor 269 is connected to a first lead of resistor 270, a second lead of which is connected to a source of potential at -5 volts. The first end of resistor 270 is also connected to an inverting input of amplifier 268 and to a first end of a resistor 271. A non-inverting input of amplifier 268 is tied to ground. A second end of resistor 271 is connected to a node 272 which provides an analog signal input ANA IN for the apparatus according to the present invention.

A/D converter 260 is connected to port B of port expander 200. This A/D has built into it its own 8 channel analog multiplexer which allows the selection of one of eight analog signals to be converted. The channel select corresponding to inputs A, B and C of converter 260 is connected to port C on bytes 5, 6 and 7.

Because A/D converter 260 operates from 0 to 5 volts, analog input at input IN should be in the range of 0 to 5 volts or an input buffer should be supplied to alter this input range. However, in keeping with general practices for safety and isolation, input IN should always be provided with an analog buffer to provide isolation for both the computer and the instrument being monitored. As illustrated, the input buffer is provided by operational amplifier 268. This amplifier converts a bipolar analog input of plus or minus 5 volts to a single unipolar input of 0 to 5 volts

- 24 -

at input IN. This analog input is used to monitor the respiration.

A/D converter 260 is set up in a free running mode such that it continuously does conversions on the analog signal. The end-of-conversion pulse at output EOC is used to generate a start pulse for the A/D so that as soon as an end of conversion occurs it a new conversion is started. This is the reason for the two gates connected between end of conversion output EOC and the start input STR. In order to prevent latchup of the device on power up, the reset line at node 187 is also used to generate a start pulse. This means that the device will always function even after being powered up. Also, in order to update A/D converter 260 as frequently as possible, the address latch enable ALE, which is used to latch in the address value for the channel to be monitored, is re-latched at every start pulse.

As illustrated in Fig. 6, digital analog (D/A) converter 300 has inputs D0, D1, D2, D3, D4, D5, D6 and D7 which are respectively connected to nodes 298, 297, 296, 295, 294, 293, 292 and 291 as illustrated in Fig. 5. Converter 300 has a write WR input connected to node 183 and has a feedback input RFB. Converter 300 also has a power supply  $V_{CC}$  input, a reference voltage  $V_{REF}$  input and an input latch enable input ILE all of which are connected to a source of potential at +5 volts by way of a node 301. Converter 300 has an analog ground AGND and a digital ground output DGND, both of which are connected by way of a node 302 to a common potential.

Converter 300 has a first output OUT1 and a second output OUT2 which are respectively connected to an inverting and a non-inverting input of an operational amplifier 303. The non-inverting input of amplifier 303 is also connected to a common potential by way of a node 305. Amplifier 303 has an input connected to a node 306

- 25 -

at a potential of +12 volts and an input connected to a node 307 at a potential of -12 volts. An output of amplifier 303 is connected to a node 308 which is connected to the RFB input of converter 300 and to a first end of a variable resistor 309. A second lead of variable resistor 309 is connected to a first lead of a variable resistor 310, a second lead of which is connected to a node 311 at a potential of +5 volts. The second lead of resistor 309 is also connected to an inverting input of operational amplifier 312 and to a first lead of a resistor 313. A non-inverting input of amplifier 312 is connected to ground. A second lead of resistor 313 is connected to an output of amplifier 312 and to a node 391 which serves as an analog output for the apparatus according to the present invention.

Port A of expander which is at location 71C, is attached to a D/A converter data bus which, includes nodes 291-298. The write latch signal for the D/A converter is provided by chip select #0. In other words, any dummy byte written to any of the addresses 700, 701, 702 or 703 hex will cause a write pulse to be sent to D/A converter 300, thereby latching the data on port A of expander 200 into the D/A converter 300 and allowing an analog signal to be generated corresponding to the digital input. The output of D/A converter 300 chip is in the form of differential currents generated at outputs OUT 1 and OUT 2. A system having two operational amplifiers is employed to convert these currents to a voltage. Amplifier 303 is a differential current to voltage converter which provides a signal from 0 to 5 volts. Amplifier 312 converts the signal to a bipolar plus or minus 5 volt signal. Feedback control for the current to voltage converter is provided in D/A converter 300 through input RFB so that in actuality three connections are made from the D/A chip to the first operational amplifier. Because the D/A converter

- 26 -

is an 8 bit device, this provides 256 voltage levels which are linearly distributed between plus and minus 5 volts. This D/A output may be used to generate calibrating signals or other control signals.

5           As illustrated in Fig. 7, a source of an ECG signal is connected by way of a node 400 to a non-inverting input of an operational amplifier 401 in an ECG trigger 60. An input of amplifier 401 is connected to a node 402 at a potential of plus 12  
10 volts. An inverting input of amplifier 401 is connected to an output of amplifier 401 and to a non-inverting input of an operational amplifier 406. A first lead of each of resistors 403a, 403b, 403c, 403d, 403e, 403f, 403g, 403h and 403i is connected to the output of  
15 amplifier 401 while the second lead of resistor 403i is permanently connected and a second lead of one other of resistors 403a through h is connected to a node 410 by a jumper. A first lead of capacitor 404 is connected to node 410 while a second lead of capacitor 404 is  
20 connected to a node 405 at a potential of minus 12 volts. An inverting input of amplifier 406 is connected to a cathode of a diode 407, an anode of which is connected to an output of amplifier 406. The cathode of diode 407 is also connected to a first lead of capacitor  
25 408 and a first lead of each of resistors 410a, 410b, 410c, 410d and 410e, the second lead of resistor 410e is permanently connected and the second lead of one other of which is connected to a node 410 (not shown) by a jumper 411g (not shown). A non-inverting input of an  
30 operational amplifier 412 is also connected to the cathode of diode 407 while an inverting input of amplifier 412 is connected to the output of amplifier 406. An input of amplifier 412 is connected to a node 413 at a potential of minus 12 volts. A first lead of  
35 resistor 414 is connected to the output of amplifier 412 while a second lead of resistor 414 is connected to a



- 27 -

cathode of a diode 415 an anode of which is connected to ground. The cathode of diode 415 is also connected to an input of a Schmitt trigger 416 an output of which is connected to a line designated IRQ 3 in PC bus 10 by way of a node 491.

ECG trigger 60 has an input buffer consisting of a non-inverting buffer of an amplifier 401 which isolates the ECG signal from the rest of the board. As illustrated in Fig. 5, the EKG trigger functions in the following manner. The R wave, which is larger than any other signal in the ECG, causes capacitor 408 to charge up to a certain value corresponding to the peak of the R wave. Any values beneath the peak of the R wave will be rejected by amplifier 403 so that no output occurs. Between R waves, the voltage on capacitor 405 decays slowly with a rate given by the RC time constant of capacitor 405 and the resistance across elements 410a-f. The voltage on the capacitor is sent to the inverting input on amplifier 403 and is used as a threshold for the R wave of the EKG. Therefore, as the electrocardiogram is being passed to the non-inverting input of amplifier 406, the only time that the operational amplifier has a positive output is when the EKG signal is larger than the voltage on capacitor 405. Whenever this occurs, capacitor 408 is immediately charged up to the value at the EKG input. In other words, the voltage on capacitor 408 is a sort of envelope on the top of the electrocardiogram, although its decay rate is limited by the RC time constant. Diode 407 insures that the envelope function which is provided by capacitor 408 is the upper envelope and not the lower envelope. The lower envelope is provided by reversing the polarity of diode 407.

The RC network of capacitor 405 and resistors 403a-i provides a low pass filtered ECG. The voltage on capacitor 405 is the baseline for the ECG, which may

- 28 -

vary. The array of jumper selected resistors 410a-e allows variation of the time constant of the RC network containing resistors 406a-e and capacitor 408. Thus, this latter network which monitors the ECG envelope is referenced to the ECG baseline present on capacitor 404 5 permitting accurate tracking of the envelope and therefore better R wave detection. As a further improvement, the jumpers may be replaced with analog switches controlled by the personal computer in order to 10 give the computer control of RC time constant selection.

An output from ECG trigger 60 is generated by connecting amplifier 412 in parallel with peak detector amplifier 406 so that the inputs are reversed. The result is that the output polarity is inverted. Because 15 the amplifiers 401, 406 and 412 are operating from a plus 12 volts to minus 12 volts supply, but the logic levels on the board are only from 0-5 volts, resistor 414 and a diode 415 are used to clamp the output value of the amplifier 412 between 0 and 12 volts. This 20 signal is then passed to a Schmitt trigger 416, which is a single conditioning device. The output of this signal conditioner is finally provided to PC bus 10 in order to drive interrupts at interrupt request 3 (IRQ3) indicating the currents of an R wave. ECG trigger 60 25 may be modified to allow selection of various decay rates for the envelope and also to provide a floating threshold for the 0 point of the EKG. The ECG triggers if the R wave passes above 0 volts. However, it can be imagined that sometimes the baseline will drift far 30 enough below 0 volts that the R wave does not cross 0 volts and in such a case this trigger would never detect the R wave. This is corrected by connecting the second leads of the charging capacitor 408 and on the selected discharging resistor of 406a-f may be connected to a low 35 pass filter consisting of a capacitor 405 and a selected one of resistors 403a-f (to choose various discharge

- 29 -

rates) which low pass filters the electrocardiograms and essentially selects out the baseline. This means that instead of measuring the R wave with respect to 0 volts, the R wave may be measured with respect to the floating baseline of the electrocardiogram. The jumper selected resistor selects an RC time constant much greater than the RR interval. So long as the baseline does not drift faster than one R wave in approximately 10 heart beats, this means that this trigger will successfully detect all R waves. Selecting one of resistors 410a-f allows variation of the RC time constant of elements 408 and 410a-f.

As illustrated in Fig. 8, in a portable calibrator 70 according to the present invention, an operational amplifier 500 has a non-inverting input connected to a first lead of each of resistors 501, 502 and 503. A second lead of resistor 501 is connected by way of a node 503a to a positive voltage source while a second lead of resistor 502 is connected by way of a node 504 to a negative voltage source. An inverting input of amplifier 500 is connected to a first lead of a capacitor 505, a second lead of which is connected by way of a node 506 to a negative voltage source. The inverting input of amplifier 505 is also connected to a first lead of a variable resistor 507 and to a first lead of a resistor 508 a second lead of which is connected to an output of amplifier 500. The output of amplifier 500 is also connected to a second lead of resistor 503. Amplifier 500 has an input connected by way of a node 509 to a positive voltage source and by way of a node 510 to a negative voltage source.

A second lead of resistor 507 is connected to a non-inverting input of an amplifier 511, an inverting input of which is connected to an output of amplifier 511 by way of a node 591 which provides an output port for a simulated respiratory frequency.

- 30 -

A first lead of a resistor 512 is connected to node 591 while a second lead of resistor 512 is connected to a first lead of a resistor 513 and to a first lead of a capacitor 514, a second lead of which is  
5 connected by way of a node 515 to a negative voltage source. A second lead of resistor 513 is connected to an output of an operational amplifier 514 and to an inverting input of amplifier 515 is connected to a first lead of a resistor 516, to a first lead of a capacitor  
10 517 and to an inverting input of an operational amplifier 518. The second lead of capacitor 517 is connected by way of a node 519 to a negative voltage source. A non-inverting input of amplifier 518 is connected to a first lead of each of resistors 520, 521  
15 and 522. A second lead of resistor 520 is connected by way of a node 523 to a positive voltage source while a second lead of resistor 521 is connected by way of a node 524 to a negative voltage source. A second lead of resistor 522 is connected to an output of amplifier 518  
20 and to a second lead of resistor 516.

An inverting input of an operational amplifier 525 is connected to the first lead of resistor 513 and to a first lead of a variable resistor 526. A non-inverting input of amplifier 525 is connected to a  
25 first lead of each of resistors 527, 528 and 529. A second lead of resistor 527 is connected to a node 530 at a positive potential while a second lead of resistor 528 is connected by way of a node 531 to a negative voltage source. A second lead of resistor 529 is  
30 connected to a second lead of resistor 526 and to an output of amplifier 525 at a node 592 which provides a square wave output simulating a modulated heart rate pulse. A first lead of a capacitor 532 is connected to node 592 while a second lead of capacitor 532 is  
35 connected by way of a node 593 to a first lead of a resistor 533, a second lead of which is connected to

- 31 -

ground. Node 593 provides an output port for a spike output simulating the R wave of an EKG.

The source of positive potential for the portable calibrator 70 may be at a voltage between about plus 5 and about plus 18 volts. Similarly, the negative voltage source for portable calibrator 70 may be at a potential of about minus 18 volts to about minus 5 volts.

Portable calibrator 70 provides test signal for the heart rate spectral analysis hardware which, although not of a truly calibrated nature, does allow one to evaluate whether or not the software and hardware is functional. Each of the output signals provided is a triangle wave which represents the respiration and a frequency modulated pulse train representing the heart rate. The modulation of the heart rate is provided at two frequencies which simulate a respiratory modulation and also a low frequency modulation.

The basic circuit of calibrator 70 for providing each pulse train consists of an oscillator having one operational amplifier as typified by the respiratory frequency modulator. A charging capacitor 505 and a variable resistor 507, provide an RC circuit which is charged by the output of the amplifier 500. It is also discharged by the amplifier 500 when the output of the amplifier 500 is low. Progressive cycles of the oscillator consist of charging and discharging the capacitor at the rate prescribed by the RC circuit. The reference level which determines whether or not one is discharging or charging is provided at the non-inverting input of the amplifier 500.

Suppose, for example, that capacitor 505 begins as being completely discharged, then the voltage at the inverting input for the operational amplifier 500 is low. The output of the operational amplifier 500 is therefore high and this means that the input at the

- 32 -

non-inverting input is  $2/3$  the voltage between the negative voltage source  $V$  and the positive voltage source  $V+$ . Thus the capacitor 505 begins to charge. When the capacitor voltage exceeds the threshold at the non-inverting input of the operational amplifier 500, the output of operational amplifier 500 changes sign and capacitor 505 begins to discharge. However, when the output of the amplifier 500 changes to the negative side, then the threshold voltage at the non-inverting input is changed and now becomes only  $1/3$  the way from the negative voltage source to the positive voltage source. This means that the voltage on the charging capacitor 505 varies between  $1/3$  and  $2/3$  the difference between the negative and the positive voltage source. This determines the range of output on capacitor 505. The voltage at capacitor 505 is buffered by a non-inverting buffer 511 and this provides the respiratory signal at node 591.

An identical oscillator is used to provide low frequency modulation. The difference in the two frequencies is obtained by adjusting the respective variable resistors, 505 and 517, which set the RC time constants. The outputs of these two modulators are fed by resistors 512 and 513 into the charging capacitor 514 for the heart rate.

The heart rate oscillator is similar in design and consists of variable resistor 526 and capacitor 532 which charges and discharges in cycles with the range of voltages on the capacitor ranging between  $1/3$  the distance from the negative voltage source to the positive voltage source to  $2/3$  the voltage between the negative voltage source and the positive voltage source. Resistors 512 and 513, which connect the outputs of the low frequency and respiratory frequency modulators to the heart rate modulator, allow a small amount of current to flow into charging capacitor 514 of

- 33 -

the heart rate modulator. This alters the charging rate of capacitor 514 and thereby affects the rate at which the heart rate oscillator oscillates. For example, on a positive cycle of the respiratory frequency modulator, the heart rate capacitor is charging more rapidly towards the plus side because more current is being supplied on the plus side of the cycle. Finally, the output of the heart rate modulator is sent through an RC filter comprising capacitor 532 and resistor 533 which converts the square wave output of the heart rate modulator into a spike output which may be sent to an R wave detector. Notice that the spike output includes both positive and negative spikes so that an R detector which depends on a high frequency filtering function may be discharging at twice the heart rate, inasmuch as it may trigger on both positive and negative spikes.

As illustrated by a block diagram in Figs. 9A and 9B, a block diagram may be constructed for the main program (designated SYNCTS19) and for sub-routine modules (SYNC7s, GWINDOW3, and FGRAPH8). This block diagram may be used in order to better interpret a complete program for heart rate fluctuation spectral analysis useful on an IBM personal computer, as illustrated in Appendix B. Although programs are provided for a Hewlett-Packard and an IBM computer herein, the software and other aspects of the present invention may be readily modified for use with other mini- and micro-computers.

In the program of Appendix B, is a routine for removing artifacts from a detected heart rate provided for by an electrocardiograph machine. This program computes histograms from the heart rate data in order to generate a tachometer waveform. The most common rate on the histogram is selected as the correct rate and other rates are interpreted in light of it. Specifically, in order to correct for a spurious extra trigger, where a

- 34 -

first and a second beat are close together while a third beat is spaced at an abnormally long interval, the second beat is discarded if the first beat to second beat interval is less than a predetermined value. The  
5 resulting interval between the first and the third beats is divided by an integer in order to provide a more normal intrabeat interval. Where a trigger has been missed, so that a first and a second beat are separated by an interval which is approximately a multiple of a  
10 normal intrabeat interval, the intrabeat interval is divided by that multiple, most commonly two, in order to provide a more correct interval length. If the slewing rate of the heartbeat samples is outside of an acceptable range of slewing rates determined as a  
15 function of a mean variance, and the problem cannot be identified as a missed trigger or as a spurious extra trigger, or if the three previous intervals have been corrected, a determined mean interval, against which all other intervals are judged, is substituted for the  
20 inappropriate interval.

The slew rate is calculated on a moving average of the heart rate waveform and corrects for triggers that fall within the parameters of 0.05 Hz (3 beats/min.) per beat and five times the maximum slew.  
25 This artifact-correcting routine never slews more than 10 percent of the heart rate waveform.

Within the software of Appendix A is a graphic routine for trending heart rate fluctuation spectral data. The parameters of LFP, RFP, LFP/RFP ratio and  
30 heart rate are plotted on a graph over time to show trends in the four parameters. These trends may then be studied in order to examine the effects of various clinical interventions. Values for the parameters heart rate, LFP/RFP ratio, LFP and RFP are stored and may be  
35 called up at any point in time through a graphing routine in order to provide a graphic depiction of the



- 35 -

course of a patient's condition. This sort of graphic depiction is illustrated for a stable patient in Fig. 10 and for an unstable patient in Fig. 11.

Also present in the program of Appendix B, a  
5 routine is provided for the segmentation of data and subsequent reanalysis. In this routine, data from the analog to digital converter 260 is collected continuously into a buffer and is dumped to a disk in blocks of 1,024 numbers (2,048 bytes equals 1,024 words  
10 and each block is referred to as a record or EPOCH). The time of heartbeat occurrence as measured by the signal provided by outputs OUT1 and OUT2 of timer 220 are collected continuously into two buffers (hb buffer 1 and hb buffer 2). These times are dumped to the disk in  
15 blocks of 1,024 pairs of numbers (1,024 from each buffer which equals 2,048 bytes or 1,024 words each). Because the heart rate is less than the sample rate of A/D converter 260 as required by signal processing, there are fewer heartbeat disk dumps.

20 In order to properly analyze data, the A/D and heartbeat data must correspond to the same time interval for the purpose of doing correlations. The correspondence may be determined from (1) the record number in a A/D file and (2) the absolute of the times  
25 stored in the heartbeat file (time differences used for intrabeat intervals). The instantaneous heart rate signal is generated backwards in time from the heartbeat corresponding to the last A/D sample in the record of interest. This means that if the heart rate signal is  
30 analyzed on a frequency scale not corresponding to the respiration data (e.g. respiration sample at 16 Hz but a heart rate analysis at 0 to 4 Hz) then the heart rate waveform extends backwards in time beyond the beginning of the present A/D record. This means that the heart  
35 rate waveform overlaps the heart rate waveform corresponding to the previous A/D records.

- 36 -

Overlapping permits lower frequency analysis than would be possible if only data corresponding to the present record were used (as in the prototype apparatus). Also, overlapping leads to the smoothing of parameters and to the subsequent reduction of fluctuating artifacts. In addition, it becomes less critical at what point analysis begins.

A calibration program providing a software driven calibrator, which may provide more realistic spectral data than the portable calibrator of Fig. 8, is contained within the program of Appendix A for a Hewlett-Packard micro-computer. Appendix C is a program which, although not tested, is believed to provide the same sort of software-driven calibration for an IBM personal computer through the data acquisition system of Figs. 4 through 7.

In general, outputs OUT0 and OUT1 of timer 220 in Fig. 5 generate a time base used via interrupt request line IRQ4 to clock data from a buffer to D/A converter 300. This buffer contains a respiratory waveform which may be a sign wave or any selected waveform as obtained by changing the contents of the buffer. Output OUT2 of timer 220 generates a heartbeat pulse as its output. In order to work properly, this pulse must be returned to the ECG trigger through node 400 or directly to interrupt request line IRQ3. If the latter course is chosen, however, node 491 must be disconnected from the output of Schmitt trigger 416. By returning the pulse to the ECG trigger, the computer is informed that the timer is through counting the present RR interval and needs a new interval to be loaded into a timer register of timer 220.

Through the use of the apparatus according to the present invention, a display of instantaneous heart rate as provided by an electrocardiograph machine, and as illustrated in Fig. 12, may be converted into an

- 37 -

instantaneous heart rate fluctuation spectrum as  
illustrated in Fig. 14. A typical spectrum for a stable  
patient is illustrated in Fig. 14 while a typical  
spectrum for an unstable patient is illustrated in Fig.

5 15.

Example I and Example II relate respectively  
to diagnosis and to treatment employing the present  
invention.

Parts suitable for use in construction of the  
10 apparatus as illustrated in Figs. 4 through 9 may  
include those as listed in Tables I, II, III and IV.

15

20

25

30

35

- 38 -

TABLE I

	<u>Element No.</u>	<u>Part No.</u>	<u>Manufacturer, Location</u>
5	100, 110, 120	74HC244	National Semiconductor Santa Clara, California
	130	74HC138	National Semiconductor Santa Clara, California
10	140, 151, 158	74HC30	National Semiconductor Santa Clara, California
15	152, 153, 154 155, 156, 157 159, 160, 161 172, 265, 266	74HC04	National Semiconductor Santa Clara, California
20	171, 173, 185 264	74HC00	National Semiconductor Santa Clara, California
	200	8255A-5	Intel Corporation Santa Clara, California
25	220	8253-5	Intel Corporation Santa Clara, California
	224		
30	240	74HC393	National Semiconductor Santa Clara, California
35	260	ADC0808	National Semiconductor Santa Clara, California

- 39 -

	268, 303, 312	LM324AN	National Semiconductor
	401, 406, 412		Santa Clara, California
	500, 511, 515		
	518, 525		
5			
	280	8286	Intel Corporation
			Santa Clara, California
	300	DAC0830	National Semiconductor
10			Santa Clara, California
	416	74HC14	National Semiconductor
			Santa Clara, California
15			
20			
25			
30			
35			

- 40 -

TABLE II

Diodes

5	<u>Element</u>	<u>Part No.</u>
	407, 415	IN4148
10		
15		
20		
25		
30		
35		

- 41 -

TABLE III

Resistors

5	<u>Element No.</u>	<u>Value (in ohms)</u>
	403i, 410	2.2k
	269	5k
10	270, 271, 409	10k
	403h	15k
	403g	27k
	403f	56k
15	313	82k
	309, 310, 501,	100k (variable)
	502, 503, 520,	
	521, 522, 527,	
	528, 529, 533,	
20	403e	
	403d, 410d	220k
	403c, 410e	560k
25	508, 516, 526,	1M(variable)
	403b, 410b	
	512, 513,	2.2M
	403a, 410a	
30		
35		

- 42 -

TABLE IV  
Capacitors

	<u>Element No.</u>	<u>Value (in microfarads)</u>
5	405	2.2
	404, 505, 517	10
10	532	0.1
	514	1

15 EXAMPLE 1

Heart rate spectral analysis was applied to the study of congestive heart failure in infants and children. Congestive heart failure is characterized by a marked alteration in cardiovascular regulation. However, many cardiovascular functions which are normally monitored in cardiac intensive care units (such as: mean heart rate; arterial blood pressure; arterial blood gases; left arterial pressure and right arterial pressure; right atrial, left atrial and pulmonary artery oxygen saturations; the peripheral pulses; peripheral perfusion; and cardiac output) may not clearly indicate a critically unstable cardiovascular condition. The usually-monitored cardiovascular function parameters may be within a normal range immediately before a major cardiovascular crisis, such as hypotension or cardiac arrest, inasmuch as the cardiovascular regulatory system maintains these parameters within a normal range up to the point of system failure.

Twenty-nine infants and children were studied in a cardiac intensive unit. Of the twenty-nine patients, twenty-six have undergone a cardiac surgical



- 43 -

procedure. The patients were studied for a minimum of three hours and a maximum of twenty-seven hours, with a mean study time of eight hours. EKG for cases were recorded and analyzed continuously in real time during the study time.

Data for a particular patient was analyzed only if the patient was in sinus rhythm. The patient's clinical course during the period of study was reviewed and, in particular, major events such as cardiac arrest, hemorrhage and profound hypotension were correlated with spectral analysis data. Administration of medication and the mode of ventilation were noted.

Real time heart rate spectral analysis was performed on a dedicated personal computer using a 6809E Motorola Microprocessor-Based System. A data acquisition system interfaced the computer with a patient monitor, available from Hewlett-Packard, Palo Alto, California, as Model No. 78341.

The heart rate power spectrum was calculated in continuous 256 second data epochs. A QRS synchronization pulse from the patient monitor was used to determine an RR interval sequence. An instantaneous heart rate signal was computed from RR interval sequence and the magnitude of the signal was set to the reciprocal of the current interbeat interval. The instantaneous heart rate signal was sampled at 4 Hz and the mean heart rate was subtracted from the resulting one thousand twenty-four point time series. A power spectrum was computed by squaring the absolute value of a Fast Fourier Transform of the one thousand twenty-four point time series. Values for low frequency power (LFP) were computed by integrating the spectrum of between 0.04 and 0.1 Hz. Respiratory frequency power (RFP) was computed by integrating the heart rate power spectrum over a 0.2 Hz-wide band centered at the mean respiratory frequency.

- 44 -

Hard copies of the heart rate time series and power spectrum were printed for each 256 second epochs. Trend graphics for the LFP, the RFP, LFP/RFP ratio, mean heart rate and respiratory rate (hereinafter referred to as the study parameters) were constructed by manually entering data in data files and analyzing the entered data by means of a computer.

Mean values for the study parameters were calculated for each period of study. The Mann-Whitney Rank Sum Test was used to determine statistically significant changes in the study parameters in individual patients and to determine differences among groups of patients. When patients were segregated into more than two groups, the Kruskal-Wallis Test, multiple comparison test, and Tukey's HSD were employed to determine statistical significance. P values of less than 0.05 were considered significant.

It was found that during each three to twenty-four hour period of study the study parameters for a given patient, the LFP, the RFP and the LFP/RFP ratio (hereinafter referred to as the spectral parameters) remain fairly stable.

Based upon the results of this study, the patients were retrospectively divided into three groups. Group I included seventeen stable patients whose median age was one month. The patients in Group I were without major post-operative complications and did not need prolonged inotropic support. The eight patients in Group II suffered cardiac arrest and died. The median age for the members of Group II was one month. In Group III, there was a total of four patients each of whom was critically ill at the time of the study but later recovered. Median age of the members of Group III was one month. Of the four members of Group III, one required re-operation, one had intermittent hypotensive episodes, and two had cardiac arrests from

- 45 -

which they were successfully resuscitated.

In order to separate all twenty-nine patients into a group of stable patients (Group A) and a group of critical patients (Group B), data from each patient in Group III was divided into the data collected during the stable period (which applied to three patients) and the data collected during the preceding critical period (which applied to four patients). When handled in this way, Group A included data for twenty patients and Group B included data for twelve patients. Typical heartrate fluctuation power spectra for Group A and B are respectively illustrated in Figs. 16 and 17.

In addition, studies were performed on three patients who had isolated coarctation of the aorta at three points in time: upon admission for congestive heart failure; during treatment; during post-operative period; and prior to discharge from an intensive care unit. An attempt was made to identify changes in cardiovascular regulatory function of each of these stages.

Patient profiles for Groups I, II and III are respectively provided in Tables V, VI and VII. These profiles include age, diagnosis and operation.

TABLE V

PATIENT PROFILE;		STABLE POST-OP		N=17
AGE	NO.	DIAGNOSIS	(NO.)	OPERATION
30	9	TGA, IVS	(3)	ARTERIAL SWITCH
		TGA, VSD, PS	(1)	L-BTS
		HLHS	(1)	STAGE 1 REPAIR
		SV	(1)	L-BTS
		SEV. COAO	(3)	SUBCL. FLAP ANGIO.

- 46 -

5	1-12 MO.	5	TGA, IVS	(1)	ARTERIAL SWITCH
			TGA, VSD, PS	(1)	BTS
			MULT. VSD'S	(1)	VSD PATCH REPAIR
			SUPRA-V. PS	(1)	PA PATCH PLASTY
			DCRV, VSD, COAO	(1)	VSD REPAIR, ANOM. B RESECTION
10	1-10 YRS.	2	PS	(1)	PULM. VALVOTOMY
			TOF	(1)	TOF REPAIR
	>10 YRS.	1	AR, MR		AVR, MVR

TABLE VI

15	PATIENT PROFILE; CRITICAL, DIED N=8				
	AGE	NO.	DIAGNOSIS	(NO.)	OPERATION
20	<30 DAYS	4	HLHS	(3)	NORWOOD PROCEDURE
			SV W/IAA	(1)	GORE-TEX GRAFT
25	1-12 MO.	3	HLHS	(1)	Fontan operation
			DORV, TAPVC,	(1)	TAPVC REPAIR, SYS.
			CCAVC		PULM. SHUNT
			HLHS	(1)	NON-OPERATIVE
30	6 1/2 YRS.	1	T OF S/P REPAIR W/ CHRONIC SEV. CARDIOMYOPATHY, S/P ARREST		NON-OPERATIVE

35

- 47 -

TABLE VII

PATIENT PROFILE: CRITICAL, RECOVERED N=4

5	AGE	NO.	DIAGNOSIS	(NO.)	OPERATION
	<30 DAYS	3	HLHS, COAO	(1)	NORWOOD PROCEDURE
			HLHS	(2)	NORWOOD PROCEDURE
10	14 YRS.	1	ACUTE MYOCARDITIS, S/P ARREST		NON-OPERATIVE

In Tables V, VI and VII: TGA is Transposition of the Great Arteries; IVS is Ventricular Septal Defect; PS is Pulmonic Stenosis; HLHS is Hypoplastic Left Heart Syndrome; SV is Single Ventricle; SEV. is severe; COAO is Coarctation of the Aorta; MULT is multiple; VSD is Ventricular Septal Defect; Supra-V. is Supravalvular; DCRV is Double Chamber Right Ventricle; TOF is Tetralogy of Fallot; AR is Aortic Regurgitation; MR is Mitral Regurgitation; W/IAA is with Interrupted Aortic Arch; DORV is Double Outlet Right Ventricle; TAPVC is Total Anomalous Pulmonary Venous Connections; CCAVC is Complete Common Atrial Ventricular Canal; S/P is Status Post; L is Left; BTS is Blalock Taussig Shunt; PA is Pulmonary Artery; ANOM. is Anomalous; B is muscle Bundle; PULM is Pulmonary; and SYS is Systemic.

Statistically significant differences were observed in the heart rates spectral parameters between the groups of patients as well as among the individual patients. However, the mean heart rate alone did not distinguish stable from critically ill patients. Both the LFP and the LFP/RFP ratio discriminated between the Group A (stable) patients and the Group B (critical) patients. The LFP/RFP ratio grew out of a statistically significant ( $p$  less than symbol 0.00001) discrimination

- 48 -

between stable and critical patients. Table VIII presents means of study parameters.

TABLE VIII

5

## MEANS OF STUDY PARAMETERS

## GROUP A, STABLE

10	PARAMETER (BEATS/MIN.)			STD. ERROR		99% CONFIDENCE	
		MEAN	STD. DEV.	OF MEAN		LOWER	UPPER
	LFP	1.77	3.35	0.75	-.37	3.91	
	RFP	0.28	0.70	0.16	-.17	0.72	
15	LFP/RFP RATIO	8.77	4.86	1.09	8.76	8.79	
	HEART RATE	139	19.60	4.38	139	139	

## 20 GROUP B, CRITICAL

	PARAMETER			STD. ERROR		99% CONFIDENCE	
		MEAN	STD.DEV.	OF MEAN		LOWER	UPPER
25	LFP	.05	.03	.01	.02	.07	
	RFP	.10	.09	.03	.01	.18	
	LFP/RFP RATIO	.83	.51	.15	.83	.83	
	HEART RATE	142	24.32	7.02	142	142	

30

The discriminate value for the LFP/RFP ratio was two. In Group A, the range of LFP/RFP ratios was 3 to 22 (arithmetic mean 8.77). The range of RFPs was 0.01 to 3.13 (arithmetic mean 0.28) and the range of LFPs was 0.09 to 13.88 (arithmetic mean 1.77). In Group

35

- 49 -

B, the range of LFP/RFP ratios was 0.17 to 1.9 (arithmetic mean 0.83), the ratio of RFPs was 0.02 to 0.32 (arithmetic mean 0.1), and the range of LFPs was 0.01 to 0.1 (arithmetic mean 0.5)

5           Although the mean value of the LFP/RFP ratio was greater than two for Group I, the ratio for the stable patients fell below two for brief periods. That which distinguishes the stable from the critical patients is the sustained value for greater than or  
10   about one hour of the LFP/RFP ratio for the critical group.

          The results are graphically depicted in Figs. 19, 17 and 18. In Figs. 16 and 17, each heavy dot A represents a geometric mean, each light line B indicates  
15   the standard error of the geometric mean and each heavy line C represents the standard deviation of the geometric mean. In Fig. 18, each heavy dot A represents an arithmetic mean, each set of slashes B1 and B2 represents the standard error of the arithmetic mean and  
20   each set of slashes C1 and C2 represents the standard deviation of the arithmetic mean.

          The significance of heart rate spectral analysis for diagnosis of cardiovascular stress and the prediction of fatality is highlighted by the fact that  
25   patients with a low LFP/RFP ratio underwent a cardiac arrest even in the presence of otherwise normal vital signs. No patient with a LFP/RFP ratio greater than two experienced a cardiac arrest.

          Infusion of pressors, alone or in combination  
30   with vasodilators, did not induce a low LFP/RFP ratio.

          Four patients in Group III had LFP/RFP ratios less than two during their critical periods. For the three of these four patients who were restudied during their recovery periods, all three had LFP/RFP ratios  
35   greater than two.

          The mean LFP for Group B [0.05 (Beats per

- 50 -

minute)<sup>2</sup>] was less than the mean LFP for Group A [1.77 beats per minute)<sup>2</sup>],  $p < 0.0001$ . There was no significant difference between the mean RFP between the groups.

5           The initial LFP/RFP ratios for the patients with isolated coarctation of the aorta ranged up to 10,000. The LFP/RFP ratios observed for this group immediately after an operation to correct the condition were within the range for Group A patients. Two  
10 patients had LFP/RFP ratios greater than 100 before discharge from the intensive care unit. These ratios were correlated with mild to moderate congestive heart failure. One of these patients died suddenly at approximately 2-1/2 months after the operation. The  
15 other two patients remained alive and well.

          Although the LFP/RFP ratio provided the sharpest discrimination between stable and critical patients in these studies, the LFP alone discriminated between Groups A and B,  $p < 0.0001$ . Neither respiratory  
20 frequency peak power nor mean heart rate distinguished between Groups A and B. On the other hand, LFP/RFP ratios and LFP levels low levels sustained for greater than or about one hour correlate with the course of the conditions of patients who experienced cardiac arrest or  
25 severe hypotensive episodes but later recovered.

          Although stable patients experienced transient depression of levels of LFP and of the LFP/RFP ratio, depression of these factors for about an hour or more never failed to predict a critical status.

30           No significant difference was observed between freely ventilating patients and mechanically ventilated patients. Eighteen out the twenty patients in Group A were mechanically ventilated and all twelve of the Group B patients were mechanically ventilated.

35           All patients in Group B received inotropic support while more than half of the patients in the



- 51 -

Group A received at least some inotropic support. The cardiac diagnoses of all of the patients in Group B and for some of the patients in Group A were known to be associated with high mortality. All of the patients in Group B underwent deep hypothermic circulatory arrest during their operations. Of the twenty patients in Group A, nine had extra cardiac surgery (i.e. not involving cardiopulmonary bypass or deep hypothermic circulatory arrest). Three of the patients in Group II did not undergo operations. Therefore, it is not believed that differences in treatment or disease specific pathology alone explained the low values LFP and the low LFP/RFP ratios in Group B patients but that the low values actually reflect a vulnerable circulatory state.

It has also been observed that the value of LFP and of the LFP/RFP ratio increase in moderate to severe heart failure but decreased to subnormal values in end stage myocardial failure. Thus, these two spectral parameters may indicate cardiovascular regulatory effectiveness (cardiovascular regulatory reserve) during the stress of heart failure.

This analysis is consistent with previous physiological studies which indicated that low frequency heart fluctuations may be mediated by both the beta-sympathetic and parasympathetic mechanisms while respiratory fluctuations are exclusively mediated by parasympathetic mechanisms. It is also consistent with this analysis that LFP has been observed to increase during conditions which elicit enhanced sympathetic activity, such as acute hypoxia, postural changes, hemorrhage and aortic constriction. In this light, the LFP/RFP ratio may represent a measure of the balance between beta adrenergic and parasympathetic modulation of cardiac function.

Thus, the increase in LFP and in the LFP/RFP

- 52 -

ratio for patients with isolated coarctation of the aorta and moderate heart failure may result from an increased activity from the sympathetic mechanism and a decreased activity of the parasympathetic mechanism. On the other hand, the decreased level of LFP and of the LFP/RFP ratio found in critical patients may be due to non-responsiveness of the sympathetic mechanism. Sympathetic non-responsiveness may be due to myocardial catecholamine depletion alone or in combination with the observed down regulation of beta receptors from cardiac tissue in the end stage of heart failure.

#### EXAMPLE 2

In patients undergoing operations, shifts in body fluid disposition during surgery may lead to changes in intervascular volume (i.e. a shift of fluid out of a circulatory tree of blood vessels). Accordingly, the availability of the method of diagnosing cardiovascular stress as described in Example 1 may be used to choose among various protocols for treatment or to justify a radical change in medical or surgical treatment.

For example, by monitoring a patient with the real time heart rate frequency spectral monitor according to the present invention during administration of anesthesia, an anesthesiologist may non-invasively monitor intravascular volume status. Upon observing an increase in the LFP or in the LFP/RFP ratio, the anesthesiologist may increase the amount of fluids administered by way of intravenous injection or may take steps to reverse effects of a particular anesthetic.

It is a particular advantage of the apparatus according to the present invention that heart rate fluctuation spectral analysis may be done in real time. This capability permits correlation of treatment administered with changes in LFP or LFP/RFP ratios.

- 53 -

Although the present invention has been described in terms of preferred embodiments, it is understood that modifications, variations and improvements will occur to those skilled in the art.

- 5 For example, it will occur to those skilled in the art to employ the present invention for monitoring cardiovascular instability in the following settings in which significant circulatory stress are commonly observed: Labor and Delivery Room; Operating Room;
- 10 Cardiac Catheterization Laboratory; Neonatal, Pediatric, Adult Medical, Adult Surgical, Cardiothoracic and Neurosurgical Intensive Care Units; Coronary Care Units; Burn Units; and Emergency Rooms.

- The present invention may also be used for
- 15 monitoring cardiovascular instability in the following patients in which adjustments in cardiovascular regulation may provide a central key to understanding the efficiency and efficacy of treatment. Ambulatory patients with known heart disease in which sudden
- 20 cardiac death is a common association, one example of which would be a patient with a congestive cardiomyopathy who is being treated with vasodialator drugs and for whom the LFP/RFP ratio has changed from a normal baseline level to decreased levels may then
- 25 subsequently be either admitted to the hospital for adjustment of medications and/or observed and monitored in the physician's office while his vasodialator drug dose is increased. A patient with renal disease (e.g. one who requires dialysis) may exhibit a marked increase
- 30 in LFP and LFP/RFP ratio secondary to the onset of incipient moderate congestive heart failure would thus be treated by dialysis to relieve a congested circulatory state; a patient with moderate to severe pulmonary disease resulting in hypoxemia and/or
- 35 hypercarbia who requires bronchodialator and/or supplementary oxygen and/or mechanical ventilation (e.g.

- 54 -

a patient who exhibits a marked decrease in LFP/RFP ratio secondary to myocardial failure due to a profound imbalance between myocardial ventricular output and oxygen demand), may be treated by adjustments in  
5 bronchodialator drugs, diuretics, and/or ventilator adjustments.

A premature infant of very low birth weight known to be at risk for intraventricular hemorrhage may, for example, develop a slow intracranial bleed  
10 associated with an abrupt increase in LFP, which may alert physicians prior to a brisk bleed thus allowing institution of appropriate changes in medical management to limit substantially known risk factors that may predispose to such an event, or may permit recognition  
15 of the presence of unsuspected circumstances that contribute to the bleed. In neurologic disease, such as one in which a patient has sustained a major intracerebral event (e.g. neurosurgical evacuation of a space occupying lesion such as a tumor or blood), a  
20 patient may, for example, exhibit a markedly attenuated LFP/RFP ratio, secondary to massively increased parasympathetic activity which would markedly increase RFP, at the expense of LFP, but which may or may not be associated with signs of increased intracranial  
25 pressure, and which may be treated by, for example, hyperventillation, rapid diuresis, or burr hole placement.

A patient with severe systemic infection may exhibit shock secondary to the infection process may,  
30 for example, exhibit an elevated LFP/RFP ratio which may then be subsequently used by the physician in managing the shock state by means of pressor agents and infusion of significant volumes of fluid, thus providing the physician an indication of how effectively he is  
35 treating the shocked state above and beyond the traditional measurements such as systemic blood pressure

- 55 -

and cardiac output. A patient with hematologic disease associated with anemia, such as Sickie Cell Anemia, exhibits an oscillation in capillary blood flow when severely anemic at the frequency associated with LFP and may exhibit large values for LFP, and for the LFP/RFP ratio may, for example, be treated by blood transfusion which may lead to an expected decrease in LFP, LFP/RFP ratio, and thus enable the physician to monitor by means of heart rate spectral analysis appropriate timing for transfusion therapy. A fetus prior to delivery, may for example, exhibit a marked attenuation in LFP associated with severe fetal distress, and may thus alert the physician to perform an emergency Caesarean section.

One skilled in the art understands that the calibrators according to the present invention may be adjusted to simulate disease states as well as normal conditions. It is also understood that the present invention is not limited to use with patients whose primary disease is of the heart but that modifications may be made for use with such patients.

Lastly, it is clear to one skilled in the art that durations and ranges for levels of LFP and LFP/RFP ratios are conservatively stated herein and that variations from these ranges and durations are contemplated within the scope of the equivalents of the present invention.

Therefore, it is intended that the methods and apparatus according to the present invention to be given the broadest scope allowable for the invention as claimed.

- 56 -

## APPENDIX A

```

5  10 Summary3:1
20      !This program takes data already collected and
      allows the data
30      !to be outputted to a printer
40      !2 MAY 1985
10  50      !
60      COM /Trends/ Mean_hr_t(60),Lfa_t(60),Rfa_
      t(60),Ratio_t(60),T_ptr,Time_now1,Mean_resp_
      t(60),Trend_dp
70      COM /Multi_param/ Start_chan,Stop_chan,Pacing_
15  bits,Pacing_rate,Num_pts,Nu
      m_xfer,Num_xfer_left,Name_len,Scr_file$(28),Scr_
      file2$(28)
80      COM /Pressure/
      Top1,Top2,Top3,Top4,Bot1,Bot2,Bot3,Bot4
20  90      COM /Editor/ Edit_msg$(80)
100     COM /Subject/ Sub_name$(25),Hos_num$(15),Id_
      age$(10),Id_wt$(10),Id_ht$(10
      ),Diag$(30),Opera$(45),Halt_pg,In_file$(6)
110     COM /Io_chart/ Io_time$(8)[10],Iv_intake(8),Fluid_
25  in(8),In_tot(8),Urine(8
      ),Chest(8),Out_tot(8),Net(8),Io_ptr
120     COM /Lab_chart/ Lab_
      time$(8)[10],Na(8),Kl(8),Cl(8),Hco3(8),Ca(8),Hct(8),G
      luc(8),Dig(8),Pt(8),Ptt(8),Creat(8),Bun(8),Lab_ptr
30  130     COM /Vent_chart/ Vent_
      time$(8)[15],Rate(8),Fio2(8),Pp(8),Peep(8),Tv(8),
      Ie_ratio$(8)[5],Airp(8),Ph(8),Po2(8),
      Pco2(8),Bgo3(8),Be(8),Vent_ptr
140     COM /Pres_chart/ Pres_time$(20)[15],Ao_s(20),Ao_
35  d(20),Ao_m(20),Pa_s(20),P
      a_d(20),Pa_m(20),La_m(20),Ra_m(20),Pres_

```

- 57 -

```
ptr,Pres_in
150  COM /Heart_index/ Heart_
      time$(15)[15],Ci(15),Pvri(15),Svri(15),Heart_ptr
160  COM /Drugs/ Drug_time$(40)[20],Drug_
5    name$(40)[40],Drug_dos$(40)[20],Drug_
      ptr
170  DIM Msg_buffer$(6400) BUFFER
180  DIM Pres_p(20),Io_p(8),Lab_p(8),Vent_p(8),Heart_
      p(5),Drug_p(40)
10  190  INPUT "enter date on which data was collected
      (ddmmyy) e.g. 22AP85",In_file$
200  Disk1$=":HP8290X,700,1"
210  INPUT "is the trend file named 'trnd'(1) or 'temp_
      trend'(2)?",Ans
15  220  IF Ans=2 THEN
230      ASSIGN @Trend_file TO "temp_
      trend"&Disk1$;FORMAT OFF
240      ASSIGN @Messages TO "messglog"&Disk1$;FORMAT
      OFF
20  250  ASSIGN @Hemo_data TO "hemo_
      data"&Disk1$;FORMAT OFF
260  ASSIGN @Io_data TO "io_data"&Disk1$;FORMAT
      OFF
270  ASSIGN @Lab_data TO "lab_data"&Disk1$;FORMAT
25  OFF
280  ASSIGN @Vent_data TO "vent_
      data"&Disk1$;FORMAT OFF
290  ASSIGN @Co_data TO "co_data"&Disk1$;FORMAT
      OFF
30  300  ASSIGN @Drug_data TO "drug_
      data"&Disk1$;FORMAT OFF
310  ASSIGN @Sub_data TO "sub_data"&Disk1$;FORMAT
      OFF
320  ON END @Trend_file GOTO Start
35  330  FOR I=0 TO 55
340      ENTER @Trend_file;Trans_t(I),Mean_hr_
```

- 58 -

```

        t(I),Lfa_t(I),Rfa_t(I),Ratio
        _t(I),Mean_resp_t(I)
350     NEXT I
360     T_ptr=I
5 370     Num_xfer=T_ptr
380     ELSE
390     ASSIGN @Trend_file TO "trnd"&In_
        file$&Disk1$;FORMAT OFF
400     ASSIGN @Messages TO "msgs"&In_
10  file$&Disk1$;FORMAT OFF
410     ASSIGN @Hemo_data TO "hemo"&In_
        file$&Disk1$;FORMAT OFF
420     ASSIGN @Io_data TO "io__"&In_
        file$&Disk1$;FORMAT OFF
15 430     ASSIGN @Lab_data TO "lab_"&In_
        file$&Disk1$;FORMAT OFF
440     ASSIGN @Vent_data TO "vent"&In_
        file$&Disk1$;FORMAT OFF
450     ASSIGN @Co_data TO "co__"&In_
20  file$&Disk1$;FORMAT OFF
460     ASSIGN @Drug_data TO "drug"&In_
        file$&Disk1$;FORMAT OFF
470     ASSIGN @Sub_data TO "sub_"&In_
        file$&Disk1$;FORMAT OFF
25 480     ENTER @Trend_file;Mean_hr_t(*),Lfa_t(*),Rfa_
        t(*),Ratio_t(*),Mean_resp
        _t(*),Trans_time(*),T_ptr
490     Num_xfer=T_ptr
500     END IF
30 510     ASSIGN @Trend_file TO *
520     ON END @Hemo_data GOTO Hemol
530     FOR I=0 TO 20
540     ENTER @Hemo_data;Pres_time$(I),Ao_s(I),Ao_
        d(I),Ao_m(I),Pa_s(I),Pa_d(I
35  ),Pa_m(I),La_m(I),Ra_m(I),Pres_p(I)
550     NEXT I

```



- 59 -

```
560 Hemol:ASSIGN @Hemo_data TO *
570   Pres_ptr=I-1
580   ON END @Io_data GOTO Iol
590   FOR I=0 TO 8
5   600       ENTER @Io_data;Io_time$(I),Iv_intake(I),Fluid_
           in(I),In_tot(I),Urine(I
           ),Chest(I),Out_tot(I),Net(I),Io_p(I)
610   NEXT I
620 Iol:ASSIGN @Io_data TO *
10  630   Io_ptr=I-1
640   ON END @Lab_data GOTO Lab1
650   FOR I=0 TO 8
660       ENTER @Lab_data;Lab_
           time$(I),Na(I),Kl(I),Cl(I),Hco3(I),Ca(I),Hct(I),G
15  luc(I),Dig(I),Pt(I),Ptt(I),Creat(I),Bun(I),Lab_p(I)
670   NEXT I
680 Lab1:ASSIGN @Lab_data TO *
690   Lab_ptr=I-1
700   ON END @Vent_data GOTO Vent1
20  710   FOR I=0 TO 8
720       ENTER @Vent_data;Vent_
           time$(I),Rate(I),Fio2(I),Pp(I),Peep(I),Tv(I),
           Ie_ratio$(I),Airp(I),Ph(I),Po2(I),
           Pco2(I),Bgo3(I),Be(I),Vent_p(I)
25  730   NEXT I
740 Vent1:ASSIGN @Vent_data TO *
750   Vent_ptr=I-1
760   ON END @Co_data GOTO Col
770   FOR I=0 TO 5
30  780       ENTER @Co_data;Heart_
           time$(I),Ci(I),Pvri(I),Svri(I),Heart_p(I)
790   NEXT I
800 Col:ASSIGN @Co_data TO *
810   Heart_ptr=I-1
35  820   ON END @Drug_data GOTO Drug1
830   FOR I=0 TO 40
```

- 60 -

```
      840      ENTER @Drug_data;Drug_time$(I),Drug_
            name$(I),Drug_dos$(I),Drug_p(I)
      850      NEXT I
      860 Drug1:ASSIGN @Drug_data TO *
5   870      Drug_ptr=I-1
      880      !
      890      !
      900      !
      910      Pacing_rate=250
10  920      Time_now1=TIMEDATE MOD 86400
      930      Out_graph=1
            !....graphics
            dump
      940      Trend_dp=2
15  950      CALL Trend_graph
      960      CALL Graph_dump(Out_graph)
      970      Trend_dp=1
      980      CALL Trend_graph
      990      CALL Graph_dump(Out_graph)
20 1000      !
      1010 Chart_dump: !
      1020      ENTER @Sub_data;Sub_name$,Hos_num$,Id_age$,Id_
            wt$,Id_ht$,Diag$,Opera$
      1030      ASSIGN @Sub_data TO *
25 1040      Out_graph=2
      1050      FOR I=1 TO 5
      1060          CALL Chart(I)
      1070      CALL Graph_dump(Out_graph)      !....chart dump
      1080      NEXT I
30 1090      !
      1100      !
      1110 Msg_dump: !
      1120      IF Ans=1 THEN
      1130          ASSIGN @Msg_file TO "msgs"&In_
35 1140          file$&Disk1$;FORMAT OFF
            ELSE
```

- 61 -

```
1150      ASSIGN @Msg_file TO "messglog"&Disk1$;FORMAT
        OFF
1160  END IF
1170  PRINTER IS 701
5  1180  ASSIGN @Msg_buffer TO BUFFER Msg_buffer$
1190  STATUS @Msg_file,3;Num_rec
1200  STATUS @Msg_file,4;Rec_len
1210  STATUS @Msg_file,7;Eof_rec
1220  STATUS @Msg_file,8;Eof_byte
10 1230  Num_bytes=(Eof_rec-1)*Rec_len+Eof_byte-1
1240 Read_msg:TRANSFER @Msg_file TO @Msg_buffer;COUNT
        Num_bytes,WAIT
1250  ASSIGN @Msg_file TO *
1260  ASSIGN @Msg_buffer TO *
15 1270  Cur_ptr=1
1280      PRINT USING Image_wt1;Sub_name$,Hos_num$,In_
        file$
1290 Image_wt1:IMAGE  "Name: ",K,XXXX,"Hosp num:
        ",K,XXXXX,K
20 1300      PRINT USING Image_wt2;Id_age$,Id_wt$,Id_
        ht$,Diag$,Opera$
1310 Image_wt2:IMAGE  "Age: ",K,XXXX,"Wt(kg):
        ",K,XXXX,"Ht(cm): ",K,XXXX,"Diag:
        ",K,XXXX,"Op: ",K
25 1320 Next_msg:!
1330  Beg_msg=POS(Msg_buffer$[4],"Time")+3
1340  IF Beg_msg=3 THEN GOTO Stopper
1350  PRINT Msg_buffer$[1,Beg_msg-1]
1360  Msg_buffer$=Msg_buffer$[Beg_msg]
30 1370  GOTO Next_msg
1380 Stopper:!PRINTER IS 1
1390  STOP
1400  END
1410      !
35 1420  !
1430  !This subroutine prints the graphics
```

- 62 -

```
1440    !
1450    !
1460    SUB Trend_graph
1470    !
5    1480    COM /Trends/ Mean_hr_t(*),Lfa_t(*),Rfa_
        t(*),Ratio_t(*),T_ptr,Time_now
        1,Meas_resp_t(*),Trend_dp,Trans_time(*)
1490    COM /Multi_param/ Start_chan,Stop_chan,Pacing_
        bits,Pacing_rate,Num_pt
10    s,Num_xfer,Num_xfer_left,Name_len,Scr_
        file$(28),Scr_
        file2$(28)
1500    COM /Pressure/
        Top1,Top2,Top3,Top4,Bot1,Bot2,Bot3,Bot4
15    1510    COM /Pres_chart/ Pres_time$(*),Ao_s(*),Ao_
        d(*),Ao_m(*),Pa_s(*),Pa_d(
        ),Pa_m(*),La_m(*),Ra_m(*),Pres_ptr,Pres_in
1520    DIM First_line(60),Sec_line(60),Third_
        line(60),Fourth_line(60)
20    1530    IF Trend_dp=1 THEN
1540        MAT First_line= Ao_m
1550        MAT Sec_line= Pa_m
1560        MAT Third_line= La_m
1570        MAT Fourth_line= Ra_m
25    1580    G_right=INT((Num_xfer*256/60)/15)
1590    Trend_ptr=Pres_ptr
1600    Top1=150
1610    Bot1=0
1620    Top2=75
30    1630    Bot2=0
1640    Top3=50
1650    Bot3=0
1660    Top4=50
1670    Bot4=0
35    1680    ELSE
1690    MAT First_line= Mean_hr_t
```

- 63 -

```
1700      MAT Sec_line= Ratio_t
1710      MAT Third_line= Lfa_t
1720      MAT Fourth_line= Rfa_t
1730      G_right=Num_xfer
5  1740      Trend_ptr=T_ptr
1750      Top1=200
1760      Bot1=0
1770      Top2=2.5
1780      Bot2=-2.5
10 1790      Top3=10
1800      Bot3=0
1810      Top4=10
1820      Bot4=0
1830      END IF
15 1840      Block_time=Pacing_rate*1.024/3600.
1850      GINIT
1860      GCLEAR
1870      PRINT CHR$(12)
1880      GRAPHICS ON
20 1890      Beg_time=Time_now1/3600-Block_time
1900      End_time=Beg_time+Num_xfer*Block_time
1910      Ibeg_time=INT(Beg_time)
1920      IF Ibeg_time<Beg_time THEN Ibeg_time=Ibeg_
time+1
25 1930 !
1940 ! label the time axes
1950 !
1960      VIEWPORT 0,128,45,50
1970      WINDOW Beg_time,End_time,0,1
30 1980      IF INT(End_time)>Beg_time THEN
1990          LDIR 0
2000          FOR T_label=Ibeg_time TO INT(End_time)
2010              MOVE T_label,.5
2020              LORG 5
35 2030              CSIZE 4
2040              LABEL T_label
```

- 64 -

```
2050          NEXT T_label
2060      END IF
2070      VIEWPORT 0,128,40,45
2080      WINDOW 0,1,0,1
5  2090      MOVE .5,0
2100      LONG 4
2110      LABEL "Time (24 hr)"
2120 !
2130 ! draw the axes
10 2140 !
2150      VIEWPORT 0,128,50,100
2160      WINDOW Beg_time,End_time,0,1
2170      AXES 1/15.,.1,Beg_time,0
2180      WINDOW 1,0,1,0
15 2190      AXES 0,.25,0,0
2200 !
2210 ! mean heart rate trends
2220 !
2230      WINDOW -1,G_right,Bot1,Top1
20 2240      MOVE 0,First_line(0)
2250      FOR I=0 TO Trend_ptr-1
2260          DRAW I,First_line(I)
2270      NEXT I
2280 !
25 2290 ! ratio trends (with a line at ratio=2)
2300 !
2310      WINDOW -1,G_right,Bot2,Top2
2320      LINE TYPE 8,5
2330      IF Trend_dp=2 THEN
30 2340          MOVE 0,LGT(Sec_line(0))
2350      ELSE
2360          MOVE 0,Sec_line(0)
2370      END IF
2380      FOR I=0 TO Trend_ptr-1
35 2390          IF Trend_dp=2 THEN
2400              DRAW I,LGT(Sec_line(I))
```

- 65 -

```
2410          ELSE
2420              DRAW I,Sec_line(I)
2430          END IF
2440      NEXT I
5  2450      IF Trend_dp=2 THEN
2460          LINE TYPE 3,5!..sparsely dotted line at
              ratio=2
2470          MOVE 0,LGT(2.)
2480          DRAW Trend_ptr-1,LGT(2.)
10 2490      END IF
2500 !
2510 ! lfa trends
2520 !
2530      WINDOW -1,G_right,Bot3,Top3
15 2540      LINE TYPE 4,5
2550      MOVE 0,Third_line(0)
2560      FOR I=0 TO Trend_ptr-1
2570          DRAW I,Third_line(I)
2580      NEXT I
20 2590 !
2600 ! rfa trends
2610 !
2620      WINDOW -1,G_right,Bot4,Top4
2630      LINE TYPE 5,5
25 2640      MOVE 0,Fourth_line(0)
2650      FOR I=0 TO Trend_ptr-1
2660          DRAW I,Fourth_line(I)
2670      NEXT I
2680 !
30 2690 ! draw a key for line types
2700 !
2710      VIEWPORT 64,128,0,50
2720      WINDOW 0,1,0,13
2730      IF Trend_dp=2 THEN
35 2740          PRINT TABXY(1,17);"trend graph"
2750          PRINT TABXY(55,15);"mean hr(0-200)"
```

- 66 -

```
2760          PRINT TABXY(55,16);"ratio(.01-100)"
2770          PRINT TABXY(55,17);"lfa      (0-10)"
2780          PRINT TABXY(55,18);"rfa      (0-10)"
2790      ELSE
5   2800          PRINT TABXY(1,17);"mean pressure graphs"
      2810          PRINT TABXY(50,15);"ao. pressure(0-150)"
      2820          PRINT TABXY(50,16);"pa pressure(0-75)"
      2830          PRINT TABXY(50,17);"la pressure(0-50)"
      2840          PRINT TABXY(50,18);"ra pressure(0-50)"
10  2850      END IF
      2860      LINE TYPE 1,5
      2870      MOVE .8,11
      2880      DRAW 1.,11
      2890      LINE TYPE 8,5
15  2900      MOVE .8,10
      2910      DRAW 1.,10
      2920      LINE TYPE 4,5
      2930      MOVE .8,9
      2940      DRAW 1.,9
20  2950      LINE TYPE 5,5
      2960      MOVE .8,8
      2970      DRAW 1.,8
      2980      SUBEND
      2990      !
25  3000      !
      3010      !This subroutine prints the charts
      3020      !
      3030      !
      3040      SUB Chart(Char_num)
30  3050          COM /Subject/ Sub_name$,Hos_num$,Id_age$,Id_
              wt$,Id_ht$,Diag$,Opera$,H
              alt_pg,In_file$
      3060          COM /Io_chart/ Io_time$(*),Iv_intake(*),Fluid_
              in(*),In_tot(*),Urine(*
35          ),Chest(*),Out_tot(*),Net(*),Io_ptr
      3070          COM /Lab_chart/ Lab_
```



- 67 -

```

time$(*),Na(*),Kl(*),Cl(*),Hco3(*),Ca(*),Hct(*),G
    luc(*),Dig(*),Pt(*),Ptt(*),Creat(*),Bun(*),Lab_
        ptr
3080      COM /Vent_chart/ Vent_
5      time$(*),Rate(*),Fio2(*),Pp(*),Peep(*),Tv(*),
    Ie_ratio$(*),Airp(*),Ph(*),Po2(*),Pco2(*),
    Bgo3(*),Be(*),Vent_ptr
3090      COM /Pres_chart/ Pres_time$(*),Ao_s(*),Ao_
    d(*),Ao_m(*),Pa_s(*),Pa_d(*
10      ),Pa_m(*),La_m(*),Ra_m(*),Pres_ptr,Pres_in
3100      COM /Pressure/
    Top1,Top2,Top3,Top4,Bot1,Bot2,Bot3,Bot4
3110      COM /Heart_index/ Heart_
    time$(*),Ci(*),Pvri(*),Svri(*),Heart_ptr
15 3120      COM /Drugs/ Drug_time$(*),Drug_name$(*),Drug_
    dos$(*),Drug_ptr
    3130      Out_graph=2
    3140      Pres_stl=0
    3150      Lab_stl=0
20 3160      Io_stl=0
    3170      Vent_stl=0
    3180      Drug_stl=0
    3190      Io_p=Io_ptr
    3200      Lab_p=Lab_ptr
25 3210      Vent_p=Vent_ptr
    3220      Pres_p=Pres_ptr
    3230      Heart_p=Heart_ptr
    3240      Drug_p=Drug_ptr
    3250      !
30 3260      ! set up identifying subject info
    3270      !
    3280      GRAPHICS OFF
    3290      PRINT CHR$(12)
    3300      PRINT TABXY(1,1);
35 3310      PRINT USING Image_wtl;Sub_name$,Hos_num$,In_
    file$

```

- 68 -

```
3320 Image_wt1:IMAGE      "Name: ",K,XXXX,"Hosp num:
                        ",K,XXXXX,K
3330      PRINT TABXY(1,2);
3340      PRINT USING Image_wt2;Id_age$,Id_wt$,Id_
5      ht$,Diag$,Opera$
3350 Image_wt2:IMAGE      "Age: ",K,XXXX,"Wt(kg):
                        ",K,XXXX,"Ht(cm): ",K,XXXX,"Diag:
                        ",K,XXXX,"Op: ",K
3360      !
10 3370      ! go to appropriate chart
3380      !
3390      ON Chart_num GOTO In_out,Lab_val,Vent_
                        val,Pres_val,Drug
3400 In_out:!                                     ....intake/output
15 3410      ! IF Io_ptr>3 THEN Io_stl=2
3420      ! IF Io_ptr>5 THEN
3430      !      DISP."do not input more Intake/Output
                        !      data; disc full"
3440      !      WAIT 3
20 3450      !      SUBEXIT
3460      ! END IF
3470      PRINT TABXY(30,3);"INTAKE/OUTPUT CHART"
3480      PRINT TABXY(1,4);"Intake (cc/hr) "
3490      PRINT TABXY(1,5);"Time"
25 3500      PRINT TABXY(4,6);"Maint. Fluid"
3510      PRINT TABXY(4,7);"Other Fluids"
3520      PRINT TABXY(1,9);"Total "
3530      PRINT TABXY(1,11);"Output (cc/hr)"
3540      PRINT TABXY(4,12);"Urine"
30 3550      PRINT TABXY(4,13);"Chest"
3560      PRINT TABXY(1,15);"Total"
3570      PRINT TABXY(1,17);"Net I/O"
3580      Start=25
3590      IF Io_ptr>3 THEN Io_p=3
35 3600 Io_dp:FOR I=Io_stl TO Io_p
3610          PRINT TABXY(Start,5);Io_time$(I)
```

- 69 -

```

3620      PRINT TABXY(Start,6);Iv_intake(I)
3630      PRINT TABXY(Start,7);Fluid_in(I)
3640      PRINT TABXY(Start,9);In_tot(I)
3650      PRINT TABXY(Start,12);Urine(I)
5  3660      PRINT TABXY(Start,13);Chest(I)
3670      PRINT TABXY(Start,15);Out_tot(I)
3680      PRINT TABXY(Start,17);Net(I)
3690      Start=Start+10
3700      NEXT I
10 3710      IF Io_ptr>Io_p THEN
3720          INPUT "more data on next page - do you
              want this dumped to printe
              r? (Y/N)",Ans$
3730          IF Ans$="Y" OR Ans$="y" THEN CALL Graph_
15      dump(Out_graph)
3740          Io_stl=4
3750          Io_p=Io_ptr
3760          Start=25
3770          FOR J=5 TO 17
20 3780          PRINT TABXY(Start,J);"      "
3790          NEXT J
3800          GOTO Io_dp
3810      END IF
3820      GOTO Finish
25 3830 !
3840 !
3850 Lab_val:!                                ...lab values
3860      !IF Lab_ptr>3 THEN Lab_stl=2
3870      !IF Lab_ptr>5 THEN
30 3880      !    DISP "do not input any more lab values;
              !    disc full"
3890      !    WAIT 3
3900      !    SUBEXIT
3910      !END IF
35 3920      PRINT TABXY(30,3);"Lab Values"
3930      PRINT TABXY(10,4);"Time"
```

- 70 -

```
3940 PRINT TABXY(1,6);"Na"
3950 PRINT TABXY(1,7);"K"
3960 PRINT TABXY(1,8);"Cl"
3970 PRINT TABXY(1,9);"HCO3"
5 3980 PRINT TABXY(1,10);"Ca"
3990 PRINT TABXY(1,11);"Hct"
4000 PRINT TABXY(1,12);"Glucose"
4010 PRINT TABXY(1,13);"Dig level"
4020 PRINT TABXY(1,14);"PT"
10 4030 PRINT TABXY(1,15);"PTT"
4040 PRINT TABXY(1,16);"Creat"
4050 PRINT TABXY(1,17);"Bun"
4060 Start=15
4070 IF Lab_ptr>3 THEN Lab_p=3
15 4080 Lab_dp:FOR I=Lab_stl TO Lab_p
4090 PRINT TABXY(Start+10,4);Lab_time$(I)
4100 PRINT TABXY(Start+10,6);Na(I)
4110 PRINT TABXY(Start+10,7);Kl(I)
4120 PRINT TABXY(Start+10,8);Cl(I)
20 4130 PRINT TABXY(Start+10,9);Hco3(I)
4140 PRINT TABXY(Start+10,10);Ca(I)
4150 PRINT TABXY(Start+10,11);Hct(I)
4160 PRINT TABXY(Start+10,12);Gluc(I)
4170 PRINT TABXY(Start+10,13);Dig(I)
25 4180 PRINT TABXY(Start+10,14);Pt(I)
4190 PRINT TABXY(Start+10,15);Ptt(I)
4200 PRINT TABXY(Start+10,16);Creat(I)
4210 PRINT TABXY(Start+10,17);Bun(I)
4220 Start=Start+10
30 4230 NEXT I
4240 IF Lab_ptr>Lab_p THEN
4250 INPUT "more data on next page - do you
want this dumped to printe
r? (Y/N)",Ans$
35 4260 IF Ans$="Y" OR Ans$="y" THEN CALL Graph_
dump(Out_graph)
```

- 71 -

```

4270          Lab_stl=4
4280          Lab_p=Lab_ptr
4290          Start=15
4300          FOR J=4 TO 17
5  4310          PRINT TABXY(Start,J);"          "
4320          NEXT J
4330          GOTO Lab_dp
4340      END IF
4350      GOTO Finish
10 4360!
4370!
4380 Vent_val:!          ....ventilation values
4390      ! IF Vent_ptr>3 THEN Vent_stl=2
4400      ! IF Vent_ptr>5 THEN Vent_stl=4
15 4410      ! IF Vent_ptr>7 THEN
4420      !     DISP "do not input any more Vent values;
         disc full"
4430      !     WAIT 3
4440      !     SUBEXIT
20 4450      ! END IF
4460      PRINT TABXY(30,3);"VENTILATION"
4470      PRINT TABXY(1,4);"Settings          Hour:"
4480      PRINT TABXY(4,5);"Rate"
4490      PRINT TABXY(4,6);"FIO2"
25 4500      PRINT TABXY(4,7);"Peak Pres"
4510      PRINT TABXY(4,8);"Peep"
4520      PRINT TABXY(4,9);"TV"
4530      PRINT TABXY(4,10);"I:E ratio"
4540      PRINT TABXY(4,11);"Mean air"
30 4550      PRINT TABXY(1,12);"Blood Gases"
4560      PRINT TABXY(4,13);"ph"
4570      PRINT TABXY(4,14);"pO2"
4580      PRINT TABXY(4,15);"pCO2"
4590      PRINT TABXY(4,16);"HCO3"
35 4600      PRINT TABXY(4,17);"BE"
4610      Start=15
```

- 72 -

```
4620      IF Vent_ptr>3 THEN Vent_p=3
4630 Vent_dp:FOR I=Vent_stl TO Vent_p
4640      PRINT TABXY(Start+10,4);Vent_time$(I)
4650      PRINT TABXY(Start+10,5);Rate(I)
5 4660      PRINT TABXY(Start+10,6);Fio2(I)
4670      PRINT TABXY(Start+10,7);Pp(I)
4680      PRINT TABXY(Start+10,8);Peep(I)
4690      PRINT TABXY(Start+10,9);Tv(I)
4700      PRINT TABXY(Start+10,10);Ie_ratio$(I)
10 4710      PRINT TABXY(Start+10,11);Airp(I)
4720      PRINT TABXY(Start+10,13);Ph(I)
4730      PRINT TABXY(Start+10,14);Po2(I)
4740      PRINT TABXY(Start+10,15);Pco2(I)
4750      PRINT TABXY(Start+10,16);Bgo3(I)
15 4760      PRINT TABXY(Start+10,17);Be(I)
4770      Start=Start+10
4780      NEXT I
4790      IF Vent_ptr>Vent_p THEN
4800          INPUT "more data on next page - do you
20          want this dumped to printe
              r? (Y/N)",Ans$
4810          IF Ans$="Y" OR Ans$="y" THEN CALL Graph_
              dump(Out_graph)
4820          Vent_stl=4
25 4830          Vent_p=Vent_ptr
4840          Start=15
4850          FOR J=4 TO 17
4860              PRINT TABXY(Start,J);"      "
4870          NEXT J
30 4880          GOTO Vent_dp
4890      END IF
4900      GOTO Finish
4910 !
4920 !
35 4930 Pres_val:!.      ....pressure values
4940      !IF Pres_ptr>12 THEN Pres_stl=5
```

- 73 -

```
4950      PRINT TABXY(9,3);"Time:"
4960      PRINT TABXY(1,4);"Systemic"
4970      PRINT TABXY(4,5);"systolic"
4980      PRINT TABXY(4,6);"diastolic"
5  4990      PRINT TABXY(4,7);"mean"
5000      PRINT TABXY(1,8);"Pulmonary"
5010      PRINT TABXY(4,9);"systolic"
5020      PRINT TABXY(4,10);"diastolic"
5030      PRINT TABXY(4,11);"mean"
10 5040      PRINT TABXY(1,12);"LA mean"
5050      PRINT TABXY(1,13);"RA mean"
5060      PRINT TABXY(9,14);"Time: "
5070      PRINT TABXY(1,15);"C.I."
5080      PRINT TABXY(1,16);"PVRI"
15 5090      PRINT TABXY(1,17);"SVRI"
5100      Start=15
5110      IF Pres_ptr>12 THEN Pres_p=12
5120 Pres_dp:FOR I=Pres_st1 TO Pres_p
5130          PRINT TABXY(Start,3);Pres_time$(I)
20 5140          PRINT TABXY(Start,5);Ao_s(I)
5150          PRINT TABXY(Start,6);Ao_d(I)
5160          PRINT TABXY(Start,7);Ao_m(I)
5170          PRINT TABXY(Start,9);Pa_s(I)
5180          PRINT TABXY(Start,10);Pa_d(I)
25 5190          PRINT TABXY(Start,11);Pa_m(I)
5200          PRINT TABXY(Start,12);La_m(I)
5210          PRINT TABXY(Start,13);Ra_m(I)
5220          Start=Start+5
5230      NEXT I
30 5240      Start=15
5250      FOR I=0 TO Heart_ptr
5260          PRINT TABXY(Start,14);Heart_time$(I)
5270          PRINT TABXY(Start,15);Ci(I)
5280          PRINT TABXY(Start,16);Pvri(I)
35 5290          PRINT TABXY(Start,17);Svri(I)
5300          Start=Start+5
```

- 74 -

```
5310      NEXT I
5320      IF Pres_ptr>Pres_p THEN
5330          INPUT "more data on next page - do you
                    want this dumped to printe
5          r? (Y/N)",Ans$
5340          IF Ans$="Y" OR Ans$="y" THEN CALL Graph_
                    dump(Out_graph)
5350          Pres_stl=13
5360          Pres_p=Pres_ptr
10 5370          Start=15
5380          FOR J=3 TO 13
5390              PRINT TABXY(Start,J);"          "
5400          NEXT J
5410          GOTO Pres_dp
15 5420      END IF
5430      GOTO Finish
5440 !
5450 !
5460 Drug:!!          ....hey man, drugs
20 5470      !IF Drug_ptr>9 THEN Drug_stl=4
5480      ! IF Drug_ptr>14 THEN Drug_stl=9
5490      !IF Drug_ptr>19 THEN Drug_stl=14
5500      !IF Drug_ptr>24 THEN Drug_stl=19
5510      !IF Drug_ptr>29 THEN Drug_stl=24
25 5520      !IF Drug_ptr>34 THEN Drug_stl=29
5530      !IF Drug_ptr>38 THEN
5540      !      DISP "do not enter more drugs; disc full"
5550      !      WAIT 3
5560      !      SUBEXIT
30 5570      ! END IF
5580      PRINT TABXY(30,4);"Drug Chart"
5590      PRINT TABXY(1,6);"Name"
5600      PRINT TABXY(30,6);"Dosage"
5610      PRINT TABXY(60,6);"Time"
35 5620      D_line=7
5630      IF Drug_ptr>9 THEN Drug_p=9
```



- 75 -

```
5640 Drug_dp:FOR I=Drug_stl TO Drug_p
5650     PRINT TABXY(1,D_line);Drug_name$(I)
5660     PRINT TABXY(30,D_line);Drug_dos$(I)
5670     PRINT TABXY(60,D_line);Drug_time$(I)
5 5680     D_line=D_line+1
5690     NEXT I
5700     IF Drug_ptr>Drug_p THEN
5710         INPUT "more data on next page - do you
want this dumped to printer? (Y/N)",Ans$
10 5720     IF Ans$="Y" OR Ans$="y" THEN CALL Graph_
dump(Out_graph)
5730     Drug_stl=Drug_stl+10
5740     Drug_p=Drug_p+10
5750     D_line=7
15 5760     FOR J=7 TO 17
5770         PRINT TABXY(1,J);"      "
5780     NEXT J
5790     GOTO Drug_dp
5800     END IF
20 5810 Finish: !
5820     SUBEND
5830     !
5840     !
5850     !
25 5860     SUB Graph_dump(A)
5870 Graph_dump:INPUT "do you want a hard copy?
<Y/N>",Ans$
5880     IF Ans$="Y" OR Ans$="y" THEN
5890         IF A=1 THEN
30 5900             DUMP GRAPHICS #701
5910             PRINTER IS 701
5911             PRINT CHR$(12)
5920             GRAPHICS OFF
5930         ELSE
35 5940             DUMP ALPHA #701
5950             PRINTER IS 701
```

- 76 -

```
5960          PRINT CHR$(12)
5970          END IF
5980          END IF
5990          PRINTER IS 1
5 6000 SUBEND
```

10

15

20

25

30

35

- 77 -

```
10 Hrsa3: !THIS IS A PROGRAM TO SET UP THE HIGH SPEED
A/D      !SYSTEM
20      ! AND CONTINUOUSLY OBTAIN INFORMATION
30      !
5  40      !
50      !.....
60      !
70      ! LAST REVISION: 30 April 1985
80      !
10 90      !.....

100      !
110      !
120      ! > FULL SET OF DECLARATIONS FOR THE HPBIB BUS
15      EXTENDED TALK ADDRESSES
130      !
140      !
150 Assignments:  !
160      ASSIGN @Multi TO 723
20 170      ASSIGN @Input_para TO 72301
180      ASSIGN @Input_intr TO 72302
190      ASSIGN @Input_ext TO 72303
200      ASSIGN @Read_format TO 72304
210      ASSIGN @Memory_input TO 72305
25 220      ASSIGN @Read_val TO 72306
230      ASSIGN @Read_status TO 72308
240      ASSIGN @Output_intr TO 72309
250      ASSIGN @Hpib_srq_status TO 72310
260      ASSIGN @Err_status_1st TO 72311
30 270      ASSIGN @Int_addr TO 72312
280      ASSIGN @Busy_instr TO 72313
290      ASSIGN @Read_clock TO 72314
300      !
310      !
35
```

- 78 -

```

320  !.....
330  !
340  ! SET UP INTERRUPT/ERROR HANDLERS
350  ! SET UP COMMON STORAGE/ARRAY STORAGE
5   360  !.....

370  !
380  !
390  COM /Intr_7/ Int_flag,Status_bytes(5)
10  400  COM /Flags/ Atod_done,Scanner_done,Memory1_
      done,Memory2_done,Timer_done,Counter_done,
      Memory3_done,Memory4_done
410  COM /Io_arrays/ Counters(3),Counters2(3),Time_
      base$(7)
15  420  COM /Multi_param/ Start_chan,Stop_chan,Pacing_
      bits,Pacing_rate,Num_pts,Num_xfer,Num_xfer_
      left,Name_len,Scr_file$(28),Scr_
      file2$(28)
430  COM /Hr_sig/ Num_pulses,Last_pulse,First_blk_
20  flg,Last_time,Num_hr_sig,Max_hr_pts,Avg_
      hr,Rollover,Hr_smooth
440  COM /Plot_par/ Plotbox,Boxcar_flg,Log_
      plot_flg,Freq_limit,Resp_search,Pct_thresh
450  COM /Graphs/
25  Hrdata(512),Hrspec(512),Respspec(512),Bpspec(512)
460  COM /Vitaldata/ Rfa,Lfa,Peakratio,Meas_resp,Next_
      time
470  COM /Messagecom/ Message$(10)[80],@Messages
480  COM /Trends/ Mean_hr_t(60),Lfa_t(60),Rfa_
30  t(60),Ratio_t(60),T_ptr,Time_now 1,Meas_resp_
      t(60),Trend_dp
490  COM /Pressure/
      Top1,Top2,Top3,Top4,Bot1,Bot2,Bot3,Bot4
500  COM /Editor/ Edit_msg$(80)
35  510  COM /Subject/ Sub_name$(25),Hos_num$(15),Id_
      age$(10),Id_wt$(10),Id_ht$(10 ),Diag$(30),

```

- 79 -

```

      Opera$(45),Halt_pg
520  COM /Io_chart/ Io_time$(8)[10],Iv_intake(8),Fluid_
      in(8),In_tot(8),Urine(8 ),Chest(8),Out_
      tot(8),Net(8),Io_ptr
5  530  COM /Lab_chart/ Lab_time$(8)[10],Na(8),Kl(8),
      Cl(8),Hco3(8),Ca(8),Hct(8),G luc(8),
      Dig(8),Pt(8),Ptt(8),Creat(8),Bun(8),Lab_ptr
540  COM /Vent_chart/ Vent_
      time$(8)[15],Rate(8),Fio2(8),Pp(8),Peep(8),Tv(8),
10  Ie_ratio$(8)[10],Airr(8),Ph(8),Po2(8),
      Pco2(8),Bgo3(8),Be(8),Vent_ptr
550  COM /Pres_chart/ Pres_time$(20)[15],Ao_s(20),Ao_
      d(20),Ao_m(20),Pa_s(20),Pa_d(20),Pa_m(20),
      La_m(20),Ra_m(20),Pres_ptr,Pres_in
15  560  COM /Heart_index/ Heart_
      time$(15)[15],Ci(15),Pvri(15),Svri(15),Heart_ptr
570  COM /Drugs/ Drug_time$(40)[20],Drug_
      name$(40)[40],Drug_dos$(40)[20],Drug_ptr
590  DIM Io$(5,15)[30],Io_msg$(5,15)[80]
20  600  DIM Msg_pad$(10)[80]
610  DIM Msg_buffer$(80) BUFFER
620  ASSIGN @Msg_buffer TO BUFFER Msg_buffer$
630  Log_plotflg=0
640  Freq_limit=1.
25  650  Resp_search=.1
660  Pct_thresh=.2
670  Scr_file$="?"
680  Halt_pg=0
690  Message$(0)="messages in "
30  700  Message$(1)="I/O chart "
710  Message$(2)="lab values"
720  Message$(3)="hemodynamics"
730  Message$(4)="Trends Display"
740  Message$(5)="messages out"
35  750  Message$(6)="STOP PROGRAM"
760  Message$(7)="ventilation"

```

- 80 -

```

770    Message$(8)="drugs"
780    Message$(9)="B.P. Display"
790    Msg_pad_ptr=0
800    P_ptr=0
5    810    !
820    ! Set up common/array storage for waveform
        analysis
830    !
840    !.....
10   850    !
860    ! Set up common/array storage for waveform
        analysis
870    !.....

15   880    !
890    COM /Directory/ Dir$(160),@Printer
900    COM /Wf1/ Printer,Plotter,String$(40]
910    COM /Wf2/ Signal(1089),Number_pnts,Type,Sampling_
        period
20   920    COM /Wf3/ Segment_size,Overlap,Num_segments,Pnts_
        used,Fft_size
930    COM /Wf5/ Refn(63),Refd(63),Refno,Refdo,Refgain
940    COM /Autoparam/ Up_down,Up_delay,Dn_delay
950    COM /Vars/ Ffthrvar,Fftrespvar
25   960    !
970    DISP "loading subroutines"
980    LOADSUB ALL FROM "multi_subs"
990    LOADSUB ALL FROM "hr_siggen8"
1000   LOADSUB ALL FROM "automaxsb2"
30   1010   LOADSUB ALL FROM "fft_anal6"
1020   DISP "load data disks and press CONTINUE"
1030   PAUSE
1040   !
1050   !.....
35   1060   ! The HP 9826/9836 flexible disk (5-1/4") has the
        ! following structure

```

- 81 -

```

1070  ! 2 sides, 33 tracks/side, 16 sectors/track, 256
      ! bytes/sector
1080  ! 1 track = 4096 bytes = 16 sectors
1090  ! 1 side = 135168 bytes = 528 sectors
5    1100 ! 1 disk = 270336 bytes = 1056 sectors
      1110 ! 1 disk = 135168 words = 132K words
      1120 !.....

      1130 !
10    1140 !
      1150 INTEGER Hpib_buffer1(2048) BUFFER
      1160 INTEGER Hpib_buffer2(2048) BUFFER
      1170 DIM Hr_signal(1024) BUFFER
      1180 Read_ptr1=0
15    1190 Read_ptr2=0
      1200 !
      1210 !
      1220 !.....
      1230 ! CLEAR MULTIPROGRAMMER
20    1240 !.....

      1250 !
      1260 !
      1270 ON INTR 7 CALL Hpib_intr
25    1280 Begin:CALL Multi_clear
      1290 !
      1300 !
      1310 !.....
      1320 ! LOAD SUPPLEMENTAL INSTRUCTION SET ("MR")
30    1330 ! usage: "MR,<card addr>,<# words>,<read
      ptr>,<mode>T"
      1340 !           <mode= 1-FIFO, 4-recirculating>
      1350 !.....

35    1360 !
      1370 !

```

- 82 -

```
1380  DISP "DOWNLOADING MR INSTRUCTION"
1390  CALL Xfer("MR")
1400  !
1410  !
5  1420  !.....
1430  ! SET UP CARDS FOR DATA COLLECTION
1440  !.....

1450  !
10 1460  !
1470 Selections:DISP "SETUP DATA COLLECTION"
1480  OUTPUT @Multi;"CY,3T"!CYCLE SCAN/PACER CARD TO
      SET DEFINITE STATE
1490  !
15 1500  !
1510  ! NOW SET UP THE SCAN CARD PARAMETERS (DEFAULT
      ! VALUES)
1520  !      START CHANNEL (3.0) - 0
1530  !      STOP CHANNEL (3.1) - 1
20 1540  !      PACING (3.2) - 40 USEC
1550  !      SEQN'L SCAN (3.3) - XXXX XXXX XXX1 ( 1)
1560  !      INTN'L PACING (3.3) - XXXX XXXX X1XX ( 4)
1570  !      MSEC TIMEBASE (3.3) - XXX1 XXXX XXXX (256)
1580  !
25 1590  CALL Get_param
1600  ASSIGN @Messages TO
      "messglog:HP8290X,700,1";FORMAT OFF
1610  ASSIGN @Temp_trend TO "temp_
      trend:HP8290X,700,1";FORMAT OFF
30 1620  ASSIGN @Hemo_data TO "hemo_
      data:HP8290X,700,1";FORMAT OFF
1630  ASSIGN @Io_data TO "io_data:HP8290X,700,1";FORMAT
      OFF
1640  ASSIGN @Lab_data TO "lab_
35  data:HP8290X,700,1";FORMAT OFF
1650  ASSIGN @Vent_data TO "vent_
```



- 83 -

```
data:HP8290X,700,1";FORMAT OFF
1660  ASSIGN @Co_data TO "co_data:HP8290X,700,1";FORMAT
OFF
1670  ASSIGN @Drug_data TO "drug_
5  data:HP8290X,700,1";FORMAT OFF
1680  IF Num_pts=0 THEN GOTO Begin
1690  Read_ptr1=0
1700  !
1710  !
10  1720  ! SET FIFO MODE AND CLEAR POINTERS IN MEMORY
1730  !
1740  !
1750  Setup_scan:DISP " NUMBER OF POINTS=";Num_pts
1760  OUTPUT @Multi;"WF,3.0",Start_chan,"3.1",Stop_
15  chan,"3.3",Pacing_bits,"3.2"
    ,Pacing_rate,"T"
1770  OUTPUT @Multi;"CC,6T"
1780  OUTPUT @Multi;"WF,5.1,1,T" ! memory set to FIFO
input mode
20  1790  OUTPUT @Multi;"AC,3,5,6T" ! cards are armed to
supply interrupts
1800  OUTPUT @Multi;"RV,6.0,6.1,6.2,6.3T" ! checking
control registers
1810  ENTER @Read_val;Counters(*)
25  1820  Read_ptr1=0
1830  Read_ptr2=0
1840  !
1850  ! setup the counter card to count
1860  !
30  1870  Setup_counter:OUTPUT @Multi;"CC,10,11,12,13T"
1880  OUTPUT @Multi;"AC,10,12,13T" !_counter not armed
1890  OUTPUT @Multi;"CY,11T"
1900  !
1910  ! setup the pacer card to generate a clock with
35  period 32 Usec
1920  ! (one half period is 16 Usec)
```

- 84 -

```
1930 !      (corresponds to 31.25KHz)
1940 !
1950 Setup_clock:OUTPUT @Multi;"WF10.2,1T"
1960   OUTPUT @Multi;"WF10,16U T"
5  1970   CALL Completer("setup completed")
1980 !
1990 !
2000 ! START THE PACERS BY CYCLING IN PARALLEL
2010 !
10 2020   OUTPUT @Multi;"GPT"
2030   CALL Init_flags
2040   ENABLE INTR 7;2
2050   OUTPUT @Multi;"CY,3,10T"
2060   OUTPUT @Multi;"GST"
15 2070   Start_pacing=TIMEDATE
2080   CALL Completer("PACING STARTED")
2090   Block_time=Pacing_rate*1.024
2100   Next_time=TIMEDATE+INT(Block_time)
2110   First_blk_flg=1
20 2120   Num_msgs=0
2130   Message_line=0
2140   Msg_dp_request=0
2150   Resp_dpflg=0
2160   Max_hr_pts=1024
25 2170   Last_time=0
2180   Trend_dp=0
2190   !Hemo_dp=0
2200   Top1=0
2210   Top2=0
30 2220   Top3=0
2230   Top4=0
2240   Bot1=0
2250   Bot2=0
2260   Bot3=0
35 2270   Bot4=0
2280 !
```

- 85 -

```
2290   Io$(1,1)="Time - hh:mm(hh=1 to 24)"
2300   Io$(1,2)="Maint. fluids"
2310   Io$(1,3)="other fluids"
2320   Io$(1,4)="urine output"
5  2330   Io$(1,5)="chest output"
2340   Io$(2,1)="Time - hh:mm"
2350   Io$(2,2)="Na"
2360   Io$(2,3)="K"
2370   Io$(2,4)="Cl"
10 2380   Io$(2,5)="HCO3"
2390   Io$(2,6)="Ca"
2400   Io$(2,7)="Hct"
2410   Io$(2,8)="Glucose"
2420   Io$(2,9)="Dig level"
15 2430   Io$(2,10)="PT"
2440   Io$(2,11)="PTT"
2450   Io$(2,12)="Creat"
2460   Io$(2,13)="Bun"
2470   Io$(3,1)="Time - hh:mm(hh=1 to 24)"
20 2480   Io$(3,2)="Resp rate"
2490   Io$(3,3)="FIO2"
2500   Io$(3,4)="Peak pres"
2510   Io$(3,5)="peep"
2520   Io$(3,6)="TV"
25 2530   Io$(3,7)="I:E"
2540   Io$(3,8)="mean airway"
2550   Io$(3,9)="ph"
2560   Io$(3,10)="pO2"
2570   Io$(3,11)="pCO2"
30 2580   Io$(3,12)="HCO3"
2590   Io$(3,13)="BE"
2600   Io$(4,1)="Time - hh:mm(hh=1 to 24)"
2610   Io$(4,2)="ao/s"
2620   Io$(4,3)="ao/d"
35 2630   Io$(4,4)="ao/m"
2640   Io$(4,5)="pa/s"
```

- 86 -

```
2650   Io$(4,6)="pa/d".
2660   Io$(4,7)="pa/m"
2670   Io$(4,8)="la/m"
2680   Io$(4,9)="ra/m"
5  2690   Io$(4,10)="Time - hh:mm(hh=1 to 24)"
2700   Io$(4,11)="C.I."
2710   Io$(4,12)="pvri"
2720   Io$(4,13)="svri"
2730   Io$(5,1)="name"
10 2740   Io$(5,2)="dosage"
2750   Io$(5,3)="Time - hh:mm:ss(hh=1 to 24)"
2760   Io_ptr=0
2770   Lab_ptr=0
2780   Vent_ptr=0
15 2790   Pres_ptr=0
2800   Heart_ptr=0
2810   Drug_ptr=0
2820   Io_in=0
2830   Lab_in=0
20 2840   Vent_in=0
2850   Pres_in=0
2860   Heart_in=0
2870   Drug_in=0
2880   Fst=1
25 2890   Fix_val=0
2900   !
2910   ! Read data continuously and write to the disk
      continuously until enough
2920   ! enough data has been obtained
30 2930   !
2940   !
2950   Reading:  !
2960   !
2970   ! set up the A/D buffers and disk files
35 2980   !
2990   ASSIGN @Memory_input TO 72305;FORMAT OFF
```

- 87 -

```
3000  ASSIGN @In_buffer TO BUFFER Hpib_buffer1(*)
3010  ASSIGN @Out_buffer TO Scr_file$;FORMAT OFF
3020  !
3030  ! set up the counter memory buffers and files
5  3040  !
3050  ASSIGN @Memory_input2 TO 72305;FORMAT OFF
3060  ASSIGN @In_buffer2 TO BUFFER Hpib_buffer2(*)
3070  ASSIGN @Out_buffer2 TO Scr_file2$;FORMAT OFF
3080  !
10 3090  Data_lockout=0
3100  !
3110  Time_now=TIMEDATE
3120  Date_now$=DATE$(TIMEDATE)
3130  Time_now1=Time_now MOD 86400
15 3140  !
3150  Blk_xfer: !
3160  CONTROL @In_buffer,3;1
      ! Reset fill pointer for buffer
3170  CONTROL @In_buffer,4;0
20      ! Reset current number of bytes in buffer
3180  CONTROL @In_buffer,5;1      ! Reset empty pointer
      for buffer
3190  !
3200  ! write an 8 byte sequence to disk as a header for
25      ! the transfer
3210  !
3220  CALL Xfheader(@Out_buffer,Num_pts,"R")
3230  !
3240  ! read A/D buffer into memory (hpib_buffer1) in 32
30      segments
3250  ! if possible
3260  !
3270  IF FRACT(Num_pts/32.)=0 THEN
3280      Num_rdseg=32
35  3290      Num_rdtype=Num_pts/32
3300  ELSE
```

- 88 -

```
3310      Num_rdseg=1
3320      Num_rdpts=Num_pts
3330      END IF
3340      !
5 3350      ! reading segments here. segmenting allows disk
      access between segments
3360      !
3370      FOR Rdseg=1 TO Num_rdseg
3380          OUTPUT @Multi;"MR,5",Num_rdpts,Read_
10      ptr1,"1T"! FIFO mode
3390          ON EOT @Memory_input GOTO Next_rdseg
3400          TRANSFER @Memory_input TO @In_buffer;COUNT
      Num_rdpts*2,CONT
3410          PRINT TABXY(1,18);
15 3420          PRINT USING Image_wtl;Num_xfer-Num_xfer_
      left+1,Num_xfer,TIMES$(Next_time),
      Rdseg,Num_rdseg
3430 Image_wtl:IMAGE      "Next xfer(",K,"/",K,"): ",K,"
      seg=",K,"/",K
20 3440 Waiter1:DISP "Now: ";TIMES$(TIMEDATE);"
      ";DATE$(TIMEDATE)
3450      IF Next_time-TIMEDATE<12 THEN
3460          OFF KEY
3470          OFF KBD
25 3480          OFF KNOB
3490          GOTO Waiter1
3500      END IF
3510      ON KEY 0 LABEL Message$(0) GOSUB Key0
3520      ON KEY 1 LABEL Message$(1) GOSUB Key1
30 3530      ON KEY 2 LABEL Message$(2) GOSUB Key2
3540      ON KEY 3 LABEL Message$(3) GOSUB Key3
3550      ON KEY 4 LABEL Message$(4) GOSUB Key4
3560      ON KEY 5 LABEL Message$(5) GOSUB Key5
3570      ON KEY 6 LABEL Message$(6) GOSUB Key6
35 3580      ON KEY 7 LABEL Message$(7) GOSUB Key7
3590      ON KEY 8 LABEL Message$(8) GOSUB Key8
```

- 89 -

```
3600      ON KEY 9 LABEL Message$(9) GOSUB Key9
3610      ON KBD GOTO Control_chars
3620      IF Msg_dp_request=2 THEN
3630          ON KNOB .05 GOSUB Move_msgs
5 3640      ELSE
3650          OFF KNOB
3660      END IF
3670      STATUS @In_buffer,10;In_xfer_stat
3680      IF In_xfer_stat<64 THEN GOTO Next_rdseg
10 3690      IF Msg_dp_request=3. THEN
3700          CALL Msg_dump(Message_chart$(*),Message_
            line,Msg_dp_request)
3710      END IF
3720      GOTO Waiter1
15 3730 Control_chars:~
3740          Kbd_hold$=KBDS$
3741          IF POS(Kbd_hold$,CHRS(6))<>0 THEN
            !..change lfa disp.range
3742              Lfa_top=Lfa_top+2.5
20 3750          IF POS(Kbd_hold$,CHRS(6))<>0 THEN
            !..change spectra disp.freq.range
3760              IF Freq_limit=1. THEN
3770                  Freq_limit=2.
3780              ELSE
25 3790                  Freq_limit=1.
3800              END IF
3810              Resp_search=.1
            !..reset resp search point each time
3820              DISP "Spectra displayed to";Freq_
30              limit;"Hz"
3830              WAIT 2
3840          END IF
3850          IF POS(Kbd_hold$,CHRS(8))<>0 THEN !..help:
            display commands
35 3860              CALL Disp_ctrls
3870          END IF
```

- 90 -

```
3880      IF POS(Kbd_hold$,CHR$(16))<>0 THEN
      !..change peak search threshold
3890      Pct_thresh=Pct_thresh+.2
3900      IF Pct_thresh>.8 THEN Pct_thresh=.2
5 3910      DISP "resp peak search threshold=";Pct_
      thresh;"%"
3920      WAIT 1
3930      END IF
3940      IF POS(Kbd_hold$,CHR$(18))<>0 THEN
10 !..display respiration time series
3950      IF Resp_dpflg=0 THEN
3960      Resp_dpflg=1
3970      DISP "resp series plot w/hr series"
3980      WAIT 2
15 3990      ELSE
4000      Resp_dpflg=0
4010      DISP "cancel resp series plot"
4020      WAIT 2
4030      END IF
20 4040      END IF
4050      IF POS(Kbd_hold$,CHR$(19))<>0 THEN
      !..change respiration peak search
4060      Resp_search=Resp_search+.1
4070      IF Resp_search>Freq_limit-.1 THEN Resp_
25 search=.1
4080      DISP "resp peak search starts at";Resp_
      search;"Hz"
4090      WAIT 1
4100      END IF
30 4110      GOTO Waiter1
4120 Next_rdseg:!
4130 !
4140 ! storing messages from soft keys if any
4150 !
35 4160      IF Msg_pad_ptr>0 THEN
4170      Num_msgs=Num_msgs+Msg_pad_ptr
```



- 91 -

```

4180          FOR I=0 TO Msg_pad_ptr-1
4190              Msg_buffer$=Msg_pad$(I)
4200              Len_message=LEN(Msg_buffer$)
4210              CONTROL @Msg_buffer,4;Len_
5              message      !...number of bytes
4220              CONTROL @Msg_buffer,5;1
              !..empty pointer to beginning
4230              TRANSFER @Msg_buffer TO
                  @Messages;COUNT Len_message,CONT
10 4240          NEXT I
4250          IF Msg_dp_request>=2 THEN
4260              DEALLOCATE Message_chart$(*)
4270              Msg_dp_request=0
4280          END IF
15 4290          OFF KNOB
4300          Msg_pad_ptr=0
4310          END IF
4320          IF Msg_dp_request=1 THEN
4330              Message_line=0
20 4340              ALLOCATE Message_chart$(17)[640]
4350              CALL Msg_dump(Message_chart$(*),Message_
                  line,Msg_dp_request)
4360              IF Msg_dp_request=0 THEN
                  !...no messages
25              yet
4370                  DEALLOCATE Message_chart$(*)
4380              END IF
4390          END IF
4400          !
30 4410          ! get read pointer for next segment
4420          !
4430          OUTPUT @Multi;"RV,6.0T"
              ! checking current read pointer
4440          ENTER @Read_val;Read_ptr1
35 4450          NEXT Rdseg
4460          !

```

- 92 -

```
4470 ! store A/D buffer on complete data file (also
      save pointers for heart rate)
4480 !
4490 !
5  4500 Resume1:OFF EOT @Memory_input
      4510 OFF KEY
      4520 OFF KBD
      4530 OFF KNOB
      4540 IF Msg_dp_request>=2 THEN
10 4550 DEALLOCATE Message_chart$(*)
      4560 Msg_dp_request=0
      4570 END IF
      4580 IF Trend_dp=1 OR Trend_dp=2 THEN DEALLOCATE
          Spectra(*)
15 4590 Next_time=Next_time+INT(Block_time)
      4600 ON EOT @Out_buffer GOTO Resume2
      4610 OUTPUT @Multi;"RV,13.0,13.1,13.2,13.3T"
          ! checking control registers
      4620 ENTER @Read_val;Counters2(*)
20 4630 Read_ptr2=Counters2(0)
      4640 Num_pulses=Counters2(1)
      4650 TRANSFER @In_buffer TO @Out_buffer;COUNT Num_
          pts*2,CONT
      4660 Waiter2:DISP TIMES$(TIMEDATE),DATES$(TIMEDATE)
25 4670 GOTO Waiter2
      4680 !
      4690 !
      4700 !
      4710 !
30 4720 Resume2:OFF EOT @Out_buffer
      4730 Num_xfer_left=Num_xfer_left-1
      4740 OUTPUT @Multi;"MR,12",Num_pulses,Read_
          ptr2,"1T" ! FIFO mode
      4750 CONTROL @In_buffer2,3;1
35 ! Reset fill pointer for buffer
      4760 CONTROL @In_buffer2,4;0
```

- 93 -

```
      ! Reset current number of bytes in buffer
4770   CONTROL @In_buffer2,5;1
      ! Reset empty pointer for buffer
4780   !
5 4790   ! write an 8 byte sequence to disk as a header for
      ! the transfer
4800   !
4810   CALL Xfheader(@Out_buffer2,Num_pulses,"H")
4820   !
10 4830   ! read multiprogrammer into computer memory (hpib_
      buffer)
4840   !
4850   ON EOT @Memory_input2 GOTO Resume4
4860   TRANSFER @Memory_input2 TO @In_buffer2;COUNT Num_
15   pulses*2,CONT
4870 Waiter4:DISP TIME$(TIMEDATE),DATE$(TIMEDATE)
4880   GOTO Waiter4
4890   !
4900   ! store computer memory on complete data file
20 4910   !
4920 Resume4:OFF EOT @Memory_input2
4930   ON EOT @Out_buffer2 GOTO Resume5
4940   TRANSFER @In_buffer2 TO @Out_buffer2;COUNT Num_
      pulses*2,CONT
25 4950 Waiter5:DISP TIME$(TIMEDATE),DATE$(TIMEDATE)
4960   GOTO Waiter5
4970   !
4980 Resume5:OFF EOT @Out_buffer2
4990   CALL Hr_sig_gen(Hpib_buffer2(*),Hr_signal(*))
30 5000   !

5010   !
5020 Resume6:!
35 5030   OUTPUT @Multi;"RV,6.0,6.1,6.2,6.3T"
      ! checking control registers
```

- 94 -

```
5040  ENTER @Read_val;Counters(*)
5050  Read_ptr1=Counters(0)
5060  IF Counters(1)=4095 THEN ! Data lockout probably
      occurred
5  5070      PRINT "DATA LOCKOUT!! TIME RECORD
      NOT      CONTINUOUS!!"
5080      PRINT "ABORTING CURRENT DATA COLLECTION."
5090      Data_lockout=1
5100      Num_xfer_left=0
10  5110  END IF
5120  OUTPUT 2;CHR$(255)&CHR$(75);
      ! Clear CRT of text
5130  GINIT
5140  PLOTTER IS 3,"INTERNAL"
15  5150  GRAPHICS ON
5160  Xscale=8
5170  Hr_max=MAX(Hr_signal(*))
5180  Hr_min=MIN(Hr_signal(*))
5190  VIEWPORT 0,64,50,100
20  5200  WINDOW 0,1,0,1
5210  AXES .1,.1,0,0
5220  CSIZE 4
5230  Hr_signal(1024)=0
5240  Hr_sigsum=SUM(Hr_signal)
25  5250  Mean_hr=INT((Hr_sigsum/1024+Avg_hr))
5260  Hr_bias=Hr_sigsum/1024
5270  LDIR 0
5280  LORG 3
5290  MOVE .2,.9
30  5300  LABEL "HR data   hr=";Mean_hr
5310  CSIZE 4
5320  MOVE .05,1
5330  LORG 3
5340  LABEL "250 bpm"
35  5350  WINDOW 1,0,1,0
5360  AXES 0,0,0,0
```

- 95 -

```
5370  IF Hr_dispflg=1 THEN
5380      WINDOW 0,1024,Hr_min,Hr_max
5390  ELSE
5400      Low_window=INT(-Avg_hr)
5  5410      High_window=Low_window+250.
5420      WINDOW 0,1024,Low_window,High_window
5430  END IF
5440  FOR I=0 TO 1023
5450      PLOT I,Hr_signal(I)
10 5460  NEXT I
5470  !
5480  ! display respirations time series also
5490  !
5500  IF Resp_dpflg=1 THEN
15 5510      Max_resp=MAX(Hpib_buffer1(*))
5520      Min_resp=MIN(Hpib_buffer1(*))
5530      IF Mean_hr>100 THEN
5540          VIEWPORT 0,64,50,65
5550      ELSE
20 5560          VIEWPORT 0,64,75,90
5570      END IF
5580      WINDOW 0,1023,Min_resp,Max_resp
5590      MOVE 0,Hpib_buffer1(0)
5600      FOR I=1 TO 1023
25 5610          PLOT I,Hpib_buffer1(I)
5620      NEXT I
5630  ELSE
5640      Resp_dpflg=0
5650  END IF
30 5660  !
5670  ! now process heart rate data with waveform
      analysis package
5680  ! make sure the hr_signal has zero mean
5690  !
35 5700  FOR I=0 TO 1023
5710      Signal(I)=Hr_signal(I)-Hr_bias
```

- 96 -

```
5720     NEXT I
5730     Plotbox=2
5740     DISP "HR fft in process"
5750     CALL Wf_analyzer(Pacing_rate)
5 5760     !
5770     ! now process respiration data with waveform
        analysis package
5780     !
5790     MAT Signal= (0)
10 5800     FOR I=0 TO 1023
5810         Signal(I)=Hpib_buffer1(I)
5820     NEXT I
5830     Signal_avg=SUM(Signal)/1024.
5840     MAT Signal= Signal-(Signal_avg)
15 5850     Plotbox=4
5860     DISP "RESP fft in process"
5870     CALL Wf_analyzer(Pacing_rate)
5880     Trend_dp=0 !..trend graph not displayed
5890     !
20 5900     ! waveform analysis completed, compile trends and
        store in temporary file
5910     !
5920     Mean_hr_t(T_ptr)=Mean_hr
5930     Lfa_t(T_ptr)=Lfa
25 5940     Rfa_t(T_ptr)=Rfa
5950     Ratio_t(T_ptr)=Peakratio
5960     Meas_resp_t(T_ptr)=Meas_resp
5961     Trans_time(T_ptr)=Xfer_time
5970     T_ptr=T_ptr+1
30 5980     OUTPUT @Temp_trend;T_ptr-1,Mean_
        hr,Lfa,Rfa,Peakratio,Meas_resp,Xfer_time
5990     IF Pres_in=1 THEN
6000         Pr=Pres_ptr-1
6010         OUTPUT @Hemo_data;Pres_time$(Pr),Ao_s(Pr),Ao_
35         d(Pr),Ao_m(Pr),Pa_s(Pr),
            Pa_d(Pr),Pa_m(Pr),La_m(Pr),Ra_m(Pr),Pr
```

- 97 -

```

6020      Pres_in=0
6030      END IF
6040      IF Io_in=1 THEN
6050          Io=Io_ptr-1
5  6060          OUTPUT @Io_data;Io_time$(Io),Iv_
            intake(Io),Fluid_in(Io),In_tot(Io),Ur
            ine(Io),Chest(Io),Out_tot(Io),Net(Io),Io
6070          Io_in=0
6080      END IF
10 6090      IF Lab_in=1 THEN
6100          L=Lab_ptr-1
6110          OUTPUT @Lab_data;Lab_
            time$(L),Na(L),Kl(L),Cl(L),Hco3(L),Ca(L),Hct(L),
            Gluc(L),Dig(L),Pt(L),Ptt(L),Creat(L),Bun(L),L
15 6120          Lab_in=0
6130      END IF
6140      IF Heart_in=1 THEN
6150          H=Heart_ptr-1
6160          OUTPUT @Co_data;Heart_
20          time$(H),Ci(H),Pvri(H),Svri(H),H
6170          Heart_in=0
6180      END IF
6190      IF Vent_in=1 THEN
6200          V=Vent_ptr-1
25 6210          OUTPUT @Vent_data;Vent_
            time$(V),Rate(V),Fio2(V),Pp(V),Peep(V),Tv(V),
            Ie_ratio$(V),Airp(V),Ph(V),Po2(V),Pco2(V),
            Bgo3(V),Be(V),V
6220          Vent_in=0
30 6230      END IF
6240      IF Drug_in=1 THEN
6250          D=Drug_ptr-1
6260          OUTPUT @Drug_data;Drug_time$(D),Drug_
            name$(D),Drug_dos$(D),D
35 6270          Drug_in=0
6280      END IF
```

- 98 -

```
6290 !
6300 ! continue with data collection
6310 !
6320 IF Num_xfer_left<=0 THEN
5 6330     Halt_pg=1
6340     GOTO Eo_blk_xfer
6350 ELSE
6360     DISP Num_xfer_left;"transfers remaining"
6370     WAIT 3
10 6380     GOTO Blk_xfer
6390 END IF
6400 Eo_blk_xfer:End_time=TIMEDATE
6410 Delta_time=End_time-Start_time
6420 !
15 6430 OUTPUT @Multi;"WF,3.2,0T"
6440 Stop_pacing=TIMEDATE
6450 !
6460 Aborter: !
6470 ASSIGN @In_buffer TO *
20 6480 ASSIGN @In_buffer2 TO *
6490 ASSIGN @Out_buffer TO *
6500 ASSIGN @Out_buffer2 TO *
6510 ASSIGN @Messages TO *
6520 ASSIGN @Temp_trend TO *
25 6530 ASSIGN @Hemo_data TO *
6540 ASSIGN @Io_data TO *
6550 ASSIGN @Lab_data TO *
6560 ASSIGN @Vent_data TO *
6570 ASSIGN @Co_data TO *
30 6580 ASSIGN @Drug_data TO *
6590 OUTPUT @Multi;"CC,3,5,6,10,11,12,13T"
6600 OUTPUT @Multi;"CC,5T"
6610 CALL Completer("READY TO RESTART")
6620 CALL Pauser
35 6630 GRAPHICS OFF
6640 CALL Get_param
```



- 99 -

```
6650  ASSIGN @Messages TO
        "messglog:HP8290X,700,1";FORMAT OFF
6660  IF Num_pts=0 THEN GOTO Begin
6670  GOTO Setup_scan
5  6680  Diag:OUTPUT 723;"RV,3.0,3.3T"
6690  ENTER 72306;C,C0
6700  PRINT "CURRENT/START CHANNEL";C,C0
6710  OUTPUT 723;"RV,6.0,6.1,6.2,6.3T"
        ! checking control registers
10  6720  ENTER 72306;Counters(*)
6730  PRINT "COUNTERS=";Counters(*)
6740  STOP
6750  Purger:!
6760  GRAPHICS OFF
15  6770  DELSUB Hpib_intr TO END
6780  PURGE "AOK:HP8290X,700,1"
6790  PURGE "hrAOK:HP8290X,700,1"
6800  PURGE "messglog:HP8290X,700,1"
6810  PURGE "temp_trend:HP8290X,700,1"
20  6820  PURGE "hemo_data:HP8290X,700,1"
6830  PURGE "co_data:HP8290X,700,1"
6840  PURGE "vent_data:HP8290X,700,1"
6850  PURGE "lab_data:HP8290X,700,1"
6860  PURGE "drug_data:HP8290X,700,1"
25  6870  PURGE "io_data:HP8290X,700,1"
6871  PURGE "sub_data:HP8290X,700,1"
6880  STOP
6890  !
6900  ! definitions for keys
30  6910  !
6920  Move_msgs:! knob is processed here
6930  IF Msg_dp_request<>2 THEN RETURN
6940  Message_line=Message_line+KNOBX
6950  IF Message_line>Num_msgs-3 THEN Message_line=Num_
35  msgs-3
6960  IF Message_line<0 THEN Message_line=0
```

- 100 -

```
6970  Msg_dp_request=3
6980  RETURN
6990  !
7000  !
5    7010  Key0:Key_id=0
      7020  Edit_msg$=""
      7030  CALL Editor
      7040  Key_msg:Msg_pad$(Msg_pad_
            ptr)="Time:"&TIME$(TIMEDATE)&" "&Edit_msg$
10   7050  Msg_pad_ptr=Msg_pad_ptr+1
      7060  DISP "only";10-Msg_pad_ptr;"more messages during
            this segment"
      7070  PRINT TABXY(1,18);"
            "
15   7080  PRINT TABXY(1,18);Edit_msg$
      7090  WAIT 3
      7100  PRINT TABXY(1,18);"
            "
      7110  PRINT TABXY(1,18);"Next transfer: ";TIME$(Next_
20   time)
      7120  GOTO Keyend
      7130  !
      7140  !
      7150  !
25   7160  Key1:Chart_num=1
            !...input/output charting
      7170  IF Next_time-TIMEDATE<45 THEN
      7180      DISP "not enough time to enter data; wait for
            next xfer"
30   7190      WAIT 2
      7200      GOTO Keyend
      7210  END IF
      7220  GRAPHICS OFF
      7230  PRINT CHR$(12)
35   7240  Num_var=5
      7250  IF Io_in=1 THEN
```

- 101 -

```

7260      DISP "data in for this xfer; chart displayed"
7270      WAIT 2
7280      Io_ptr=Io_ptr-1
7290      CALL Chart(Char_num)
5 7300      Io_ptr=Io_ptr+1
7310      GOTO Keyend
7320  ELSE
7330      INPUT "Input values=1 or display
      chart=2?",Inp
10 7340      IF Inp=1 THEN
7350          IF Io_ptr>5 THEN
7360              DISP "Do not enter more I/O data;
              disc full"
7370              WAIT 3
15 7380              GOTO Keyend
7390          ELSE
7400              GOTO I_o
7410          END IF
7420      ELSE
20 7430          CALL Chart(Char_num)
7440          GOTO Keyend
7450      END IF
7460  END IF
7470  Data1:
25 7480      Io_time$(Io_ptr)=Io_msg$(Char_num,1)
7490      Iv_intake(Io_ptr)=FNLval(Io_msg$(Char_num,2))
7500      IF Iv_intake(Io_ptr)=9999.999 THEN
7510          Ionum=2
7520          Fix_val=1
30 7530          GOTO Data_edit
7540      END IF
7550      Fluid_in(Io_ptr)=FNLval(Io_msg$(Char_num,3))
7560      IF Fluid_in(Io_ptr)=9999.999 THEN
7570          Ionum=3
35 7580          Fix_val=1
7590          GOTO Data_edit
```

- 102 -

```
7600     END IF
7610     Urine(Io_ptr)=FNLval(Io_msg$(Chart_num,4))
7620     IF Urine(Io_ptr)=9999.999 THEN
7630         Ionum=4
5 7640         Fix_val=1
7650         GOTO Data_edit
7660     END IF
7670     Chest(Io_ptr)=FNLval(Io_msg$(Chart_num,5))
7680     IF Chest(Io_ptr)=9999.999 THEN
10 7690         Ionum=5
7700         Fix_val=1
7710         GOTO Data_edit
7720     END IF
7730     In_tot(Io_ptr)=Iv_intake(Io_ptr)+Fluid_in(Io_ptr)
15 7740     Out_tot(Io_ptr)=Urine(Io_ptr)+Chest(Io_ptr)
7750     Net(Io_ptr)=In_tot(Io_ptr)-Out_tot(Io_ptr)
7760     CALL Chart(Chart_num)
7770     Io_ptr=Io_ptr+1
7780     Io_in=1
20 7790     Fix_val=0
7800     GOTO Keyend
7810     !
7820     !
7830     Key2:Chart_num=2
25     !...ventilation charting
7840     GRAPHICS OFF
7850     PRINT CHR$(12)
7860     IF Next_time-TIMEDATE<45 THEN
7870         DISP "not enough time to enter data; wait for
30         next xfer"
7880         WAIT 2
7890         GOTO Keyend
7900     END IF
7910     Num_var=13
35 7920     IF Lab_in=1 THEN
7930         DISP "data in for this xfer; chart displayed"
```

- 103 -

```
7940      WAIT 2
7950      Lab_ptr=Lab_ptr-1
7960      CALL Chart(Char_num)
7970      Lab_ptr=Lab_ptr+1
5 7980      GOTO Keyend
7990 -ELSE
8000      INPUT "Input values=1 or display
      chart=2?",Inp
8010      IF Inp=1 THEN
10 8020          IF Lab_ptr>7 THEN
8030              DISP "Do not enter more Lab data;
              disc full"
8040              WAIT 3
8050              GOTO Keyend
15 8060          ELSE
8070              GOTO I_o
8080          END IF
8090      ELSE
8100          CALL Chart(Char_num)
20 8110          GOTO Keyend
8120      END IF
8130  END IF
8140 Data2:
8150      Lab_time$(Lab_ptr)=Io_msg$(Char_num,1)
25 8160      Na(Lab_ptr)=FNLval(Io_msg$(Char_num,2))
8170      IF Na(Lab_ptr)=9999.999 THEN
8180          Ionum=2
8190          Fix_val=1
8200          GOTO Data_edit
30 8210      END IF
8220      K1(Lab_ptr)=FNLval(Io_msg$(Char_num,3))
8230      IF K1(Lab_ptr)=9999.999 THEN
8240          Ionum=3
8250          Fix_val=1
35 8260          GOTO Data_edit
8270      END IF
```

- 104 -

```

      8280      Cl(Lab_ptr)=FNLval(Io_msg$(Chart_num,4))
      8290      IF Cl(Lab_ptr)=9999.999 THEN
      8300          Ionum=4
      8310          Fix_val=1
5      8320          GOTO Data_edit
      8330      END IF
      8340      Hco3(Lab_ptr)=FNLval(Io_msg$(Chart_num,5))
      8350      IF Hco3(Lab_ptr)=9999.999 THEN
      8360          Ionum=5
10     8370          Fix_val=1
      8380          GOTO Data_edit
      8390      END IF
      8400      Ca(Lab_ptr)=FNLval(Io_msg$(Chart_num,6))
      8410      IF Ca(Lab_ptr)=9999.999 THEN
15     8420          Ionum=6
      8430          Fix_val=1
      8440          GOTO Data_edit
      8450      END IF
      8460      Hct(Lab_ptr)=FNLval(Io_msg$(Chart_num,7))
20     8470      IF Hct(Lab_ptr)=9999.999 THEN
      8480          Ionum=7
      8490          Fix_val=1
      8500          GOTO Data_edit
      8510      END IF
25     8520      Gluc(Lab_ptr)=FNLval(Io_msg$(Chart_num,8))
      8530      IF Gluc(Lab_ptr)=9999.999 THEN
      8540          Ionum=8
      8550          Fix_val=1
      8560          GOTO Data_edit
30     8570      END IF
      8580      Dig(Lab_ptr)=FNLval(Io_msg$(Chart_num,9))
      8590      IF Dig(Lab_ptr)=9999.999 THEN
      8600          Ionum=9
      8610          Fix_val=1
35     8620          GOTO Data_edit
      8630      END IF
```

- 105 -

```
8640 Pt(Lab_ptr)=FNLval(Io_msg$(Chart_num,10))
8650 IF Pt(Lab_ptr)=9999.999 THEN
8660     Ionum=10
8670     Fix_val=1
5 8680     GOTO Data_edit
8690 END IF
8700 Ptt(Lab_ptr)=FNLval(Io_msg$(Chart_num,11))
8710 IF Ptt(Lab_ptr)=9999.999 THEN
8720     Ionum=11
10 8730     Fix_val=1
8740     GOTO Data_edit
8750 END IF
8760 Creat(Lab_ptr)=FNLval(Io_msg$(Chart_num,12))
8770 IF Creat(Lab_ptr)=9999.999 THEN
15 8780     Ionum=12
8790     Fix_val=1
8800     GOTO Data_edit
8810 END IF
8820 Bun(Lab_ptr)=FNLval(Io_msg$(Chart_num,13))
20 8830 IF Bun(Lab_ptr)=9999.999 THEN
8840     Ionum=13
8850     Fix_val=1
8860     GOTO Data_edit
8870 END IF
25 8880 CALL Chart(Chart_num)
8890 Lab_ptr=Lab_ptr+1
8900 Lab_in=1
8910 Fix_val=0
8920 GOTO Keyend
30 8930 !
8940 !
8950 Key3:Chart_num=4
      !...hemodynamic graphics
8960 IF Next_time-TIMEDATE<45 THEN
35 8970     DISP "not enough time to enter data; wait for
      next xfer"
```

- 106 -

```
      8980      WAIT 2
      8990      GOTO Keyend
      9000      END IF

5
      9010      GRAPHICS OFF
      9020      PRINT CHR$(12)
      9030      INPUT "Blood pressures(1) or cardiac
      indices(2)?",Bp
10     9040      IF Bp=1 THEN
      9050          Num_var=9
      9060      ELSE
      9070          Fst=10
      9080          Num_var=13
15     9090      END IF
      9100      IF Pres_in=1 AND Bp=1 THEN
      9110          DISP "data in for this xfer; chart displayed"
      9120          WAIT 2
      9130          Pres_ptr=Pres_ptr-1
20     9140          IF Heart_in=1 THEN Heart_ptr=Heart_ptr-1
      9150          CALL Chart(Char_num)
      9160          IF Heart_in=1 THEN Heart_ptr=Heart_ptr+1
      9170          Pres_ptr=Pres_ptr+1
      9180          GOTO Keyend
25     9190      ELSE
      9200          IF Heart_in=1 AND Bp=2 THEN
      9210              DISP "data in for this xfer; chart
              displayed"
      9220              WAIT 2
30     9230              IF Pres_in=1 THEN Pres_ptr=Pres_ptr-1
      9240              Heart_ptr=Heart_ptr-1
      9250              CALL Chart(Char_num)
      9260              Heart_ptr=Heart_ptr+1
      9270              IF Pres_in=1 THEN Pres_ptr=Pres_ptr-1
35     9280              GOTO Keyend
      9290      ELSE
```



- 107 -

```

9300      INPUT "Input values=1 or display
          chart=2?",Inp
9310      IF Inp=1 THEN
9320          IF Bp=1 AND Pres_ptr>17 THEN
5 9330          DISP "Do not enter more Pressure
          data; disc full"
9340          WAIT 3
9350          GOTO Keyend
9360      ELSE
10 9370          GOTO I_o
9380      END IF
9390      ELSE
9400          IF Heart_in=1 THEN Heart_ptr=Heart_
          ptr-1
15 9410          IF Pres_in=1 THEN Pres_ptr=Pres_ptr-1
9420          CALL Chart(Char_num)
9430          IF Heart_in=1 THEN Heart_ptr=Heart_
          ptr+1
9440          IF Pres_in=1 THEN Pres_ptr=Pres_ptr+1
20 9450          GOTO Keyend
9460      END IF
9470      END IF
9480      END IF
9490      Data4: !
25 9500      IF Bp=1 THEN
9510          Pres_time$(Pres_ptr)=Io_msg$(Chart_num,1)
9520          Ao_s(Pres_ptr)=FNLval(Io_msg$(Chart_num,2))
9530          IF Ao_s(Pres_ptr)=9999.999 THEN
9540              Ionum=2
30 9550              Fix_val=1
9560              GOTO Data_edit
9570          END IF
9580          Ao_d(Pres_ptr)=FNLval(Io_msg$(Chart_num,3))
9590          IF Ao_d(Pres_ptr)=9999.999 THEN
35 9600              Ionum=3
9610              Fix_val=1
```

- 108 -

```

9620          GOTO Data_edit
9630      END IF
9640      Ao_m(Pres_ptr)=FNLval(Io_msg$(Chart_num,4))
9650      IF Ao_m(Pres_ptr)=9999.999 THEN
5  9660          Ionum=4
9670          Fix_val=1
9680          GOTO Data_edit
9690      END IF
9700      Pa_s(Pres_ptr)=FNLval(Io_msg$(Chart_num,5))
10 9710      IF Pa_s(Pres_ptr)=9999.999 THEN
9720          Ionum=5
9730          Fix_val=1
9740          GOTO Data_edit
9750      END IF
15 9760      Pa_d(Pres_ptr)=FNLval(Io_msg$(Chart_num,6))
9770      IF Pa_d(Pres_ptr)=9999.999 THEN
9780          Ionum=6
9790          Fix_val=1
9800          GOTO Data_edit
20 9810      END IF
9820      Pa_m(Pres_ptr)=FNLval(Io_msg$(Chart_num,7))
9830      IF Pa_m(Pres_ptr)=9999.999 THEN
9840          Ionum=7
9850          Fix_val=1
25 9860          GOTO Data_edit
9870      END IF
9880      La_m(Pres_ptr)=FNLval(Io_msg$(Chart_num,8))
9890      IF La_m(Pres_ptr)=9999.999 THEN
9900          Ionum=8
30 9910          Fix_val=1
9920          GOTO Data_edit
9930      END IF
9940      Ra_m(Pres_ptr)=FNLval(Io_msg$(Chart_num,9))
9950      IF Ra_m(Pres_ptr)=9999.999 THEN
35 9960          Ionum=9
9970          Fix_val=1
```

- 109 -

```

9980          GOTO Data_edit
9990      END IF
10000      IF Heart_in=1 THEN Heart_ptr=Heart_ptr-1
10010      CALL Chart(Chart_num)
5  10020      IF Heart_in=1 THEN Heart_ptr=Heart_ptr+1
10030      Pres_ptr=Pres_ptr+1
10040      Pres_in=1
10050      Fix_val=0
10060      GOTO Keyend
10  10070  ELSE
10080      Heart_time$(Heart_ptr)=Io_msg$(Chart_num,10)
10090      Ci(Heart_ptr)=FNLval(Io_msg$(Chart_num,11))
10100      IF Ci(Heart_ptr)=9999.999 THEN
10110          Ionum=11
15  10120          Fix_val=1
10130          GOTO Data_edit
10140      END IF
10150      Pvri(Heart_ptr)=FNLval(Io_msg$(Chart_num,12))
10160      IF Pvri(Heart_ptr)=9999.999 THEN
20  10170          Ionum=12
10180          Fix_val=1
10190          GOTO Data_edit
10200      END IF
10210      Svri(Heart_ptr)=FNLval(Io_msg$(Chart_num,13))
25  10220      IF Svri(Heart_ptr)=9999.999 THEN
10230          Ionum=13
10240          Fix_val=1
10250          GOTO Data_edit
10260      END IF
30  10270      IF Pres_in=1 THEN Pres_ptr=Pres_ptr-1
10280      CALL Chart(Chart_num)
10290      IF Pres_in=1 THEN Pres_ptr=Pres_ptr+1
10300      Heart_ptr=Heart_ptr+1
10310      Heart_in=1
35  10320      Fst=1
10330      Fix_val=0
```

- 110 -

```
10340  END IF
10350  GOTO Keyend
10360  !
10370  !
5  10380  Key4:Key_id=4
10390  IF Trend_dp=0 THEN
10400      ALLOCATE INTEGER Spectra(7499)
10410      GSTORE Spectra(*)
10420      Trend_dp=2
10  10430      Top1=200
10440      Top2=2.5
10450      Bot2=-2.5
10460      Top3=10
10470      Top4=10
15  10480      CALL Trend_graph
10490  ELSE
10500      IF Trend_dp=2 THEN
10510          GRAPHICS ON
10520          GLOAD Spectra(*)
20  10530          DEALLOCATE Spectra(*)
10540          CALL Offgraph
10550          Trend_dp=0
10560      ELSE
10570          Trend_dp=2
25  10580          Top1=200
10590          Top2=2.5
10600          Bot2=-2.5
10610          Top3=10
10620          Top4=10
30  10630          CALL Trend_graph
10640      END IF
10650  END IF
10660  GOTO Keyend
10670  !
35  10680  !
10690  Key5:Key_id=5
```

- 111 -

```
      !...display message file
10700 IF Msg_dp_request<2 THEN
10710     DISP "messages will be recalled soon"
10720     Msg_dp_request=1
5  10730     WAIT 1
10740 ELSE
10750     Msg_dp_request=3
10760 END IF
10770 GOTO Keyend
10 10780 !
10790 !
10800 Key6:Key_id=6                !..premature program
                                   termination
10810 DISP "To halt program hit KEY 6 again (within 10
15      sec)"
10820 ON TIME (TIMEDATE+10) MOD 86400,4 GOTO Keyend
10830 ON KEY 6,3 GOTO Halter
10840 Cancel_wait:GOTO Cancel_wait
10850 Halter:Num_xfer_left=1
20 10860 Halt_pg=1
10870 GOTO Key_msg
10880 !
10890 !
10900 Key7:Chart_num=3
25 10910 IF Next_time-TIMEDATE<45 THEN
10920     DISP "not enough time to enter data; wait for
           next xfer"
10930     WAIT 2
10940     GOTO Keyend
30 10950 END IF
10960 GRAPHICS OFF
10970 PRINT CHR$(12)
10980 Num_var=13
10990 IF Vent_in=1 THEN
35 11000     DISP "data in for this xfer; chart displayed"
11010     WAIT 2
```

- 112 -

```
11020      Vent_ptr=Vent_ptr-1
11030      CALL Chart(Char_num)
11040      Vent_ptr=Vent_ptr+1
11050      GOTO Keyend
5 11060  ELSE
11070      INPUT "Input values=1 or display
chart=2?",Inp
11080      IF Inp=1 THEN
11090          IF Vent_ptr>7 THEN
10 11100              DISP "Do not enter more Vent data;
disc full"
11110              WAIT 3
11120              GOTO Keyend
11130          ELSE
15 11140              GOTO I_o
11150          END IF
11160      ELSE
11170          CALL Chart(Char_num)
11180          GOTO Keyend
20 11190      END IF
11200  END IF
11210 Data3:
11220  Vent_time$(Vent_ptr)=Io_msg$(Chart_num,1)
11230  Rate(Vent_ptr)=FNLval(Io_msg$(Chart_num,2))
25 11240  IF Rate(Vent_ptr)=9999.999 THEN
11250      Ionum=2
11260      Fix_val=1
11270      GOTO Data_edit
11280  END IF
30 11290  Fio2(Vent_ptr)=FNLval(Io_msg$(Chart_num,3))
11300  IF Fio2(Vent_ptr)=9999.999 THEN
11310      Ionum=3
11320      Fix_val=1
11330      GOTO Data_edit
35 11340  END IF
11350  Pp(Vent_ptr)=FNLval(Io_msg$(Chart_num,4))
```

- 113 -

```
11360 IF Pp(Vent_ptr)=9999.999 THEN
11370     Ionum=4
11380     Fix_val=1
11390     GOTO Data_edit
5 11400 END IF
11410 Peep(Vent_ptr)=FNLval(Io_msg$(Chart_num,5))
11420 IF Peep(Vent_ptr)=9999.999 THEN
11430     Ionum=5
11440     Fix_val=1
10 11450     GOTO Data_edit
11460 END IF
11470 Tv(Vent_ptr)=FNLval(Io_msg$(Chart_num,6))
11480 IF Tv(Vent_ptr)=9999.999 THEN
11490     Ionum=6
15 11500     Fix_val=1
11510     GOTO Data_edit
11520 END IF
11530 Ie_ratio$(Vent_ptr)=Io_msg$(Chart_num,7)
11540 Airp(Vent_ptr)=FNLval(Io_msg$(Chart_num,8))
20 11550 IF Airp(Vent_ptr)=9999.999 THEN
11560     Ionum=8
11570     Fix_val=1
11580     GOTO Data_edit
11590 END IF
25 11600 Ph(Vent_ptr)=FNLval(Io_msg$(Chart_num,9))
11610 IF Ph(Vent_ptr)=9999.999 THEN
11620     Ionum=9
11630     Fix_val=1
11640     GOTO Data_edit
30 11650 END IF
11660 Po2(Vent_ptr)=FNLval(Io_msg$(Chart_num,10))
11670 IF Po2(Vent_ptr)=9999.999 THEN
11680     Ionum=10
11690     Fix_val=1
35 11700     GOTO Data_edit
11710 END IF
```

- 114 -

```
11720 Pco2(Vent_ptr)=FNLval(Io_msg$(Chart_num,11))
11730 IF Pco2(Vent_ptr)=9999.999 THEN
11740     Ionum=11
11750     Fix_val=1
5 11760     GOTO Data_edit
11770 END IF
11780 Bgo3(Vent_ptr)=FNLval(Io_msg$(Chart_num,12))
11790 IF Bgo3(Vent_ptr)=9999.999 THEN
11800     Ionum=12
10 11810     Fix_val=1
11820     GOTO Data_edit
11830 END IF
11840 Be(Vent_ptr)=FNLval(Io_msg$(Chart_num,13))
11850 IF Be(Vent_ptr)=9999.999 THEN
15 11860     Ionum=13
11870     Fix_val=1
11880     GOTO Data_edit
11890 END IF
11900 CALL Chart(Chart_num)
20 11910 Vent_ptr=Vent_ptr+1
11920 Vent_in=1
11930 Fix_val=0
11940 GOTO Keyend
11950 !
25 11960 !
11970 Key8:Chart_num=5
11980 IF Next_time-TIMEDATE<45 THEN
11990     DISP "not enough time to enter data; wait for
        next xfer"
30 12000     WAIT 2
12010     GOTO Keyend
12020 END IF
12030 GRAPHICS OFF
12040 PRINT CHR$(12)
35 12050 Num_var=3
12060 IF Drug_in=1 THEN
```



- 115 -

```
12070      DISP "data in for this xfer; chart displayed"
12080      WAIT 2
12090      Drug_ptr=Drug_ptr-1
12100      CALL Chart(Char_num)
5 12110      Drug_ptr=Drug_ptr+1
12120      GOTO Keyend
12130  ELSE
12140      INPUT "Input values=1 or display
      chart=2?",Inp
10 12150      IF Inp=1 THEN
12160          IF Drug_ptr>38 THEN
12170              DISP "Do not enter more Drug data;
              disc full"
12180              WAIT 3
15 12190              GOTO Keyend
12200          ELSE
12210              GOTO I_o
12220          END IF
12230      ELSE
20 12240          CALL Chart(Char_num)
12250          GOTO Keyend
12260      END IF
12270  END IF
12280  Data5: !
25 12290  Drug_time$(Drug_ptr)=Io_msg$(Char_num,3)
12300  Drug_name$(Drug_ptr)=Io_msg$(Char_num,1)
12310  Drug_dos$(Drug_ptr)=Io_msg$(Char_num,2)
12320  CALL Chart(Char_num)
12330  Drug_ptr=Drug_ptr+1
30 12340  Drug_in=1
12350  GOTO Keyend
12360  !
12370  !
12380  Key9:Key_id=9
35 12390  Bp_graph: !
12400  IF Next_time-TIMEDATE<12 THEN GOTO Waiter1
```

- 116 -

```
12410 IF Trend_dp=0 THEN
12420     Trend_dp=1
12430     Top1=150
12440     Top2=75
5 12450     Bot2=0
12460     Top3=50
12470     Top4=50
12480     ALLOCATE INTEGER Spectra(7499)
12490     GSTORE Spectra(*)
10 12500     CALL Trend_graph
12510 ELSE
12520     IF Trend_dp=1 THEN
12530         GRAPHICS ON
12540         GLOAD Spectra(*)
15 12550         DEALLOCATE Spectra(*)
12560         CALL Offgraph
12570         Trend_dp=0
12580     ELSE
12590         Trend_dp=1
20 12600         Top1=150
12610         Top2=75
12620         Bot2=0
12630         Top3=50
12640         Top4=50
25 12650         CALL Trend_graph
12660     END IF
12670 END IF
12680 GOTO Keyend
12690 !
30 12700 !
12710 I_o:1
12720 IF TIMEDATE>Next_time-20 THEN
12730     DISP "not enough time to enter data; wait for
        next xfer"
35 12740     WAIT 2
12750     GOTO Keyend
```

- 117 -

```
12760 END IF
12770 PRINT TABXY(1,1);"enter values"
12780 FOR I=Fst TO Num_var
12790     PRINT TABXY(1,17);"          "
5  12800     PRINT TABXY(1,17);Io$(Chart_num,I)
12810     Edit_msg$=""
12820     CALL Editor
12830     Io_msg$(Chart_num,I)=Edit_msg$
12840     PRINT TABXY(1,I+2);Io$(Chart_num,I);"=";Io_
10      msg$(Chart_num,I)
12850 NEXT I
12860 PRINT TABXY(1,17);"          "
12870 PRINT TABXY(1,18);"          "
12880!
15  12890!....editting the data
12900!
12910 Io_fix:DISP "Do you want to edit I/O
values?              (Y/N)"
12920 ENTER 2;Ans$
20  12930 DISP "          "
12940 IF Ans$="Y" OR Ans$="y" THEN
12950     IF TIMEDATE>Next_time-15 THEN
12960         DISP "not enough time; data not stored;
retry next xfer"
25  12970         GOTO Keyend
12980     END IF
12990     ON Chart_num GOTO Value,Lab,Vent,Pres,Drug
13000 Value:DISP "which value? 1=time, 2=maint. fluid,
3=other fluids, 4=urine, 5=chest"
30  13010     ENTER 2;Ionum
13020     IF Ionum<1 OR Ionum>5 THEN GOTO Value
13030     GOTO Data_edit
13040 Lab: DISP "which value?
1=time, 2=Na, 3=K, 4=Cl, 5=HCO3, 6=Ca, 7=Hct, 8=Gluc, 9=Di
35      g, 10=PT, 11=PTT, 12=Creat, 13=Bun"
13050     ENTER 2;Ionum
```

- 118 -

```
13060      IF Ionum<1 OR Ionum>13 THEN GOTO Lab
13070      GOTO Data_edit
13080 Vent:PRINT TABXY(1,17);"which value?
           1=time,2=rate,3=FIO2,4=PP,5=peep,6=TV,
5          7=I:E,8=airway"
13090      PRINT TABXY(1,18);
           "9=ph,10=pO2,11=pCO2,12=HCO3,13=Be"
13100      ENTER 2;Ionum
13110      IF Ionum<1 OR Ionum>13 THEN GOTO Vent
10 13120      GOTO Data_edit
13130 Pres:IF Bp=1 THEN
13140          PRINT TABXY(1,17);"which value? 1=pres
           time,2=ao/s,3=ao/d,4=ao/m,
           5=pa/s,6=pa/d,7=pa/m,8=1a,9=ra"
15 13150      ELSE
13160          PRINT TABXY(1,18);"which value? 10=heart
           time,11=c.i.,12=pvri,13=svri"
13170      END IF
13180      ENTER 2;Ionum
20 13190      IF Ionum<1 OR Ionum>13 THEN GOTO Pres
13200      GOTO Data_edit
13210 Drug:DISP "which value? 1=name,2=dosage,3=time"
13220      ENTER 2;Ionum
13230      IF Ionum<1 OR Ionum>10 THEN GOTO Drug
25 13240      GOTO Data_edit
13250 Data_edit:
13260      IF TIMEDATE>Next_time-15 THEN
13270          DISP "not enough time; data not stored;
           retry next xfer"
30 13280      WAIT 2
13290      GOTO Keyend
13300      END IF
13310      C_num=Chart_num
13320      R_num=2
35 13330      IF Fix_val=1 THEN
13340          PRINT TABXY(1,17);"Error on input; enter
```

- 119 -

```

                                value again"
13350      PRINT TABXY(1,18);Io$(C_num,Ionum)
13360      END IF
13370      PRINT TABXY(1,18);Io_msg$(C_num,Ionum)
5  13380      Edit_msg$=Io_msg$(C_num,Ionum)
13390      CALL Editor
13400      Io_msg$(C_num,Ionum)=Edit_msg$
13410      PRINT TABXY(1,Ionum+R_num);"      "
13420      PRINT TABXY(1,Ionum+R_num);Io$(C_
10      num,Ionum);"=";Edit_msg$
13430      PRINT TABXY(1,17);"      "
13440      PRINT TABXY(1,18);"      "
13460      GOTO Io_fix
13470  ELSE
15  13480      ON Chart_num GOTO
          Data1,Data2,Data3,Data4,Data5
13490  END IF
13500  Keyend:OFF TIME
13510  OFF KBD
20  13520  RETURN
13530  END
13540  !
13550  !
13560  !
25  13570  !
13580  !
13590  SUB Pauser
13600      DISP "press CONTINUE to continue"
13610      PAUSE
30  13620      DISP
13630  SUBEND
13640  !
13650  !
13660  !
35  13670  !
13680  !
```

- 120 -

```

13690 SUB Get_param
13700 COM /Multi_param/ Start_chan,Stop_chan,Pacing_
bits,Pacing_rate,Num_pt
s,Num_xfer,Num_xfer_left,Name_len,Scr_
5 file$(28),Scr_
file2$(28)
13710 COM /Messagecom/ Message$(10)[80],@Messages
13720 COM /Trends/ Mean_hr_t(*),Lfa_t(*),Rfa_
t(*),Ratio_t(*),T_ptr,Time_now
10 l,Meas_resp_t(*),Trend_dp
13730 COM /Vitaldata/ Rfa,Lfa,Peakratio,Meas_
resp,Next_time
13740 COM /Pressure/
Top1,Top2,Top3,Top4,Bot1,Bot2,Bot3,Bot4
15 13750 COM /Pres_chart/ Pres_time$(*),Ao_s(*),Ao_
d(*),Ao_m(*),Pa_s(*),Pa_d(*
),Pa_m(*),La_m(*),Ra_m(*),Pres_ptr,Pres_in
13760 COM /Subject/ Sub_name$(25),Hos_num$(15),Id_
age$(10),Id_wt$(10),Id_ht
20 $(10),Diag$(30),Opera$(45),Halt_pg
13770 COM /Io_chart/ Io_time$(*),Iv_intake(*),Fluid_
in(*),In_tot(*),Urine(*
),Chest(*),Out_tot(*),Net(*),Io_ptr
13780 COM /Lab_chart/ Lab_
25 time$(*),Na(*),Kl(*),Cl(*),Hco3(*),
Ca(*),Hct(*),Gluc(*),Dig(*),Pt(*),
Ptt(*),Creat(*),Bun(*),Lab_ptr
13790 COM /Vent_chart/ Vent_
time$(*),Rate(*),Fio2(*),Pp(*),Peep(*),Tv(*),Ie
30 _ratio$(*),Airlp(*),Ph(*),Po2(*),Pco2(*),
Bgo3(*),Be(*),Vent_ptr
13800 COM /Heart_index/ Heart_
time$(*),Ci(*),Pvri(*),Svri(*),Heart_ptr
13810 COM /Drugs/ Drug_time$(*),Drug_name$(*),Drug_
35 dos$(*),Drug_ptr
13820 DIM Mo$(24)

```

- 121 -

```
13830      Mo$="JAFBMRAPMYJNJLAUSPOCNODC"
13840      !   INTEGER Id_buffer(255) BUFFER
13850      Disk_name$=":HP8290X,700,1"
13860      IF Halt_pg=1 THEN GOTO Purger_get!.....quit
5         program
13870      !
13880      ! change soft key messages
13890      !
13900      Oldmsg:PRINT CHR$(12)
10      13910      PRINT "These are the current soft key
           messages:"
13920      FOR I=0 TO 9
13930          PRINT "KEY";I;":";Message$(I)
13940      NEXT I
15      14100      DISP "Press cont when ready to continue"
14110      PAUSE
14120!
14130      INPUT "Enter subject name, 10 chars (Doe if
           unknown)",Sub_name$
20      14140      Sub_name$=Sub_name$[1,10]
14150      INPUT "Enter hospital number, 8 chars (00 if
           unknown):",Hos_num$
14160      Hos_num$=Hos_num$[1,8]
14170      INPUT "Enter subject age(00 if unknown):",Id_
25      age$
14180      INPUT "Enter subject weight,kg (00 if
           unknown):",Id_wt$
14190      INPUT "Enter subject height,cm (00 if
           unknown):",Id_ht$
30      14200      INPUT "Enter diagnosis, 10 chars (Unk if
           unknown):",Diag$
14210      Diag$=Diag$[1,10]
14220      INPUT "Enter operation, 15 chars (Unk if
           unknown):",Opera$
35      14230      Opera$=Opera$[1,15]
14240!
```

- 122 -

```
14250 Ch_sel:!  
14260      Start_chan=0  
14270      Stop_chan=0  
14280 !  
5 14290      Pacing_bits=0  
14300 Pacing_sel:!  
14310      Base$="M"  
14320      Pacing_bits=261  
14330 !  
10 14340      Base$=Base$&"SEC"  
14350 !  
14360 !  
14370 ! FINDOUT BLOCKSIZE FOR DATA TRANSFER  
14380 !  
15 14390      Num_xfer=55  
14400!  
14410! since new data is to be taken, zero the trend  
      graphs (120 pts=8hrs)  
14420!  
20 14430      MAT Mean_hr_t= (0)  
14440      MAT Rfa_t= (0)  
14450      MAT Lfa_t= (0)  
14460      MAT Ratio_t= (0)  
14470      MAT Meas_resp_t= (0)  
25 14471      MAT Trans_time= (0)  
14480      T_ptr=0  
14490      MAT Pres_time$= ("")  
14500      MAT Ao_s= (0)  
14510      MAT Ao_d= (0)  
30 14520      MAT Ao_m= (0)  
14530      MAT Pa_s= (0)  
14540      MAT Pa_d= (0)  
14550      MAT Pa_m= (0)  
14560      MAT La_m= (0)  
35 14570      MAT Ra_m= (0)  
14580      MAT Io_time$= ("")
```



- 123 -

	14590	MAT Iv_intake= (0)
	14600	MAT Fluid_in= (0)
	14610	MAT In_tot= (0)
	14620	MAT Urine= (0)
5	14630	MAT Chest= (0)
	14640	MAT Out_tot= (0)
	14650	MAT Net= (0)
	14660	MAT Lab_time\$= ("")
	14670	MAT Na= (0)
10	14680	MAT K1= (0)
	14690	MAT Cl= (0)
	14700	MAT Hco3= (0)
	14710	MAT Ca= (0)
	14720	MAT Hct= (0)
15	14730	MAT Gluc= (0)
	14740	MAT Dig= (0)
	14750	MAT Pt= (0)
	14760	MAT Ptt= (0)
	14770	MAT Creat= (0)
20	14780	MAT Bun= (0)
	14790	MAT Vent_time\$= ("")
	14800	MAT Rate= (0)
	14810	MAT Fio2= (0)
	14820	MAT Pp= (0)
25	14830	MAT Peep= (0)
	14840	MAT Tv= (0)
	14850	MAT Ie_ratio\$= ("")
	14860	MAT Airp= (0)
	14870	MAT Ph= (0)
30	14880	MAT Po2= (0)
	14890	MAT Pco2= (0)
	14900	MAT Bgo3= (0)
	14910	MAT Be= (0)
	14920	MAT Heart_time\$= ("")
35	14930	MAT Ci= (0)
	14940	MAT Pvri= (0)

- 124 -

```
14950      MAT Svri= (0)
14960      MAT Drug_time$= ("")
14970      MAT Drug_name$= ("")
14980      MAT Drug_dos$= ("")
5  14990      Pres_ptr=0
15000      Trend_ptr=0

15010      Ratio_t(0)=1 !..prevent trend graph errors on
10      startup
15020      Rfa=0
15030      Lfa=0
15040      Meas_resp=0
15050      Peakratio=1
15 15060      !
15070      !
15080      Pacing_rate=250
15090      Num_pts=1024*Num_xfer
15100      Num_header=256+8*Num_xfer
20 15110      IF Scr_file$="?" THEN GOTO Skipl
15120 Purger_get:DISP "PURGE FILE?"
15130      ENTER 2;Resp$
15140      IF Resp$="Y" OR Resp$="YES" THEN
15150          PURGE Scr_file$
25 15160          PURGE Scr_file2$
15170          PURGE "messglog:HP8290X,700,1"
15180          PURGE "temp_trend:HP8290X,700,1"
15190          PURGE "hemo_data:HP8290X,700,1"
15200          PURGE "io_data:HP8290X,700,1"
30 15210          PURGE "drug_data:HP8290X,700,1"
15220          PURGE "lab_data:HP8290X,700,1"
15230          PURGE "co_data:HP8290X,700,1"
15231          PURGE "sub_data:HP8290X,700,1"
15240      ELSE
35 15250!
15260! the data files are named according to the date
```

- 125 -

```

15270! in the following format:
15280!      xxxxxmddy
15290! where
15300!      xxxx - resp,hr__,msgs,errs,trnd
5  15310!      dd  - day
    15320!      mm  - month
                        (JA,FB,MR,AP,MY,JN,JL,AU,SP,OC,NO,DC)
15330!      yy   - year
15340      Date_now$=DATE$(TIMEDATE)
10  15350      Month_now=FNMonth(Date_now$)*2-1
    15360      Mm$=Mo$(Month_now;2]
    15370      Id_field$=Date_now$[1;2]&Mm$&Date_
                        now$[10;2]
15380! new name for respiratory file: respddmmyy
15  15390      RENAME Scr_file$ TO "resp"&Id_field$&Disk_
                        name$
15400! new name for heart rate file: hr__ddmmyy
15410      RENAME Scr_file2$ TO "hr__"&Id_
                        field$&Disk_name$
20  15420! new name for message log: msgsdmmyy
    15430      RENAME "messglog:HP8290X,700,1" TO
                        "msgs"&Id_field$&Disk_name$
15440! new name for hemo data: dataddmmyy
15450      RENAME "hemo_data:HP8290X,700,1" TO
25  "hemo"&Id_field$&Disk_name$
15460! new name for io data
15470      RENAME "io_data:HP8290X,700,1" TO "io__
                        "&Id_field$&Disk_name$
15480! new name for lab data
30  15490      RENAME "lab_data:HP8290X,700,1" TO "lab_
                        "&Id_field$&Disk_name$
15500! new name for vent data
15510      RENAME "vent_data:HP8290X,700,1" TO
                        "vent"&Id_field$&Disk_name$
35  15520! new name for co data
    15530      RENAME "co_data:HP8290X,700,1" TO "co__

```

- 126 -

```

                                "&Id_field$&Disk_name$
15540! new name for drug data
15550                RENAME "drug_data:HP8290X,700,1" TO
                                "drug"&Id_field$&Disk_name$
5  15551! new name for subject data
15552                RENAME "sub_data:HP8290X,700,1" TO "sub_
                                "&Id_field$&Disk_name$
15560! name for trend summary file: trndddmmyy
15570                PURGE "temp_trend:HP8290X,700,1"
10 15580                CREATE BDAT "trnd"&Id_field$&Disk_
                                name$,19,256
15590                ASSIGN @Trend_file TO "trnd"&Id_
                                field$&Disk_name$;FORMAT OFF
15600                OUTPUT @Trend_file;Mean_hr_t(*),Lfa_
15  t(*),Rfa_t(*),Ratio_t(*),Meas
                                _resp_t(*),Trans_time(*),T_ptr
15610                ASSIGN @Trend_file TO *
15620                END IF
15630                IF Halt_pg=1 THEN      !..terminate program
20 15640                DISP "PROGRAM COMPLETED"
15650                STOP
15660                END IF
15670 Skipl:DISP
15680                Scr_file$="AOK"&Disk_name$
25 15690                Num_rec=-INT(-(Num_pts+Num_header)/128.)
15700                Scr_file2$="hr"&Scr_file$
15710                CREATE BDAT Scr_file$,Num_rec,256
15720                CREATE BDAT Scr_file2$,Num_rec,256
15730                CREATE BDAT "messglog:HP8290X,700,1",20,640
30 15740                CREATE BDAT "temp_trend"&Disk_name$,19,256
15750                CREATE BDAT "hemo_data"&Disk_name$,10,256
15760                CREATE BDAT "io_data"&Disk_name$,10,256
15770                CREATE BDAT "lab_data"&Disk_name$,10,256
15780                CREATE BDAT "vent_data"&Disk_name$,10,256
35 15790                CREATE BDAT "co_data"&Disk_name$,10,256
15800                CREATE BDAT "drug_data"&Disk_name$,10,256

```

- 127 -

```
15801      CREATE BDAT "sub_data"&Disk_name$,1,256
15802      ASSIGN @Sub_data TO "sub_data"&Disk_
           name$;FORMAT OFF
15803      OUTPUT @Sub_data;Sub_name$,Hos_num$,Id_
5         age$,Id_wt$,Id_ht$,Diag$,Opera$
15804      ASSIGN @Sub_data TO *
15810      Halt_pg=0
15820      Num_pts=1024
15830      PRINT Num_pts*Num_xfer;"points will be
10         transferred in";Num_xfer;"bloc
           ks of";Num_pts;"points"
15840      !
15850      Num_xfer_left=Num_xfer
15860      SUBEND
15870      !
15880      !
15890      !
15900      !
15910      DEF FNMonth(Date_now$)
20 15920      Month$=Date_now$[4;3]
15930      Month=0
15940      IF Month$="Jan" THEN Month=1
15950      IF Month$="Feb" THEN Month=2
15960      IF Month$="Mar" THEN Month=3
25 15970      IF Month$="Apr" THEN Month=4
15980      IF Month$="May" THEN Month=5
15990      IF Month$="Jun" THEN Month=6
16000      IF Month$="Jul" THEN Month=7
16010      IF Month$="Aug" THEN Month=8
30 16020      IF Month$="Sep" THEN Month=9
16030      IF Month$="Oct" THEN Month=10
16040      IF Month$="Nov" THEN Month=11
16050      IF Month$="Dec" THEN Month=12
16060      RETURN Month
35 16070      FNEND
16080!
```

- 128 -

```
16090!
16100!
16110!
16120!
5 16130 SUB Xfheader(@Disk,Num_bytes,File_id$)
16140     INTEGER Xheader(7) BUFFER
16150     Xheader(0)=(TIMEDATE MOD 86400)/60
16160     Xheader(1)=Num_bytes
16170     Xheader(2)=NUM(File_id$[1;1])
10 16180     Xheader(3)=0
16190     Xheader(4)=0
16200     Xheader(5)=0
16210     Xheader(6)=0
16220     Xheader(7)=0
15 16230     ASSIGN @Xheader TO BUFFER Xheader(*)
16240     CONTROL @Xheader,5;1 ! Reset empty pointer
    for buffer
16250     CONTROL @Xheader,4;16 ! Reset current number
    of bytes in buffer
20 16260     TRANSFER @Xheader TO @Disk;COUNT 16,WAIT
16270     ASSIGN @Xheader TO *
16280 SUBEND
16290!
16300!
25 16310!
16320!
16330!
16340!
16350 SUB Trend_graph
30 16360!
16370     COM /Trends/ Mean_hr_t(*),Lfa_t(*),Rfa_
    t(*),Ratio_t(*),T_ptr,Time_now
    1,Meas_resp_t(*),Trend_dp,Trans_time(*),Lfa_
    top,Rfa_top
35 16380     COM /Multi_param/ Start_chan,Stop_chan,Pacing_
    bits,Pacing_rate,Num_pt
```

- 129 -

```

s,Num_xfer,Num_xfer_left,Name_len,Scr_
  file$(28),Scr_
  file2$(28]
16390 COM /Pressure/
5   Top1,Top2,Top3,Top4,Bot1,Bot2,Bot3,Bot4
16400 COM /Pres_chart/ Pres_time$(*),Ao_s(*),Ao_
    d(*),Ao_m(*),Pa_s(*),Pa_d(*
      ),Pa_m(*),La_m(*),Ra_m(*),Pres_ptr,Pres_in
16410 DIM First_line(60),Sec_line(60),Third_
10  line(60),Fourth_line(60)
16420 IF Trend_dp=1 THEN
16430   MAT First_line= Ao_m
16440   MAT Sec_line= Pa_m
16450   MAT Third_line= La_m
15  16460   MAT Fourth_line= Ra_m
16470   G_right=INT((Num_xfer*256/60)/15)
16480   ! IF Pres_in=0 THEN ! Trend_ptr=Pres_
      ptr+1
16490   ! Trend_ptr=Pres_ptr+1
20  16500   ! ELSE
16510   Trend_ptr=Pres_ptr
16520   ! END IF
16530 ELSE
16540   MAT First_line= Mean_hr_t
25  16550   MAT Sec_line= Ratio_t
16560   MAT Third_line= Lfa_t
16570   MAT Fourth_line= Rfa_t
16580   G_right=Num_xfer
16590   Trend_ptr=T_ptr
30  16600 END IF
16610 Block_time=Pacing_rate*1.024/3600.
16620 GINIT
16630 GCLEAR
16640 PRINT CHR$(12)
35  16650 GRAPHICS ON
16660 Beg_time=Time_now1/3600-Block_time

```

- 130 -

```
16670      End_time=Beg_time+Num_xfer*Block_time
16680      Ibeg_time=INT(Beg_time)
16690      IF Ibeg_time<Beg_time THEN Ibeg_time=Ibeg_
        time+1
5  16700!
16710! label the time axes
16720!
16730      VIEWPORT 0,128,45,50
16740      WINDOW Beg_time,End_time,0,1
10 16750      IF INT(End_time)>Beg_time THEN
16760          LDIR 0
16770          FOR T_label=Ibeg_time TO INT(End_time)
16780              MOVE T_label,.5
16790              LONG 5
15 16800              CSIZE 4
16810              LABEL T_label
16820          NEXT T_label
16830      END IF
16840      VIEWPORT 0,128,40,45
20 16850      WINDOW 0,1,0,1
16860      MOVE .5,0
16870      LONG 4
16880      LABEL "Time (24 hr)"
16890!
25 16900! draw the axes
16910!
16920      VIEWPORT 0,128,50,100
16930      WINDOW Beg_time,End_time,0,1
16940      AXES 1/15.,.1,Beg_time,0
30 16950      WINDOW 1,0,1,0
16960      AXES 0,.25,0,0
16970!
16980! mean heart rate trends
16990!
35 17000      WINDOW -1,G_right,Bot1,Top1
17010      MOVE 0,First_line(0)
```



- 131 -

```
17020      FOR I=0 TO Trend_ptr-1
17030          DRAW I,First_line(I)
17040      NEXT I
17050!
5 17060! ratio trends (with a line at ratio=2)
17070!
17080      WINDOW -1,G_right,Bot2,Top2
17090      LINE TYPE 8,5
17100      IF Trend_dp=2 THEN
10 17110          MOVE 0,LGT(Sec_line(0))
17120      ELSE
17130          MOVE 0,Sec_line(0)
17140      END IF
17150      FOR I=0 TO Trend_ptr-1
15 17160          IF Trend_dp=2 THEN
17170              DRAW I,LGT(Sec_line(I))
17180          ELSE
17190              DRAW I,Sec_line(I)
17200          END IF
20 17210      NEXT I
17220      IF Trend_dp=2 THEN
17230          LINE TYPE 3,5!..sparsely dotted line at
              ratio=2
17240          MOVE 0,LGT(2.)
25 17250          DRAW Trend_ptr-1,LGT(2.)
17260      END IF
17270!
17280! lfa trends
17290!
30 17300      WINDOW -1,G_right,Bot3,Top3
17310      LINE TYPE 4,5
17320      MOVE 0,Third_line(0)
17330      FOR I=0 TO Trend_ptr-1
17340          DRAW I,Third_line(I)
35 17350      NEXT I
17360!
```

- 132 -

```
17370! rfa trends
17380!
17390      WINDOW -1,G_right,Bot4,Top4
17400      LINE TYPE 5,5
5  17410      MOVE 0,Fourth_line(0)
17420      FOR I=0 TO Trend_ptr-1
17430          DRAW I,Fourth_line(I)
17440      NEXT I
17450!
10 17460! draw a key for line types
17470!
17480      VIEWPORT 64,128,0,50
17490      WINDOW 0,1,0,13
17500      IF Trend_dp=2 THEN
15 17510          PRINT TABXY(1,17);"trend graph"
17520          PRINT TABXY(55,15);"mean hr(0-200)"
17530          PRINT TABXY(55,16);"ratio(.01-100)"
17540          PRINT TABXY(55,17);"lfa      (0-10)"
17550          PRINT TABXY(55,18);"rfa      (0-10)"
20 17560      ELSE
17570          PRINT TABXY(1,17);"mean pressure graphs"
17580          PRINT TABXY(50,15);"ao pressure(0-150)"
17590          PRINT TABXY(50,16);"pa pressure(0-75)"
17600          PRINT TABXY(50,17);"la pressure(0-50)"
25 17610          PRINT TABXY(50,18);"ra pressure(0-50)"
17620      END IF
17630      LINE TYPE 1,5
17640      MOVE .8,11
17650      DRAW 1.,11
30 17660      LINE TYPE 8,5
17670      MOVE .8,10
17680      DRAW 1.,10
17690      LINE TYPE 4,5
17700      MOVE .8,9
35 17710      DRAW 1.,9
17720      LINE TYPE 5,5
```

- 133 -

```
17730      MOVE .8,8
17740      DRAW 1.,8
17750  SUBEND
17760!
5  17770!
17780!
17790!
17800!
17810  SUB Msg_dump(Message_chart$(*),Message_line,Flg)
10 17820      COM /Messagecom/ Message$(10)[80],@Messages
17830      DIM Msg_buffer$(1280) BUFFER
17840      IF Flg>=2 THEN GOTO Chart_filled
17850      ASSIGN @Msg_buffer TO BUFFER Msg_
buffer$;FORMAT OFF
15 17860      STATUS @Messages,3;Num_rec
17870      STATUS @Messages,4;Rec_len
17880      STATUS @Messages,5;Cur_rec
17890      STATUS @Messages,6;Cur_byte
17900      IF Cur_rec<=1 AND Cur_byte<=1 THEN !.. no
20      messages yet
17910          Flg=0
17920          DISP "no messages yet"
17930          WAIT 2
17940          SUBEXIT
25 17950      END IF
17960      Flg=2
17970      CONTROL @Messages,5;1
17980      CONTROL @Messages,6;1
17990      FOR Rec=1 TO Cur_rec-1
30 18000  Read_msg:TRANSFER @Messages TO @Msg_buffer;COUNT
Rec_len,WAIT
18010      Message_chart$(Rec-1)=Msg_buffer$(1;Rec_
len]
18020      CONTROL @Msg_buffer,4;0
35 18030      CONTROL @Msg_buffer,5;1
18040      NEXT Rec
```

- 134 -

```
18050      IF Cur_byte>1 THEN
18060          TRANSFER @Messages TO @Msg_buffer;COUNT
              Cur_byte-1,WAIT
18070          Message_chart$(Cur_rec-1)=Msg_
5          buffer$[1;Cur_byte-1]
18080      END IF
18090      ASSIGN @Msg_buffer TO *
18100 Reset_msg_file:~
18110      CONTROL @Messages,5;Cur_rec
10 18120      CONTROL @Messages,6;Cur_byte
18130 Chart_filled:~
18140      STATUS @Messages,5;Cur_rec
18150      STATUS @Messages,6;Cur_byte
18160      Flg=2
15 18170      Cur_msg_ptr=0
18180      Chart_line=1
18190      Msg_buffer$=Message_chart$(0)
18200      Last_msg=Message_line+17
18210      Clear$=CHR$(255)&CHR$(75)
20 18220      OUTPUT 2;Clear$
18230      GRAPHICS OFF
18240 Next_msg:~
18250      Beg_msg=POS(Msg_buffer$[4],"Time")+3
18260      IF Beg_msg=3 THEN GOTO Next_chart_line
25 18270      Cur_msg_ptr=Cur_msg_ptr+1
18280      IF Cur_msg_ptr>Message_line THEN
18290          Tab_line=Cur_msg_ptr-Message_line
18300          PRINT TABXY(1,Tab_line);"      "
18310          PRINT TABXY(1,Tab_line);Msg_buffer$[1,Beg_
30          msg-1]
18320      END IF
18330      Msg_buffer$=Msg_buffer$[Beg_msg]
18340      IF Cur_msg_ptr=Last_msg THEN Subend_msg
18350      GOTO Next_msg
35 18360 Next_chart_line:IF Chart_line<Cur_rec THEN
18370          Msg_buffer$=Msg_buffer$&Message_
```

- 135 -

```

                                chart$(Chart_line)
18380      Chart_line=Chart_line+1
18390      GOTO Next_msg
18400      END IF
5  18410 Stopper:PRINT Msg_buffer$
18420 Subend_msg:PRINT
18430 SUBEND
18440 !
18450 !
10 18460 !
18470 !
18480 !
18490 SUB Disp_ctrls
18500     DISP "f̂ - freq range adjust (1 or 2 Hz)"
15 18510     WAIT 2
18520     DISP "ĥ - help: display these controls"
18530     WAIT 2
18540     DISP "p̂ - peak threshold adjust (+20%)"
18550     WAIT 2
20 18560     DISP "r̂ - resp time series display"
18570     WAIT 2
18580     DISP "ŝ - search for resp peak (+.1 Hz)"
18590     WAIT 2
18600 SUBEND
25 18610 !
18620 !
18630 !
18640 SUB Offgraph
18650     COM /Vitaldata/ Rfa,Lfa,Peakratio,Meas_
30      resp,Next_time
18660     PRINT CHR$(12)
18670     PRINT TABXY(1,14);"RR=";PROUND(Meas_resp,-
      2);"Hz"
18680     PRINT TABXY(1,15);"lfa=";Lfa
35 18690     PRINT TABXY(1,16);"rfa=";Rfa
18700     PRINT TABXY(1,17);"ratio=";Peakratio
```

- 136 -

```
18710      PRINT TABXY(1,18);"next transfer:
           ";TIME$(Next_time)
18720  SUBEND
18730      !
5  18740      !
18750      ! This subroutine edits the data
18760      !
18770      !
18780  SUB Editor
10  18790      COM /Editor/ Edit_msg$(80)
18800      COM /Vitaldata/ Rfa,Lfa,Peakratio,Meas_
           resp,Next_time
18810  Key_in: !
18820      PRINT TABXY(1,18);"
15  18830      PRINT TABXY(1,18);Edit_msg$
18840      IF TIMEDATE>Next_time-15 THEN GOTO Keyend
18850      ON TIME (TIMEDATE+10) MOD 86400,3 GOTO Keyend
18860      DISP "type message"
18870      GRAPHICS OFF
20  18880      ON KBD,2 GOTO Next_char
18890  Key_wait:GOTO Key_wait
18900  Next_char:Key$=KBD$
18910      ON TIME (TIMEDATE+10) MOD 86400,3 GOTO Keyend
18920      IF NUM(Key$)=255 THEN
25  18930          IF NUM(Key$[2])=69 THEN GOTO End_key
18940          IF NUM(Key$[2])=66 THEN !..backspacing
18950              New_msg_len=LEN(Edit_msg$)-1
18960              IF New_msg_len<=0 THEN New_msg_len=0
18970              Edit_msg$=Edit_msg$[1;New_msg_len]
30  18980          END IF
18990          IF NUM(Key$[2])=35 THEN !..clear line
19000              Edit_msg$=""
19010          END IF
19020      ELSE
35  19030          IF LEN(Edit_msg$)<66 THEN !..can add
           ! characters
```

- 137 -

```

19040          Edit_msg$=Edit_msg$&Key$
19050          ELSE
19060          BEEP
19070          END IF
5 19080      END IF
19090      PRINT TABXY(1,18);"
19100      PRINT TABXY(1,18);Edit_msg$
19110      GOTO Key_wait
19120 Keyend: !
10 19130 End_key:OFF KBD
19140      OFF TIME
19150 SUBEND
19160 !
19170 !
15 19180 !
19190 SUB Chart(Char_num)
19200      COM /Subject/ Sub_name$,Hos_num$,Id_age$,Id_
          wt$,Id_ht$,Diag$,Opera$,Halt_pg
19210      COM /Io_chart/ Io_time$(*),Iv_intake(*),Fluid_
20      in(*),In_tot(*),Urine(*),Chest(*),Out_
          tot(*),Net(*),Io_ptr
19220      COM /Lab_chart/ Lab_
          time$(*),Na(*),Kl(*),Cl(*),Hco3(*),Ca(*),Hct(*),G
          luc(*),Dig(*),Pt(*),Ptt(*),Creat(*),Bun(*),Lab_
25      ptr
19230      COM /Vent_chart/ Vent_
          time$(*),Rate(*),Fio2(*),Pp(*),Peep(*),Tv(*),
          Ie_ratio$(*),Airp(*),Ph(*),Po2(*),Pco2(*),
          Bgo3(*),Be(*),Vent_ptr
30 19240      COM /Pres_chart/ Pres_time$(*),Ao_s(*),Ao_
          d(*),Ao_m(*),Pa_s(*),Pa_d(
          ),Pa_m(*),La_m(*),Ra_m(*),Pres_ptr,Pres_in
19250      COM /Pressure/
          Top1,Top2,Top3,Top4,Bot1,Bot2,Bot3,Bot4
35 19260      COM /Heart_index/ Heart_
          time$(*),Ci(*),Pvri(*),Svri(*),Heart_ptr

```

- 138 -

```
19270      COM /Drugs/ Drug_time$(*),Drug_name$(*),Drug_
          dos$(*),Drug_ptr
19280      Pres_stl=0
19290      Lab_stl=0
5  19300      Io_stl=0
19310      Vent_stl=0
19320      Drug_stl=0
19330      !
19340      ! set up identifying subject info
10 19350      !
19360      PRINT CHR$(12)
19370      PRINT TABXY(1,1);
19380      PRINT USING Image_wt1;Sub_name$,Hos_
          num$,TIMES$(TIMEDATE),DATES$(TIMEDATE)
15 19390 Image_wt1:IMAGE  "Name: ",K,XXXX,"Hosp num:
          ",K,XXXXX,K,XXXXX,K
19400      PRINT TABXY(1,2);
19410      PRINT USING Image_wt2;Id_age$,Id_wt$,Id_
          ht$,Diag$,Opera$
20 19420 Image_wt2:IMAGE  "Age: ",K,XXXX,"Wt(kg):
          ",K,XXXX,"Ht(cm): ",K,XXXX,"Diag
          : ",K,XXXX,"Op: ",K
19430      !
19440      ! go to appropriate chart
25 19450      !
19460      ON Chart_num GOTO In_out,Lab_val,Vent_
          val,Pres_val,Drug
19470 In_out:!.intake/output
19480      IF Io_ptr>3 THEN Io_stl=2
30 19490      IF Io_ptr>5 THEN
19500          DISP "do not input more Intake/Output
          data; disc full"
19510          WAIT 3
19520          SUBEXIT
35 19530      END IF
19540      PRINT TABXY(30,3);"INTAKE/OUTPUT CHART"
```



- 139 -

```
19550      PRINT TABXY(1,4);"Intake (cc/hr) "
19560      PRINT TABXY(1,5);"Time"
19570      PRINT TABXY(4,6);"Maint. Fluid"
19580      PRINT TABXY(4,7);"Other Fluids"
5  19590      PRINT TABXY(1,9);"Total "
19600      PRINT TABXY(1,11);"Output (cc/hr)"
19610      PRINT TABXY(4,12);"Urine"
19620      PRINT TABXY(4,13);"Chest"
19630      PRINT TABXY(1,15);"Total"
10 19640      PRINT TABXY(1,17);"Net I/O"
19650      Start=25
19660      FOR I=Io_st1 TO Io_ptr
19670          PRINT TABXY(Start,5);Io_time$(I)
19680          PRINT TABXY(Start,6);Iv_intake(I)
15 19690          PRINT TABXY(Start,7);Fluid_in(I)
19700          PRINT TABXY(Start,9);In_tot(I)
19710          PRINT TABXY(Start,12);Urine(I)
19720          PRINT TABXY(Start,13);Chest(I)
19730          PRINT TABXY(Start,15);Out_tot(I)
20 19740          PRINT TABXY(Start,17);Net(I)
19750          Start=Start+10
19760      NEXT I
19770      GOTO Finish
19780!
25 19790!
19800 Lab_val:!      ...lab values
19810      IF Lab_ptr>3 THEN Lab_st1=2
19820      IF Lab_ptr>7 THEN
19830          DISP "do not input any more lab values;
30          disc full"
19840          WAIT 3
19850          SUBEXIT
19860      END IF
19870      PRINT TABXY(30,3);"Lab Values"
35 19880      PRINT TABXY(10,4);"Time"
19890      PRINT TABXY(1,6);"Na"
```

- 140 -

```
19900      PRINT TABXY(1,7);"K"
19910      PRINT TABXY(1,8);"Cl"
19920      PRINT TABXY(1,9);"HCO3"
19930      PRINT TABXY(1,10);"Ca"
5  19940      PRINT TABXY(1,11);"Hct"
19950      PRINT TABXY(1,12);"Glucose"
19960      PRINT TABXY(1,13);"Dig level"
19970      PRINT TABXY(1,14);"PT"
19980      PRINT TABXY(1,15);"PTT"
10 19990      PRINT TABXY(1,16);"Creat"
20000      PRINT TABXY(1,17);"Bun"

20010      Start=15
15 20020      FOR I=Lab_stl TO Lab_ptr
20030          PRINT TABXY(Start+10,4);Lab_time$(I)
20040          PRINT TABXY(Start+10,6);Na(I)
20050          PRINT TABXY(Start+10,7);Kl(I)
20060          PRINT TABXY(Start+10,8);Cl(I)
20 20070          PRINT TABXY(Start+10,9);Hco3(I)
20080          PRINT TABXY(Start+10,10);Ca(I)
20090          PRINT TABXY(Start+10,11);Hct(I)
20100          PRINT TABXY(Start+10,12);Gluc(I)
20110          PRINT TABXY(Start+10,13);Dig(I)
25 20120          PRINT TABXY(Start+10,14);Pt(I)
20130          PRINT TABXY(Start+10,15);Ptt(I)
20140          PRINT TABXY(Start+10,16);Creat(I)
20150          PRINT TABXY(Start+10,17);Bun(I)
20160          Start=Start+10
30 20170      NEXT I
20180      GOTO Finish
20190!
20200!
20210 Vent_val:!!          ....ventilation values
35 20220      IF Vent_ptr>3 THEN Vent_stl=2
20230      IF Vent_ptr>5 THEN Vent_stl=4
```

- 141 -

```
20240 IF Vent_ptr>7 THEN
20250     DISP "do not input any more Vent values;
           disc full"
20260     WAIT 3
5 20270     SUBEXIT
20280 END IF
20290 PRINT TABXY(30,3);"VENTILATION"
20300 PRINT TABXY(1,4);"Settings           Hour:".
20310 PRINT TABXY(4,5);"Rate"
10 20320 PRINT TABXY(4,6);"FIO2"
20330 PRINT TABXY(4,7);"Peak Pres"
20340 PRINT TABXY(4,8);"Peep"
20350 PRINT TABXY(4,9);"TV"
20360 PRINT TABXY(4,10);"I:E ratio"
15 20370 PRINT TABXY(4,11);"Mean air"
20380 PRINT TABXY(1,12);"Blood Gases"
20390 PRINT TABXY(4,13);"ph"
20400 PRINT TABXY(4,14);"pO2"
20410 PRINT TABXY(4,15);"pCO2"
20 20420 PRINT TABXY(4,16);"HCO3"
20430 PRINT TABXY(4,17);"BE"
20440 Start=15
20450 FOR I=Vent_stl TO Vent_ptr
20460     PRINT TABXY(Start+10,4);Vent_time$(I)
25 20470     PRINT TABXY(Start+10,5);Rate(I)
20480     PRINT TABXY(Start+10,6);Fio2(I)
20490     PRINT TABXY(Start+10,7);Pp(I)
20500     PRINT TABXY(Start+10,8);Peep(I)
20510     PRINT TABXY(Start+10,9);Tv(I)
30 20520     PRINT TABXY(Start+10,10);Ie_ratio$(I)
20530     PRINT TABXY(Start+10,11);Airp(I)
20540     PRINT TABXY(Start+10,13);Ph(I)
20550     PRINT TABXY(Start+10,14);Po2(I)
20560     PRINT TABXY(Start+10,15);Pco2(I)
35 20570     PRINT TABXY(Start+10,16);Bgo3(I)
20580     PRINT TABXY(Start+10,17);Be(I)
```

- 142 -

```
20590          Start=Start+10
20600      NEXT I
20610      GOTO Finish
20620!
5  20630!
20640 Pres_val:!!          ....pressure values
20650      IF Pres_ptr>12 THEN Pres_stl=5
20660      IF Pres_ptr>17 THEN
20670          DISP "Do not input any more pressures;
10          disc full"
20680          WAIT 3
20690          SUBEXIT
20700      END IF
20710      PRINT TABXY(9,3);"Time:"
15  20720      PRINT TABXY(1,4);"Systemic"
20730      PRINT TABXY(4,5);"systolic"
20740      PRINT TABXY(4,6);"diastolic"
20750      PRINT TABXY(4,7);"mean"
20760      PRINT TABXY(1,8);"Pulmonary"
20  20770      PRINT TABXY(4,9);"systolic"
20780      PRINT TABXY(4,10);"diastolic"
20790      PRINT TABXY(4,11);"mean"
20800      PRINT TABXY(1,12);"LA mean"
20810      PRINT TABXY(1,13);"RA mean"
25  20820      PRINT TABXY(9,14);"Time: ".
20830      PRINT TABXY(1,15);"C.I."
20840      PRINT TABXY(1,16);"PVRI"
20850      PRINT TABXY(1,17);"SVRI"
20860      Start=15
30  20870      FOR I=Pres_stl TO Pres_ptr
20880          PRINT TABXY(Start,3);Pres_time$(I)
20890          PRINT TABXY(Start,5);Ao_s(I)
20900          PRINT TABXY(Start,6);Ao_d(I)
20910          PRINT TABXY(Start,7);Ao_m(I)
35  20920          PRINT TABXY(Start,9);Pa_s(I)
20930          PRINT TABXY(Start,10);Pa_d(I)
```

- 143

```
20940      PRINT TABXY(Start,11);Pa_m(I)
20950      PRINT TABXY(Start,12);La_m(I)
20960      PRINT TABXY(Start,13);Ra_m(I)
20970      Start=Start+5
5  20980      NEXT I
20990      Start=15
21000      FOR I=0 TO Heart_ptr
21010          PRINT TABXY(Start,14);Heart_time$(I)
21020          PRINT TABXY(Start,15);Ci(I)
10 21030          PRINT TABXY(Start,16);Pvri(I)
21040          PRINT TABXY(Start,17);Svri(I)
21050          Start=Start+5
21060      NEXT I
21070      GOTO Finish
15 21080!
21090!
21100 Drug:!                      ....hey man, drugs
21110      IF Drug_ptr>9 THEN Drug_stl=4
21120      IF Drug_ptr>14 THEN Drug_stl=9
20 21130      IF Drug_ptr>19 THEN Drug_stl=14
21140      IF Drug_ptr>24 THEN Drug_stl=19
21150      IF Drug_ptr>29 THEN Drug_stl=24
21160      IF Drug_ptr>34 THEN Drug_stl=29
21170      IF Drug_ptr>38 THEN
25 21180          DISP "do not enter more drugs; disc full"
21190          WAIT 3
21200          SUBEXIT
21210      END IF
21220      PRINT TABXY(30,4);"Drug Chart"
30 21230      PRINT TABXY(1,6);"Name"
21240      PRINT TABXY(30,6);"Dosage"
21250      PRINT TABXY(60,6);"Time"
21260      D_line=7
21270      FOR I=Drug_stl TO Drug_ptr
35 21280          PRINT TABXY(1,D_line);Drug_name$(I)
21290          PRINT TABXY(30,D_line);Drug_dos$(I)
```

- 144 -

```
21300      PRINT TABXY(60,D_line);Drug_time$(I)
21310      D_line=D_line+1
21320      NEXT I
21330 Finish: !
5  21340  SUBEND
21350  !
21360  !
21370  DEF FNLval(Lnum$)
21380      Numval=VAL("9"&Lnum$)
10 21390      If Num val=9 THEN
21400      Rval=9999.999
21410      RETURN Rval
21420      ELSE
21430      Numval=VAL(Lnum$)
15 21440      RETURN Numval
21450      END IF
21460      FNEND
```

20

25

30

35

- 145 -

```

10 Teaser7: ! This program reviews data taken by sgrape
20      ! and allows all the graphs to be printed (when
30      ! its done)
40      !
5  50      ! .....
60      !
70      ! LAST REVISION: 1 May 1985
80      ! .....

10  90      !
100     !
110     ! .....
120     !
130     ! SET UP ERROR HANDLERS
15  140     ! SET UP COMMON STORAGE/ARRAY STORAGE
150     ! .....

160     !
170     !
20  171     COM /Vars/ Ffthrvar, Fftrespvar
180     COM /Intr_7/ Int_flag, Status_bytes(5)
190     COM /Flags/ Atod_done, Scanner_done, Memory1_
done, Memory2_done, Timer_done, Counter_done,
Memory3_done, Memory4_done
25  200     COM /Io_arrays/ Counters(3), Counters2(3), Time_
base$(7)
210     COM /Multi_param/ Start_chan, Stop_chan, Pacing_
bits, Pacing_rate, Num_pts, Num_xfer,
Num_xfer_left, Name_len, Scr_file$(28), Scr_
30  file2$(28)
220     COM /Hr_sig/ Num_pulses, Last_pulse, First_blk_
flg, Last_time, Num_hr_sig, Max
_hr_pts, Avg_hr, Rollover, Hr_smooth
230     COM /Hr_stats/ Hr_histo(128), Histo_min, Histo_
35  max, Num_fudge, Num_histo_pnts
,@Err_log

```

- 146 -

```
240    COM /Plot_par/ Plotbox,Boxcar_flg,Log_
      plotflg,Freq_limit,Resp_search,Pct_thresh
250    COM /Graphs/
      Hrdata(512),Hrspec(512),Respspec(512),Bpspec(512)
5    260    COM /Vitaldata/ Rfa,Lfa,Peakratio,Meas_resp,Next_
      time
      270    COM /Idfield/ Id_field$(18)
      280    COM /Messagecom/ Message$(10)[80],@Messages
      290    COM /Trends/ Mean_hr_t(60),Lfa_t(60),Rfa_
10    t(60),Ratio_t(60),T_ptr,Time_now
      1,Meas_resp_t(60)
      300    DIM Msg_pad$(20)[80],Edit_msg$(80)
      310    DIM Msg_buffer$(80) BUFFER
      320    ASSIGN @Msg_buffer TO BUFFER Msg_buffer$
15    330    Log_plotflg=0
      340    Freq_limit=1.
      350    Resp_search=.1
      360    Pct_thresh=.2
      370    Scr_file$="?"
20    380    !
      390    ! Set up common/array storage for waveform
      analysis
      400    !
      410    !.....
25    420    !
      430    ! Set up common/array storage for waveform
      ! analysis
      440    !.....

30    450    !
      460    COM /Directory/ Dir$(160),@Printer
      470    COM /Wf1/ Printer,Plotter,String$(40)
      480    COM /Wf2/ Signal(8257),Number_pnts,Type,Sampling_
      period
35    490    COM /Wf3/ Segment_size,Overlap,Num_segments,Pnts_
      used,Fft_size
```



- 147 -

```

500    COM /Wf5/ Refn(63),Refd(63),Refno,Refdo,Refgain
510    COM /Autoparam/ Up_down,Up_delay,Dn_delay
520    COM /Fftcom/ INTEGER Bitrev(512),Sincos(512)
530    !
5  540    DISP "loading subroutines"
550    LOADSUB ALL FROM "hr_siggen8"
560    LOADSUB ALL FROM "automaxsb2"
570    LOADSUB ALL FROM "fft_anal6"
580    DISP "load data disks and press CONTINUE"
10 590    PAUSE
600    !
610    !.....
620    ! The HP 9826/9836 flexible disk (5-1/4") has the
        following structure
15 630    ! 2 sides, 33 tracks/side, 16 sectors/track, 256
        bytes/sector
640    ! 1 track = 4096 bytes = 16 sectors
650    ! 1 side = 135168 bytes = 528 sectors
660    ! 1 disk = 270336 bytes = 1056 sectors
20 670    ! 1 disk = 135168 words = 132K words
680    !.....

690    !
700    !
25 710    INTEGER Hpib_buffer1(2048) BUFFER
720    INTEGER Hpib_buffer2(2048) BUFFER
730    DIM Hr_signal(1024) BUFFER
740    Read_ptr1=0
750    Read_ptr2=0
30 760    Begin: !
770    Selections: !
780    !
790    !
800    ! NOW SET UP THE SCAN CARD PARAMETERS (DEFAULT
35    ! VALUES)
810    ! START CHANNEL (3.0) - 0

```

- 148 -

```
820  !      STOP CHANNEL (3.1) - 1
830  !      PACING (3.2) - 40 USEC
840  !      SEQN'L SCAN (3.3) - XXXX XXXX XXX1 ( 1)
850  !      INTN'L PACING (3.3) - XXXX XXXX X1XX ( 4)
5   860  !      MSEC TIMEBASE (3.3) - XXX1 XXXX XXXX (256)
870  !
880  CALL Get_param
890  !
900  ! set up the bit reverse index
10  910  !
920  Npair=Num_pts/2
930  K=0
940  FOR J=1 TO Npair-1
950      I=2
15  960      Ndivi=Npair/I
970      IF K<Ndivi THEN 1010
980      K=K-Ndivi
990      I=I+I
1000     GOTO 960
20  1010     K=K+Ndivi
1020     Bitrev(J+1)=K+1
1030  NEXT J
1040  !
1050  ! set up the sin/cosine table
25  1060  !
1070  Angl=ATN(1)*8/Npair
1080  FOR J=0 TO Npair-1
1090      Sincos(J)=SIN(Angl*J)
1100  NEXT J
30  1110  !
1120  ! set up other data paths
1130  !
1140  ! ASSIGN @Err_log TO "errs"&Id_
      field$&":HP8290X,700,1";FORMAT OFF
35  1150  ! ASSIGN @Messages TO "msgs"&Id_
      field$&":HP8290X,700,1";FORMAT OFF
```

- 149 -

```
1160 ! ASSIGN @Temp_trend TO "trnd"&Id_
      field$&":HP8290X,700,1";FORMAT OFF
1170 IF Num_pts=0 THEN GOTO Begin
1180 Read_ptr1=0
5 1190 Setup_scan:DISP " NUMBER OF POINTS=";Num_pts
1200 Read_ptr1=0
1210 Read_ptr2=0
1220 Setup_counter:~
1230 Setup_clock:~
10 1240 Block_time=Pacing_rate*1.024
1250 First_blk_flg=1
1260 Num_msgs=0
1270 Message_line=0
1280 Msg_dp_request=0
15 1290 Resp_dp_flg=0
1300 Max_hr_pts=1024
1310 Last_time=0
1320 !
1330 ! setup control parameters
20 1340 !
1350 Defaultset:~
1360 INPUT "use default settings?",Resp$
1370 IF Resp$="N" THEN Frqlimset
1380 Freq_limit=2.
25 1390 Pct_thresh=.2
1400 Resp_dp_flg=1
1410 Resp_search=.2
1420 Hcdopy_flg=0
1430 PRINT "Spectra displayed to";Freq_limit;"Hz"
30 1440 PRINT "resp peak search threshold=";Pct_thresh
1450 PRINT "resp series plot w/hr series"
1460 PRINT "resp peak search starts at";Resp_
      search;"Hz"
1470 PRINT "no hard copy will be printed"
35 1480 INPUT "is this ok?",Resp$
1490 IF Resp$<>"Y" THEN Defaultset
```

- 150 -

```
1500  GOTO Skipset
1510  Frqlimset: !
1520  INPUT "frequency limit?", Freq_limit  !..change
      spectra disp.freq.range
5  1530  IF Freq_limit<>1. THEN Freq_limit=2.
1540  PRINT "Spectra displayed to"; Freq_limit; "Hz"
1550  INPUT "is this ok?", Resp$
1560  IF Resp$<>"Y" THEN Frqlimset
1570  Searchset: !
10 1580  INPUT "resp peak threshold?", Pct_thresh !..change
      peak search threshold
1590  IF Pct_thresh>.8 THEN Pct_thresh=.2
1600  PRINT "resp peak search threshold="; Pct_thresh
1610  INPUT "is this ok?", Resp$
15 1620  IF Resp$<>"Y" THEN Searchset
1630  Respdpset: !
1640  INPUT "display resp time series?", Resp$
      !..display respiration time series
1650  IF Resp$<>"N" THEN
20 1660      Resp_dpflg=1
1670      PRINT "resp series plot w/hr series"
1680  ELSE
1690      Resp_dpflg=0
1700      PRINT "cancel resp series plot"
25 1710  END IF
1720  INPUT "is this ok?", Resp$
1730  IF Resp$<>"Y" THEN Respdpset
1740  Resppkset: !
1750  INPUT "start for resp peak search?", Resp_
30      search  !..change respiration
      peak search
1760  IF Resp_search>Freq_limit-.1 THEN Resp_search=.1
1770  PRINT "resp peak search starts at"; Resp_
      search; "Hz"
35 1780  INPUT "is this ok?", Resp$
1790  IF Resp$<>"Y" THEN Resppkset
```

- 151 -

```
1800 Hdcopyset: !
1810 INPUT "print hardcopy?",Resp$
1820 IF Resp$="N" THEN
1830     Hdcopyflg=0
5 1840     PRINT "no hard copy will be printed"
1850 ELSE
1860     Hdcopyflg=1
1870     PRINT "hard copy will be printed"
1880 END IF
10 1890 INPUT "is this ok?",Resp$
1900 IF Resp$<>"Y" THEN Hdcopyset
1910 Skipset: !
1920 !
1930 ! Read data continuously
15 1940 !
1950 ! Set up the memory buffers and disk files
1960 !
1970 Reading: !
1980 ASSIGN @In_buffer TO BUFFER Hpib_buffer1(*)
20 1990 ASSIGN @Diskbuffer TO Scr_file$;FORMAT OFF
2000 ASSIGN @In_buffer2 TO BUFFER Hpib_buffer2(*)
2010 ASSIGN @Diskbuffer2 TO Scr_file2$;FORMAT OFF
2020 !
2030 Data_lockout=0
25 2040 !
2050 ! generate id fields to identify data files
2060 !.....
2070 ! the first 256 bytes of the file are reserved for
    identification
30 2080 !
2090 ! the reserved data are:
2100 !     byte 1 - 72 ("H") or 82 ("R"): hr or resp_
    !     file
2110 !     byte 2 - year (at beginnig of expt.)
35 2120 !     byte 3 - month
2130 !     byte 4 - day
```

- 152 -

```

2140 !      byte    5 - hour
2150 !      byte    6 - minute
2160 !      byte    7 - collecting program date (0-365)
2170 !      byte    8 - collecting program year (1984-?)
5  2180 !      byte   9-16: unused
2190 !      byte   17 - pacing rate (0-32768)
2200 !      byte   18 - pacing rate units(77 ="M" or 85
      !      ="U")
2210 !      byte   19 - number of transfers
10 2220 !      byte   20 - number of point/transfer (=1024)
2230 !      byte   21 - number of A/D channels used (=1)
2240 !      byte   22-256 : unassigned
2250 !
2260 !      the remainder of the file is data
15 2270 !      each transfer is preceded by an identifying
      !      string of 8 bytes
2280 !      byte 1 - time of day (timedate mod 86400)/60
2290 !      byte 2 - number of points in next transfer
2300 !      byte 3 - H/R (check to make sure this is the
20      right file)
2310 !.....
2320 !
2330 ! INTEGER Id_buffer(255) BUFFER
2340   Time_now=TIMEDATE
25 2350 ! Id_buffer(0)=72           !..Heart rate file
2360   Date_now$=DATE$(TIMEDATE)
2370 ! Day_now=VAL(Date_now$)
2380 ! Year_now=VAL(Date_now$[8;4])
2390 ! Month_now=FNMonth(Date_now$)
30 2400 ! Id_buffer(1)=Year_now           !..year
2410 ! Id_buffer(2)=Month_now           !..month
2420 ! Id_buffer(3)=Day_now             !..day
2430   Time_now1=Time_now MOD 86400
2440 ! Id_buffer(4)=Time_now1/3600       !..hour
35 2450 ! Id_buffer(5)=(Time_now1 MOD 3600)/60 !..min
2460 ! Id_buffer(6)=348                 !..pgm date

```

- 153 -

```
2470 ! Id_buffer(7)=1984                !..pgm year
2480 ! Id_buffer(16)=Pacing_rate
2490 ! Id_buffer(17)=77                  !..MSEC
2500 ! Id_buffer(18)=Num_xfer
5  2510 ! Id_buffer(19)=1024             !..num_pts
2520 ! Id_buffer(20)=1                  !..# channels
2530 !
2540 !
2550 ! read id field for heart rate file
10 2560 !
2570 ! ASSIGN @Id_buffer TO BUFFER Id_buffer(*)
2580 ! TRANSFER @Diskbuffer2 TO @Id_buffer;COUNT
      ! 256,WAIT
2590 ! ASSIGN @Id_buffer TO *
15 2600 !
2610 ! read id field for respiratory file
2620 !
2630 ! Id_buffer(0)=82                  !..Resp file
2640 ! ASSIGN @Id_buffer TO BUFFER Id_buffer(*)
20 2650 ! TRANSFER @Diskbuffer TO @Id_buffer;COUNT 256,WAIT
2660 ! ASSIGN @Id_buffer TO *
2670 !
2680 !
2690 !
25 2700 ! begin transferring data from the A/D buffer
2710 !
2720 Blk_xfer:
2730   CONTROL @In_buffer,3;1
      ! Reset fill pointer for buffer
30 2740   CONTROL @In_buffer,4;0
      ! Reset current number of bytes in buffer
2750   CONTROL @In_buffer,5;1
      ! Reset empty pointer for buffer
2760 !
35 2770 ! read an 8 byte sequence to disk as a header for
      ! the transfer
```

- 154 -

```
2780  !
2790  CALL Rdheader(@Diskbuffer,Num_pts,"R")
2800  !
2810  Num_rdpts=Num_pts
5  2820  TRANSFER @Diskbuffer TO @In_buffer;COUNT Num_
      rdpts*2,CONT
2830  PRINT TABXY(1,18);
2840  PRINT USING Image_wtl;Num_xfer-Num_xfer_
      left+1,Num_xfer,TIMES$(Next_time),
10  Rdseg,Num_rdseg
2850  image_wtl:IMAGE      "Next xfer(",K,"/",K,"): ",K,"
                        seg=",K,"/",K
2860  !
2870  ! store A/D buffer on complete data file (also
15      save pointers for heart rate)
2880  !
2890  !
2900  Resume1:!
2910  Next_time=Next_time+INT(Block_time)
20  2920  !
2930  !
2940  !
2950  Resume2:!
2960  Num_xfer_left=Num_xfer_left-1
25  2970  CONTROL @In_buffer2,3;1
      ! Reset fill pointer for buffer
2980  CONTROL @In_buffer2,4;0
      ! Reset current number of bytes in buffer
2990  CONTROL @In_buffer2,5;1
30      ! Reset empty pointer for buffer
3000  !
3010  ! read an 8 byte sequence to disk as a header for
      ! the transfer
3020  !
35  3030  CALL Rdheader(@Diskbuffer2,Num_pulses,"H")
3040  TRANSFER @Diskbuffer2 TO @In_buffer2;COUNT Num_
```



- 155 -

```
        pulses*2, WAIT
3050  !
3060 Resume5: !
3070  Histo_max=8000
5  3080  Histo_min=-8000
3090  CALL Hr_sig_gen(Hpib_buffer2(*), Hr_signal(*))
3100  !
3110  !
3120 Resume6: !
10 3130  OUTPUT 2; CHR$(255)&CHR$(75);
    ! Clear CRT of text
3140  GINIT
3150  PLOTTER IS 3, "INTERNAL"
3160  GRAPHICS ON
15 3170  Xscale=8
3180  Hr_max=MAX(Hr_signal(*))
3190  Hr_min=MIN(Hr_signal(*))
3200  VIEWPORT 0,64,50,100
3210  WINDOW 0,1,0,1
20 3220  AXES .1,.1,0,0
3230  CSIZE 4
3240  Hr_signal(1024)=0
3250  Hr_sigsum=SUM(Hr_signal)
3260  Mean_hr=INT((Hr_sigsum/1024+Avg_hr))
25 3270  LDIR 0
3280  LORG 3
3290  MOVE .2,.9
3300  LABEL "HR data   hr="; Mean_hr
3310  CSIZE 4
30 3320  MOVE .05,1
3330  LORG 3
3340  LABEL "250 bpm"
3350  WINDOW 1,0,1,0
3360  AXES 0,0,0,0
35 3370  IF Hr_dispflg=1 THEN
3380      WINDOW 0,1024,Hr_min,Hr_max
```

- 156 -

```
3390 ELSE
3400     Low_window=INT(-Avg_hr)
3410     High_window=Low_window+250.
3420     WINDOW 0,1024,Low_window,High_window.
5  3430 END IF
3440     FOR I=0 TO 1023
3450         PLOT I,Hr_signal(I)
3460     NEXT I
3470 !CALL Pauser
10 3480 IF Fftskpflg=1 THEN GOTO Skip_fft
3490 !
3500 ! display respirations time series also
3510 !
3520 IF Resp_dpflg=1 THEN
15 3530     Max_resp=MAX(Hpib_buffer1(*))
3540     Min_resp=MIN(Hpib_buffer1(*))
3550     IF Mean_hr>100 THEN
3560         VIEWPORT 0,64,50,65
3570     ELSE
20 3580         VIEWPORT 0,64,75,90
3590     END IF
3600     WINDOW 0,1023,Min_resp,Max_resp
3610     MOVE 0,Hpib_buffer1(0)
3620     FOR I=1 TO 1023
25 3630         PLOT I,Hpib_buffer1(I)
3640     NEXT I
3650 ELSE
3660     Resp_dpflg=0
3670 END IF
30 3680 !
3690 ! now process heart rate data with waveform
      analyse package
3700 ! make sure the hr_signal has zero mean
3710 !
35 3711 MAT Signal= (0)
3720     Hr_bias=Hr_sigsum/1024
```

- 157 -

```
3730  FOR I=0 TO 1023
3740      Signal(I)=Hr_signal(I)-Hr_bias
3750  NEXT I
3751  Hr_var=DOT(Signal,Signal)/1024
5   3760  Plotbox=2
3770  DISP "HR fft in process"
3780  CALL Wf_analyzer(Pacing_rate)
3790  !
3800  ! now process respiration data with waveform
10  analysis package
3810  !
3820  MAT Signal= (0)
3830  FOR I=0 TO 1023
3840      Signal(I)=Hpib_buffer1(I)
15  3850  NEXT I
3860  Signal_avg=SUM(Signal)/1024.
3870  MAT Signal= Signal-(Signal_avg)
3880  Plotbox=4
3881  Respvar=DOT(Signal,Signal)/1024
20  3890  DISP "RESP fft in process"
3900  CALL Wf_analyzer(Pacing_rate)
3901  PRINT "hr_var,respvar";Hr_var;Respvar
3902  PRINT "fft vars: ";Ffthrvar,Fftrespvar
3910  Trend_dp=0 !..trend graph not displayed
25  3920  !
3930  ! waveform analysis completed, compile trends and
      store in temporary file
3940  !
3950  Mean_hr_t(T_ptr)=Mean_hr
30  3960  Lfa_t(T_ptr)=Lfa
3970  Rfa_t(T_ptr)=Rfa
3980  Ratio_t(T_ptr)=Peakratio
3990  Meas_resp_t(T_ptr)=Meas_resp
4000  T_ptr=T_ptr+1
35  4010  IF Hdcopyflg=1 THEN
4011      DUMP DEVICE IS 701
```

- 158 -

```
4020      DUMP GRAPHICS
4030      PRINTER IS 701
4040      PRINT "hr=";Mean_hr
4050      PRINT "lfa=";Lfa
5 4060      PRINT "rfa=";Rfa
4070      PRINT "ratio";Peakratio
4080      PRINT "RR";Meas_resp
4090      PRINT "transfer#";T_ptr
4091      PRINT "hr_var,respvar";Hr_var;Respvar
10 4092      PRINT "fft vars: ";Ffthrvar,Fftrespvar
4100      PRINTER IS 1
4110      END IF
4120      !
4130      ! continue with data collection
15 4140      !
4150      Skip_fft: !
4160      IF Num_xfer_left<=0 THEN
4170          GOTO Eo_blk_xfer
4180      ELSE
20 4190          DISP Num_xfer_left;"transfers remaining"
4200          WAIT 3
4210          GOTO Blk_xfer
4220      END IF
4230      Eo_blk_xfer:End_time=TIMEDATE
25 4240      Delta_time=End_time-Start_time
4250      !
4260      Stop_pacing=TIMEDATE
4270      !
4280      Aborter: !
30 4290      ASSIGN @In_buffer TO *
4300      ASSIGN @In_buffer2 TO *
4310      ASSIGN @Diskbuffer TO *
4320      ASSIGN @Diskbuffer2 TO *
4330      ! ASSIGN @Err_log TO *
35 4340      ! ASSIGN @Messages TO *
4350      ! ASSIGN @Temp_trend TO *
```

- 159 -

```
4360 CALL Pauser
4370 GRAPHICS OFF
4380 CALL Get_param
4390 ! ASSIGN @Err_log TO "errs"&Id_
5      field$&":HP8290X,700,1";FORMAT OFF
4400 ! ASSIGN @Messages TO "msgs"&Id_
      field$&":HP8290X,700,1";FORMAT OFF
4410 IF Num_pts=0 THEN GOTO Begin
4420 GOTO Setup_scan
10 4430 END
4440 !
4450 !
4460 !
4470 !
15 4480 !
4490 SUB Pauser
4500     DISP "press CONTINUE to continue"
4510     PAUSE
4520     DISP
20 4530 SUBEND
4540 !
4550 !
4560 !
4570 !
25 4580 !
4590 SUB Get_param
4600     COM /Multi_param/ Start_chan,Stop_chan,Pacing_
        bits,Pacing_rate,Num_pt
        s,Num_xfer,Num_xfer_left,Name_len,Scr_
30      file$[28],Scr_
        file2$[28]
4610     COM /Trends/ Mean_hr_t(*),Lfa_t(*),Rfa_
        t(*),Ratio_t(*),T_ptr,Time_now
        l,Meas_resp_t(*)
35 4620     COM /Vitaldata/ Rfa,Lfa,Peakratio,Meas_
        resp,Next_time
```

- 160 -

```
4630      COM /Idfield/ Id_field$
4640      DIM Mo$(24)
4650      Mo$="JAFBMRAPMYJNJLAUSPOCNODC"
4660      INTEGER Id_buffer(255) BUFFER
5 4670      Disk_name$=":HP8290X,700,1"
4680 Oldmsg:PRINT CHR$(12)
4690 !
4700 !
4710 Ch_sel:~
10 4720      Start_chan=0
4730      Stop_chan=0
4740 !
4750      Pacing_bits=0
4760 Pacing_sel:~
15 4770      Base$="M"
4780      Pacing_bits=261
4790 !
4800      Base$=Base$&"SEC"
4810 !
20 4820 !
4830 ! FINDOUT BLOCKSIZE FOR DATA TRANSFER
4840 !
4850 Get_xfer:DISP "Enter number of transfers: (0 -
      change scan, <0 - quit)"
25 4860      OUTPUT 2;55;
4870      ENTER 2;Num_xfer
4880      IF Num_xfer<0 THEN      !..terminate program
4890      INPUT "to lose trend data type
      'lose'",Response$
30 4900      IF Response$<>"lose" THEN
4910      CREATE BDAT
      "teasertrnd:HP8290X,700,1",19,256
4920      ASSIGN @Trndfile TO
      "teasertrnd:HP8290X,700,1";FORMAT OFF
35 4930      OUTPUT @Trndfile;Mean_hr_t(*),Lfa_
      t(*),Rfa_t(*),Ratio_t(*),Me
```

- 161 -

```

                                as_resp_t(*),T_ptr
4940      ASSIGN @Trndfile TO *
4950      END IF
4960      DISP "PROGRAM COMPLETED"
5 4970      STOP
4980      END IF
4990      IF Num_xfer=0 THEN
5000          Num_pts=0
5010          SUBEXIT
10 5020      END IF
5030 !
5040 ! since new data is to be taken, zero the trend
      graphs (120 pts=8hrs)
5050 !
15 5060      MAT Mean_hr_t= (0)
5070      MAT Rfa_t= (0)
5080      MAT Lfa_t= (0)
5090      MAT Ratio_t= (0)
5100      MAT Meas_resp_t= (0)
20 5110      T_ptr=0
5120      Ratio_t(0)=1 !..prevent trend graph errors on
      startup
5130      Rfa=0
5140      Lfa=0
25 5150      Meas_resp=0
5160      Peakratio=1
5170 !
5180 Intvl_sel:DISP "ENTER PACING RATE (IN
      ";Base$[1,4];"):"
30 5190      OUTPUT 2;250;
5200      ENTER 2;Pacing_rate
5210      IF Pacing_rate<0 OR Pacing_rate>65535 THEN
      GOTO Intvl_sel
5220 !
35 5230      Num_pts=1024*Num_xfer
5240      Num_header=256+8*Num_xfer
```

- 162 -

```
5250      INPUT "type in date on which data was
          taken",Datdate$
5251      INPUT "is trend file named 'trnd' (1) or
          'temp_trend' (2)?",File_nm
5  5260      Datdate$=DATE$(DATE(Datdate$))
          5270 !
          5280 ! the data files are named according to the date
          5290 ! in the following format:
          5300 !      xxxxxmmddyy
10 5310 ! where
          5320 !      xxxx - resp,hr__,msgs,errs,trnd
          5330 !      dd   - day
          5340 !      mm   - month
          (JA,FB,MR,AP,MY,JN,JL,AU,SP,OC,NO,DC)
15 5350 !      yy    - year
          5360      Month_now=FNMonth(Datdate$)*2-1
          5370      Mm$=Mo$(Month_now;2]
          5380      Id_field$=Datdate$[1;2]&Mm$&Datdate$[10;2]
          5390 ! new name for respiratory file: respddmmyy
20 5391 IF File_nm=1 THEN
          5400      Scr_file$="resp"&Id_field$&Disk_name$
          5410 ! new name for heart rate file: hr__ddmmyy
          5420      Scr_file2$="hr__"&Id_field$&Disk_name$
          5421      ELSE
25 5422      Scr_file$="AOK"&Disk_name$
          5423      Scr_file2$="hrAOK"&Disk_name$
          5424      END IF
          5430 ! new name for errorlog: errsddmmyy
          5440 ! new name for message log: msgsdmmyy
30 5450 ! name for trend summary file: trndddmmyy
          5460      Num_rec=-INT(-(Num_pts+Num_header)/128.)
          5470      Num_pts=1024
          5480      PRINT Num_pts*Num_xfer;"points were
          transferred in";Num_xfer;"blocks
35      of";Num_pts;"points"
          5490      !
```



- 163 -

```
5500      Num_xfer_left=Num_xfer
5510  SUBEND
5520  !
5530  !
5  5540  !
5550  !
5560  DEF FNMonth(Date_now$)
5570      Month$=Date_now$[4;3]
5580      Month=0
10  5590      IF Month$="Jan" THEN Month=1
5600      IF Month$="Feb" THEN Month=2
5610      IF Month$="Mar" THEN Month=3
5620      IF Month$="Apr" THEN Month=4
5630      IF Month$="May" THEN Month=5
15  5640      IF Month$="Jun" THEN Month=6
5650      IF Month$="Jul" THEN Month=7
5660      IF Month$="Aug" THEN Month=8
5670      IF Month$="Sep" THEN Month=9
5680      IF Month$="Oct" THEN Month=10
20  5690      IF Month$="Nov" THEN Month=11
5700      IF Month$="Dec" THEN Month=12
5710      RETURN Month
5720  FNEND
5730  !
25  5740  !
5750  !
5760  !
5770  !
5780  SUB Rdheader(@Disk,Num_bytes,File_id$)
30  5790      INTEGER Xheader(7) BUFFER
5800      ASSIGN @Xheader TO BUFFER Xheader(*)
5810      TRANSFER @Disk TO @Xheader;COUNT 16,WAIT
5820      ASSIGN @Xheader TO *
5830      Num_bytes=Xheader(1)
35  5840      File_id$=CHR$(Xheader(2))
5850  SUBEND
```

- 164 -

```
5860 !
5870 !
5880 !
5890 !
5 5900 !
5910 !
5920 SUB Trend_graph
5930 !
5940 COM /Trends/ Mean_hr_t(*),Lfa_t(*),Rfa_
10 t(*),Ratio_t(*),T_ptr,Time_now
    1,Meas_resp_t(*)
5950 COM /Multi_param/ Start_chan,Stop_chan,Pacing_
    bits,Pacing_rate,Num_pt
    s,Num_xfer,Num_xfer_left,Name_len,Scr_
15 file$(28),Scr_
    file2$(28)
5960 Block_time=Pacing_rate*1.024/3600.
5970 GINIT
5980 GCLEAR
20 5990 PRINT CHR$(12)
6000 GRAPHICS ON
6010 PRINT TABXY(1,18);"trend graph"
6020 Beg_time=Time_now1/3600-Block_time
6030 End_time=Beg_time+Num_xfer*Block_time
25 6040 Ibeg_time=INT(Beg_time)
6050 IF Ibeg_time<Beg_time THEN Ibeg_time=Ibeg_
    time+1
6060 !
6070 ! label the time axes
30 6080 !
6090 VIEWPORT 0,128,45,50
6100 WINDOW Beg_time,End_time,0,1
6110 IF INT(End_time)>Beg_time THEN
6120 LDIR 0
35 6130 FOR T_label=Ibeg_time TO INT(End_time)
6140 MOVE T_label,.5
```

- 165 -

```

        6150          LORG 5
        6160          CSIZE 4
        6170          LABEL T_label
        6180          NEXT T_label
5   6190          END IF
        6200          VIEWPORT 0,128,40,45
        6210          WINDOW 0,1,0,1
        6220          MOVE .5,0
        6230          LORG 4
10  6240          LABEL "Time (24 hr)"
        6250 !
        6260 ! draw the axes
        6270 !
        6280          VIEWPORT 0,128,50,100
15  6290          WINDOW Beg_time,End_time,0,1
        6300          AXES 1/15.,.1,Beg_time,0
        6310          WINDOW 1,0,1,0
        6320          AXES 0,.25,0,0
        6330 !
20  6340 ! mean heart rate trends
        6350 !
        6360          WINDOW -1,Num_xfer,0,200.
        6370          MOVE 0,Mean_hr_t(0)
        6380          FOR I=0 TO T_ptr-1
25  6390          DRAW I,Mean_hr_t(I)
        6400          NEXT I
        6410 !
        6420 ! lfa trends
        6430 !
30  6440          WINDOW -1,Num_xfer,0,10.
        6450          LINE TYPE 4,5
        6460          MOVE 0,Lfa_t(0)
        6470          FOR I=0 TO T_ptr-1
        6480          DRAW I,Lfa_t(I)
35  6490          NEXT I
        6500 !
```

- 166 -

```
6510 ! rfa trends
6520 !
6530     WINDOW -1,Num_xfer,0,10.
6540     LINE TYPE 5,5
5 6550     MOVE 0,Rfa_t(0)
6560     FOR I=0 TO T_ptr-1
6570         DRAW I,Rfa_t(I)
6580     NEXT I
6590 !
10 6600 ! ratio trends (with a line at ratio=2)
6610 !
6620     WINDOW -1,Num_xfer,-2.5,2.5
6630     LINE TYPE 8,5
6640     MOVE 0,LGT(Ratio_t(0))
15 6650     FOR I=0 TO T_ptr-1
6660         DRAW I,LGT(Ratio_t(I))
6670     NEXT I
6680     LINE TYPE 3,5 !..sparsely dotted line at
        ratio=2
20 6690     MOVE 0,LGT(2.)
6700     DRAW T_ptr-1,LGT(2.)
6710 !
6720 ! respiration trends
6730 !
25 6740     WINDOW -1,Num_xfer,0,200
6750     LINE TYPE 5,10
6760     MOVE 0,Meas_resp_t(0)
6770     FOR I=0 TO T_ptr-1
6780         DRAW I,Meas_resp_t(I)
30 6790     NEXT I
6800 !
6810 ! draw a key for line types
6820 !
6830     VIEWPORT 64,128,0,50
35 6840     WINDOW 0,1,0,13
6850     PRINT TABXY(55,15);"mean hr(0-200)"
```

- 167 -

```
6860      PRINT TABXY(55,16);"lfa      (0-10)"
6870      PRINT TABXY(55,17);"rfa      (0-10)"
6880      PRINT TABXY(55,18);"ratio(.01-100)"
6890      LINE TYPE 1,5
5  6900      MOVE .8,11
6910      DRAW 1.,11
6920      LINE TYPE 4,5
6930      MOVE .8,10
6940      DRAW 1.,10
10 6950      LINE TYPE 5,5
6960      MOVE .8,9
6970      DRAW 1.,9
6980      LINE TYPE 8,5
6990      MOVE .8,8
15 7000      DRAW 1.,8
7010      LINE TYPE 1,5
7020      SUBEND
```

20

25

30

35

- 168 -

' CALIB - program to calibrate instruments using  
board#1

' last revision: 4 April 1985

5

defint a-y ' only z denotes a real number

dim buffer(12800)

10

hrbpm=0

zfqglow=0.

zfqgres=0.

zlfa=0.

zrfa=0.

15

cls

'define ports on 8253

timer0=&h720

20

timer1=&h721

timer2=&h722

con8253=&h723

25

' set timer modes to 16 bit square wave rate  
generators

out con8253,&h36

out con8253,&h76

out con8253,&hB6

30

'for testing set timer 0 to 100Hz timebase

'2.38MHz/23864: 23864=93\*256+56

'set timer 0 to 1280Hz timebase

'(2.38MHz/1864) (1864=7\*256+72)

35

'set timer 1 as a 1Hz clock at startup

'(gives a heart rate signal at

- 169 -

```
'60bpm) 'set timer 2 as a flip flop
      out timer0,56
      out timer0,93
      out timer0,72
5      out timer0,7
      out timer1,0
      out timer1,5
      hrbpm=60
      out timer2,2
10     out timer2,0

      ' turn the gates on using the 8255 at bits 0,1,2
on portc
15     porta=&H700
      portb=&H708
      portc=&H710
      con8255=&H718

20     ' first set all 8255 ports to output, then set
portc to 0FFH
      out con8255,128
      out portc,&H0FF

25

      ' first print out the present value of the
interrupt vectors
      locate 4,1
30     gosub 10000

      ' install the interrupt with a dummy buffer and
      print vectors
35     reseter=256
      call wrbuffer(reseter)
```

- 170 -

```

        reseter=128
        call wrbuffer(reseter)
        call instinct
        locate 5,1
5         gosub 10000

        ' now go through required startup subroutines
        gosub 90          ' set up breathing

10    signal
        gosub 70          ' set up heart rate
        variations
        gosub 50          ' put some information
        on screen
15        gosub 80          ' turn D/A on
        locate 1,1
        print "commands: h(rvar),i(nt
on),q(uit),r(beats),b(reath),c(ounts)"

20
        ' wait until user hits a key
        savekey$=""
40        while
len(savekey$)=0:savekey$=savekey$+inkey$:wend
25        if savekey$="r" then gosub 50  'print heart
        beats
        if savekey$="q" then goto 9996  'quit
        if savekey$="c" then gosub 60  'print timers
        if savekey$="h" then gosub 70  'set up heart
30    rate
        ' variations
        if savekey$="i" then gosub 80  'unmask
        interrupts
        if savekey$="b" then gosub 90  'set up
35    breathing signal
        savekey$=""
```



- 171 -

```
        goto 40

        'print present value of heartbeats

5    50    locate 7,1
        call rdbeat(n)
        print "present heart beats are: ";n;time$
        return

10

        ' print present value of counters
        60    out control,0          'latch timer0
            tlow0=inp(timer0)
            thigh0=inp(timer0)
15    out control,&h40          'latch timer1
            tlow1=inp(timer1)
            thigh1=inp(timer1)
            out control,&h80          'latch timer2
            tlow2=inp(timer2)
20    thigh2=inp(timer2)
        locate 8,1
        print "timer0: ";tlow0+thigh0*16;tab(20);"
timer1:
        ";tlow1+thigh1*16;
25    print tab(40);"timer2: ";tlow2+thigh2*16
        return

30    ' set up the heart rate variations
        '      respiratory frequency is given by
        1280Hz/buffer
        '      length
        '      low frequency is 1280Hz/low frequency
35    divider
        ,
```

- 172 -

```
70      if numval<=0 then beep:print "setup analog
buffer
          first":return
71      locate 17,1
5        print "present lfa,rfa(bpm)= ";zlf,zrfa,"at
freqs(Hz):
          ";zfqlow,zfqres
          input "lfa,rfa,low freq: ",zlfan,zrfan,zfqlown
          if zlfan>30. then beep:goto 71 else zlf=zlfan
10         if zrfa>30. then beep:goto 71 else zrfa=zrfan
          if zfqlown<.02 or zfrlown>zfqres then beep:goto
71 else
          zfqlow=zfqlown
          locate 21,1
15         print "mean heart rate(bpm)= ";hrbpm
72         locate 22,1
          input "new mean heart rate(bpm): .",newhrbpm
          if newhrbpm>150 or newhrbpm<30 then beep:goto 72
else
20         hrbpm=newhrbpm
          'clear screen after input
          locate 17,1
          print space$(72)
          print space$(72)
25         print space$(72)
          print space$(72)
          print space$(72)
30         ' now compute values for hrsetup subroutine
          meandiv=76800#/hrbpm      '1280*60 ticks/min gives
          ticks/beat
          rfascal=76800#/(hrbpm-zrfa)-76800#/(hrbpm+zrfa)
          ' rfascal is the total excursion
35      of
          ' respiration
```

- 173 -

```
lfascal=76800#/(hrbpm-zlfa)-76800#/(hrbpm+zlfa)
      ' lfascal is the total excursion
      of low frequency
lowdiv=meandiv-(rfascal+lfascal)/2#
5
tbaserst=1280#/zfqlow
locate 17,1
print "tbaserst,rfascal,lfascal,lowdiv:
      ";tbaserst;rfascal;lfascal;
10 print lowdiv
call hrsetup(tbaserst,rfascal,lfascal,lowdiv)

return

15

      ' print out interrupt controller parameters
80 locate 10,1
mask=inp(&h21)
20 if (mask mod 16)<8 then mask=mask+8 else
mask=mask-8
      out &h21,mask
      mask=inp(&h21)
      print "8259 IMR(interrupt mask regisiter)=
25 ";mask;"
      =";hex$(mask)
      return

30

      ' this subroutine will change the analog buffer
90 locate 12,1
input "enter breathing rate (bpm): ",brate
if brate>75 or brate<7 then beep:goto 90
35 zfgres=brate/60#
numval=76800#/brate
```

- 174 -

```
        ztincr=8*ATN(1#)/numval
        locate 12,40
        color 31:print "calculating respiratory
signal...":color 7
5         call exstint          ' turn off interrupts
                                   while resetting buffer

        reseter=256
        call wrbuffer(reseter)
        for itime=0 to numval
10         ztnow=ztnow+ztincr
            analogval=127*(1#+SIN(ztnow))
            call wrbuffer(analogval)
        next itime
        call instint
15         locate 12,40
        print "respiratory signal active now      "
        return

20
        ' exstall the interrupt and print vector
9996      cls
            locate 4,1
            gosub 10000
25         call exstint
            locate 5,1
            gosub 10000
            locate 21,1
9999      stop

30
        ' subroutine to print out the interrupt vectors

10000     def seg=0
35         print "IRQ3 @0B*4H: ";hex$(peek(&h2C));"
            ";hex$(peek(&h2D));" ";
```

- 175 -

```
        print hex$(peek(&h2E));"  
";hex$(peek(&h2F));tab(40);  
        print "IRQ4 @0C*4H: ";hex$(peek(&h30));"  
        ";hex$(peek(&h31));" "  
5      print hex$(peek(&h32));" ";hex$(peek(&h33))  
        return  
  
        end  
  
10  
  
15  
  
20  
  
25  
  
30  
  
35
```

- 176 -

```

        page      66,80
; bdzint.asm - an assembler routine to handle interrupts
;               from IRQ3
; Last.revision: 1 April 1985
5      ;
;
;-----;
; 8088 interrupt location      ;
;-----;

10      abs0      segment at 0      ;absolute memory segment
                                   ;allows placement of
                                   ;interrupt address
                                   ;future timebase
15                                   ; interrupt handler
                                   ; resides at int 0B

      IRQ3_int    dw      2 dup(?);offset value is a word

                                   org      0CH*4      ;heart beat interrupt
                                   ;handler resides at int
20                                   ; 0C

      IRQ4_int    dw      2 dup(?);offset value is a word

      abs0      ends                ;

25

;-----;
; int_buffer: area to save DOS ;
30 ;      dummy interrupt ptr      ;
;-----;

      int_buffer  segment            ;data segment containing
35                                   ;user interrupt buffer

      save_int    dw      4 dup(?);offset for two DOS

```

- 177 -

```

;interrupts saved
;to be restored using
;exstint

5  int_buffer      ends      ;

;-----;
10 ; working storage for      ;
; time base interrupts      ;
;-----;

15 dseg_tbase      segment    ;data segment for timebase
; interrupt
    heartbeats      dw      ?  ;keep track of heart beats
; here (for debugging)
    base_rate       dw      ?  ;lowest divisor for heart
20 ; rate
    lfa_scal        db      ?  ;low frequency modulation
    rfa_scal        db      ?  ;high frequency modulation
    tbase_ctr       dw      ?  ;counter for timebase
; interrupt
25 ;(use for low frequency
; generation)
    tbase_rst       dw      ?  ;reset value for tbase_ctr
; used to set low frequency
    tbase_ptr       dw      ?  ;pointer to present analog
30 ; value
    tbase_len       dw      ?  ;length of analog data buffer
    tbase_buffer    db      2800dup(?) ;buffer for A/D values
    dseg_tbase      ends      ;

35

```

- 178 -

```

;-----;
; setup structures to allow access to;
; arguments passed by BASIC      ;
;-----;

5

; subroutine rdbeat(BASIC_beats)
frame_rd      struc      ;define the stack
                        ;structure for passing
10                        ;arguments to BASIC
savebp1        dw        ?      ;caller's base pointer
saveret1       dd        ?      ;return offset and
                        ;segment pushed by BASIC
BASIC_beats    dw        ?      ;place to return heart
15                        ;beats to BASIC
frame_rd      ends

;subroutine wrbuffer (analog)
frame_wr      struc      ;define the stack structure
20                        ; for passing
                        ;arguments from BASIC to
                        ; analog buffer
savebp2        dw        ?      ;caller's base pointer
saveret2       dd        ?      ;return offset and segment
25                        ; pushed by BASIC
analog         dw        ?      ;place to receive analog value
                        ; from BASIC
frame_wr      ends

30
;subroutine hrsetup(B_lreset,
; Brfa_scal,Blfa_scal,Bbase_
; rate)
frame_hr      struc      ;define the stack structure for
                        ; passing
35                        ;arguments from BASIC to heart
                        ; rate controls

```



- 179 -

```

    savebp3    dw    ?    ;caller's base pointer
    saveret3   dd    ?    ;return offset and segment pushed
                        ;    by BASIC
    Bbase_rate dw    ?    ;BASIC's lowest divider for heart
5      ;    rate
    Blfa_scal  dw    ?    ;BASIC's low frequency scaler
                        ;    (amplitude)
    Brfa_scal  dw    ?    ;BASIC's high frequency scaler
                        ;    (amplitude)
10     B_lreset dw    ?    ;BASIC's low frequency timer
                        ;    reset value

    frame_hr   ends

                        ;.....code segment begins here
15
    cseg_calibs    segment 'code'
    basic_dgroup   group  data,stack,const,heap,memory
                        ;defining link to BASIC
    porta          equ    0700H    ;port definitions for
20                ;8255 port expander
    portb          equ    0708H    ;these addresses are
                        ;decoded on the homemade
    portc          equ    0710H    ;board
    control        equ    0718H    ;control word in the
25                ;8255
    timer0         equ    0720H    ;8253 timer0 register
    timer1         equ    0721H    ;8253 timer1 register
    timer2         equ    0722H    ;8253 timer2 register
    con8253        equ    0723H    ;8253 control register
30

                        ;-----;
                        ; timebase interrupt handler (not accessible to;
                        ; BASIC)                                ;
35                ;-----;
                        ;this routine reads the A/D every timer0

```

- 180 -

```

;tick
;with the next point in the analog
;buffer

5
tbase_int    proc    far    ;this procedure is not
                    ;made public
                    assume  cs:cseg_sync,ds:dseg_
                        base,es:nothing,ss:nothing
10            push    ax    ;save registers used
                        ;during interrupt
                    push    bx    ;
                    push    dx    ;
                    push    ds    ;
15
                    mov     ax,dseg_base    ;set up segment
                        ;register for data area
                    mov     ds,ax    ;
20
                    ;.....increment counter used for
                        ;low frequency generation
                    dec     tbase_ctr    ;decrement
25                    ;interrupt counter
                    jnz     ctr_ok    ;if not zero then
                        ;continue
                    mov     ax,tbase_rst ;else reload reset
                        ;value
30                    mov     tbase_ctr,ax ;
                    ctr_ok:
                        ;.....get analog value from
                        ;buffer and send to DAC
35
                    mov     bx,tbase_ptr    ;get pointer to
```

- 181 -

```
                                ;analog data
dec    bx                      ;
mov     al,tbase_buffer[bx]    ;get analog
                                ;value

5
    mov     dx,porta           ;send analog value
                                ;to DAC
    out     dx,al              ;

10
    mov     dx,control         ;toggle the write
                                ;latch for the DAC
    mov     al,6               ;by using direct
                                ;bit reset
    out     dx,al              ;and
15
    inc     al                  ;reset commands
    out     dx,al              ;

    dec     tbase_ptr          ;point to next
                                ;value
20
    jnz     tbase_eoi          ;if zero, reset
                                ;pointer
    mov     ax,tbase_len       ;reset with buffer
                                ;length
    mov     tbase_ptr,ax       ;

25
    ;.....acknowledge interrupt to
    ;      8259A
tbase_eoi: mov     al,20H       ;send EOI to 8259A
    out     20H,al            ;

30
    pop     ds                 ;restore registers which
                                ;were used
    pop     dx                 ;
    pop     bx                 ;
35
    pop     ax                 ;
    iret                      ;return to place where
```



- 183 -

```

                                push    ax      ;save registers during
                                ;interrupt
                                push    bx      ;
                                push    cx      ;
5                                push    dx      ;
                                push    ds      ;

                                mov     ax,dseg_tbase ;set up segment
                                ;register for data area
10                               mov     ds,ax      ;

                                inc     heartbeats ;increment heart
                                ; beat counter

15                               ;.....calculate low frequency modulation
                                ; (the tbase buffer is used as a trig
                                ; table here)
                                mov     ax,tbase_ctr ;get number of 1280Hz
                                ;pulses
20                               dec     ax      ;
                                mul     tbase_len ;scale by length of
                                ; respiratory
                                ; buffer
                                div     tbase_rst ;divided by reset
25                               ;value to get
                                ; pointer
                                mov     bx,ax    ;to low frequency
                                ; modulation
                                mov     al,tbase_buffer[bx] ;get sinusoidal
30                               ; modulation
                                mul     lfa_scal ;and scale
                                ; appropriately
                                mov     cx,ax    ;cx accumulate
                                ;divider for 1280Hz
35                               clock
```

- 184 -

```

;.....calculate respiratory modulation
mov     bx,tbase_ptr      ;get present
                                ;respiration signal
mov     al,tbase_buffer[bx] ;from buffer
5      mul     rfa_scal      ;scale with rfa scaler
add     cx,ax              ;and add to cx

add     cx,base_rate       ;finally add base rate
                                ;to get
10      ; value for
                                ;timer1 (heart rate
                                ;generator on
                                ; 8253)

15      ;.....send new divider to 8253 timer
mov     al,76H              ;set timer 1 to square
                                ; wave
                                ; generator
mov     dx,con8253          ;
20      out     dx,al        ;

mov     dx,timer1           ;send divider to
                                ;timer1
mov     al,cl               ;low byte first
25      out     dx,al        ;
mov     al,ch               ;high byte next
out     dx,al              ;

;.....acknowledge interrupt to
30      ; 8259A
mov     al,20H              ;send EOI to 8259A
out     20H,al             ;

35      pop     ds           ;restore registers and
pop     dx                  ;
pop     cx                  ;
```

- 185 -

```

                                pop    bx      ;
                                pop    ax      ;
                                ired         ;return to place where
                                           ;interrupt occurred
5
                                debugmsg2      db      'this is the end of the heart
                                           beat interrupt'

10  hbeat_int      endp

                                ;-----;
15  ; subroutine instint (install_interrupts) ;
                                ;-----;

                                instint      proc    far
                                           public  instint
20  ;public symbol allows external references
                                ;es,ds used to access interrupt and must
                                ; be restored movsw
                                ;uses (ds:si)(es:di) addr
                                           assume  cs:cseg_calibs,ss:basic_
25  ;                                dgroup,ds:basic_dgroup
                                           assume  es:int_buffer

                                ;.....save registers
                                push    ds      ;save ds register on the
30  ;                                ; stack
                                push    es      ;save es register on the
                                ; stack

                                push    bp      ;save BASIC base pointer
35  ;                                ; for return to BASIC
                                mov     bp,sp    ;point stack pointer at

```

- 186 -

```
                                ;frame reference to
                                ;address of BASIC analog
                                ;data buffer

5          push    ax          ;save additional
                                ;registers
          push    si          ;
          push    di          ;

10         ;set up the segment registers as assumed

          mov     ax,int_buffer ;
          ;es points to buffer area to save
          ;DOS dummy interrupt vector
15         mov     es,ax        ;
          mov     ax,0          ;ds points to
                                ;abs0 (interrupt table)
          mov     ds,ax        ;
          assume  ds:abs0      ;

20         ;setup access to interrupt vectors
          lea     di,save_int   ;load offset of
                                ;save_int in es,di
          lea     si,IRQ3_int   ;load offset of
25         ;IRQ3_int in ds,si
          movsw                   ;save DOS dummy
                                ;interrupt vectors to be
          movsw                   ;restored later
          movsw                   ;now saving IRQ4
30         movsw                   ;

          ;install the DAC timebase (IRQ3)
          mov     IRQ3_int+2,cseg_calibs
35         mov     IRQ3_int,offset tbase_int;
          ;interrupt handler now
```



- 187 -

```
;install the heart beat (IRQ4) interrupt handler now
      mov     IRQ4_int+2,cseg_calibs;
      mov     IRQ4_int,offset hbeat_int;

5
      ;.....return to BASIC

      pop     di      ;restore additional
                      registers
10     pop     si      ;
      pop     ax      ;

      pop     bp      ;restore BASIC's base
                      ;pointer and
15     pop     es      ;segment registers
                      before returning
      pop     ds      ;
      ret     0        ;delete 0 parameters (0
                      ;bytes) from the stack
20                      ;and return to the
                      ;calling routine

      debugmsg3      db     'this is the end of the
                      interrupt installation'
25

      instint        endp

30
      ;-----;
      ; subroutine exstint (exstall_      ;
      ; interrupts)                      ;
      ;-----;
```

- 188 -

```

exstint      proc    far
              public  exstint ;public symbol allows
                                ;external references
              assume  cs:cseg_calibs,ss:basic_dgroup
5             assume  ds:int_buffer,es:abs0
              ;es,ds used to access interrupt
              ;vectors and must be restored
              ;movsw uses (ds:si)(es:di) addr

10            ;.....save registers

              push    ds        ;save ds register on the
                                ; stack
              push    es        ;save es register on the
15             ; stack
              push    bp        ;save BASIC base pointer
                                ; for return to BASIC
              mov     bp,sp      ;point stack pointer at
                                ; frame reference to
20             ;access arguments passed
                                ; by BASIC (none here)

              push    ax        ;save additional
                                ;registers
25             push    si        ;
              push    di        ;
                                ;set up the segment
                                ; registers as assumed
              mov     ax,0       ;es points to
30             ;abs0 (interrupt table)
              mov     es,ax      ;
              mov     ax,int_buffer ;ds points to
                                ;buffer area to save
              mov     ds,ax      ;DOS dummy
35             ;interrupt vector

```

- 189 -

```

                                ;setup access to interrupt vectors
lea    di,IRQ3_int             ;load offset of
                                ;IRQ3_int in es,di
lea    si,save_int             ;load offset of
5      ;save_int in ds,si
movsw                                     ;restore DOS
                                ;dummy interrupt vectors
movsw                                     ;for IRQ3
movsw                                     ;and IRQ4
10    movsw                                     ;

                                ;.....return to BASIC

15    pop    di                 ;restore additional
                                ; registers
pop    si                       ;
pop    ax                       ;

20    pop    bp                 ;restore BASIC's base
pop    es                       ;pointer and segment
pop    ds                       ;registers before
                                ;returning
ret    0                       ;delete 0 parameters (0
25    ;bytes) from the stack
                                ;and return to the
                                ;calling routine

debugmsg4    db    'this is the end of the
30            interrupt exstallation'

exstint      endp

35
```

- 190 -

```

;-----;
; subroutine rdbeat (read_heart_beats ;
;-----;
5
rdbeat      proc    far
              public rdbeat    ;public symbol allows
                                ;external references
10          assume cs:cseg_calibs,es:dseg_tbase
              assume ds:basic_dgroup,ss:basic_dgroup

              ;.....save registers
15
              push    bp        ;save BASIC base pointer
                                ;for return to BASIC
              mov     bp,sp      ;point stack pointer at
                                ;frame reference to
20                                ;access arguments passed
                                ;by BASIC (one here)

              push    ax        ;save additional
                                ;registers
25          push    es          ;
              push    di        ;

              mov     ax,dseg_tbase    ;set up segment
                                ;register for data area
30          mov     es,ax        ;

              mov     ax,heartbeats    ;get
                                ;beats from local memory
35          mov     di,[bp].BASIC_beats ;
              mov     [di],ax          ;send
```

- 191 -

;beats to BASIC

;.....return to BASIC

```
5
    pop    di        ;restore additional
                        registers
    pop    es        ;
    pop    ax        ;
10
    pop    bp        ;restore BASIC's base
                        ;pointer,
    ret     2         ;delete 2 parameters (4
                        ;bytes) from the stack
15
                        ;and return to the
                        ;calling routine

    debugmsg5        db    'this is the end of the heart
                        beat read routine'
20
    rdbeat    endp
```

```

;-----;
25  ; subroutine wrbuffer(analog)      ;
;-----;
```

```
    wrbuffer    proc    far
    public wrbuffer    ;public symbol allows
30
                        ;external references
    assume cs:cseg_calibs,es:dseg_tbase
    assume ds:basic_dgroup,ss:basic_dgroup
```

```
35
;.....save registers
```

- 192 -

```

                    push    bp      ;save BASIC base pointer
                                   ;for return to BASIC
                    mov     bp,sp   ;point stack pointer at
                                   ;frame reference to
5                    ;access arguments passed
                                   ;by BASIC (one here)

                    push    ax      ;save additional
                                   ;registers

10                   push    bx      ;
                    push    es      ;
                    push    si      ;
                    mov     ax,dseg_tbase ;set up segment
                                   ;register for data area

15                   mov     es,ax    ;

                    mov     si,[bp].analog ;get analog value
                                   ;from BASIC

                    mov     ax,[si]   ;
20                   test    ah,OFFH   ;if upper byte is
                                   ;zero
                    jz      new_buff  ;then install a
                                   ; new point in
                                   ; the buffer

25                   mov     tbase_len,0 ;otherwise reset
                                   ;the buffer

                    mov     tbase_ptr,1 ;
                    jmp     wr_ret     ;

30                   mov     bx,tbase_len ;get present
                                   ;pointer and
                                   ;use it
                    mov     tbase_buffer[bx],al ;to store
                                   ; buffer value

35                   inc     tbase_len ;point to next
                                   ;buffer value

```

- 193 -

```

;.....return to BASIC

5          pop     si      ;restore additional
                                ;registers
wr_ret:    pop     es      ;wr_ret:
          pop     bx      ;
          pop     ax      ;
10
          pop     bp      ;restore BASIC's base
                                ;pointer,
          ret     2        ;delete 1 parameters (2
                                ;bytes) from the stack
15                                ;and return to the
                                ;calling routine

debugmsg6  db        'this is the end of the buffer
                                write routine'
20
wrbuffer   endp

;-----;
25 ; subroutine hrsetup(B_lreset,Brfa_scal,Blfa_scal,
; Bbase_rate)
;-----;

proc      far
30 public hrsetup ;public symbol allows
                                external references
assume cs:cseg_calibs,es:dseg_tbase
assume ds:basic_dgroup,ss:basic_dgroup

35
;.....save registers

```

- 194 -

```

5      push    bp                ;save BASIC base
                                   ;pointer for return
                                   ;to BASIC
      mov     bp,sp            ;point stack pointer
                                   ;at frame
                                   ;reference to
                                   ;access arguments
                                   ;passed by BASIC
10     ;(one here)

      push    ax                ;save additional
                                   ;registers
      push    es                ;
15     push    si                ;

      mov     ax,dseg_tbase     ;set up segment
                                   ;register for
                                   ;data area
20     mov     es,ax            ;

      mov     si,[bp].Bbase_rate ;get lowest
                                   ;divisor for heart
      mov     ax,[si]           ;rate from BASIC
25     mov     base_rate,ax      ;and save in local
                                   ;      data
                                   ;      segment

      mov     si,[bp],Blfa_sacl ;get low freq
30     ;                        ;      modulation
                                   ;      scale
      mov     ax,[si]           ;      from BASIC
      mov     lfa_scal,al       ;and save LSbyte in
                                   ;local data
35     ;                        ;      segment

```



- 195 -

```

      mov     si,[bp].Brfa_scal  ;get high freq
                                ; modulation scale
      mov     ax,[si]            ;from BASIC
      mov     rfa_scal,al        ;and save
5      ;LSbyte in local data
      ;segment
      mov     si,[bp].B_lreset   ;get low freq
                                ; timer reset value
      mov     ax,[si]            ;from BASIC
10     mov     tbase_rst,ax       ;and save in
                                ; local data segment

      ;.....return to BASIC

15     pop     si                ;restore additional
                                ;registers
      pop     es                ;
      pop     ax                ;

20     pop     bp               ;restore BASIC's base
                                ;pointer,
      ret     8                 ;delete 4 parameters (8
                                ; bytes) from the stack
                                ;and return to the
25     ; calling routine

      debugmsg 7 db 'this is the end of the heart rate
                        setup routine'

30     hrsetup endp

      cseg_calibs ends

      end

35
```

- 196 -

5

10

15

20

25

30

35

- 197 -

## APPENDIX B

1985 - Makoto R. Arai  
Laura E. McAlpine, and  
Daivd Gordon

5

10

15

20

25

30

35

```
' SYNCTS19 - program to test synchronous data
'           acquisition and also
'           test asynchronous processing using
'           board#2
'           addition: asynchronous data
'           archiving (poll driven)
'           reviewing old data
' last revision: 15 May 1985
'
' REQUIRED SUBROUTINES: <MODULE>
'
'           instint(fdbuf1ptr,fdbuf2ptr,fdbuf3ptr)
'           <SYNC7S>
'           exstint                                <SYNC7S>
'           rdbeat(heart,sync)                    <SYNC7S>
'           rdbuf(dataptr,bufferno)                <SYNC7S>
'           rdptrs(adrd,hbrd,adflag,hbflag) <SYNC7S>
'
'           swindow(xmins,xmaxs,ymins,ymaxs)
'           <GWINDOW3>
'           dwindow(xmind,xmaxd,ymind,ymaxd)
'           <GWINDOW3>
'           clrwindw                                <GWINDOW3>
'           axes                                    <GWINDOW3>
'           scaler(dataptr,gdataptr,numval)
'           <GWINDOW3>
'
'           fgraph(dataptr,numval,xnow,linemask)
'           <FGRAPH8>
'           [for scaled graphs, use
```

- 198 -

```

'          xnow=xmins,
'          numval=numvalg=xmaxs-xmins+1,
'          and gdataptr]
'          dumpgr   [to dump graphics]      <DUMPGR>
5
'

defint a-y          ' only z denotes a real number
defdbl z
10  dim zreal(514),zrimag(514),zdata(1025)
    dim ydata(1025),ydatag(1025)
    dim hb1(1025),hb2(1025),zhr(1025)
    dim zspec.hb.real(512),zspec.hb.imag(512)
    dim sresetval(5),resprstval(5)
15  dim linetype(3),histogram(100)
    def fnzmag(z1,z2)=z1*z1+z2*z2
    def fnzcoher(zr1,zil,zr2,zi2)=fnzmag
        (zr1*zr2+zil*zi2,zil*zr2-zr1*zi2)

20
    ' initialize timer reset values
1    sval=27 : for i=1 to 5 sresetval(i)=sval :
        sval=sval+sval : next i
2    sval=1381 : for i=0 to 3 : resprstval(i)=sval :
25    sval=sval+sval : next i
3    resprstval(4)=sval

    ' define fft parameters
30  4    fftsize=1024 : npair=fftsize/2 :
        znpair=cdbl(npair) : lpower=9

5    for i=0 to 514 : zreal(i)=0# : zrimag(i)=0# :
        next i
35
    datacycle=0

```

- 199 -

```
        ' flag for automatic fft: when non-zero,
        ' marks stage of data
        ' processing (semi asynchronous)
cyclewait=0
5      ' define linetype for plots
      linetype(0)=&HFFFF
      linetype(1)=&HAAAA
      linetype(2)=&HCCCC
      linetype(3)=&HFAFA
10     req.cls=0
      sounder=1

      'define ports on 8253
15     timer0=&h704
      timer1=&h705
      timer2=&h706
      con8253=&h707

20     'define ports on 8255
      porta=&H71C
      portb=&H71D
      portc=&H71E
25     con8255=&H71F

      ' set up sampling rate for heart rate timer and
      ' respirations
30     gosub 100

      ' first set 8255 ports A,C to output, port B to
35     ' input
      ' turn the gates on using the 8255 at bits 0,1,2
```

- 200 -

```

'   on portc
'   by setting portc to 1FH
'   this also selects channel 0 for the A/D
out con8255,130
5 out portc,&H1F

' now go through required startup subroutines to
10 '   set up data archives
    open "R",1,"resp.dat",2048
    open "R",2,"hb1.dat",2048
    open "R",3,"hb2.dat",2048
    open "R",10,"trends.dat",128
15
31 field #1,2048 as analog$
    field #2,2048 as fdhb1$
    field #3,2048 as fdhb2$
    field #10,128 as trends$
20
    fdflag=0
    fdrecord=1
    record1no=0 : record2no=0 : record3no=0 :
        record10no=0
25 adflag1st=0 : hbflag1st=0
    fdbuf1ptr=varptr(#1)+188      ' set up
        'pointers to disk buffers
    fdbuf2ptr=varptr(#2)+188
    fdbuf3ptr=varptr(#3)+188
30
        '.....field definitions for
            trend data file
    field #10,8 as hr$,8 as rr$,8 as rcf$,8
        as lfa$,8 as rfa$,8 as coher$
35 field #10,48 as dummy1$,8 as ratio$,8
        as cratio$,8 as hrintegral$

```

- 201 -

```
field #10,72 as dummy2$,8 as respintegral$,8
  as timestamp$
field #10,120 as dummy3$,2 as hbrecord$,2
  as adrecord$
5 field #10,124 as dummy4$,2 as hbeat$,2
  as samprate$

' first print out the present value of the
10 ' interrupt vectors
locate 23,1 : gosub 20000
gosub 19000

15 ' make sure interrupts are off before installing
' handlers
mask=inp(&h21) : mask=mask or 24 : out &h21,mask

' install the interrupts
20 call instint(fdbuf1ptr,fdbuf2ptr,fdbuf3ptr)
locate 24,1 : gosub 20000
gosub 19000

25 ' turn interrupts back on
mask=inp(&h21) : mask=mask and &h0e7 : out
&h21,mask

30 40 locate 1,1 : gosub 20000
print "commands: c(ounts), f(ft), g(raph),
i(in on), q(uit), r(beats)";
print "s(tore), x(cls), #(samples);

35 ' wait until user hits a key
41 savekey$=""
```

- 202 -

```
while len(savekey$)=0 and datacycle<=0
    savekey$=savekey$+inkey$:gosub 30000:locate
    24,70:print time$;wend

5      while datacycle=1
        fdrecord=recordlno : fdflag=1
        'set up future A/D analysis
        analrec.ad=recordlno : analrec.hr=record2no+1
        if req.cls=1 then cls : req.cls=0
10      'clear screen if needed

        gosub 950 '.....analyze heart rate
42      hrspecsum#=zspectsum*2#

15      gosub 900 '.....analyze A/D data (from floppy
43      respspecsum#=zspectsum*2#

        gosub 15000
        ' calculate spectral amplitudes
20      gosub 16000

        save trend data

        datacycle=cyclewait : wend
        'end auto data analysis cycle
25

49

30      if savekey$="c" then gosub 60
        ' print timer counts
        if savekey$="f" then gosub 900
        ' fft A/D buffer contents
        if savekey$="F" then gosub 950
35      ' fft heart rate buffer contents
        if savekey$="g" then gosub 12700
```



- 203 -

```
' graph current A/D buffer
if savekey$="G" then gosub 12710
' graph current heart rate buffer
if savekey$="h" then gosub 90
5  ' (no) plot histogram
   if savekey$="p" then gosub 91
   ' (no) print trends
   if savekey$="i" then gosub 80
   ' unmask interrupt 3
10  if savekey$="I" then gosub 81
   'unmask interrupt 4
   if savekey$="q" then goto 9996
   ' quit
   if savekey$="r" then gosub 50
15  ' print heart beats
   if savekey$="S" then gosub 800
   ' analyze data in disk file (set fdflag)
   if savekey$="t" then gosub 16500
   ' print out the trends
20  if savekey$="x" then cls          'clear screen
   if savekey$="#" then gosub 100
   ' reset sampling rate
   if savekey$="?" then gosub 700  'help
   savekey$=""
25
   goto 41

'print present value of heartbeats
30
50  locate 24,1 : gosub 20000
   call rdbeat(heart, sync)
   print "heart beats: ";heart,"sync pulses:
      ";sync;time$;
35  return
```

- 204 -

```
' print present value of counters
60 out con8253,0          'latch timer0
   tlow0=inp(timer0)
5   thigh0=inp(timer0)
   out con8253,&h40        'latch timer1
   tlow1=inp(timer1)
   thigh1=inp(timer1)
   out con8253,&h80        'latch timer2
10  tlow2=inp(timer2)
   thigh2=inp(timer2)
   locate 24,1 : gosub 20000
   print "timer0: ";tlow0+thigh0*256;tab(20);"
     timer1: ";tlow1+thigh1*256;
15 61 print tab(40);"timer2: ";tlow2+thigh2*256#;
   return

20 ' print out interrupt controller parameters:
   ' entry point for IRQ3
80 mask=inp(&h21) : mask=mask xor 8 : out &h21,mask
   goto 82
   ' entry point for IRQ4
25 81 mask=inp(&h21) : mask=mask xor 16 : out
     &h21,mask
82 mask=inp(&h21)
   locate 24,1 : gosub 20000
   print "8259 IMR(interrupt mask regsiter)=
30     ";mask;" =";hex$(mask);
   return

35 ' (re)set sampling rates
   ' set timer0 to 16 bit square wave rate
```

- 205 -

```
' generator mode
' set timers 1,2 to 16 bit rate generator mode
100 out con8253,&h36
    out con8253,&h74
5    out con8253,&hB4

    '.....set real time multiplier
105 locate 23,1 : gosub 20000
    input "real time multiplier: ",rt.mult
10    rt.multqual=0
    if rt.mult=1 then rt.multqual=1
    if rt.mult=2 then rt.multqual=2
    if rt.mult=4 then rt.multqual=3
    if rt.mult=8 then rt.multqual=4
15    if rt.multqual<>0 then goto 110
    beep : goto 105

    ' get heart rate resolution desired to reset
20    ' timer0 reset value
110 locate 1,1 : gosub 20000
    input "heart rate resolution: (11,23,45,91,181
        usec) ",hrresol

25    ' check heart rate resolution validity
    hrqual=0
    if hrresol=11 then hrqual=1
    if hrresol=23 then hrqual=2
    if hrresol=45 then hrqual=3
30    if hrresol=91 then hrqual=4
    if hrresol=181 then hrqual=5
    if hrqual<>0 then sreset=sresetval(hrqual) :
        goto 120
    beep : goto 110
35    ' invalid heart rate resolution
```

- 206 -

```

5      '      sreset=27      'set timer 0 to 88384Hz
      '      'timebase (11.3 usec res
      '      '(2.38MHz/27)(max resp
      '      'samples then 64Hz)

10     '      sreset=54      'set timer 0 to 44192Hz
      '      'timebase (22.6 usec res
      '      '(2.38MHz/54)(max resp
      '      'samples then 32Hz)

15     '      sreset=108     'set timer 0 to 22096Hz
      '      'timebase (45.3 usec res
      '      '(2.38MHz/108)(max resp
      '      'samples then 16Hz)

20     '      sreset=216     'set timer 0 to 11048Hz
      '      'timebase (90.5 usec res
      '      '(2.38MHz/216)(max resp
      '      'samples then 8Hz)

25     '      sreset=432     'set timer 0 to 5524Hz
      '      'timebase (181 usec res
      '      '(2.38MHz/432)(max resp
      '      'samples then 4Hz)

120    '.....set respiratory sampling rate
30    locate 2,1 : gosub 20000
      print "respiratory sampling rate: ( 4";
      twopwr=4 : for i=hrqual+rt.multqual to 5 :
      twopwr=twopwr+twopwr
      print using ",##";twopwr; : next i : print "
      Hz) ";
35    input respsampl

```

- 207 -

```
' check respiratory sampling rate validity
respqual=0 : respsampl.eff=respsampl*rt.mult
if respsampl=4 then respqual=1
if respsampl=8 then respqual=2
5   if respsampl=16 then respqual=3
    if respsampl=32 then respqual=4
    if respsampl=64 then respqual=5
    if respqual=0 or respqual+hrqual+rt.multqual>7
      then beep : goto 120
10
    resprst=resprstval(7-hrqual-respqual-
      rt.multqual)

      '.....set cycle delay time between
15   '      analyses
    130 locate 3,1 : gosub 20000
        input "waiting time between cycles: ",dropcycle
        if dropcycle<0 or dropcycle>5 then beep : goto
          130
20   cyclewait=0-dropcycle

    out timer0,(sreset mod 256)
    ' system timebase generated here out
    '   timer0,(sreset\256)
25
    out timer1,(resprst mod 256)
    ' timer 1 counts timebase and outputs out
    '   timer1,(resprst\256)
    ' the respiratory sampling rate
30
    out timer2,0
    ' set timer 2 as an overflow counter for the
    '   out timer2,0
    ' number of overflows (65536 counts)
35
```

- 208 -

```

200      timer2over#=65536#
          ' overflow value for timer2
201      zlover=resprst          ' reset count for timer1
202      zlfreq=14318180#/6#/sreset
5        ' timer1 input clock frequency
203      zhrsampler=zlover/zlfreq
          ' timer1 output=sampling interval
204      segment.time=fftsize*zhrsampler
205      zlfreq.real=zlfreq/rt.mult
10       ' real time used to calculate HR
206      zhrsampler.real=zlover/zlfreq.real

          '.....respiratory peak search
          '                parameters
15      210      minrespfrq#=.2#
          ' start at frequency (in pixels)
211      minresp=minrespfrq#/respsampl*1024
212      combwidth#=.032#
          use comb tooth width (in pixels)
20      213      combpix=combwidth#/respsampl*1024
214      if combpix<=0 then combpix=0

          '.....low frequency
          '                peak/integration parameters
25      220      pixel.04=cint(40.96#/fftsampl)+1
          ' pixel for .04Hz
221      pixel.10=cint(102.4#/fftsampl)+1
          ' pixel for .10Hz
222      fft.expansion=respsampl/fftsampl
30
          if datacycle=0 then datacycle=-1
          if recordlno=0 then return
          ' on startup don't delay
          ' exclude the current data segment
35      ' from analysis since changes in
          ' sampling rate will introduce glitches

```

- 209 -

```

    return

5
    ' set floppy disk flag (fdflag) to analyze data
    ' stored on floppy (resp)
800    fdflag=1
    locate 23,1 : gosub 20000 : input "record
10    number: ",fdrecord
    if fdrecord>=1 and fdrecord<=recordlno then
        gosub 12700 : return
    locate 24,1 : gosub 20000 : beep : print
        "invalid record number";
15    return

    ' set up data for fft here
20    ' get analog data from the A/D
900    gosub 12700      ' get analog data and plot
901    for i=1 to fftsize : zdata(i)=ydata(i) : next i
902    locate 23,1 : gosub 20000 : print "A/D buffer is
        transformed";

25

    xmin=330 : xmax=630 : ymin=102 : ymax=167
    call swindow(xmin,xmax,ymin,ymax)

    glabel=3      ' plot label is "resp spect"
30    gosub 10000  ' fft
    return

    ' get heart rate data for fft
950    locate 23,1 : gosub 20000 : print "heart rate is
35    transformed";
951    gosub 12710  ' get hr function and plot it
```

- 210 -

```
952      for i=1 to fftsize : zdata(i)=zhr(i) : next i

953      xmins=330 : xmaxs=630 : ymins=28 : ymaxs=93
954      call swindow(xmins,xmaxs,ymins,ymaxs)

5
955      glabel=4          ' plot label is "hr spect"
956      gosub 10000        ' fft

          ' save spectrum in spec.hr buffers
10 960      for i=0 to 512
961          zspec.hb.real(i)= zreal(i) :
          zspec.hb.imag(i)=zrimag(i)
962      next i

15      return

20
          ' exstall the interrupt and print vector
9996      cls
          ' make sure interrupts are off before
          removing handlers
25      mask=inp(&h21) : mask=mask or 24 : out &h21,mask

          ' remove interrupt handlers
          screen 0

30      locate 4,1
          gosub 19000
          call exstint
          locate 5,1
          gosub 19000
35      locate 21,1
```



- 211 -

```
' close files after storing last bit of data
bufferno=0
call rdbuf(fdbuf1ptr,bufferno)
5 put #1,record1no+1
bufferno=1
call rdbuf(fdbuf2ptr,bufferno)
put #2,record2no+1
bufferno=2
10 call rdbuf(fdbuf3ptr,bufferno)
put #3,record3no+1

close #1,#2,#3,#10

15 ' and quit
9999 stop

20

' FFT ROUTINE
'
25 ' set up the data
'
10000 zreal(0)=0#
10001 zrimag(0)=0#
10002 zreal(npair+1)=0#
30 10003 zrimag(npair+1)=0#

' compute mean value of array
10004 zmean=0#
10005 for i=1 to fftsize : zmean=zmean+zdata(i)
35 : next i
10006 zmean=zmean/1024#
```

- 212 -

```

10007   for k=1 to npair : j=k+k-1 : zreal(k)=zdata
        (j)-zmean
5  10008   zrimag(k)=zdata(j+1)-zmean : next k

10009   ' locate 24,1 : gosub 20000
10010   ' print "arrays initialized at
        ' ";time$;space$(20);
10
        '
        ' fft routine <fftandift> begins here
        '
10011   ' locate 24,1 : print "entering fft routine at
15      ' ";time$;space$(20);
10012   k=0
10013   for j=1 to npair-1 : i=2
10014       ndivi=npair/i
10015       if k<ndivi then 10017
20  10016       k=k-ndivi : i=i+i : goto 10014
10017       k=k+ndivi
10018       if k<=j then 10025
10019       za=zreal(j+1)
10020       zreal(j+1)=zreal(k+1)
25  10021       zreal(k+1)=za
10022       za=zrimag(j+1)
10023       zrimag(j+1)=zrimag(k+1)
10024       zrimag(k+1)=za
10025   next j
30  10026   ' locate 24,1:print "bit reversal completed at
        ' ";time$;space$(20);

10030   g=1 : zp=1#
10031   for i=1 to lpower : gosub 30000
35      'check if disk requires service
10032   'locate 24,1:print "entering stage ";g;" at

```

- 213 -

```

      time ";time$;space$(20);
10033   if i=1 then zsign=-1# else zsign=1#
10034   zc=1# : ze=0#
10035   zq2=(1#-zp)/2# : if zq2<=0# then zq=0# : else
5       zq=sqr(zq2)
10036   zp2=(1#+zp)/2# : if zp2<=0# then zp=0# : else
       zp=zsign*sqr(zp2)
10037   itwog=g+g

10 10040   for r=1 to g
      10041   for j=r to npair step itwog
            k=j+g : if k>npair then print "k j g
            over>> ";k;j;g
10042       za=zc*zreal(k)+ze*zrimag(k)
15 10043       zb=ze*zreal(k)-zc*zrimag(k)
10044       zreal(k) =zreal(j) -za
10045       zrimag(k)=zrimag(j)+zb
10046       zreal(j) =zreal(j) +za
10047       zrimag(j)=zrimag(j)-zb
20 10048       next j
10049       za=ze*zp+zc*zq
10050       zc=zc*zp-ze*zq
10051       ze=za
10052       next r
25 10053       g=itwog
10054   next i
10055   'locate 24,1:print "entering final stage at
      ";time$;space$(20);
10056   gosub 30000
30   ' check if disk requires service

10060   za=4#*atn(1#)/znpair
10061   zp=cos(za)
10062   zq=sin(za)
35 10063   za=zreal(1)
10064   zreal(1)=za+zrimag(1)

```

- 214 -

```
10065  zrimag(1)=za-zrimag(1)
10066  zreal(1)=zreal(1)/2#
10067  zrimag(1)=zrimag(1)/2#
10068  zc=1# : ze=0#

5
10070  j=2
10071  while j<npair/2
10072      za=ze*zp+zc*zq
10073      zc=zc*zp-ze*zq
10  10074      ze=za
10075      k=npair-j+2
10076      za=zreal(j)+zreal(k)
10077      zb=(zrimag(j)+zrimag(k))*zc-(zreal(j)-
          zreal(k))*ze
15  10078      zu=zrimag(j)-zrimag(k)
10079      zv=(zrimag(j)+zrimag(k))*ze+(zreal(j)-
          zreal(k))*zc
10080      zreal(j)=(za+zb)/2#
10081      zrimag(j)=(zu-zv)/2#
20  10082      zreal(k)=(za-zb)/2#
10083      zrimag(k)=- (zu+zv)/2#
10084      j=j+1 : wend
10085      zrimag(npair/2+1)=-zrimag(npair/2+1)

25  10090  for j=2 to npair
10091      zreal(j)=zreal(j)/znpair/2#
10092      zrimag(j)=zrimag(j)/znpair/2#
10093  next j
10094  zreal(1)=zreal(1)/znpair
30  10095  zrimag(1)=zrimag(1)/znpair

,
' fft routine now completed
,

35  10100  locate 24,1:print "fft completed
          ";time$;space$(20);
```

- 215 -

```

      '...integrate spectrum
      '      sum up the spectrum noting that only the
      '      first npair elements of
5      '      the fft are valid
      '      (npair+1 to fftsize are complex conjugates
      '      of 1 to npair and are
      '      not calculated)
10101  zspectsum=0#
10      10102  zsummax=0#
10103  ipeak=-1
10104  for i=1 to npair
10105      zadd=fnzmag(zreal(i),zrimag(i))
10106      zspectsum=zspectsum+zadd
15      10107      if zadd<=zsummax then 10110
10108          zsummax=zadd
10109          ipeak=i
10110  next i

20
      '
      ' graphing routine for fft spectra
      '
10111  'locate 1,1 : gosub 20000
25      10113  'print "total spectral weight
      '      <variance>:";zspectsum*2#;
10114  'locate 2,1 : gosub 20000
10115  'print "peak weight : ";zsummax;" peak
      '      frequency= ";
30      10116  'print (ipeak-1#)/fftsize*respsampl;

      10117  gosub 12730
      ' fgraph of spectrum
10118  return
35
```

- 216 -

```
5      '-----'
      ' UTILITY ROUTINES HERE '
      '-----'

      ' graphing routine: gets data from A/D buffer
      ' and displays graph
12700 glabel=1
10    numpts=fftsize
      indata=0
      ' local flag
      ' indicating data is read while indata=0 and
      ' fdflag=0
15    dataptr=varptr(ydata(1))
      bufferno=0      'read A/D buffer
      call rdbuf(dataptr,bufferno)
      indata=1
      wend

20    while indata=0 and fdflag=1
      gosub 30000
      ' check file buffer to see if service is
      ' required
25    get #1,fdrecord
      for i=1 to 1024 :
      ydata(i)=cvi(mid$(analog$,i+i-1,2)) : next i
      indata=1
      wend

30

      xmin=10 : xmax=310 : ymin=102 : ymax=167
      call swindow(xmin,xmax,ymin,ymax)

35    xmind=0 : xmaxd=300 : ymind=0 : ymaxd=255
      call dwindow(xmind,xmaxd,ymind,ymaxd)
```

- 217 -

```
      ' max A/D value is 255
      call clrwindw
      call axes
      goto 12770

5

      ' entry point for plot of heart rate function
12710 screen 2          ' get heart rate function
12711 glabel=2
10 12712 numpts=fftsize
12713 gosub 13000
12714 ibeg=adrd+2
12715 for i=1 to fftsize : if ibeg=i then
      ibeg=ibeg+fftsize
15 12716 ydata(i)=cint(zhr(ibeg-i)) : next i

      xmins=10 : xmaxs=310 : ymins=28 : ymaxs=93
      call swindow(xmins,xmaxs,ymins,ymaxs)

20      xmind=0 : xmaxd=300 : ymind=0 : ymaxd=250
      call dwindow(xmind,xmaxd,ymind,ymaxd)
      ' max hr is 250 bpm

      goto 12770

25

      ' entry point for plotting spectra (screen
      ' windows already setup)
12730 zgain=250#/zsummax
30 12731 for i=1 to npair
12732 ydata(i)=cint(zgain*fnzmag
      (zreal(i),zrimag(i))) +1
12733 next i
12734 numpts=npair

35

      ' max spectral element (scaled to 250)
```

- 218 -

```
xmind=0 : xmaxd=300 : ymind=0 : ymaxd=255
call dwindow(xmind,xmaxd,ymind,ymaxd)

12770  call clrwindw
5      call axes

12780  dataptr=varptr(ydata(1))
      gdataptr=varptr(ydatag(1))
      call scaler(dataptr,gdataptr,numpts)
10     'correctly selects screen width

      ' entry point for plot of ydatag(i)
12790  x=xmins
      numvalg=xmaxs-xmins+1
15     linemask=&hffff
      gdataptr=varptr(ydatag(1))
      call fgraph(gdataptr,numvalg,x,linemask)

      ' graph labels printed here
20     on glabel goto 12800,12810,12820,12830
      return 'invalid label

      ' respirations in time domain
12800  if fdflag=1 then locate 14,30 : print
25     "rec#";fdrecord : fdflag=0
      return

      ' heart rate in time domain
12810  locate 5,3
30     print using "HR= ### bpm";cint(zavghr)
      return

      ' respiratory spectrum
12820  locate 14,63 :      print " Resp Spect ";
35     locate 15,63 : print using " (0-
      ##Hz)";respsampl\2
```



- 219 -

```
gosub 14000
' respiratory rate from spectrum by comb method
locate 14,3
' print respiratory rate with time tracing
5 print using "RR=### bpm
   (rcf=#.###)";cint(respfreq#*60),respcombfrac#
return

' heart rate spectrum
10 12830 locate 4,63 :      print " HR Spect ";
      locate 5,63 : print using " (0-##Hz)";fftsampl\2
      return

15
      'heart rate functions:
      ' read times from memory
      ' convert to heart rate function
      ' FFT resulting buffer
20      ' display the spectral amplitudes

13000 call rdptrs(adrd,hbrd,adflag,hbflag)
13002 if record2no=0 then startup=1 else startup=0
      startup is special

25
13003 hbptr1=varptr(hb1(1))
13004 bufferno=1
      'read heart beat buffer 1 (least sig. cts
13005 call rdbuf(hbptr1,bufferno)
30 13006 locate 24,1 : gosub 20000
13007 print "hbrd= ";hbrd; : anal.beat=hbrd

13008 hbptr2=varptr(hb2(1))
13009 bufferno=2
35      'read heart beat buffer 2 (most sig. cts
13010 call rdbuf(hbptr2,bufferno)
```

- 220 -

```
13011  for i=0 to 100 : histogram(i)=0 : next i
      'initialize histogram for deglitching (.4-40Hz)
13012  histomax#=zlfreq.real*2.5#
5  13013  histoscal#=zlfreq.real/40#

      ' compute time differences for entire hb array
      ' and save in zdata
      ' from the top down
10      ' zdata will contain the latest hr intervals,
      ' with the latest in
      ' (hbrd) and older intervals for decreasing
      ' array index
      ' since the timers are decrementing,
15      ' lstbeat<thisbeat
      ' (lstbeat is later, therefore smaller)
      ' this relation fails whenever there is a carry
      ' over (timer overflow)
      ' note: timer1 overflows exactly fftsize times
20      ' during one data segment
13020  lstbeat#=hbl(hbrd) : lstover#=hb2(hbrd)
13022  hbnow=hbrd-1
13023  if hbnow<=0 then hbnow=fftsize
13024  if startup=1 and hbnow=fftsize then return ' no
25      data yet

13025  numint=1
      ' valid intervals only (1 less than
      ' buffer size
30 13026  while numint<fftsize
13027      thisbeat#=hbl(hbnow)
      ' check for overflow of overflow counter
13028      thisover#=hb2(hbnow)
13029      if hb2(hbnow)<cint(lstover#) then
35          lstover#=lstover#-timer2over#
13030      hbnow=hbnow-1
```

- 221 -

```

13031      if hbnow=0 then hbnow=fftsize
13032      if hbnow=fftsize and startup=1 then goto
           13048
13033      zdatnow=thisbeat#-lstbeat#+overdif#*zlover
5  13034      if zdatnow>=0 then goto 13047 '?error

13040      if zdatnow>histomax# then goto 13044
13041      index=cint(zdatnow/histoscal#)
13042      histogram(index)=histogram(index)+1
10  13043      goto 13045
           'keep histogram of intervals (.2-20Hz:
           '  give 10% resolution @2Hz) extended
           '    data lapses
13044      histogram(100)=histogram(100)+1
15      ' extended data lapses

13045      zdata(numint)=zdatnow : numint=numint+1
13046      lstbeat#=thisbeat# : lstover#=thisover#
13047      wend
20  13048      numint=numint-1

           '.....find the interval
           '      corresponding to mean heart rate
           '          1) find largest peak in
25      '              .5-4Hz (2 pixels wide)
           '          2) calculate corrected
           '              mean interval
           '          3) calculate corrected
           '              interval variance
30      '          4) set slewing
           '              parameters for HR
           '              generation

13050      lstint=histogram(4) : hpeak=0 : hpeak.ht=0
35  13051      for i=3 to 40 : thisint=histogram(i)
13052      if (thisint+lstint)>hpeak.ht then

```

- 222 -

```

        hpeak.ht=thisint+1stint : hpeak=i
13053      1stint=thisint : next i
13054      approx.avg#=(hpeak-0.5#)*histoscal#

5  13060      zhistsum=0# : zhistsum2=0#
    13061      for i=1 to numint :
        index=cint(zdata(i)/approx.avg#)
13062          if index<=0 then index=1
13063          zhistsum=zhistsum+zdata(i)/index : next i
10 13064      avgint#=zhistsum/numint

    13070      for i=1 to numint : index=cint(zdata(i)/avgint#)
13071          if index<=0 then index=1
13072          zdif=zdata(i)/index-avgint# :
15      zhistsum2=zhistsum2+zdif*zdif
    13073      next i
13074      histvar#=zhistsum2/numint

        ' calculate deglitching parameters
20 13081      varslew#=31.4#*sqr(histvar#)/respsampl
        '5x max slew (1Hz rfa) slew at least .05Hz
        ' (3bpm)/beat infslew has infimum of slew
        '      maxima
    13082      min.maxslew#=.05
25 13083      infslew#=1#/(1#/avgint#-
        min.maxslew#/zlfreq.real)-avgint#
    13084      if maxslew#<infslew# then maxslew#=infslew#
    13085      supslew#=avgint#/5#
        'never slew more than 20% HR
30 13086      if maxslew#>supslew# then maxslew#=supslew#
    13087      locate 1,1 : gosub 20000 ' : print "maxslew:
        ";maxslew#

        ' compute heart rate waveform next
35 13100      ztime=0#
        ' time for present heart rate signal

```

- 223 -

```
      ' pointer in zdata to present beat number
      ' of beats accepted
13101  intnow=1
13102  beatno=1 :
5  13103  while zdata(intnow)<=0
13104      intnow=intnow+1 : if intnow>numint then goto
          13140 : wend
13105  zintlst=avgint# : zdropper=avgint# :
          zintnow=zdata(intnow)
10  13106  znext=zintnow/zlfreq.real
      ' time of previous heart beat deglitch first
      ' beat present heart rate keep statistics for
      ' deglitching sampling rate determined by
      ' timers
15  13107  avgnow#=avgint# : gosub 13500
13108  zhrnow=60#*zlfreq.real/zintnow
13109  zsum=zhrnow
13110  zsum2=zhrnow*zhrnow
13111  zincr=zhrsampler.real
20  13120  numsig=1
      ' point to heart rate function

13121  while numsig<=fftsize and ztime<=znext
13122      zhr(numsig)=zhrnow : numsig=numsig+1 :
25  13123      ztime=ztime+zincr
          wend:zintlst=zintnow
13124  if numsig=fftsize+1 then goto 13142
13125      intnow=intnow+1 : if intnow>numint then goto
          13140
30  13126      zintnow=zdata(intnow) : if zintnow<=0 then
          goto 13125
13127      znext=znext+zintnow/zlfreq.real : gosub
          13500      ' deglitcher
13128      zhrnow=60#*zlfreq.real/zintnow
35  13129      zsum=zsum+zhrnow : zsum2=zsum2+zhrnow*zhrnow
          : beatno=beatno+1
```

- 224 -

```
13130      goto 13121

13140      zavghr=zsum/beatno
          ' averaged over number of beats
5  13141      while numsig<=fftsize : zhr(numsig)=zavghr :
          numsig=numsig+1 : wend
13142      zavghr=zsum/beatno

13400      locate 24,13 : print "  avg hr(bpm): ";zavghr;
10
          ' zhr now has heart rate function
13401      print " ...heart rate function computed";

          return
15
          ' deglitching of three types employed here:
          '      correction of premature triggers (not
          '      yet)
          '      correction of dropped beats (not yet)
20          '      slew rate limiting of final output (a
          '      crude bandlimiter)

13500      if abs(zintnow-zintlst)<maxslew# then return
          'check for dropped beats
13501      numdrop=cint(zintnow/avgnow#) : if numdrop<=0
25          then goto 13510
13502      if abs(zintlst-zintnow/numdrop)>maxslew# then
          1350#
13503      zintnow=zintnow/numdrop : sound 1200,sounder :
          return          'dropped beat
30 13504      if numdrop>1 then goto 13520 else goto 13510

          ' check for premature trigger (note:
          '      premature trigger assump-
          '      -tion remains in effect
35          '      only for glitched time
          '      (if added portion is an
```

- 225 -

```

        ' acceptable beat,
        ' (that's how it's used;
        ' otherwise slew rate
        ' (limiter extends
5
        ' assumption to added portion
13510 if abs(zintnow+zdata(intnow+1)_
        zintlstr)>maxslew# then 13520
13511 zintnow=zintnow+zdata(intnow+1)
10 ' assume premature trigger here
13512 sound 1400,sounder : return

        ' slew rate limiter
13520 sound 600,sounder : zintnow=zintlstr
15
        return

        ' calculating the respiratory rate using the
20 ' comb method
        ' [spectrum in ydata(*)]
        ' start at frequency : minrespfrq#
        ' (in pixels): minresp
        ' use comb tooth width: combwidth#
25 ' (in pixels): combpix

14000 maxcomb#=0# : respcomb=0 : combstep=combpix\2+1
        ' for loop shifts comb beginning to different
30 ' frequencies
14001 for comb=minresp to npair step combstep
14002 curcomb#=0# : harmbeg=comb-combstep+2
14003 lastbeg=harmbeg+9*comb : if lastbeg>npair
        then lastbeg=npair
35
        ' while loop adds up 10 teeth

```

- 226 -

```

      ' (harmonics) in the comb
14004   while harmbeg<=lastbeg
14005       toothptr=harmbeg
14006       lstooth=harmbeg+combpix : if
5         lstooth>npair then lstooth=npair

      ' this while loop adds one tooth's
      ' contribution to comb
14007   while toothptr<=lstooth
10 14008       curcomb#=curcomb#+ydata(toothptr)
14009       toothptr=toothptr+1
14010   wend
14011       harmbeg=harmbeg+comb
14012   wend
15 14013   if curcomb#>maxcomb# then maxcomb#=curcomb# :
       respcomb=comb
14014   next comb

14050   locate 3,1 : gosub 20000 : print "respiratory
20       comb fraction: ";
14051   curcomb#=0# : for i=1 to npair :
       curcomb#=curcomb#+ydata(i) : next i
14052   respcombfrac#=maxcomb#/curcomb# : print using
       "#.###";respcombfrac#;

25
      ' respcomb now has respiratory frequency or a
      ' subharmonic
      ' to decide which is the first harmonic look at
      ' weight in each tooth
30      ' of the comb; a higher harmonic comb must
      ' contribute at least double
      ' amplitude to be designated as the fundamental
      ' (4xspectral weight)

35 14100   maxtooth#=0 : resptooth=0 : harmbeg=respcomb+1-
          combpix

```



- 227 -

```
14101  lastbeg=harmbeg+9*respcomb : if lastbeg>npair
      then lastbeg=npair

14102      while harmbeg<=lastbeg
5  14103          toothptr=harmbeg : curtooth#=0#
      14104          lstooth=harmbeg+combpix+combpix
      14105          if lstooth>npair then lstooth=npair

          ' add up one widened tooth
10  14110      while toothptr<=lstooth
      14111          curtooth#=curtooth#+ydata(toothptr)
          : toothptr=toothptr+1
      14112      wend

          ' compare to previous teeth
15  14120      if curtooth#<4*maxtooth# then goto 14130
      14121          maxtooth#=curtooth# :
          resptooth=harmbeg

      14130      harmbeg=harmbeg+respcomb
20  14131      wend

          ' compute respiratory frequency as peak
          ' average

14200      toothptr=resptooth : respfreq#=0#
25  14201      lstooth=toothptr+combpix+combpix
      14202      if lstooth>npair then lstooth=npair

          ' average frequency over fundamental
          ' peak

30  14210      while toothptr<=lstooth
      14211          respfreq#=respfreq#+ydata(toothptr)
          *cdbl(toothptr-1)
      14212          toothptr=toothptr+1
      14213      wend
35  14214      respfreq#=respfreq#/maxtooth#/1024#*respsampl
```

- 228 -

```

14220  resp.lopixel=cint((respfreq#-
      .06#)/respsampl*1024#)+1
      ' integration limits
14221  resp.hipixel=cint((respfreq#+.06#)
5      /respsampl*1024#)+1

      return

10      ' spectral amplitude calculations
15000  lfa#=0# : rfa#=0# : coherence#=0#
15001  for i=pixel.04 to pixel.10
15002      lfa#=lfa#+fnzmag(zspec.hb.real(i),
15      zspec.hb.imag(i))
15003      next i
15004  lfa#=lfa#+lfa#

15010  for i=resp.lopixel to resp.hipixel
20 15011      rfa#=rfa#+fnzmag(zspec.hb.real(i),
      zspec.hb.imag(i))
15012      next i
15013  rfa#=rfa#+rfa#

25 15020  for i=1 to 512
15021      coherence#=coherence#+fnzcoher
      (zreal(i),zrimag(i),_
      zspec.hb.real(i),
      zspec.hb.imag(i))
30 15022      next i
15023  coherence#=coherence#/zspectsum

15030  ratio#=lfa#/rfa#
15031  cratio#=lfa#/coherence#

35 15040  locate 6,60 : print using  "lfa: ##.###";lfa#;

```

- 229 -

```
15041  locate 7,60 : print using "rfa: ##.###
                                     (##.###)";rfa#,coherence#;
15042  locate 8,58 : print using "ratio: ##.###
                                     (##.###)";ratio#,cratio#;
5      return

      ' storing trend data on floppy disk (file #10)
10 16000 lset hr$=mkd$(zavghr)
    16001 lset rr$=mkd$(respfreq#)
    16002 lset rcf$=mkd$(respcombfrac#)
    16003 lset lfa$=mkd$(lfa#)
    16004 lset rfa$=mkd$(rfa#)
15 16005 lset coher$=mkd$(coherence#)
    16006 lset ratio$=mkd$(ratio#)
    16007 lset cratio$=mkd$(cratio#)
    16008 lset hrintegral$=mkd$(hrspectsum#)
    16009 lset respintegral$=mkd$(respspectsum#)
20 16010 lset timestamp$=time$
    16011 lset hbrecord$=mki$(analrec.hr)
    16012 lset adrecord$=mki$(analrec.ad)
    16013 lset hbeat$=mki$(anal.beat)
    16014 lset samprate$=mki$(respsampl)
25
    recordl0no=recordl0no+1 : put #10,recordl0no

    return

30

      ' reading trend data from floppy disk (file #10)
16500 if recordl0no<=1 then return
16501 cls
35
16510 xmin=10 : xmax=310 : ymin=2 : ymax=197 :
```

- 230 -

```

        numvalg=xmaxs-xmins+1
16511  call swindow(xmins,xmaxs,ymins,ymaxs)

16512  call clrwindw
5  16513  call axes

16520  numpts=recordl0no
16521  lfa.beg=recordl0no
16522  rfa.beg=2*recordl0no
10 16523  ratio.beg=3*recordl0no
16524  lastydata=4*recordl0no
16525  ln10#=log(10#)
16526  xscale#=numvalg/recordl0no

15      ' get trend information from the disk file
16530  for temprec=1 to recordl0no
16531      get #10,temprec
16532      ydata(temprec)=197-.78#*cvd(hr$)
16533      ydata(temprec+lfa.beg)=197-19.5*cvd(lfa$)
20 16534      ydata(temprec+rfa.beg)=197-19.5*cvd(rfa$)
16535      ydata(temprec+ratio.beg)=100-
          log(cvd(ratio$))/ln10#*45#
16536      next temprec

25 16537  for i=1 to lastydata : if ydata(i)<ymins then
          ydata(i)=ymins
16538      if ydata(i)>ymaxs then ydata(i)=ymaxs :
          next i

30      ' plot trends here
16540  for trend=0 to 3 : trendoff=trend*recordl0no
16542      gctr=1 : ydata1st=ydata(1) :
          ydatag(1)=ydata1st
16543      for temprec=2 to recordl0no :
35          gctrmax=temprec*xscale#
16544          gdif=gctrmax-gctr : if gdif<=0 then goto
```

- 231 -

```

16550
16545      ydatadif=ydata(temprec+trendoff)-
          ydatalst : part=0
16546      while gctr<gctrmax : gctr=gctr+1 :
5          part=part+1
16547      ydatag(gctr)=ydatalst+
          (part/gdif)*ydatadif : wend
16548      ydatalst=ydata(temprec+trendoff)
16550      next temprec
10 16551      linemask=linetype(trend) : x=xmins
16552      gdataptr=varptr(ydatag(1)) : numvalg=xmaxs-
          xmins+1
16553      call fgraph(gdataptr,numvalg,x,linemask)
16554      next trend
15
16560      locate 2,42 : print "HR (0-250 bpm)";
16561      locate 3,42 : print "lfa (0-10 bpm^2)";
16562      locate 4,42 : print "rfa (0-10 bpm^2)";
16563      locate 5,42 : print "ratio (.01-100)";
20
16600      req.cls=1

      return

25

      ' subroutine to print out the interrupt vectors

30 19000      def seg=0
          print "IRQ3 @0B*4H: ";hex$(peek(&h2C));
          '   ";hex$(peek(&h2D));" ";
          print hex$(peek(&h2E));
          '   ";hex$(peek(&h2F));tab(40);
35      print "IRQ4 @0C*4H: ";hex$(peek(&h30));
          '   ";hex$(peek(&h31));" ";

```

- 232 -

```
        print hex$(peek(&h32));" ";hex$(peek(&h33));
        return

5          ' routine to clear the present line
20000      csnow=csrlin:locate csnow,1:print
          return

10          ' check pointers to see if any disk files need
          ' to be written
30000      call rdptrs(adwr,hbwr,adflag,hbflag)
30001      if adflag=adflaglst and hbflag=hbflaglst then
          return

15          30010      while adflag>recordlno+1 : beep : locate 23,1 :
          print "data #1 loss";
          30011          recordlno=adflag-1 : wend
          30020      while hbflag>record2no+1 : beep : locate 23,1 :
20          print "data #2 loss";
          30021          record2no=hbflag-1 : wend
          30030      while hbflag>record3no+1 : beep : locate 23,1 :
          print "data #3 loss";
          30031          record3no=hbflag-1 : wend

25          30040      if adflag<recordlno+1 then goto 30050
          30041          recordlno=adflag : put #1,adflag
          30042          if datacycle<=0 then datacycle=datacycle+1
          'if not processing, begin

30          30050      if hbflag=record2no+1 then record2no=hbflag :
          put #2,hbflag
          30060      if hbflag=record3no+1 then record3no=hbflag :
          put #3,hbflag

35          locate 3,1 : gosub 20000 : print "current file
          records: ";adflag; print " (#1) ";hbflag;"
```

- 233 -

```
(#2)";  
adflaglst=adflag : hbflaglst=hbflag
```

```
return
```

5

```
end
```

10

15

20

25

30

35

- 234 -

```

                                page    66,80
; sync7s.asm - an assembler routine to handle interrupts
;
5 ;      from IRQ4 and collect
;      synchronous data from the A/D (board 2
;      configuration assumed)
;      The routine checks A/D readings for
;      output validity
;      Data is loaded by interrupts into both a
10 ;      processing buffer and
;      a disk file I/O buffer to allow quick
;      archival; an overflow
;      flag signals when a disk file buffer
;      should be stored and
15 ;      also indicates whether the disk buffer
;      was corrupted.
;      To acknowledge storage of a disk buffer
;      one must reset the
;      overflow flag using <ackfdio>
20 ; Last revision:  3 May 1985
;
;
;-----;
; 8088 interrupt location      ;
25 ;-----;

abs0      segment at 0      ;absolute memory segment
                                ;allows placement of
                                ;interrupt address
30      org      0BH*4      ;future heart beat
                                interrupt handler resides
IRQ3_int   dw      2 dup(?);at int 0B

                                org      0CH*4      ;8253 timebase interrupt
                                ;handler resides
35      IRQ4_int   dw      2 dup(?);at int 0C

```



- 235 -

```

abs0          ends          ;

5

;-----;
; int_buffer: area to save DOS ;
;   dummy interrupt ptrs      ;
;-----;

10

int_buffer    segment        ;data segment containing
                           ;user interrupt buffer

15  save_int   dw      4 dup(?);offset for two DOS
                           ;interrupts saved
                           ;to be restored using
                           ;exstint

20  int_buffer ends          ;

;-----;
25  ; working storage for      ;
; interrupts                  ;
;-----;

dseg_sync     segment        ;data segment for
30                                     ;interrupts

;.....declare all variables public
;           for use by other
;           assembly level routines
35  public ad_buffer,ad_rd,ad_wr,sync_ctr
      public hb_buffer1,hb_buffer2,hb_rd,hb_

```

- 236 -

wr,heartbeats

```

;.....timebase local storage and
buffer

```

```

5      ad_buffer      db      1024 dup(?) ;buffer for A/D
                                   values
      ad_rd           dw      ?          ;read indicator for A/D
                                   ;disk buffer
10     ad_wr           dw      ?          ;write pointer for A/D
                                   ;buffer (incrementing)
      sync_ctr        dw      ?          ;counter for timebase
                                   ;interrupt (overflows)

15
;.....heart beat local storage and
;      buffer
;      note:for main clock
;      14.318 180 MHz (osc)
20     ;      system clock
;      4.772 727 MHz (clock)
;      8253 clock
;      2.386 363 MHz (ck8253)
;      (ck8253 / 432)
25     ;      5.524 KHz      (hb.clk)
;      (ck8253 /596592) 4 Hz
;      (respck)
;      hb.clk = 1381*respck
;      sync.ctr overflow =
30     ;      16384 sec (4:33:04)

      hb_buffer1      dw      1024 dup(?) ;heart beat time
                                   stamps for previous 1024
      hb_buffer2      dw      1024 dup(?) ;beats (2 words:
35     hb.clk,sync.ctr)

```

- 237 -

```

      hb_rd      dw      ?      ;read indicator for
                                ;heart beat disk buffers
      hb_wr      dw      ?      ;write pointer
                                ;(incrementing) for hb_buffer
5
      heartbeats dw      ?      ;keep track of number of
                                ;beats processed

                                ;.....pointers to disk file buffers
10
      fd1ptr     label   dword   ;pointer to floppy disk
                                file #1 buffer
      fd1ptroff  dw      ?      ; (offset)
      fd1ptrseg  dw      ?      ; (segment)
15
      fd2ptr     label   dword   ;pointer to floppy disk
                                file #2 buffer
      fd2ptroff  dw      ?      ; (offset)
      fd2ptrseg  dw      ?      ; (segment)
20
      fd3ptr     label   dword   ;pointer to floppy disk
                                file #3 buffer
      fd3ptroff  dw      ?      ; (offset)
      fd3ptrseg  dw      ?      ; (segment)
25
      dseg_sync  ends          ;

30
                                ;-----;
                                ; setup structures to allow access to;
                                ; arguments pased by BASIC          ;
                                ;-----;
35
                                ; subroutine

```



- 239 -

```

whichbuff      dw      ?      ;place to select which
                                ;buffer to read
BASIC_ptr      dw      ?      ;place to get pointer to
                                ;BASIC data array
5  frame_rdbuf  ends

                                ; subroutine rdptrs
                                ;(adwr,hbwr,adflag,hbflag)
frame_rdptrs    struc          ;define the stack
10                                ;structure for passing
                                ;arguments to BASIC
savebp3         dw      ?      ;caller's base pointer
saveret3        dd      ?      ;return offset and
                                ;segment pushed by BASIC
15 hbflag       dw      ?      ;flag indicating disk
                                ;file #1,#2 buffers full
adflag          dw      ?      ;flag indicating disk
                                ;file #1 buffer is full
BASIC_hbwr      dw      ?      ;write pointer for heart
20                                ;beat buffer
BASIC_adwr      dw      ?      ;write pointer for ad
                                ;buffer
frame_rdptrs    ends

25

                                ;.....code segment begins here

cseg_sync       segment 'code'
30 basic_dgroup  group  data,stack,const,heap,memory
                                ;defining link to BASIC
porta          equ      071CH   ;port definitions for
                                ;8255 port expander
portb          equ      071DH   ;these addresses are
35                                ;decoded on the homemade
portc          equ      071EH   ;board

```

- 240 -

```

control      equ      071FH      ;control word in the
                                   ;8255
timer0       equ      0704H      ;8253 timer0 register
timer1       equ      0705H      ;8253 timer1 register
5  timer2     equ      0706H      ;8253 timer2 register
con8253      equ      0707H      ;8253 control register

;-----;
10  ; timebase interrupt handler (not accessible to;
; BASIC)                                     ;
;-----;
      ;this routine reads the A/D every timer1
      ;tick
15  ;and stores the point in the analog
      ;buffer

tbase_int    proc      far      ;this procedure is not
20                                     ;made public
      assume  cs:cseg_sync,ds:dseg_
      sync,es:nothing,ss:nothing
      push    ax      ;save registers used
                                   ;during interrupt
25  push      bx      ;
      push      cx      ;
      push      dx      ;
      push      si      ;
      push      di      ;
30  push      ds      ;
      push      es      ;

      mov      ax,dseg_sync      ;set up segment
                                   ;register for data area
35  mov      ds,ax      ;

```

- 241 -

```

;.....increment counters/ decrement
;           pointers
inc      sync_ctr      ;increment
                        ;interrupt counter
5      mov      cx,20      ;allow up to 20
                        ;rereads of A/D

;.....get analog value from A/D and
;           send to buffer
10     mov      dx,portb      ;get analog
                        ;value from A/D
      in       al,dx      ;

      mov      bx,ad_wr      ;and put analog
15     ;data pointer in bx
      mov      ad_buffer[bx],al
      ;save analog value in ad_buffer

      chk_adc:      in      al,dx      ;reread adc and
20     ;check if previous
      cmp      ad_buffer[bx],al      ;value agrees
      je       adc_ok      ;if value is the
                        ;same we're done
      loop     retry      ;retry if retry
25     ;counter is not depleted
                        ;failure returns
                        ;last value read

30     adc_ok:      inc      ad_wr      ;increment write
                        ;pointer
      cmp      ad_wr,1023      ;see if write
                        ;pointer<=1023
      jle      tbase_eoi      ;if pointer is
35     ;in range then finish

      int

```

- 242 -

```

;.....reset local ptr and load disk
;      buffer for file #1

5      xor     ah,ah           ;zero ah as
                                ;upper byte of A/D reading
mov     cx,1024               ;load counter
                                ;for 1024 repetitions
10     lea     si,ad_buffer    ;load local
                                ;buffer address
les     di,fdlptr             ;load pointer to
                                ;disk file #1 buffer
fdllp: lodsb                   ;repeat moves
                                ;1024 times (ds:si->es:di)
15     stosw                   ;converting
                                ;bytes to words
loop    fdllp                 ;
mov     ad_wr,cx              ;reset write
                                ;pointer (wrap around)
20     inc     ad_rd           ;increment read
                                ;request for disk

;.....acknowledge interrupt to
;      8259A
25     tbase_eoi: mov     al,20H ;send EOI to 8259A
out     20H,al ;

pop     es                    ;restore registers which
                                ;were used
30     pop     ds ;
pop     di ;
pop     si ;
pop     dx ;
pop     cx ;
35     pop     bx ;
pop     ax ;

```



- 243 -

```

                                ired                ;return to place where
                                                ;interrupt occurred

5  debugmsg1      db      'this is the end of the time
                        base interrupt'

tbase_int      endp

10

                                ;-----;
                                ; heart beat interrupt handler (not accessible ;
                                ; to BASIC)                                     ;
15                                ;-----;
                                ;this routine reads the local system
                                ;timers
                                ;every heart beat and stores the time in
                                ;the heart beat buffer for use in
20                                ;spectral analysis
                                ;

hbeat_int      proc      far      ;this procedure is not
25                                ;made public
                                assume cs:cseg_sync,ds:dseg_sync
                                assume es:nothing,ss:nothing

                                push     ax      ;save registers during
30                                ;interrupt
                                push     bx      ;
                                push     cx      ;
                                push     dx      ;
                                push     si      ;
35                                push     di      ;
                                push     ds      ;
```

- 244 -

```
push    es    ;

mov     ax,dseg_sync    ;set up segment
                        ;register for data area
5      mov     ds,ax    ;

inc     heartbeats    ;increment heart
                        ; beat counter

10      ;.....read counters and store
        ;          result in hb_buffer
mov     dx,con8253    ;prepare to read
                        ;hbl.clk from timer1
mov     al,40H        ;by latching
15      ;counts in timer1
out     dx,al        ;

mov     dx,timer1    ;prepare to read
                        ;the latched value
20      in     al,dx    ;from the timer
                        ;(low byte first)
mov     ah,al        ;save low byte
                        ;in ah
in     al,dx        ;(high byte
25      ;last)
xchg    al,ah        ;get the bytes'
                        ;order right

mov     bx,hb_wr    ;get write
30      ;pointer for hb_buffer
add     bx,bx        ;double to
                        ;point to a word
mov     hb_buffer1[bx],ax    ;and store
                        ;hbl.clk counts
35      ;.....read overflow counter from
```

- 245 -

```

; timer2
mov dx,con8253 ;prepare to read
;hb2.clk from timer2
mov al,80H ;by latching
5
;counts in timer2
out dx,al ;
mov dx,timer2 ;prepare to read
10 ;the latched value
in al,dx ;from the timer
; (low byte first)
mov ah,al ;save low byte
;in ah
15 in al,dx ;(high byte
;last)
xchg al,ah ;get the bytes'
;order right in ax
20 mov hb_buffer2[bx],ax ;store
result in hb2.clk buffer
;.....increment write pointer and
; check for buffer overflow
25 inc hb_wr ;increment write
;pointer
cmp hb_wr,1023 ;if hb_wr<=1023
jle hb_eoi ;then finish up
30 ;.....reset local ptr/load disk
; buffers for files #2,#3
; (routine takes about 15-20
; msec to fill disk buffer)
mov cx,1024 ;load counter
35 ;for 1024 repetitions
lea si,hb_buffer1 ;load local

```

- 246 -

```

;buffer address
les     di,fd2ptr      ;load pointer to
;disk file #2 buffer
fd2lp:  movsw          ;repeat moves
5         ;1024 times (ds:si->es:di)
loop    fd2lp          ;
mov     cx,1024        ;load counter
;for 1024 repetitions
lea     si,hb_buffer2  ;load local
10         ;buffer address
les     di,fd3ptr      ;load pointer to
;disk file #3 buffer
fd3lp:  movsw          ;repeat moves
;1024 times (ds:si->es:di)
15         loop    fd3lp          ;
mov     hb_wr,cx        ;reset write
;pointer (wrap around)
inc     hb_rd           ;increment read
;request
20
;.....acknowledge interrupt to
;      8259A
hb_eoi: mov     al,20H  ;send EOI to 8259A
out     20H,al        ;
25
pop     es             ;restore registers and
pop     ds             ;
pop     di             ;
pop     si             ;
30         pop     dx             ;
pop     cx             ;
pop     bx             ;
pop     ax             ;
iret                  ;return to place where
35         ;interrupt occurred

```

- 247 -

```

debugmsg2      db      'this is the end of the heart
                    beat interrupt'

5  hbeat_int      endp

;-----;
10  ; subroutine instint [install_interrupts]      ;
; (fillptr,fil2ptr,fil3ptr)                      ;
;-----;

instint        proc    far
15      public  instint
;public symbol allows external references
;es,ds vectors and must be restored movsw
;uses (ds:si)(es:di) addr
      assume  cs:cseg_sync,ss:basic_
20      dgroup,ds:basic_dgroup
      assume  es:basic_dgroup
      used to access interrupt

;.....save registers

25      push   bp      ;save BASIC base pointer
                    ;   for return to BASIC
      mov     bp,sp    ;point stack pointer at
                    ;frame reference to
30      ;address of BASIC analog
                    ;data buffer

      push    ax      ;save additional
                    ;registers
35      push    si      ;
      push    di      ;

```

- 248 -

```

push    ds      ;
push    es      ;
pushf                    ;and flags

5                                ;set up the segment
                                ;registers
mov     ax,dseg_sync    ;set up access
                                ;to floppy disk data ptrs
mov     es,ax          ;
10    assume es:dseg_sync    ;

                                ;.....put disk file pointers into
                                ;      local memory
15    mov     di,[bp].B_fillptr    ;get
                                pointers from BASIC
mov     ax,[di]          ;and
                                ;      save in dseg_sync areas
mov     fdlptroff,ax     ;

20

mov     di,[bp].B_fil2ptr    ;
mov     ax,[di]            ;
mov     fd2ptroff,ax        ;

25    mov     di,[bp].B_fil3ptr    ;
mov     ax,[di]            ;
mov     fd3ptroff,ax        ;

mov     ax,ds            ;put segment
30                                ;registers into
mov     fdlptrseg,ax      ;pointers
mov     fd2ptrseg,ax      ;
mov     fd3ptrseg,ax      ;

35                                ;set up the segment
```

- 249 -

```

;registers
mov     ax,int_buffer
;es points to buffer area to save
mov     es,ax           ;DOS dummy
5
;interrupt vector
assume  es:int_buffer  ;
mov     ax,0           ;ds points to
;abs0 (interrupt table)
mov     ds,ax          ;
10
assume  ds:abs0        ;

;setup access to
;interrupt vectors
lea     di,save_int    ;load offset of
15
;save_int in es,di
lea     si,IRQ3_int    ;load offset of
;IRQ3_int in ds,si
cld
;clear direction
;flag to increment ptrs
20
movsw
;save DOS dummy
;interrupt vectors to be
movsw
;restored later
movsw
;now saving IRQ4
movsw
;
25

mov     IRQ3_int+2,cseg_sync ;install
;the heart beat (IRQ3)
mov     IRQ3_int,offset hbeat_int
;interrupt handler now
30
mov     IRQ4_int+2,cseg_sync ;install
;the DAC timebase (IRQ4)
mov     IRQ4_int,offset tbase_int
;interrupt handler now
35

```

- 250 -

;.....initialization of buffer  
control variables

```
5      mov     ax,dseg_sync      ;setup data  
      ;segment for initialization  
      mov     ds,ax             ;  
      assume  ds:dseg_sync      ;ds segment  
      ;register now redefined  
  
10     xor     ax,ax             ;zero ax  
      ;register to initialize  
      mov     heartbeats,ax     ;counters  
      mov     sync_ctr,ax       ;  
15     mov     ad_wr,ax          ;initialize  
      ;read/write pointers to top  
      mov     hb_wr,ax          ;of buffer  
      mov     ad_rd,ax          ;  
      mov     hb_rd,ax          ;  
20  
  
      ;.....return to BASIC  
  
25     popf                    ;restore flags  
      pop     es                ;restore additional  
      registers  
      pop     ds                ;  
      pop     di                ;  
      pop     si                ;  
30     pop     ax                ;  
  
      pop     bp                ;restore BASIC's base  
      ;pointer and  
35     ret     6                ;delete 3 parameters (6  
      ;bytes) from the stack  
      ;and return to the
```



- 251 -

```
                                ;calling routine

debugmsg3      db      'this is the end of the
                                interrupt installation'

5

instint        endp

10

                                ;-----;
                                ; subroutine exstint (exstall_      ;
                                ; interrupts)                        ;
                                ;-----;

15

exstint        proc      far
                public   exstint ;public symbol allows
                                ;external references

20              assume   cs:cseg_sync,ss:basic_dgroup
                assume   ds:int_buffer,es:abs0
                ;es,ds used to access interrupt
                ;vectors and must be restored
                ;movsw uses (ds:si)(es:di) addr

25

                ;.....save registers

                push     bp      ;save BASIC base pointer
                                ;    for return to BASIC

30              mov      bp,sp   ;point stack pointer at
                                ;    frame reference to
                                ;access arguments passed
                                ;    by BASIC (none here)

35              push     ax      ;save additional
                                ;registers
```

- 252 -

```

                    push    si        ;
                    push    di        ;
                    push    ds        ;
                    push    es        ;
5                   pushf          ;and flags

                                ;set up the segment
                                ;registers as assumed
                    mov     ax,0      ;es points to
10                  ;abs0 (interrupt table)
                    mov     es,ax     ;
                    mov     ax,int_buffer ;ds points to
                                ;buffer area to save
                    mov     ds,ax     ;DOS dummy
15                  ;interrupt vector

                                ;setup access to
                                ;interrupt vectors
                    lea     di,IRQ3_int ;load offset of
20                  ;IRQ3_int in es,di
                    lea     si,save_int ;load offset of
                                ;save_int in ds,si
                    cld              ;clear direction
                                ;flag to increment ptrs
25                  movsw           ;restore DOS
                                ;dummy interrupt vectors
                    movsw           ;for IRQ3
                    movsw           ;and IRQ4
                    movsw           ;
30

                                ;.....return to BASIC

                    popf            ;restore flags
35                  pop     es      ;restore additional
                                ;registers

```

- 253 -

```

    pop    ds      ;
    pop    di      ;
    pop    si      ;
    pop    ax      ;
5
    pop    bp      ;restore BASIC's base
                        ;pointer and
    ret     0      ;delete 0 parameters (0
                        ;bytes) from the stack
10
                        ;and return to the
                        ;calling routine

    debugmsg4      db      'this is the end of the
                        interrupt exstallation'
15
    exstint        endp

20

                        ;-----;
                        ; subroutine rdbeat (heartbeats, sync_ ;
                        ; pulses)                                ;
25
                        ;-----;

    rdbeat        proc    far
                        public rdbeat    ;public symbol allows
30
                        external references
                        assume cs:cseg_sync,es:dseg_sync
                        assume ds:basic_dgroup,ss:basic_dgroup

35
                        ;.....save registers
```

- 254 -

```

                    push    bp        ;save BASIC base poin
                                   ;ter for return to BASIC
                    mov     bp,sp     ;point stack pointer at
                                   ;frame reference to
5
                                   ;access arguments passed
                                   ;by BASIC (one here)

                    push    ax        ;save additional
                                   registers
10
                    push    di        ;
                    push    es        ;

                    mov     ax,dseg_sync ;set up segment
                                   register for data area
15
                    mov     es,ax      ;

                    mov     ax,heartbeats ;get
                                   ;beats from local memory
20
                    mov     di,[bp].BASIC_beats ;
                    mov     [di],ax      ;send
                                   ;beats to BASIC

                    mov     ax,sync_ctr ;get
                                   ;sync pulses from local
25
                    mov     di,[bp].BASIC_syncs ;memory
                    mov     [di],ax      ;send
                                   ;sync pulses to BASIC

30
                    ;.....return to BASIC

                    pop     es        ;restore additional
                                   registers
35
                    pop     di        ;
                    pop     ax        ;
```

- 255 -

```

5      pop    bp      ;restore BASIC's base
                        ;pointer,
      ret     4        ;delete 2 parameters (4
                        ;bytes) from the stack
                        ;and return to the
                        ;calling routine

      debugmsg5      db      'this is the end of the heart
10                        beat read routine'

      rdbeat  endp

15      ;-----;
      ; subroutine rdbuf (BASIC_      ;
      ; ptr,whichbuff)                ;
      ;-----;

      ;this routine dumps a buffer
20      ;from the
      ;assembly routine data area to a
      ;BASIC array
      ;pointed to by BASIC_ptr;
      ;whichbuff selects
25      ;the assembler buffer to be
      ;dumped.
      ;choices of buffer are:
      ; 0 - ad_buffer      (bytes)
      ; 1 - hb_buffer1     (words)
30      ; 2 - hb_buffer2     (words)

      rdbuf      proc      far
                        public rdbuf      ;public symbol allows
                        ;external references

35      assume cs:cseg_sync,es:basic_dgroup
      assume ds:basic_dgroup,ss:basic_dgroup

```

- 256 -

```
;.....save registers

5      push    bp      ;save BASIC base pointer
                        ;for return to BASIC
      mov     bp,sp    ;point stack pointer at
                        ;frame reference to
                        ;access arguments passed
10     ;by BASIC (one here)

      push    ax      ;save additional
                        ;registers

      push    cx      ;
15     push    si      ;
      push    di      ;
      push    ds      ;
      push    es      ;
      pushf                ;and flags
20

;.....get pointers from BASIC
      mov     di,[bp].whichbuff      ;get
                        ;buffer choice from BASIC
25     mov     ax,[di]                ;

      mov     di,[bp].BASIC_ptr
      ;get pointer to BASIC's data area
      mov     di,[di]                ;and put pointer
30     ;into di

;.....set up extra segment register
;      and counter
35     mov     cx,dseg_sync      ;set up segment
                        register for data area
```

- 257 -

```

mov     ds,cx          ;
assume  ds:dseg_sync
mov     cx,1024        ;load counter
                        ;with number of objects

5

;.....select buffer here and place
;      pointer in si
or      ax,ax          ;compare
10                        ;selector with 0

jz      rd_adbuf
;if zero (select =0) read ad_buffer
dec     ax              ;decrement to
                        ;see if select was 1

15      jz      rd_hbbuf1
;if zero (select =1) read hb_buffer1
dec     ax              ;decrement to
                        ;see if select was 2

20      jz      rd_hbbuf2
;if zero (select =2) read hb_buffer2
jmp     rdbuf_end
;not a valid buffer, so return to BASIC

rd_adbuf:  lea     si,ad_buffer    ;point source
25                        ;index to ad_buffer

jmp     move_dta_byte    ;

rd_hbbuf1: lea     si,hb_buffer1   ;point source
                        ;index to hb_buffer1

30      jmp     move_dta_word    ;

rd_hbbuf2: lea     si,hb_buffer2   ;point source
                        ;index to hb_buffer2

35      jmp     move_dta_word    ;

```

- 258 -

```

;.....move byte data from local
;      storage to BASIC array
move_dta_byte: xor     ah,ah    ;zero upper byte of ax

5          cld                ;clear direction flag to
                                ;increment si,di by 2
byt_lp:    lodsb              ;move data bytes from
                                ;local storage (ds:si)
                                ;and store as a word in
10         stosw              ;BASIC's area (es:di)
                                ;
                                loop    byt_lp  ;
                                jmp      rdbuf_end ;finished

15         ;.....move word data from local
;      storage to BASIC array
move_dta_word: cld           ;clear direction flag to
                                ;increment si,di by 2
wd_lp:      movsw            ;get data word from
                                ;local storage (ds:si)
20         loop    wd_lp      ;and store as a word in
                                ;BASIC's area (es:di)

25         ;.....return to BASIC

rdbuf_end: popf              ;restore flags
            pop     es        ;restore additional
                                ;registers
30         pop     ds        ;
            pop     di        ;
            pop     si        ;
            pop     cx        ;
            pop     ax        ;
35         pop     bp        ;restore BASIC's base
```





- 260 -

```
assume ds:basic_dgroup,ss:basic_dgroup

;.....save registers

5      push    bp      ;save BASIC base pointer
                        ;for return to BASIC
      mov     bp,sp    ;point stack pointer at
                        ;frame reference to
10     ;access arguments passed
                        ;by BASIC (one here)

      push    ax      ;save additional
                        ;registers
15     push    di      ;
      push    es      ;

      mov     ax,dseg_sync ;set up segment
                        ;register for data area
20     mov     es,ax    ;

      mov     ax,ad_wr   ;get write
                        ;pointer for A/D buffer
25     mov     di,[bp].BASIC_adwr ;and send
                        ;to BASIC
      mov     [di],ax    ;

      mov     ax,hb_wr   ;get
                        ;write pointer for heart
30     mov     di,[bp].BASIC_hbwr ;beat
                        ;buffer and send to BASIC
      mov     [di],ax    ;

35     mov     ax,ad_rd   ;get
                        ;disk file flag for A/D
```

- 261 -

```

                                mov     di,[bp].adflag      ;buffer
                                ;and send to BASIC
                                mov     [di],ax           ;
5                                mov     ax,hb_rd           ;get
                                ;disk file flag for heart
                                mov     di,[bp].hbflag     ;beat
                                ;buffers and send to BASIC
                                mov     [di],ax           ;
10                               ;.....return to BASIC

                                pop     es               ;restore additional
                                ;registers
15                               pop     di               ;
                                pop     ax               ;

                                pop     bp               ;restore BASIC's base
                                ;pointer,
20                               ret     8               ;delete 4 parameters (8
                                ;bytes) from the stack
                                ;and return to the
                                ;calling routine

25  debugmsg7                   db     'this is the end of the pointer
                                read routine'

                                rdptrs   endp

30                               cseg_sync               ends
                                ;
                                ; module gwindow1.asm - a collection of routines useful
                                ; for preparing data
35                               ; for the fast graphics routine.
                                ;
```

- 262 -

```
;      subroutines:
;
;      dwindow(xmin,xmax,ymin,ymax) - establish
;      data value limits corresponding to
5 ;      screen window.
;
;      swindow(xmin,xmax,ymin,ymax) - establish
;      screen boundaries for data to be
;      plotted.
10 ;
;      clrwindw - clear contents of present
;      window
;
;      axes - prepare axes for current window
;      (no tick marks yet)
15 ;      (first version: only draws a box
;      around window)
;
;      scaler(indata_ptr,outdata_ptr,numval) -
20 ;      scale data to fit into window requires
;      correct initialization
;      using dwindow
;      and swindow
;      (first version: only scales y-
25 ;      coordinate with dwindow)
;      (
;      scaled by numval)
;      (
;      maximum y-value
30 ;      is plotted)
;
;-----
;
;      arguments passed by BASIC
35 ;
;
```

- 263 -

```

;      indata_ptr      - offset of BASIC array
;                        containing y-coordinates of
;                        points to be plotted
;      outdata_ptr     - offset of BASIC array
5  ;      containing scaled y-coordinates
;      numval          - number of values to plot
;
;-----
10

;.....screen memory definition

screen_memory segment at 0B800H
15 even_pixels db      8000 dup(?)      ;pixels with
;even y-coordinates
org      2000H      ;beginning of
;high screen memory
odd_pixels db      8000 dup(?)      ;pixels with odd
20 ;y-coordinates
screen_memory ends

;.....local memory definitions
25

dseg_wind segment      ;valid default values
;present at startup

xmin_s dw      0      ;minimum screen ordinate
30 ;for window
xmax_s dw      639    ;maximum screen ordinate
;for window
ymin_s dw      0      ;minimum screen abscissa
;for window
35 ymax_s dw      199   ;maximum screen abscissa
;for window

```

- 264 -

	xmin_d	dw	0	;minimum data ordinate ;for window
	xmax_d	dw	16384	;maximum data ordinate ;for window
5	ymin_d	dw	0	;minimum data abscissa ;for window
	ymax_d	dw	16384	;maximum data abscissa ;for window
10	ulh_cor	dw	0	;offset for upper left ;hand corner of screen
	urh_cor	dw	79	;offset for upper right ;hand corner of screen
15	llh_cor	dw	3EF0H	;offset for lower left ;hand corner of screen
	lrh_cor	dw	3F3FH	;offset for lower right ;hand corner of screen
20	outptr	dw	?	;pointer to output array ;in BASIC (must be ;at least as large as ;input array)
	rndoff	dw	?	;roundoff correction (if ;fraction>.5 round up)
25	numvalt	dw	?	;save number of points ;in input array for xpass
	bx_last	dw	?	;save pointer during x- ;scaling to allow
30				;use of largest y per x ;pixel
	dseg_wind	ends		

35

- 265 -

```

;-----;
; define structures for passing arguments from ;
; BASIC ;
;-----;

5
; subroutines
; dwindow/swindow(xmin,xmax,ymin,ymax)

frame_lim      struc      ;define structure
savebp1        dw         ?   ;caller's base pointer
10 saveret1     dd         ?   ;return offset and
;segment pushed by BASIC
ymax          dw         ?   ;maximum abscissa
; (screen or data coordinate)
ymin          dw         ?   ;minimum abscissa
15 ; (screen or data coordinate)
xmax          dw         ?   ;maximum ordinate
; (screen or data coordinate)
xmin          dw         ?   ;minimum ordinate
; (screen or data coordinate)
20 frame_lim    ends

; subroutine scaler(indata_ptr,outdata_
; ptr,numval)

25 frame_scl    struc      ;define structure
savebp2        dw         ?   ;caller's base pointer
saveret2       dd         ?   ;return offset and
;segment pushed by BASIC
numval         dw         ?   ;number of values in
30 ;BASIC's data array
outdata_ptr     dw         ?   ;scaled values are
;passed to a BASIC
;array pointed to by
;this pointer(for fgraph)
35 indata_ptr   dw         ?   ;values to be graphed
;are passed from a BASIC

```

- 266 -

```

;array pointed to by
;this pointer.

```

```

frame_scl      ends

```

5

```

;.....subroutines' code begins here

```

```

10  cseg_gr segment 'code'
    dgroup group data,stack,const,heap,memory
        ;defining link to BASIC

```

15

```

;-----;
; subroutine dwindow(xmin,xmax,ymin,ymax) ;
;-----;
;subroutine to establish data value
;limits
;corresponding to screen window.

```

20

```

dwindow      proc    far
              public  dwindow
25            ;public symbols allow external references
              assume  cs:cseg_gr,ds:dgroup
              ;BASIC defines regs
              assume  ss:dgroup,es:dseg_wind

```

30

```

              push    bp      ;save base pointer for the
                              ;return to BASIC
              mov     bp,sp    ;point stack pointer at frame
                              ;structure

```

35



- 267 -

```

;.....save additional registers and
;          set up extra data seg
push    ax      ;
push    di      ;
5      push    es      ;

mov      ax,dseg_wind    ;set up extra data
                                ;segment as assumed
mov      es,ax          ;

10

;.....get specifications for window from
;    BASIC and store locally

15      mov      di,[bp].ymax    ;
mov      ax,[di]          ;
mov      ymax_d,ax        ;

mov      di,[bp].ymin    ;
20      mov      ax,[di]          ;
mov      ymin_d,ax        ;

mov      di,[bp].xmax    ;
mov      ax,[di]          ;
25      mov      xmax_d,ax        ;

mov      di,[bp].xmin    ;
mov      ax,[di]          ;
30      mov      xmin_d,ax        ;

;.....restore all registers which
;          were corrupted
pop      es              ;
35      pop      di              ;
pop      ax              ;
```

- 268 -

```

                pop    bp            ;restore BASIC base
                                ;pointer before returning
5                ret     8          ;delete 4 parameter
                                ;addresses (8 bytes) from
                                ;stack and return to
                                ;calling routine

dwindow endp
10

                ;-----;
                ; subroutine swindow(xmin,xmax,ymin,ymax) ;
                ;-----;
15                ;subroutine to establish absolute screen
                ;coordinate limits
                ;corresponding to screen window.

swindow        proc    far
20                public  swindow ;public symbols allow
                external references
                assume  cs:cseg_gr,ss:dgroup
                ;BASIC defines regs
                assume  ds:dseg_wind,es:dgroup
25

                push    bp          ;save base pointer for the
                                ;return to BASIC
                mov     bp,sp       ;point stack pointer at frame
30                ;structure

                                ;.....save additional registers and
                                ;          set up extra data seg
35                push    ax        ;
                push    cx        ;

```

- 269 -

```

        push    dx        ;
        push    di        ;
        push    ds        ;

5         mov     ax,dseg_wind    ;set up extra data
                                   ;segment as assumed
        mov     ds,ax        ;

10        ;.....get specifications for window from
        ;    BASIC and store locally
        ;.....first y coordinate ranges
        mov     di,es:[bp].ymax ;
        mov     ax,es:[di]      ;
15        cmp     ax,199        ;make sure ymax_s <=199
        jg      y_bad         ;use default value if
                                   ;value sent is bad
        mov     ymax_s,ax      ;

20        mov     di,es:[bp].ymin ;
        mov     ax,es:[di]      ;
        mov     ymin_s,ax      ;

        ;.....y range limits examined
25        add     ax,8          ;make sure that ymax
                                   ;exceeds ymin by at least 8
        cmp     ax,ymax_s      ;
        jng     y_ok          ;if ymax_s <= ymin_s+8
y_bad:    mov     ax,199        ;then set ymax_s,ymin_s
30        ;to default values
        mov     ymax_s,ax      ;ymax_s default=199
        xor     ax,ax          ;ymin_s default=0
        mov     ymin_s,ax      ;

35        ;.....x coordinate ranges set up

```

- 270 -

```

y_ok:  mov     di,es:[bp].xmax ;
        mov     ax,es:[di]      ;
        cmp     ax,639          ;make sure xmax_s <=639
        jg      x_bad           ;use default value if
5                                     ;value sent is bad
        mov     xmax_s,ax       ;

        mov     di,es:[bp].xmin ;
        mov     ax,es:[di]      ;
10      mov     xmin_s,ax       ;

        ;.....x range limits examined
        cmp     ax,xmax_s       ;make sure that xmax
        ;exceeds xmin
15      jnge    x_ok            ;if xmax_s < xmin_s
        x_bad:  mov     ax,639   ;then set xmax_s,xmin_s
        ;to default values
        mov     xmax_s,ax       ;xmax_s default=199
        xor     ax,ax           ;xmin_s default=0
20      mov     xmin_s,ax       ;

        ;.....set up the pointers to the
25      ;          four screen corners

        ; --ymin
        x_ok:   xor     dx,dx    ;put lowest screen
        ;memory location (=0) into dx
30      mov     ax,ymin_s       ;first calculate y
        ;contribution to offset of
        shr     ax,1           ;upper corners by
        ;multiplying (ymin/2) by 80.
        jnc     y0_even        ;if ymin was not even
35      mov     dx,2000H       ;then the upper corners
        ;are odd pixels (2000H)

```

- 271 -

```

y0_even:mov    cl,80          ;[promised
                                ;multiplication by 80]
                                ;
mul            cl              ;
add            dx,ax           ;y contribution to
                                ;offset is here
5
mov            ulh_cor,dx      ;save partial result
mov            urh_cor,dx      ;

                                ; --ymax
10 xor          dx,dx          ;put lowest screen
                                ;memory location (=0) into dx
mov            ax,ymax_s       ;first calculate y
                                ;contribution to offset of
shr            ax,1            ;lower corners by
15                                ;multiplying (ymax/2) by 80.
jnc            yl_even         ;if ymax was not even
mov            dx,2000H         ;then the upper corners
                                ;are odd pixels (2000H)
yl_even:mov     cl,80          ;[promised
20                                ;multiplication by 80]
mul            cl              ;
add            dx,ax           ;y contribution to
                                ;offset is here
mov            llh_cor,dx      ;save partial result
25 mov            lrh_cor,dx      ;

mov            ax,xmin_s        ;x contribution is
                                ;xmin/8
mov            cl,3            ;calculated by shifting
30                                ;right 3 bits
shr            ax,cl           ;and
add            ulh_cor,ax       ;adding the result to
                                ;the stored partial result
add            llh_cor,ax       ;
35
mov            ax,xmax_s        ;x contribution is

```

- 272 -

```

                                xmin/8
                                ;calculated by shifting
                                ;right 3 bits
                                ;and
5      mov     cl,3
      shr     ax,cl
      add     urh_cor,ax        ;adding the result to
                                ;the stored partial result
      add     lrh_cor,ax        ;

10      ;.....restore all registers which
      ;      were corrupted
      pop     ds                ;
      pop     di                ;
      pop     dx                ;
15      pop     cx                ;
      pop     ax                ;

      pop     bp                ;restore BASIC base
20      ret     8                ;pointer before returning
                                ;delete 4 parameter
                                ;addresses (8 bytes) from
                                ;stack and return to
                                ;calling routine

25  swindow endp

                                ;-----;
                                ; subroutine clrwindw                ;
30      ;-----;
                                ;subroutine to clear
                                ;the screen window.

      clrwindw   proc     far
35      public   clrwindw ;public symbols allow
                                ;external references

```

- 273 -

```

        assume  cs:cseg_gr,ss:dgroup
        ;BASIC defines regs
        assume  ds:dseg_wind,es:screen_memory

5
        push    bp        ;save base pointer for the
                           ;return to BASIC
        mov     bp,sp      ;point stack pointer at frame
                           ;structure

10
                           ;.....save additional registers and
                           ;          set up data segments
        push    ax         ;
15        push    bx         ;
        push    cx         ;
        push    dx         ;
        push    si         ;
        push    di         ;
20        push    ds         ;
        push    es         ;

                           ;.....set up data segments as
25        ;          assumed
        mov     ax,dseg_wind ;
        mov     ds,ax         ;
        mov     ax,screen_memory;
        mov     es,ax         ;

30
                           ;.....clear screen by zeroing out
                           ;          graphics memory
                           ;          register usage:
35        ;          ax - marker for
                           ;          rightmost column
```

- 274 -

```

;          bh - # x bytes
;          bl - pixel mask
;          cx - y
;          coordinate counter
5          dx - # y lines
;          si - offset of
;          top of column
;          di - offset of
;          present byte
10          ;.....first clear leftmost part of window
mov        dx,ymax_s      ;compute number of
                          ;vertical lines
sub        dx,ymin_s      ;
inc        dx              ;and save in dx
15          mov        ax,urh_cor      ;compute number of
                          ;horizontal bytes
sub        ax,ulh_cor      ;(a number 1-79)
mov        bh,al           ;and save in bh
20          xor        ax,ax           ;clear ax register to
                          ;indicate clearing of all
                          ;columns except the
                          ;rightmost one

25          ;.....set up to blank leftmost
;          column
mov        cx,xmin_s      ;compute mask for
                          ;blanking leftmost column
30          call       mask0          ;

lea        di,even_pixels ;get offset of
add        di,ulh_cor      ;upper left hand corner
                          ;of window
35          mov        si,di          ;save location in si

```



- 275 -

```

;.....blank all columns except
;      rightmost
nxt_col:call  clr_col      ;
5          xor    bl,bl      ;subsequent columns
                        ;blank all bits (bl mask=0)
          inc    si          ;compute offset of
                        ;present column
          mov    di,si        ;and load into di
10         dec    bh          ;see if there are any
                        ;columns left
          jnz    nxt_col      ;

;.....blank rightmost column
15         mov    cx,xmax_s    ;compute mask for
                        ;rightmost column
          inc    cx            ;include rightmost pixel
          and    cl,7          ;using cx mod 8
          mov    bl,0FFH      ;put mask in bl
20         jz     mask_r       ;if cx mod 8 <>0 then
          shr    bl,cl         ;shift mask
                        ;appropriately
          jmp     lst_clr      ;
mask_r: xor    bl,bl          ;set bl mask to blank
25         ;all bits
lst_clr:call   clr_col        ;clear rightmost column

;.....restore all registers which
30         ;      were corrupted
          pop    es            ;
          pop    ds            ;
          pop    di            ;
          pop    si            ;
35         pop    dx            ;
          pop    cx            ;

```

- 276 -

```

        pop    bx          ;
        pop    ax          ;

5         pop    bp          ;restore BASIC base
                                ;pointer before returning
        ret     0          ;delete 0 parameter
                                ;addresses (0 bytes) from
                                ;stack and return to
10        ;calling routine
        clrwindw endp

15        ;-----;
        ; subroutine axes          ;
        ;-----;
                                ;subroutine to draw a box
                                ;enclosing the screen window.

20        axes          proc    far
                                public axes      ;public symbols allow
                                                ;external references
                                assume cs:cseg_gr,ss:dgroup
25        ;BASIC defines regs
                                assume ds:dseg_wind,es:screen_memory

        push    bp          ;save base pointer for the
30        ;return to BASIC
        mov     bp,sp      ;point stack pointer at frame
                                ;structure

35        ;.....save additional registers and
        ;          set up data segments

```

- 277 -

```

    push    ax      ;
    push    bx      ;
    push    cx      ;
    push    dx      ;
5   push    si      ;
    push    di      ;
    push    ds      ;
    push    es      ;

10
    ;.....set up data segments as
    ;          assumed
    mov     ax,dseg_wind    ;
    mov     ds,ax           ;
15   mov     ax,screen_memory;
    mov     es,ax           ;

    ;.....draw box screen by setting
20   ;          appropriate bits
    ;          register usage:
    ;          ax - marker for
    ;          rightmost column
    ;          bh - # x bytes
25   ;          bl - pixel mask
    ;          cx - y
    ;          coordinate counter
    ;          dx - # y lines
    ;          si - offset of
30   ;          top of column
    ;          di - offset of
    ;          present byte
    ;.....first calculate number of
    ;          vertical, horizontal counts
35   mov     dx,ymax_s      ;compute number of
    ;          ;vertical lines
```

- 278 -

```

sub    dx,ymin_s    ;
inc    dx            ;and save in dx

mov    ax,urh_cor    ;compute number of
5      ;horizontal bytes

sub    ax,ulh_cor    ;(a number 1-79)
mov    bh,al         ;and save in bh

10      ;.....left edge of box
lea    di,even_pixels ;get offset of
add    di,ulh_cor     ;upper left hand corner
                        ;of window

15      mov    cx,xmin_s    ;compute mask to draw
                        ;left end of top line
call   mask0          ;[mask0 gives pixels to
                        ;left of x coordinate]
xor    bl,0FFH        ;[requiring
20      ;complementation here]
or     es:[di],bl     ;

mov    cx,xmin_s      ;compute mask for
                        ;setting leftmost box edge
25      call   mask1        ;
call   drw_ln         ;draw the left most
                        ;border of the box

lea    di,even_pixels ;get offset of
30      add    di,llh_cor    ;lower left hand corner
                        ; of window
mov    cx,xmin_s      ;compute mask to draw
                        ;left end of bottom line
call   mask0          ;[mask0 gives pixels to
35      ;left of x coordinate]
xor    bl,0FFH        ;[requiring

```

- 279 -

```

                                ;complementation here]
or      es:[di],bl      ;

5      ;.....bottom edge of box
mov     bl,bh           ;save number of
                                ;horizontal bytes in bl
call    hbar           ;draw horizontal bar

10     ;.....top edge of box
mov     bh,bl           ;get number of
                                ;horizontal bytes from bl
lea     di,even_pixels  ;get offset of
15     add     di,ulh_cor ;upper left hand corner
                                ;of window
call    hbar           ;draw horizontal bar

20     ;.....right edge of box
lea     di,even_pixels  ;get offset of
add     di,urh_cor      ;upper left hand corner
                                ;of window

25     mov     cx,xmax_s  ;compute mask to draw
                                ;right end of top line
call    mask0           ;
or      es:[di],bl      ;

30     mov     cx,xmax_s  ;compute mask for
                                ;setting rightmost box edge
call    mask1           ;
call    drw_ln          ;set rightmost box edge

35     lea     di,even_pixels ;get offset of
add     di,lrh_cor      ;lower right hand corner
                                ;of window

```

- 280 -

```

    mov     cx,xmax_s      ;compute mask to draw
                           ;right end of bottom line
    call    mask0          ;
    or      es:[di],bl     ;

5
                           ;.....restore all registers which
                           ;           were corrupted
    pop     es             ;
10    pop     ds           ;
    pop     di             ;
    pop     si             ;
    pop     dx             ;
    pop     cx             ;
15    pop     bx           ;
    pop     ax             ;

    pop     bp             ;restore BASIC base
20    ret     0            ;pointer before returning
                           ;delete 0 parameter
                           ;addresses (0 bytes) from
                           ;stack and return to
                           ;calling routine

25    axes endp

    ;-----;
30    ; subroutine scaler(indata_ptr,outdata_      ;
    ; ptr,numval)                                  ;
    ;-----;
                           ;subroutine to scale data values within
                           ;limits
35    ;corresponding to data window. As a
    ;convenience,

```

- 281 -

```

;the data is inverted so ymax_d is at
;top of
;the window (screen values increase
;towards
5      ;bottom of the screen)
;
;scaling occurs in two passes: first y
;is scaled, then x
scaler proc far
10      public scaler ;public symbols allow
                    external references
        assume cs:cseg_gr,es:dgroup
        ;BASIC defines regs
        assume ss:dgroup,ds:dseg_wind
15
        push bp      ;save base pointer for the
                    ;return to BASIC
        mov bp,sp     ;point stack pointer at frame
20
                    ;structure

        ;.....save additional registers and
        ;          set up extra data seg
25      push ax      ;
        push bx      ;
        push cx      ;
        push dx      ;
        push si      ;
30      push di      ;
        push ds      ;

        mov ax,dseg_wind ;set up extra data
                    ;segment as assumed
35      mov ds,ax      ;

```

- 282 -

```

;....get data from BASIC point by point
;           and scale according to
;           data window. (use di,bx as
5           holding registers)

mov     si,es:[bp].outdata_ptr
;get pointer for scaled data output
mov     si,es:[si]           ;pointer is now in si
10      mov     outptr,si      ;save output pointer

mov     si,es:[bp].numval
;get number of points to scale into cx
mov     cx,es:[si]           ;
15      mov     numvalt,cx     ;save value for second
                                ;pass

mov     si,es:[bp].indata_ptr
;get pointer to BASIC's array of data
20      mov     si,es:[si]     ;pointer for
                                ;input is now in si
mov     di,outptr            ;pointer for
                                ;output is now in di

25      mov     bx,ymax_s      ;put screen scale into
                                ;bx
sub     bx,ymin_s            ;

mov     ax,bx                ;use half screen scale
30      ;as a roundoff correction
shr     ax,1                 ;
mov     rndoff,ax            ;

mov     bp,ymax_d            ;put data scale into bp
35      sub     bp,ymin_d      ;
```







- 285 -

```

      xchg    ax,bx          ;save source ptr in ax,
                              ;using bx to point to
                              ;offset of destination
                              ;(which is a word)
5      cmp     bx,bx_last    ;see if we are on the
                              ;same x-coordinate
      jne     y_save        ;if not put a valid
                              ;abscissa at this coordinate
10     cmp     es:[di][bx],si ;compare yscaled value
                              ;to last yscaled value
      jle     y_more        ;stored. if y was
                              ;greater or equal then keep it
      y_save: mov    es:[di][bx],si ;else store yscaled
15                              ;value in output array
      mov     bx_last,bx    ;save current
                              ;destination pointer

      y_more: xchg    bx,ax    ;restore bx register
20
      inc     bx            ;point to next input
                              ;point
      inc     bx            ;
      loop    get_ysc       ;continue scaling x
25                              ;until counter cx is zero

                              ;.....restore all registers which
30      ;           were corrupted
      pop     ds            ;
      pop     di            ;
      pop     si            ;
      pop     dx            ;
35      pop     cx          ;
      pop     bx            ;

```

- 286 -

```

                pop    ax                ;
                pop    bp                ;restore BASIC base
                ;pointer before returning
5              ret     6                ;delete 3 parameter
                ;addresses (6 bytes) from
                ;stack and return to
                ;calling routine

scaler endp
10

                ;-----;
                ; utility routines local to the window ;
15              ; module                      ;
                ;-----;

                ;.....utility procedure for fast
                ;          clearing of vertical cols
20  clr_col proc near

                mov     cx,dx            ;set up counter for
                ;clearing first column
                clr_lp: and     es:[di],bl    ;clear a graphics byte
25              ;using mask
                xor     di,2000H          ;switch even/odd pixel
                test    di,2000H          ;if odd pixel go to
                loop    ;          statement
                jnz     go_clr            ;
30              add     di,80              ;go to next even/odd
                ;pair
                go_clr: loop   clr_lp      ;continue clearing this
                ;column
                ret                     ;
35
clr_col endp

```

- 287 -

```

;.....utility procedure for fast
;      drawing of vertical lines
5      drw_ln  proc  near

          mov     cx,dx      ;set up counter for
                              ;clearing first column
10      drw_lp: or     es:[di],bl  ;set a graphics bit
                              ;using mask
          xor     di,2000H    ;switch even/odd pixel
          test    di,2000H    ;if odd pixel go to loop
                              ;statement
15          jnz    go_drw     ;
          add     di,80      ;go to next even/odd
                              ;pair
          go_drw: loop    drw_lp  ;continue clearing this
                              ;column
20          ret              ;

      drw_ln  endp

25

;.....utility for fast drawing of
;      horizontal lines
      hbar    proc  near  ;requires di to have byte before
                              ;first byte of line
30          ;bh is used as a decrementing
                              ;byte counter for number
                              ;of bytes drawn

          dec     bh          ;check to make sure at
35          jz     hbar_ok    ;least one byte to plot
                              ;if bh=0 then done

```

- 288 -

```

    hbar_lp:inc    di            ;go to next byte
                mov    byte ptr es:[di],0FFH    ;set byte
                dec    bh            ;decrement number of
                ;bytes remaining
5                jnz    hbar_lp        ;continue if more bytes
                ;need to be drawn

    hbar_ok:ret    ;

10    hbar    endp

                ;.....utility procedure for
15    mask0    proc    near    ;computing bit mask for clears
                ;uses value in cx to compute bit
                ;mask in bl

                and    cl,7        ;using cx mod 8
20                mov    bl,0FFH    ;put mask in bl
                jz    mask0_ok    ;if cx mod 8 <>0 then
                shr    bl,cl        ;shift mask
                ;appropriately
    mask0_ok:xor    bl,0FFH        ;complement mask to set
25                ;bits to be retained
                ret

    mask0    endp

30

                ;.....utility procedure for
                ;computing bit mask for drawing
    mask1    proc    near    ;uses value in cx to compute bit
35                ;mask in bl

```

- 289 -

```

        and    cl,7            ;using cx mod 8
*        mov    bl,80H         ;put mask in bl
        jz     mask1_ok       ;if cx mod 8 <>0 then
*        shr    bl,cl          ;shift mask
5        ;appropriately
mask1_ok:ret

mask1    endp

10

cseg_gr ends
end

15    ; subroutine fgraph (data_ptr,numval,x_coord,line_type)
    ;        called from BASIC this routine graphs an array
    ;        on the screen
    ;        this routine is designed to allow rapid access
    ;        to the screen to allow
20    ;        real time graph generation.
    ;
    ;-----
    ;
25    ;        arguments passed by BASIC
    ;
    ;        data_ptr        - offset of BASIC array
    ;                          containing y-coordinates of
    ;                          points to be plotted
30    ;
    ;        numval          - number of values to plot
    ;        x_coord         - absolute (screen) x coordinate
    ;                          of first point
    ;                          succeeding values are plotted
    ;                          at succeeding pixels
35    ;
    ;        line_type       - if 0 then just plot points

```

- 290 -

```

;           if not zero this byte value
;           gives the line mask for
;           plotting various lines
;           (eg. 55H interpolates a line
5 ;           between adjacent
;           points with every other point
;           on the interpolation
;           line; in other words, a fine
;           dotted line)
10 ;
;-----

;.....screen memory definition

15 screen_memory segment at 0B800H
even_pixels db 8000 dup(?) ;pixels with
;even y-coordinates
org 2000H ;beginning of
;high screen memory
20 odd_pixels db 8000 dup(?) ;pixels with odd
;y-coordinates
screen_memory ends

25
frame struc ;define structure
savebp dw ? ;caller's base pointer
save_es dw ? ;save es on stack for
;return to BASIC
30 saveret dd ? ;return offset and
;segment pushed by BASIC
line_type dw ? ;mask for plotting
;various line types
35 x_coord dw ? ;x_coordinate of first
;point to be plotted

```



- 291 -

```

numval      dw      ?      ;number of values in
                                ;graph_data(*) array
data_ptr    dw      ?      ;values to be graphed
                                ;are passed in an array
5           ;graph_data(*) pointed
                                ;to by this pointer.

frame      ends

10
cseg      segment 'code'
dgroup    group  data,stack,const,heap,memory
                                ;defining link to BASIC
                                assume cs:cseg,ds:dgroup,ss:dgroup
15         ;BASIC defines regs
                                assume es:screen_memory      ;use extra data
                                                ;segment to access the
                                                ;screen memory

20  fgraph  proc      far
                                public  fgraph      ;public symbols allow
                                                ;external references

25         push    es      ;save BASIC's es
                                ;register
                                push    bp      ;save base pointer for
                                                ;the return to BASIC
                                mov      bp,sp    ;point stack pointer at
30         ;frame structure

                                ;.....save additional registers
                                push    ax      ;
35         push    bx      ;
                                push    cx      ;

```



- 293 -

```

                                ;to generate various
                                ;dotted/dashed lines
pixel_mask    db      ?      ;pixel mask is used to
                                ;set one pixel in the
5              ;screen memory (using an
                                ;OR instruction)

10  setup:  mov      last_di,0ffffH ;initialize last_di to
                                ;ffff
              mov      numval_t,ax  ;save number of points
                                ;to plot
              mov      si,[bp].line_type ;get line type mask
                                ;from BASIC
15          mov      ax,[si]      ;
              mov      line_mask,al ;and store lower byte in
                                ;local storage

              mov      si,[bp].x_coord ;get x coordinate of
20          ;first point from BASIC
              mov      ax,[si]      ;
              mov      bx,numval_t  ;get number of points in
                                ;order
              dec      bx           ;to compute
25          add      bx,ax          ;the last x-coordinate
              cmp      bx,640       ;x-coordinate is modulo
                                ;640
              jle      lst_x        ;if less than 640 store
                                ;value
30          sub      bx,640         ;else make less than 640
lst_x:  mov      last_x,bx          ;store last_x value for
                                ;return to BASIC

              mov      bx,seg even_pixels ;set up screen
35          ;memory as extra segment
              mov      es,bx        ; (note: cannot move an

```

- 294 -

```

;immediate direct to es)

mov     cl,al           ;get low byte of x_
                        ;coordinate
5       and     cl,7     ;modulo 8
mov     pixel_mask,80H  ;initialize pixel mask
                        ;to first bit
        jz      mask_ok  ;if x_coord mod 8 is
                        ;zero, the mask is ok
10      shr     pixel_mask,cl ;rotate mask bit to
                        ;correct position

mask_ok:mov     cl,3     ;x_coord/8 is byte
                        ;offset for pixel
15      shr     ax,cl     ;this result is termed x_
                        ;now
mov     x_now,ax        ;

mov     di,[bp].data_ptr
20      ;use [si] with offset in bx to access y
mov     si,[di]         ;coordinates in BASIC
                        ;data(*) array
mov     bx,0           ;initialize to first
                        ;element of array
25      mov     dx,[si][bx] ;get first y-coordinate
                        ;from BASIC
mov     last_y,dx      ;and initialize last_y

get_y:  mov     dx,[si][bx] ;get y-coordinate from
30      ;BASIC
mov     ax,dx          ;ax is used to calculate
                        ;screen memory offset
shr     ax,1           ;divide by two to get
                        ;rid of lsb
35      mov     cl,80    ;80 bytes per line (lsb
                        ;gives interlace)

```

- 295 -

```

        mul    cl                ;ax is offset for y-
                                   ;coord in screen memory
        add    ax,x_now          ;add offset for x-
                                   ;coordinate to y offset in ax
5        mov    di,ax            ;and put x,y offset into
                                   ;di
        test   dx,1             ;if y_coordinate was
                                   ;even
        jz     ln_beg           ;then we are ready to
10        add    di,2000H        ;plot a point or a line
                                   ;odd pixels require the
                                   ;interlace offset

ln_beg:  cmp    last_di,0ffffH   ;if last_di is not ffff
15        jne    lst_di         ;(first point)
                                   ;then go to set next
                                   ;pixel
        mov     last_di,di       ;else initialize di
                                   ;properly
20  lst_di: cmp    line_mask,0    ;if line mask is not 0
        jne     draw_line       ;then draw the
                                   ;appropriate line
        set_px: mov    al,pixel_mask ;else set pixel using OR
                                   ;with mask
25        or     even_pixels[di],al
        jmp     more            ;and go to next point

                                   ;.....drawing the required line
30
        draw_line:xchg di,last_di ;get old screen memory
                                   ;location to start
        mov     cx,last_y       ;cx will be the y
                                   ;distance to current pixel
35        sub    cx,dx           ;dx still has current y-
                                   ;coord.

```

- 296 -

```

        jcxz    ln_done        ;if cx is zero then plot
                                ;only one point
        jg      nxt_pxu        ;if last_y>y-coord then
                                ;draw up on screen
5          ;since lowest y is at
                                ;top of screen

                                ;.....draw a line down on screen
10         ;          (increasing y)

        neg     cx             ;cx was negative
        jmp     nxt_pix        ;only plot one point per
                                ;y-coord if possible
15  dn_lp:  shl   line_mask,1   ;set up line mask for
                                ;next pixel
        jnc     nxt_pix        ;if no bits are shifted
                                ;out then no pixel here
        or      line_mask,1     ;is msb was shifted out,
20          ;now set lsb
        mov     al,pixel_mask   ;load pixel mask and
        or      even_pixels[di],al ;set pixel using
                                ;OR with mask

25         ;.....now find next pixel position
        ;          for line
        nxt_pix:xor    di,2000H ;change from high to low
                                ;memory (or vice versa)
        test    di,2000H        ;if in high screen
30          ;memory
        jnz     dn_di          ;then di points to next
                                ;pixel
        add     di,80           ;else go to next line in
                                ;lower memory
35  dn_di:  loop   dn_lp        ;do another pixel in
                                ;this line

```

- 297 -

```

        jmp     ln_done           ;plot last pixel when
                                   ;done

                                   ;.....draw a line up on screen
5      ;          (decreasing y)

        up_lp:  shl     line_mask,1    ;set up line mask for
                                   ;next pixel
                                   jnc     nxt_pxu      ;if no bits are shifted
10      ;out then no pixel here
        or      line_mask,1    ;is msb was shifted out,
                                   ;now set lsb
        mov     al,pixel_mask    ;load pixel mask and
        or      even_pixels[di],al    ;set pixel using
15      ;OR with mask

                                   ;.....now find next pixel position
                                   ;          for line
        nxt_pxu: xor     di,2000H      ;change from high to low
20      ;memory (or vice versa)
        test    di,2000H      ;if in low screen memory
        jz      up_di         ;then di points to next
                                   ;pixel
        sub     di,80          ;else go to next line in
25      ;upper memory
        up_di:  loop     up_lp        ;do another pixel in
                                   ;this line
        ;      jmp     ln_done      ;plot last pixel when
                                   ;done(statement not needed
30      ;          here)

                                   ;.....finish up with line by
35      ;          storing current data

```

- 298 -

```

ln_done:shl    line_mask,1    ;set up line mask for
                                ;next pixel
                                ;if no bits are shifted
                                ;out then no pixel here
                                ;is msb was shifted out,
5          or    line_mask,1    ;now set lsb
                                ;load pixel mask and
                                ;set pixel using
                                ;OR with mask
                                ;save present y-
                                ;coordinate
10   end_pix:mov  last_y,dx
                                ;save present
                                ;pixel byte pointer
                                ;.....prepare for next point if
                                ;there is one
15
                                ;one less point left now
more:   dec    numval_t
20      jz     finish          ;finished if none left
                                ;if not done increment
                                ;base index by 2 to point
                                ;to next y-coord in
                                ;BASIC array
                                ;move pixel mask to next
                                ;x-coord
25      shr    pixel_mask,1
                                ;if mask points to some
                                ;pixel get the y-coord
                                ;otherwise set up mask
                                ;for next 8 x-coordinates
30      mov    pixel_mask,80H
                                ;x_now points to next
                                ;byte (for next 8 pts)
                                ;fix last di to point to
                                ;present column
                                ;there are only 80 bytes
                                ;per line, so
35      cmp    x_now,80

```



- 299 -

```

        jl      go_gety      ;if x_now<80 then x_now
                               ;is ok to get next y
        mov     x_now,0      ;otherwise wrap around
                               ;to x_now=0
5         sub    last_di,80   ;also reset di to first
                               ;column
        go_gety: jmp     get_y      ;

10        ;.....finish up and send present
        ;           pointers,mask to BASIC

        finish: mov     al,line_mask    ;get present line mask
        xor     ah,ah      ;zero upper byte
15         mov     si,[bp].line_type ;and
        mov     [si],ax     ;send to BASIC
        mov     ax,last_x   ;get last x-coordinate
        mov     si,[bp].x_coord ;and send to BASIC
        mov     [si],ax     ;

20        ;.....restore all registers which
        ;           were corrupted

        pop     di          ;
25         pop     si          ;
        pop     dx          ;
        pop     cx          ;
        pop     bx          ;
        pop     ax          ;

30        ;

        pop     es          ;restore the es register
                               ;and
        pop     bp          ;restore BASIC base
35         ;pointer before returning
        ret     8           ;delete 4 parameter
```

- 300 -

;addresses (8 bytes) from  
;stack and return to  
;calling routine

5      fgraph    endp  
      cseg      ends  
                end

10

15

20

25

30

35

- 301 -

## APPENDIX C

```
5      ' CALIB - program to calibrate instruments using
      '   board#1
      ' last revision: 1985

10

      defint a-y
      ' only z denotes a real number
      dim buffer(12800)
      hrbpm=0
15      zfqlow=0.
      zfqres=0.
      zlfa=0.
      zrfa=0.
      cls

20

      'define ports on 8253
      timer0=&h704
      timer1=&h705
25      timer2=&h706
      con8253=&h707

      ' set timer modes to 16 bit square wave rate
30      '   generators
      out con8253,&h36
      out con8253,&h76
      out con8253,&hB6

35      'for testing set timer 0 to 100Hz timebase
      out timer1,164
```

- 302 -

```
out timer1,3

out timer2,0
'set timer 0 to 1280Hz timebase
5 out timer2,5
' (2.38MHz/1864) (1864=2*256+104)
                                'set timer 2 as a 1Hz
                                ' clock at
                                'startup
10 hrbpm=60                      '(gives a heart rate
                                ' signal at
                                '60bpm)
out timer0,1                    'set timer 0 as a flip
                                ' flop
15 out timer0,0                  '

' turn the gates on using the 8255 at bits 0,1,2
' on portc
20 porta=&H70C
portb=&H70D
portc=&H71E
con8255=&H71F
' port A output port B input port C output
25 ' first set all 8255 ports to output, then set
' portc to
' 0FFH
out con8255,130
out portc,&H0FF

30

' first print out the present value of the
' interrupt
35 ' vectors
locate 4,1
```

- 303 -

```
gosub 10000

' install the interrupt with a dummy buffer and
5  ' print
   ' vectors
   reseter=256
   call wrbuffer(reseter)
   reseter=128
10  call wrbuffer(reseter)
   call instint
   locate 5,1
   gosub 10000

15  ' now go through required startup subroutines
   gosub 90
   ' set up breathing signal
   gosub 70
20  ' set up heart rate variations
   gosub 50
   ' put some information on screen
   gosub 80           ' turn D/A on
   locate 1,1
25  print "commands: h(rvar),i(nt
      on),q(uit),r(beats),b(reath),c(ounts)"

   ' wait until user hits a key
30  savekey$=""
   40  while
len(savekey$)=0:savekey$=savekey$+inkey$:wend
   if savekey$="r" then gosub 50
   'print heart beats
35  if savekey$="q" then goto 9996 'quit
   if savekey$="c" then gosub 60  'print timers
```

- 304 -

```

    if savekey$="h" then gosub 70
    'set up heart rate variations unmask interrupts
    if savekey$="i" then gosub 80
    if savekey$="b" then gosub 90
5   'set up breathing signal
    savekey$=""
    goto 40

    'print present value of heartbeats
10
    50   locate 7,1
        call rdbeat(n)
        print "present heart beats are: ";n;time$
        return

15

    ' print present value of counters
    60   out control,0           'latch timer0
        tlow0=inp(timer0)
    20   thigh0=inp(timer0)
        out control,&h40       'latch timer1
        tlow1=inp(timer1)
        thigh1=inp(timer1)
        out control,&h80       'latch timer2
    25   tlow2=inp(timer2)
        thigh2=inp(timer2)
        locate 8,1
        print "timer0:
30   ";tlow0+thigh0*16;tab(20);"           timer1:
        ";tlow1+thigh1*16;
        print tab(40);"timer2: ";tlow2+thigh2*16
        return

35

    ' set up the heart rate variations
```

- 305 -

```

      '      respiratory frequency is given by
      '      1280Hz/buffer
      '      length
      '      low frequency is 1280Hz/low frequency
5      '      divider
      '
70      if numval<=0 then beep:print "setup analog
          buffer first":return
71      locate 17,1
10      print "present lfa,rfa(bpm)= ";zlf,zrfa,"at
          freqs(Hz):";zfqlow,zfqres
          input "lfa,rfa,low freq: ",zlfan,zrfan,zfqlown
          if zlfan>30. then beep:goto 71 else zlf=zlfan
          if zrfan>30. then beep:goto 71 else zrfa=zrfan
15      if zfqlown<.02 or zfrlown>zfqres then beep:goto
          71 else
              zfqlow=zfqlown
              locate 21,1
              print "mean heart rate(bpm)= ";hrbpm
20      72      locate 22,1
              input "new mean heart rate(bpm): ",newhrbpm
              if newhrbpm>150 or newhrbpm<30 then beep:goto 72
          else
              hrbpm=newhrbpm
25      'clear screen after input
              locate 17,1
              print space$(72)
              print space$(72)
              print space$(72)
30      print space$(72)
              print space$(72)

      ' now compute values for hrsetup subroutine
35      meandiv=76800#/hrbpm      '1280*60 ticks/min gives
          '      ticks/beat

```

- 306 -

```

rfascal=76800#/(hrbpm-zrfa)-76800#/(hrbpm+zrfa)
' rfascal is the total excursion
'   of respiration
lfascal=76800#/(hrbpm-zlfa)-76800#/(hrbpm+zlfa)
5      ' lfascal is the total excursion
      '   of low frequency
lowdiv=meandiv-(rfascal+lfascal)/2#

tbaserst=1280#/zfqlow
10    locate 17,1
      print "tbaserst,rfascal,lfascal,lowdiv:
          ";tbaserst;rfascal;lfascal;
      print lowdiv
      call hrsetup(tbaserst,rfascal,lfascal,lowdiv)
15
      return

20    ' print out interrupt controller parameters
80    locate 10,1
      mask=inp(&h21)
      mask=maskx or 24
      out &h21,mask
25    mask=inp(&h21)
      print "8259 IMR(interrupt mask regsiter)=
";mask;"
          =";hex$(mask)
      return

30

90    ' this subroutine will change the analog buffer
      locate 12,1
35    input "enter breathing rate (bpm): ",brate
      if brate>75 or brate<7 then beep:goto 90

```



- 307 -

```

        zfqres=brate/60#
        numval=76800#/brate
        ztincr=8*ATN(1#)/numval
        locate 12,40
5         color 31:print "calculating respiratory
signal...":color
        7
        call exstint          ' turn off interrupts
                                ' while
10         resetting buffer
        reseter=256
        call wrbuffer(reseter)
        for itime=0 to numval
            ztnow=ztnow+ztincr
15         analogval=127*(1#+SIN(ztnow))
            call wrbuffer(analogval)
        next itime
        call instint
        locate 12,40
20         print "respiratory signal active now    "
        return

25         ' exstall the interrupt and print vector
9996        cls
            locate 4,1
            gosub 10000
            call exstint
30         locate 5,1
            gosub 10000
            locate 21,1
9999        stop

35         ' subroutine to print out the interrupt vectors
```

- 308 -

```
10000  def seg=0
        print "IRQ3 @0B*4H: ";hex$(peek(&h2C));"
        ";hex$(peek(&h2D));" ";
5      print hex$(peek(&h2E));"
        ";hex$(peek(&h2F));tab(40);
        print "IRQ4 @0C*4H: ";hex$(peek(&h30));"
        ";hex$(peek(&h31));" ";
        print hex$(peek(&h32));" ";hex$(peek(&h33))
10      return

        end
```

15

20

25

30

35

- 309 -

```

                                page    66,80
; bdzint.asm - an assembler routine to handle interrupts
;               from IRQ3
; Last revision: 1 April 1985
5  ;
;
;-----;
; 8088 interrupt location      ;
;-----;

10 abs0          segment at 0      ;absolute memory segment
                                ;allows placement of
                                ;interrupt address
                                ;future timebase
15                                ; interrupt handler
                                ; resides at int 0B
IRQ3_int        dw      2 dup(?);offset value is a word

                                org      0CH*4  ;heart beat interrupt
                                ;handler resides at int
20                                ; 0C
IRQ4_int        dw      2 dup(?);offset value is a word

abs0            ends              ;

25

;-----;
; int_buffer: area to save DOS ;
;       dummy interrupt ptr    ;
;-----;

30

int_buffer      segment          ;data segment containing
35                                ;user interrupt buffer
save_int        dw      4 dup(?);offset for two DOS

```

- 310 -

```

;interrupts saved
;to be restored using
;exstint

5  int_buffer      ends      ;

;-----;
10 ; working storage for      ;
; time base interrupts      ;
;-----;

15 dseg_tbase      segment    ;data segment for timebase
; interrupt
heartbeats         dw         ? ;keep track of heart beats
; here (for debugging)
base_rate          dw         ? ;lowest divisor for heart
20 ; rate
lfa_scal           db         ? ;low frequency modulation
rfa_scal           db         ? ;high frequency modulation
tbase_ctr          dw         ? ;counter for timebase
; interrupt
25 ;(use for low frequency
; generation)
tbase_rst          dw         ? ;reset value for tbase_ctr
; used to set low frequency
tbase_ptr          dw         ? ;pointer to present analog
30 ; value
tbase_len          dw         ? ;length of analog data buffer
tbase_buffer       db 2800dup(?) ;buffer for A/D values
dseg_tbase         ends      ;

35

```

- 311 -

```

;-----;
; setup structures to allow access to;
; arguments passed by BASIC      ;
;-----;

5

; subroutine rdbeat(BASIC_beats)
frame_rd      struc      ;define the stack
                        ;structure for passing
10                        ;arguments to BASIC
savebp1        dw        ?      ;caller's base pointer
saveret1       dd        ?      ;return offset and
                        ;segment pushed by BASIC
BASIC_beats    dw        ?      ;place to return heart
15                        ;beats to BASIC
frame_rd      ends

;subroutine wrbuffer (analog)
frame_wr      struc      ;define the stack structure
20                        ; for passing
                        ;arguments from BASIC to
                        ; analog buffer
savebp2        dw        ?      ;caller's base pointer
saveret2       dd        ?      ;return offset and segment
25                        ; pushed by BASIC
analog         dw        ?      ;place to receive analog value
                        ; from BASIC
frame_wr      ends

30
;subroutine hrsetup(B_lreset,
; Brfa_scal,Blfa_scal,Bbase_
; rate)
frame_hr      struc      ;define the stack structure for
                        ; passing
35                        ;arguments from BASIC to heart
                        ; rate controls

```

- 312 -

```

savebp3    dw    ?    ;caller's base pointer
saveret3    dd    ?    ;return offset and segment pushed
                    ;    by BASIC
Bbase_rate dw    ?    ;BASIC's lowest divider for heart
5           ;    rate
Blfa_scal   dw    ?    ;BASIC's low frequency scaler
                    ;    (amplitude)
Brfa_scal   dw    ?    ;BASIC's high frequency scaler
                    ;    (amplitude)
10          B_lreset dw    ?    ;BASIC's low frequency timer
                    ;    reset value

frame_hr    ends

                    ;.....code segment begins here
15
cseg_calibs    segment 'code'
basic_dgroup   group  data,stack,const,heap,memory
                    ;defining link to BASIC
porta          equ    0700H    ;port definitions for
20             ;8255 port expander
portb          equ    0708H    ;these addresses are
                    ;decoded on the homemade
portc          equ    0710H    ;board
control        equ    0718H    ;control word in the
25             ;8255
timer0         equ    0720H    ;8253 timer0 register
timer1         equ    0721H    ;8253 timer1 register
timer2         equ    0722H    ;8253 timer2 register
con8253        equ    0723H    ;8253 control register
30

                    ;-----;
                    ; timebase interrupt handler (not accessible to;
                    ; BASIC) ;
35             ;-----;
                    ;this routine reads the A/D every timer0

```

- 313 -

```

;tick
;with the next point in the analog
;buffer

5
    tbase_int    proc    far        ;this procedure is not
                                ;made public
                                assume cs:cseg_sync,ds:dseg_
                                base,es:nothing,ss:nothing
10    push    ax        ;save registers used
                                ;during interrupt
    push    bx        ;
    push    dx        ;
    push    ds        ;
15
    mov     ax,dseg_base    ;set up segment
                                ;register for data area
    mov     ds,ax        ;

20

    ;.....increment counter used for
                                ;low frequency generation
    dec     tbase_ctr    ;decrement
25                                ; interrupt counter
    jnz     ctr_ok        ;if not zero then
                                ; continue
    mov     ax,tbase_rst ;else reload reset
                                ;value
30    mov     tbase_ctr,ax ;
    ctr_ok:
    ;.....get analog value from
    ;buffer and send to DAC

35
    mov     bx,tbase_ptr ;get pointer to
```

- 314 -

```

                                ;analog data
dec    bx                      ;
mov    al,tbase_buffer[bx]    ;get analog
                                ;      value
5
mov    dx,porta                ;send analog value
                                ;to DAC
out    dx,al                   ;

10   mov    dx,control          ;toggle the write
                                ;latch for the DAC
mov    al,6                    ;by using direct bit
                                ;reset
out    dx,al                   ;and
15   inc    al                  ;reset commands
out    dx,al                   ;

dec    tbase_ptr               ;point to next
                                ;value
20   jnz    tbase_eoi           ;if zero, reset
                                ;pointer
mov    ax,tbase_len            ;reset with buffer
                                ;length
mov    tbase_ptr,ax            ;

25
    ;.....acknowledge interrupt to
    ;      8259A
tbase_eoi: mov    al,20H        ;send EOI to 8259A
out    20H,al                  ;

30   pop    ds                  ;restore registers which
                                ;were used
pop    dx                      ;
pop    bx                      ;
35   pop    ax                  ;
iret                           ;return to place where

```



- 315 -

```
;interrupt occurred
```

```

5      debugmsg1      db      'this is the end of the time
                                base interrupt'

      tbase_int      endp

10      ;-----;
      ; heart beat interrupt handler (not accessible ;
      ; to BASIC) ;
      ;-----;

15      ;this routine updates the timer1 rate generator
      ;every heart beat with the divider necessary to
      ;generate the next heart beat
      ;
20      ;the respiratory modulation is given by a scaler
      ;      (0-255)
      ;times the present value of the respiratory
      ;      signal.
      ;the low frequency modulation is given by scaler
25      ;      (0-255)
      ;times a value selected from the respiratory
      ;      buffer.
      ;the value selected is the
      ;      (tbase_ctr/tbase_rst)*buffer_length
30      ;element

      hbeat_int      proc      far      ;this procedure is not
                                ;made public
                                assume cs:cseg_calibs,ds:dseg_tbase
35      assume es:nothing,ss:nothing

```

- 316 -

```

                    push    ax      ;save registers during
                                   ;interrupt
                    push    bx      ;
                    push    cx      ;
5                   push    dx      ;
                    push    ds      ;

                    mov     ax,dseg_tbase    ;set up segment
                                   ;register for data area
10                  mov     ds,ax      ;

                    inc     heartbeats    ;increment heart
                                   ; beat counter

15                  ;.....calculate low frequency
                    ;          modulation
                    ; (the tbase buffer is used as a trig
                    ; table here)
                    mov     ax,tbase_ctr    ;get number of
                                   ;1280Hz pulses
20                  dec     ax          ;
                    mul     tbase_len    ;scale by length
                                   ;of respiratory
                                   ; buffer
25                  div     tbase_rst    ;divided by reset
                                   ;value to get
                                   ; pointer
                    mov     bx,ax        ;to low frequency
                                   ; modulation
30                  mov     al,tbase_buffer[bx] ;get
                                   ; sinusoidal
                                   ; modulation
                    mul     lfa_scal    ;and scale
                                   ; appropriately
35                  mov     cx,ax        ;cx accumulate
                                   ;divider for
```

- 317 -

```

; 1280Hz clock

;.....calculate respiratory
;      modulation
5      mov     bx,tbase_ptr      ;get present
                                   ;respiration
                                   ;signal
      mov     al,tbase_buffer[bx] ;from
                                   ;buffer
10     mul     rfa_scal          ;scale with rfa
                                   ;scaler
      add     cx,ax              ;and add to cx

      add     cx,base_rate       ;finally add base
15     ;rate to get
                                   ; value for
                                   ;timer1 (heart
                                   ;rate
generator                                ; on
20     ; 8253)

;..... send new divider to 8253
;      timer
      mov     al,76H             ;set timer 1 to
25     ;square wave
                                   ; generator

      mov     dx,con8253         ;
      out     dx,al              ;

30     mov     dx,timer1         ;send divider to
                                   ;timer1
      mov     al,cl              ;low byte first
      out     dx,al              ;
      mov     al,ch              ;high byte next
35     out     dx,al              ;

```

- 318 -

```

;.....acknowledge interrupt to
;          8259A
mov     al,20H ;send EOI to 8259A
out     20H,al ;

5
    pop     ds      ;restore registers and
    pop     dx      ;
    pop     cx      ;
    pop     bx      ;
10    pop     ax      ;
    iret          ;return to place where
                  ;interrupt occurred

15    debugmsg2    db      'this is the end of the heart
                        beat interrupt'

    hbeat_int      endp

20

;-----;
; subroutine instint (install_interrupts) ;
;-----;

25    instint      proc     far
                        public instint
;public symbol allows external references
;es,ds used to access interrupt and must
30    ; be restored movsw
;uses (ds:si)(es:di) addr
    assume cs:cseg_calibs,ss:basic_
        dgroup,ds:basic_dgroup
    assume es:int_buffer

35
        ;.....save registers

```

- 319 -

```

    push    ds        ;save ds register on the
                        ; stack
    push    es        ;save es register on the
                        ; stack
5
    push    bp        ;save BASIC base pointer
                        ; for return to BASIC
    mov     bp,sp      ;point stack pointer at
                        ;frame reference to
10
                        ;address of BASIC analog
                        ;data buffer

    push    ax        ;save additional
                        ;registers
15
    push    si        ;
    push    di        ;

                        ;set up the segment registers as assumed

20
    mov     ax,int_buffer ;
    ;es points to buffer area to save
    ;DOS dummy interrupt vector
    mov     es,ax      ;
    mov     ax,0        ;ds points to
25
                        ;abs0 (interrupt table)
    mov     ds,ax      ;
    assume  ds:abs0    ;

                        ;setup access to interrupt vectors
30
    lea     di,save_int ;load offset of
                        ;save_int in es,di
    lea     si,IRQ3_int ;load offset of
                        ;IRQ3_int in ds,si
    movsw                   ;save DOS dummy
35
                        ;interrupt vectors to be
    movsw                   ;restored later
```

- 320 -

```
movsw                                ;now saving IRQ4
movsw                                ;

5  ;install the DAC timebase (IRQ3)
    mov    IRQ3_int+2,cseg_calibs
    mov    IRQ3_int,offset tbase_int;
        ;interrupt handler now
;install the heart beat (IRQ4) interrupt handler now
10    mov    IRQ4_int+2,cseg_calibs;
    mov    IRQ4_int,offset hbeat_int;

        ;.....return to BASIC
15
    pop     di        ;restore additional
                        registers
    pop     si        ;
    pop     ax        ;

20
    pop     bp        ;restore BASIC's base
                        ;pointer and
    pop     es        ;segment registers
                        before returning
25    pop     ds        ;
    ret     0         ;delete 0 parameters (0
                        ;bytes) from the stack
                        ;and return to the
                        ;calling routine
30
    debugmsg3    db    'this is the end of the
                        interrupt installation'

35    instint     endp
```

- 321 -

```

;-----;
; subroutine exstint (exstall_      ;
5      ; interrupts)                ;
;-----;

exstint      proc      far
10      public exstint ;public symbol allows
                        ;external references
      assume cs:cseg_calibs,ss:basic_dgroup
      assume ds:int_buffer,es:abs0
      ;es,ds used to access interrupt
15      ;vectors and must be restored
      ;movsw uses (ds:si)(es:di) addr

      ;.....save registers

20      push      ds      ;save ds register on the
                        ; stack
      push      es      ;save es register on the
                        ; stack
      push      bp      ;save BASIC base pointer
25      ; for return to BASIC
      mov      bp,sp    ;point stack pointer at
                        ; frame reference to
                        ;access arguments passed
                        ; by BASIC (none here)

30      push      ax      ;save additional
                        ;registers
      push      si      ;
      push      di      ;

35      ;set up the segment
      ; registers as assumed

```

- 322 -

```

mov     ax,0                ;es points to
                                ;abs0 (interrupt table)
mov     es,ax               ;
mov     ax,int_buffer       ;ds points to
5      ;buffer area to save
mov     ds,ax               ;DOS dummy
                                ;interrupt vector

                                ;setup access to interrupt vectors
10      lea     di,IRQ3_int   ;load offset of
                                ;IRQ3_int in es,di
        lea     si,save_int   ;load offset of
                                ;save_int in ds,si
movsw                   ;restore DOS
15      ;dummy interrupt vectors
movsw                   ;for IRQ3
movsw                   ;and IRQ4
movsw                   ;

20      ;.....return to BASIC

        pop     di           ;restore additional
                                registers
25      pop     si           ;
        pop     ax           ;

        pop     bp           ;restore BASIC's base
        pop     es           ;pointer and segment
30      pop     ds           ;registers before
                                ;returning
        ret     0            ;delete 0 parameters (0
                                ;bytes) from the stack
                                ;and return to the
35      ;calling routine

```



- 323 -

```

    debugmsg4      db      'this is the end of the
                        interrupt exstallation'

    exstint         endp

5

10                ;-----;
                ; subroutine rdbeat (read_heart_beats ;
                ;-----;

15  rdbeat         proc    far
                public rdbeat ;public symbol allows
                        ;external references
                assume cs:cseg_calibs,es:dseg_tbase
                assume ds:basic_dgroup,ss:basic_dgroup

20

                ;.....save registers

25                push    bp      ;save BASIC base pointer
                        ;for return to BASIC
                mov     bp,sp     ;point stack pointer at
                        ;frame reference to
                        ;access arguments passed
                        ;by BASIC (one here)

30                push    ax      ;save additional
                        ;registers
                push    es      ;
                push    di      ;

35                mov     ax,dseg_tbase ;set up segment
```

- 324 -

```
                                ;register for data area
mov     es,ax                    ;
5      mov     ax,heartbeats      ;get
                                ;beats from local memory
mov     di,[bp].BASIC_beats      ;
mov     [di],ax                  ;send
                                ;beats to BASIC
10
                                ;.....return to BASIC

15      pop     di                ;restore additional
                                registers
      pop     es                ;
      pop     ax                ;

20      pop     bp              ;restore BASIC's base
                                ;pointer,
      ret     2                 ;delete 2 parameters (4
                                ;bytes) from the stack
                                ;and return to the
                                ;calling routine

25      debugmsg5      db      'this is the end of the heart
                                beat read routine'

      rdbeat  endp

30
                                ;-----;
                                ; subroutine wrbuffer(analog)      ;
                                ;-----;

35      wrbuffer      proc     far
```

- 325 -

```
public wrbuffer    ;public symbol allows
                   ;external references
assume cs:cseg_calibs,es:dseg_tbase
assume ds:basic_dgroup,ss:basic_dgroup

5

;.....save registers

push    bp        ;save BASIC base pointer
10                                ;for return to BASIC
mov     bp,sp      ;point stack pointer at
                                ;frame reference to
                                ;access arguments passed
                                ;by BASIC (one here)

15

push    ax        ;save additional
                                ;registers
push    bx        ;
push    es        ;
20  push    si        ;
mov     ax,dseg_tbase    ;set up segment
                                ;register for data area
mov     es,ax      ;

25  mov     si,[bp].analog    ;get analog value
                                ;from BASIC
mov     ax,[si]        ;
test    ah,OFFH        ;if upper byte is
                                ;zero
30  jz      new_buff      ;then install a
                                ;new point in the
                                ;buffer
mov     tbase_len,0    ;otherwise reset
                                ;the buffer
35  mov     tbase_ptr,1    ;
jmp     wr_ret        ;
```

- 326 -

```

                                mov    bx,tbase_len    ;get present
                                                ;pointer and use
                                                ;it
5      mov    tbase_buffer[bx],al    ;to store
                                                ;buffer value
      inc    tbase_len                ;point to next
                                                ;buffer value

10      ;.....return to BASIC

                                pop     si            ;restore additional
                                                ;registers
15  wr_ret:    pop     es                ;wr_ret:
                                pop     bx            ;
                                pop     ax            ;

                                pop     bp            ;restore BASIC's base
20      ret     2                      ;pointer,
                                                ;delete 1 parameters (2
                                                ;bytes) from the stack
                                                ;and return to the
                                                ;calling routine

25  debugmsg6    db        'this is the end of the buffer
                                write routine'

                                wrbuffer    endp

30      ;-----;
      ; subroutine hrsetup(B_lreset,Brfa_scal,Blfa_scal,
      ; Bbase_rate)
35      ;-----;
```

- 327 -

```
proc    far
public hrsetup    ;public symbol allows
                  external references
assume cs:cseg_calibs,es:dseg_tbase
5         assume ds:basic_dgroup,ss:basic_dgroup

;.....save registers

10        push    bp        ;save BASIC base
                  ;pointer for return
                  ;to BASIC
        mov     bp,sp        ;point stack pointer
                  ;at frame
15        ;reference to
                  ;access arguments
                  ;passed by BASIC
                  ;(one here)

20        push    ax        ;save additional
                  ;registers
        push    es        ;
        push    si        ;

25        mov     ax,dseg_tbase    ;set up segment
                  ;register for
                  ;data area
        mov     es,ax        ;

30        mov     si,[bp].Bbase_rate ;get lowest
                  ;divisor for heart
        mov     ax,[si]        ;rate from BASIC
        mov     base_rate,ax    ;and save in local
                  ;    data
35        ;    segment
```

- 328 -

```

5      mov     si,[bp],Blfa_sacl    ;get low freq
                                   ;      modulation
                                   ;      scale
      mov     ax,[si]              ;      from BASIC
      mov     lfa_scal,al ;and save LSbyte in
                                   ;local data
                                   ;  segment

10     mov     si,[bp].Brfa_scal    ;get high freq
                                   ;      modulation scale
      mov     ax,[si]              ;from BASIC
      mov     rfa_scal,al          ;and save
                                   ;LSbyte in local data
                                   ;segment

15     mov     si,[bp].B_lreset     ;get low freq
                                   ;      timer reset value
      mov     ax,[si]              ;from BASIC
      mov     tbase_rst,ax         ;and save in
                                   ;  local data segment

20     ;.....return to BASIC

      pop     si                   ;restore additional
                                   ;registers

25     pop     es                   ;
      pop     ax                   ;

      pop     bp                   ;restore BASIC's base
                                   ;pointer,

30     ret     8                   ;delete 4 parameters (8
                                   ;  bytes) from the stack
                                   ;and return to the
                                   ;  calling routine

35     debugmsg 7      db    'this is the end of the heart rate
                           setup routine'
```

- 329 -

hrsetup           endp

cseg\_calibs       ends

5

end

10

15

20

25

30

35

- 330 -

## WHAT IS CLAIMED IS:

1. An apparatus for correcting artifacts in a series of heartbeats comprising:

5 means for collecting a series of heartbeat samples;

means coupled to said means for collecting, for selecting an appropriate interval between heartbeats;

10 means for identifying a mean variance among the intervals between heartbeat samples coupled to said means for determining;

means, coupled to said means for identifying, for establishing an acceptable range of slewing rates as a function of the mean variance;

15 means, coupled to said means for determining, for particularizing the absolute value of the slewing rate of a heartbeat sample relative to the mean interval; and

20 means, coupled to said means for particularizing, for substituting the appropriate interval between heartbeats for all heartbeat interval samples having an absolute value outside the range of acceptable slewing rates.

25

2. The apparatus as recited in claim 1 wherein said means for selecting an appropriate interval comprises means for dividing intervals having a length equal to a multiple of the appropriate interval by the multiple.

35



- 331 -

3. The apparatus as recited in claim 1 wherein said means for selecting an appropriate interval comprises means for discarding interval shorter than a predetermined length.

5

4. The apparatus as recited in claim 1 wherein said means for selecting an appropriate interval comprises means for determining a mean interval and means for substituting a mean interval for intervals having preceded by a preselected number of intervals having an absolute value outside the range of acceptable slewing rates and having an absolute value outside of the range of acceptable slewing rates.

15

5. A method for correcting artifacts in a series of heartbeats comprising the steps of:

collecting a series of heartbeat samples;

selecting an appropriate interval between heartbeats;

20

identifying variances in the intervals between heartbeats;

establishing an acceptable range of slewing rates as a function of a mean variance;

particularizing the absolute value of the slewing rate of a heartbeat sample relative to the mean interval; and

25

substituting the selected interval for all heartbeat interval samples having an absolute value outside the range of acceptable slewing rates.

30

6. The method as recited in claim 5 wherein said selecting step comprises the step of dividing intervals having a length equal to a multiple of the appropriate interval by the multiple.

35

- 332 -

7. The method as recited in claim 5 wherein said selecting step comprises the steps of determining a mean interval and substituting a mean interval for intervals having preceded by a preselected number of intervals having an absolute value outside the range of acceptable slewing rates and having an absolute value outside of the range of acceptable slewing rates.

10

8. The method as recited in claim 5 wherein said selecting step comprises the step of discarding interval shorter than a predetermined length.

15

9. Apparatus for calibrating a heart rate power spectrum monitor comprising:

means for supplying a signal simulating a heart rate;

means for generating a signal simulating a respiratory frequency fluctuation in heart rate;

means for providing a signal simulating a low frequency fluctuation in heart rate; and

means for applying signals from said means for supplying, said means for generating and said means for providing to a power spectrum monitor.

30

35

- 333 -

10. Apparatus for heart rate fluctuation  
power spectral analysis comprising:

5           means for providing an electrocardiogram  
          signal;

          means for supplying an electroplethysmogram  
          signal;

10           means, coupled to said means for providing and  
to said means for supplying, for obtaining a heart rate  
fluctuation power spectrum from an electrocardiogram  
signal and an electroplethysmogram signal; and

          relative means, coupled to said means for  
obtaining, for displaying a heart rate fluctuation power  
15   spectrum.

20

25

30

35

- 334 -

11. Apparatus for trending heart rate  
fluctuation power spectral data comprising:

5           means for providing an electrocardiogram  
signal;

          means for supplying an electroplethysmogram  
signal;

          means, coupled to said means for providing and  
10 to said means for supplying, for obtaining a heart rate  
fluctuation power spectrum from an electrocardiogram  
signal and from an electroplethysmogram signal; and

          means, coupled to said means for obtaining,  
for storing heart rate fluctuation power spectral data;

15           addressable means, coupled to said means for  
storing, for transmitting stored heart rate fluctuation  
power spectral data;

          means, coupled to said addressable means for  
transmitting, for converting heart rate fluctuation  
20 power spectral data into graphic form; and

          real time means, coupled to said means for  
converting, for displaying heart rate fluctuation power  
spectra.

25

30

35

- 335 -

12. The apparatus according to claim 11 further comprising:

means, coupled between said means for obtaining and said means for storing, for segmenting  
5 data into overlapping samples.

13. A method for treatment of a condition related to malfunctions of the cardiovascular control system in a patient comprising the steps of:

10 monitoring a power spectrum of heart rate fluctuations in the patient;

identifying a level below about 0.1 (beats/min.)<sup>2</sup> in the power spectrum of heart rate fluctuations at a frequency between about 0.04 and about  
15 0.10 Hz as indicative of cardiovascular instability; and

applying a procedure to treat the condition and thereby to increase the level of heart rate fluctuations between about 0.04 and about 0.10 Hz.

20

14. A method for treatment of a condition related to malfunctions of the cardiovascular control system in a patient comprising the steps of:

25 monitoring a power spectrum of heart rate fluctuations in the patient; and

identifying a marked increase to above about 10 (beats/min.)<sup>2</sup> in heart rate fluctuations at a frequency between about 0.04 to about 0.10 Hz as  
30 indicative of cardiovascular stress; and

applying a procedure to treat the condition and thereby to decrease the level of heart rate fluctuations between about 0.04 and about 0.10 Hz.

35

- 336 -

15. A method for treatment of a condition related to cardiovascular control system in a patient comprising the steps of:

- 5 monitoring a power spectrum of heart rate fluctuations in the patient; and
- identifying a ratio of the area under a heart rate fluctuation power spectrum of a peak at a frequency between about 0.04 and about 0.1 Hz to the area under a peak in the heart rate fluctuation power spectrum
- 10 centered at the mean respiratory rate about 0.1 Hz as having an absolute value less than 2.0 as indicative of cardiac instability; and
- applying a procedure to treat the condition and thereby to increase the ratio.

15


16. A method for treatment of a condition related to cardiovascular control system in a patient comprising the steps of:

- 20 monitoring a power spectrum of heart rate fluctuations in the patient; and
- identifying a ratio of the area under a heart rate fluctuation power spectrum of a peak at a frequency between about 0.04 and about 0.1 Hz to the area under a peak in the heart rate fluctuation power spectrum
- 25 centered at the mean respiratory rate about 0.1 Hz as having an absolute value greater than or about 50 as
- 30 indicative of cardiac instability; and
- applying a procedure to treat the condition and thereby to increase the ratio.

35

1 / 17

FIG. 1  
PRIOR ART

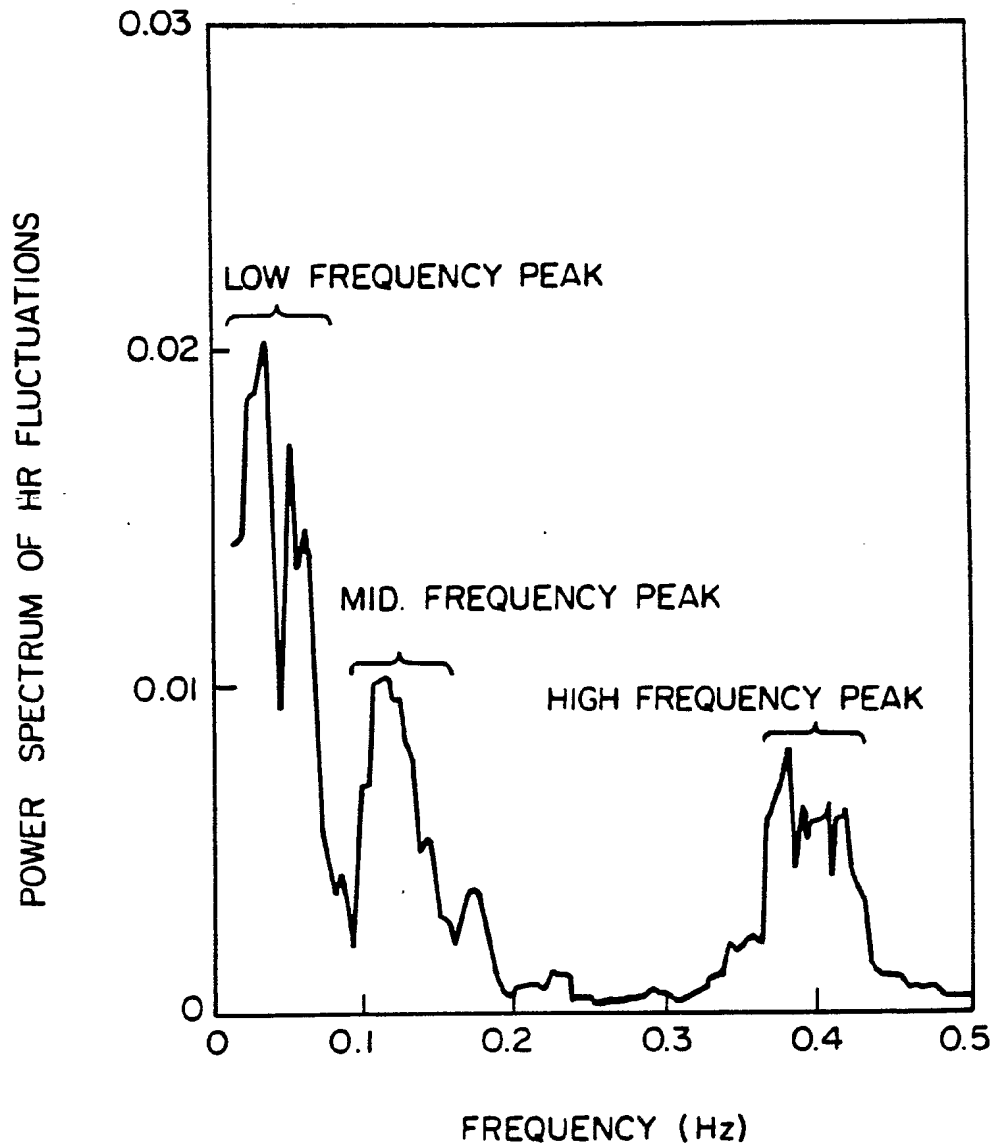

2 / 17


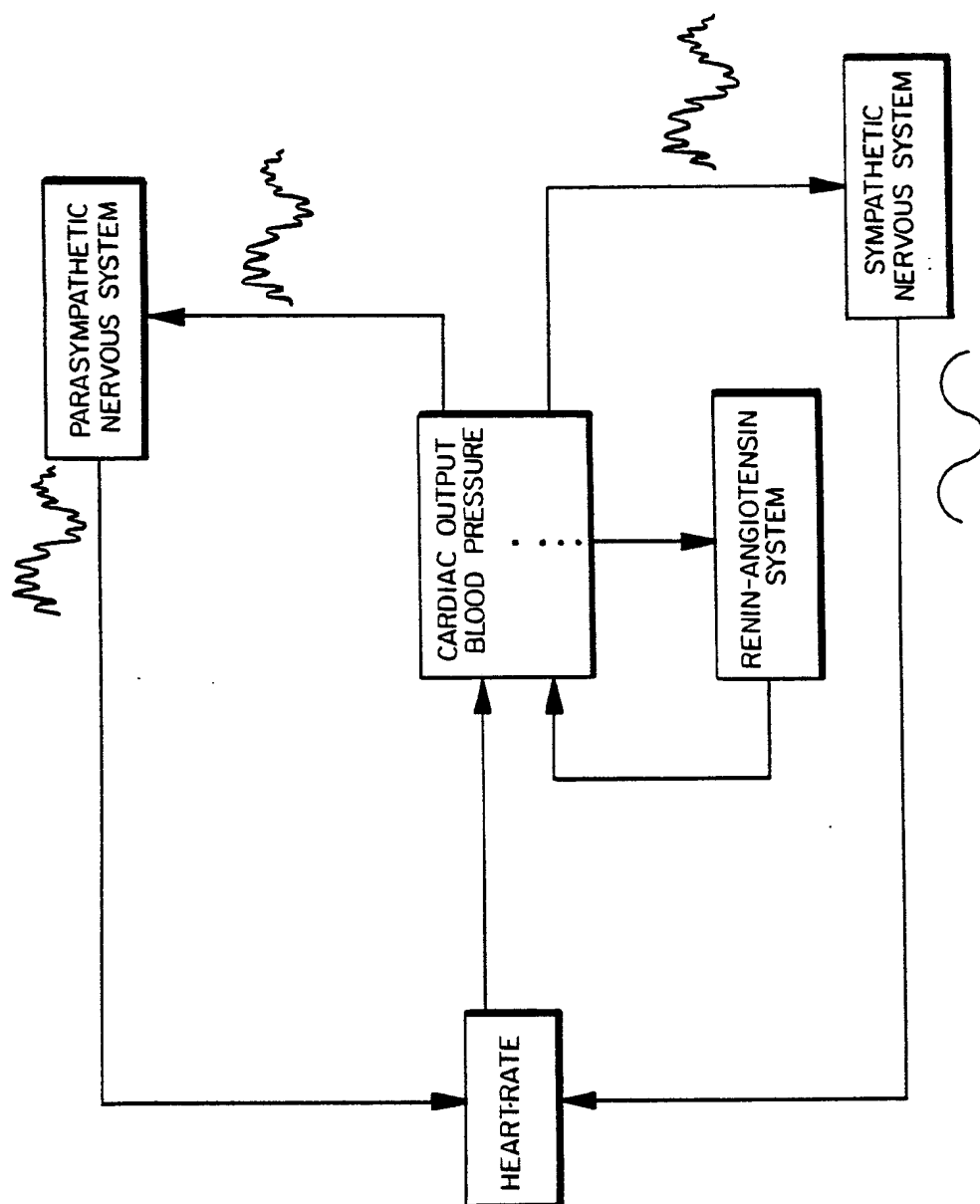


FIG. 2  
PRIOR ART



3 / 17

FIG. 3

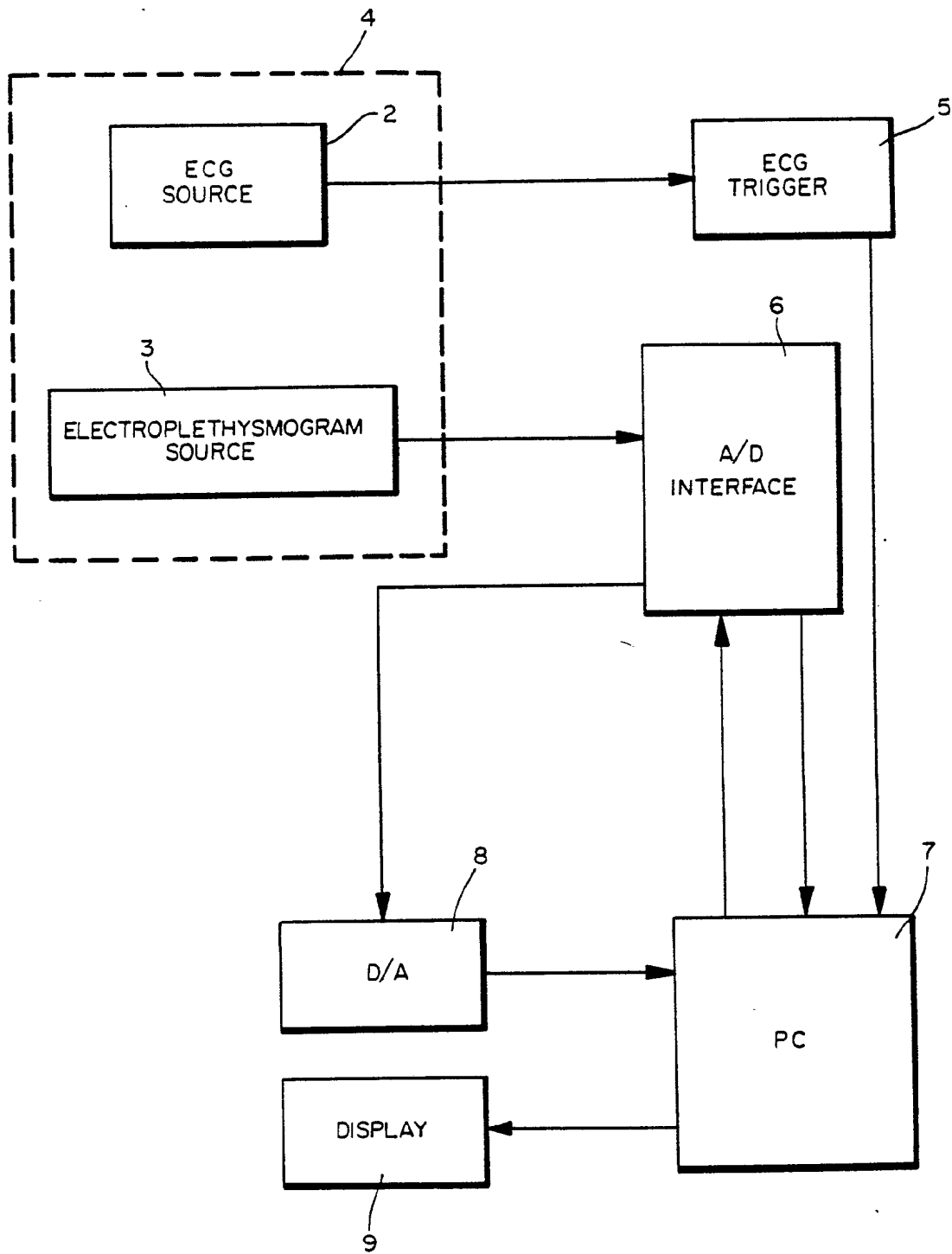

4 / 17

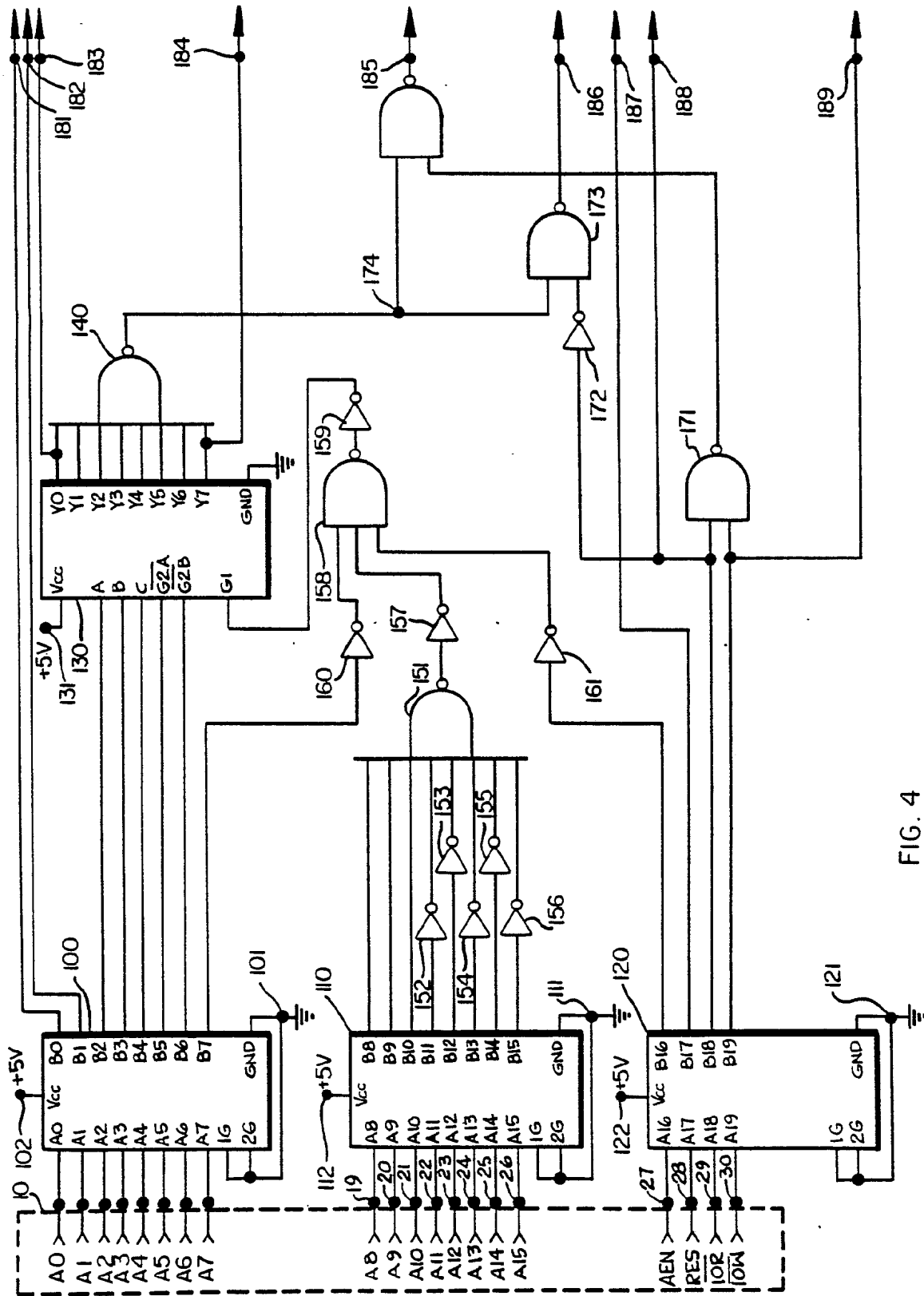


FIG. 4

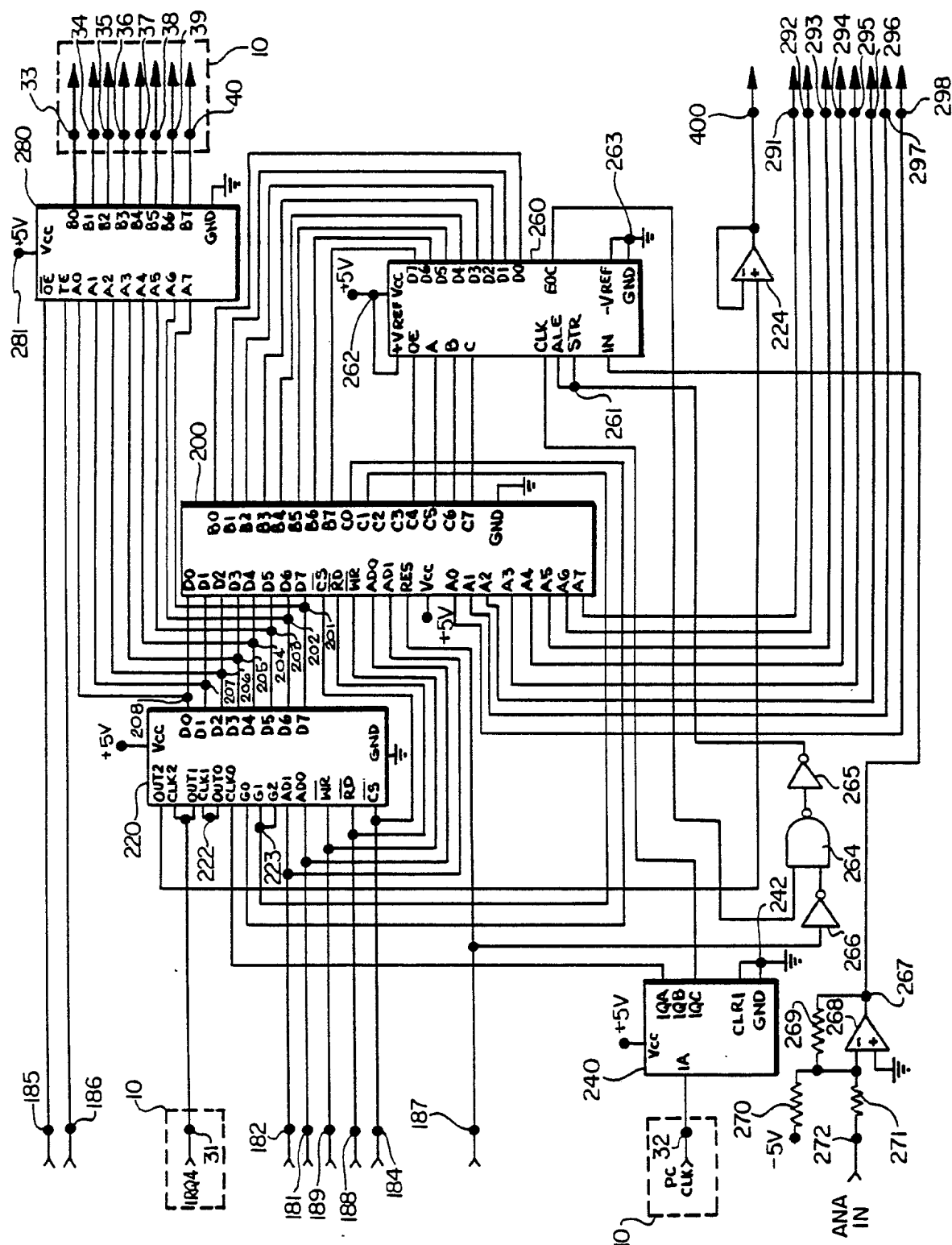


FIG. 5

6 / 17

FIG. 6



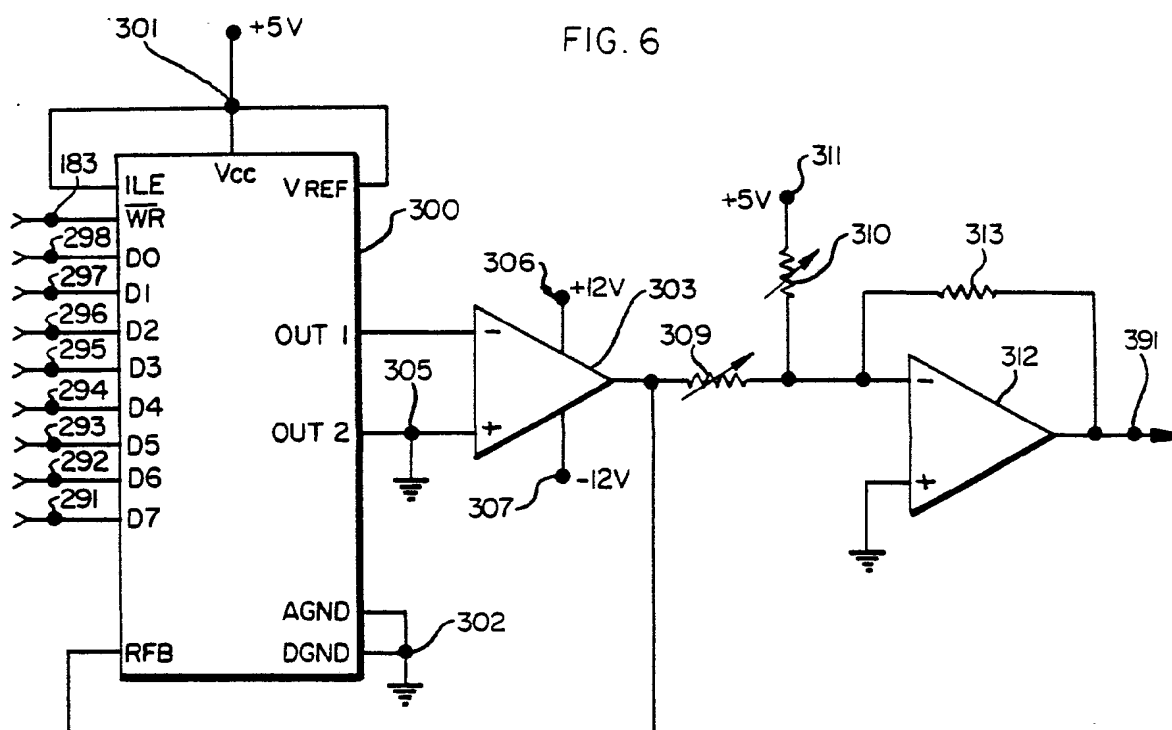


FIG. 7

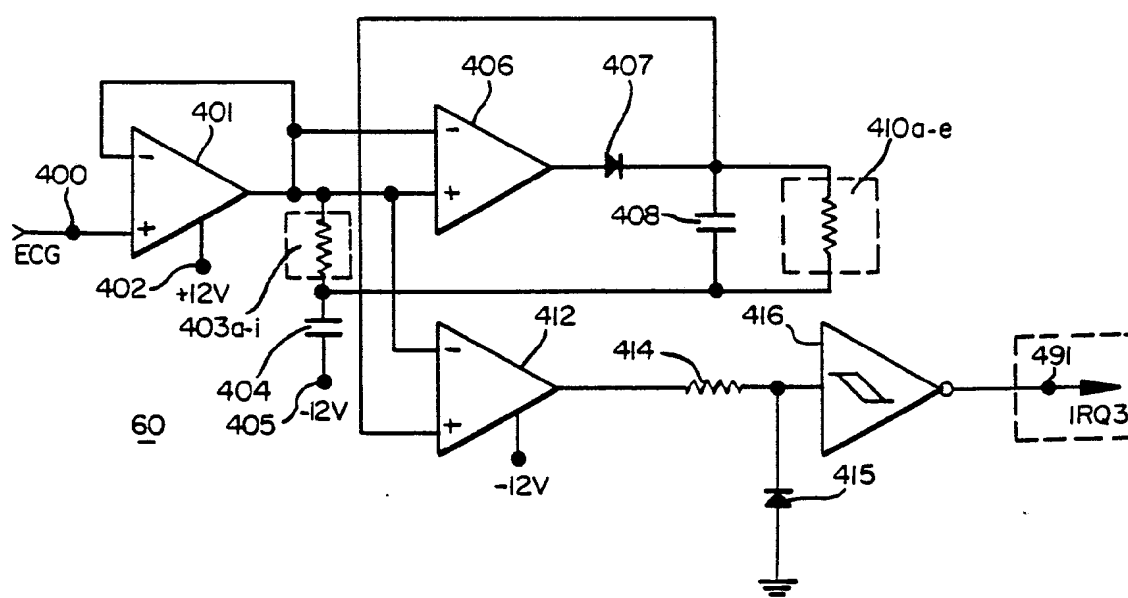


7/17

FIG. 8

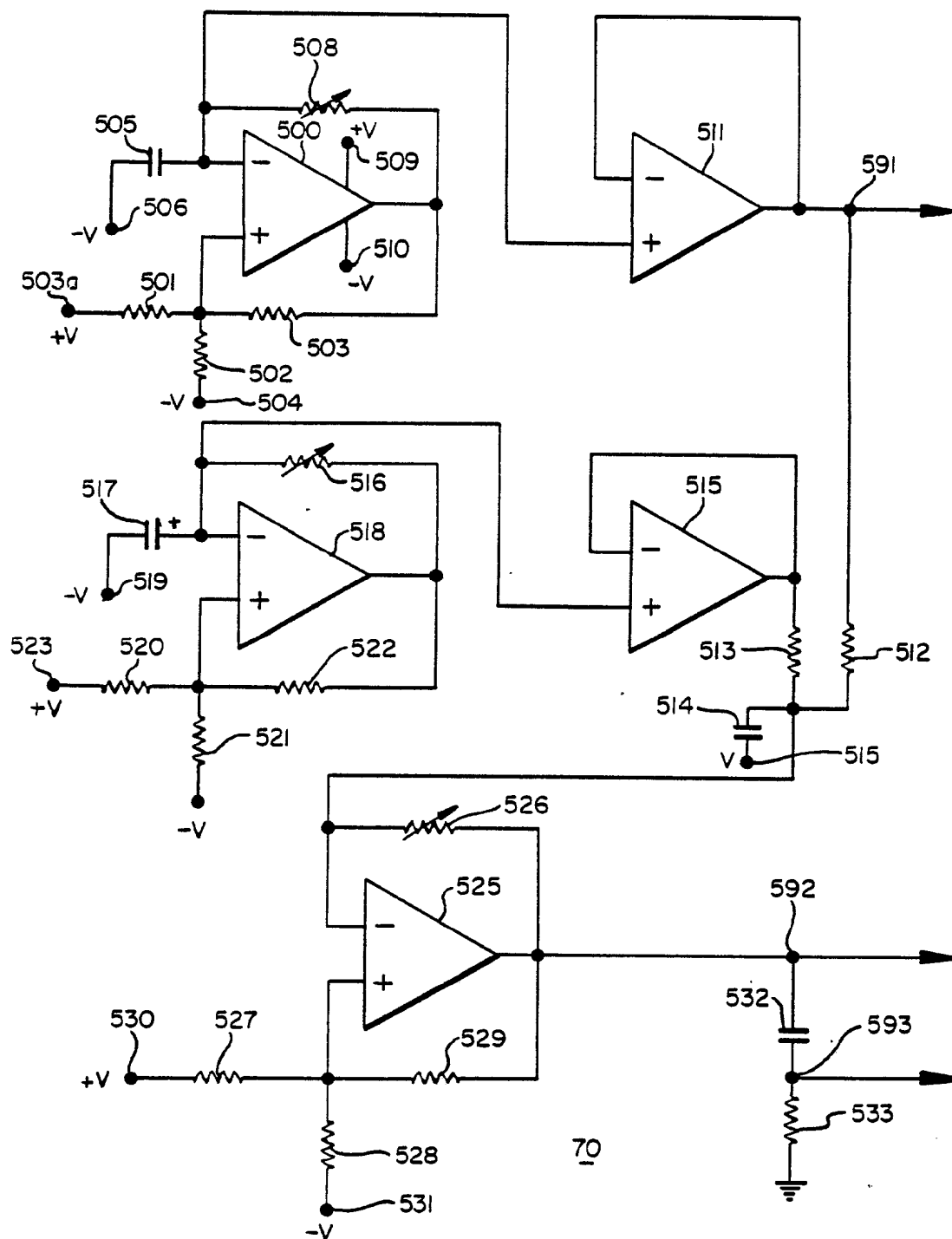

8/17

FIG. 9A

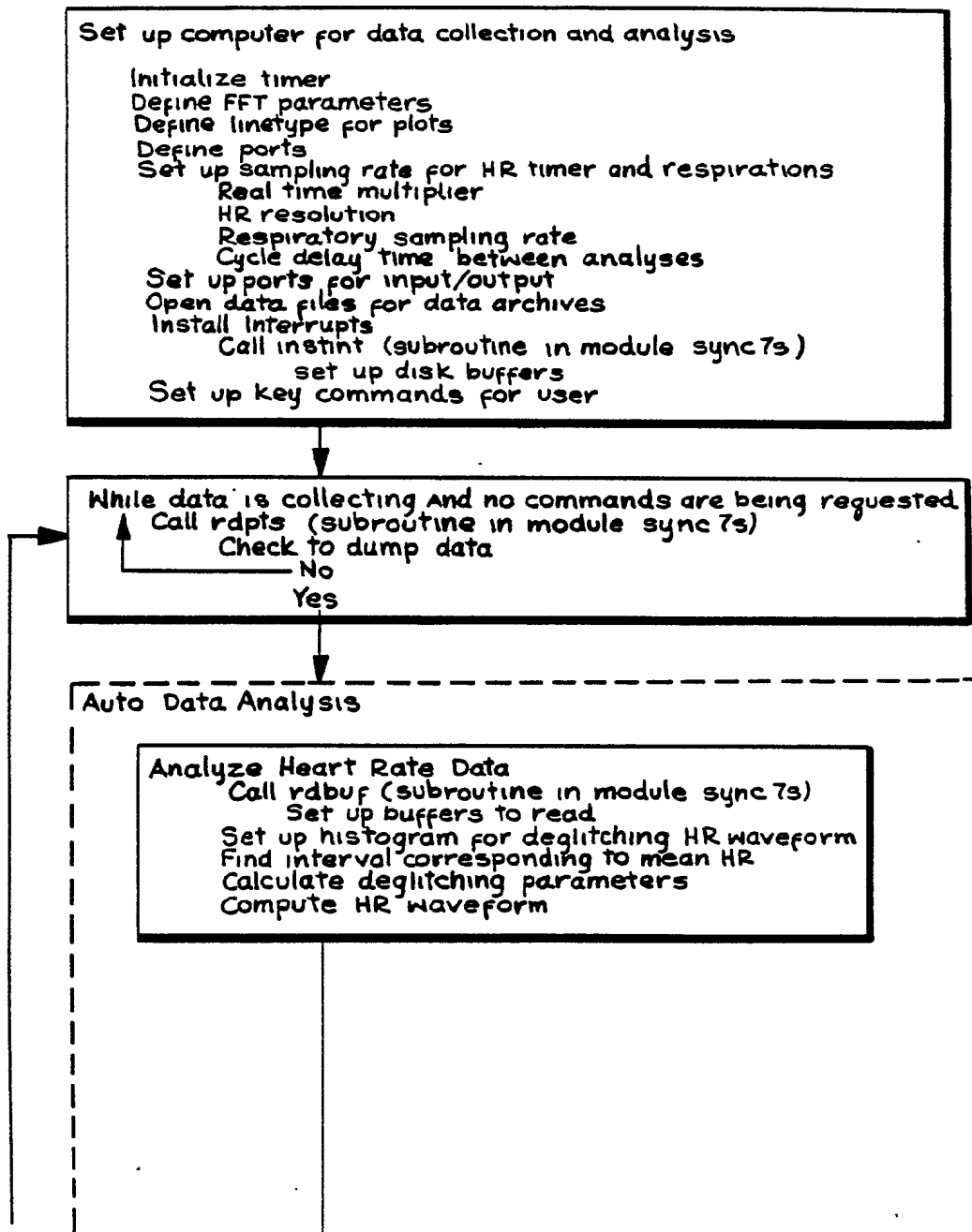

9 / 17


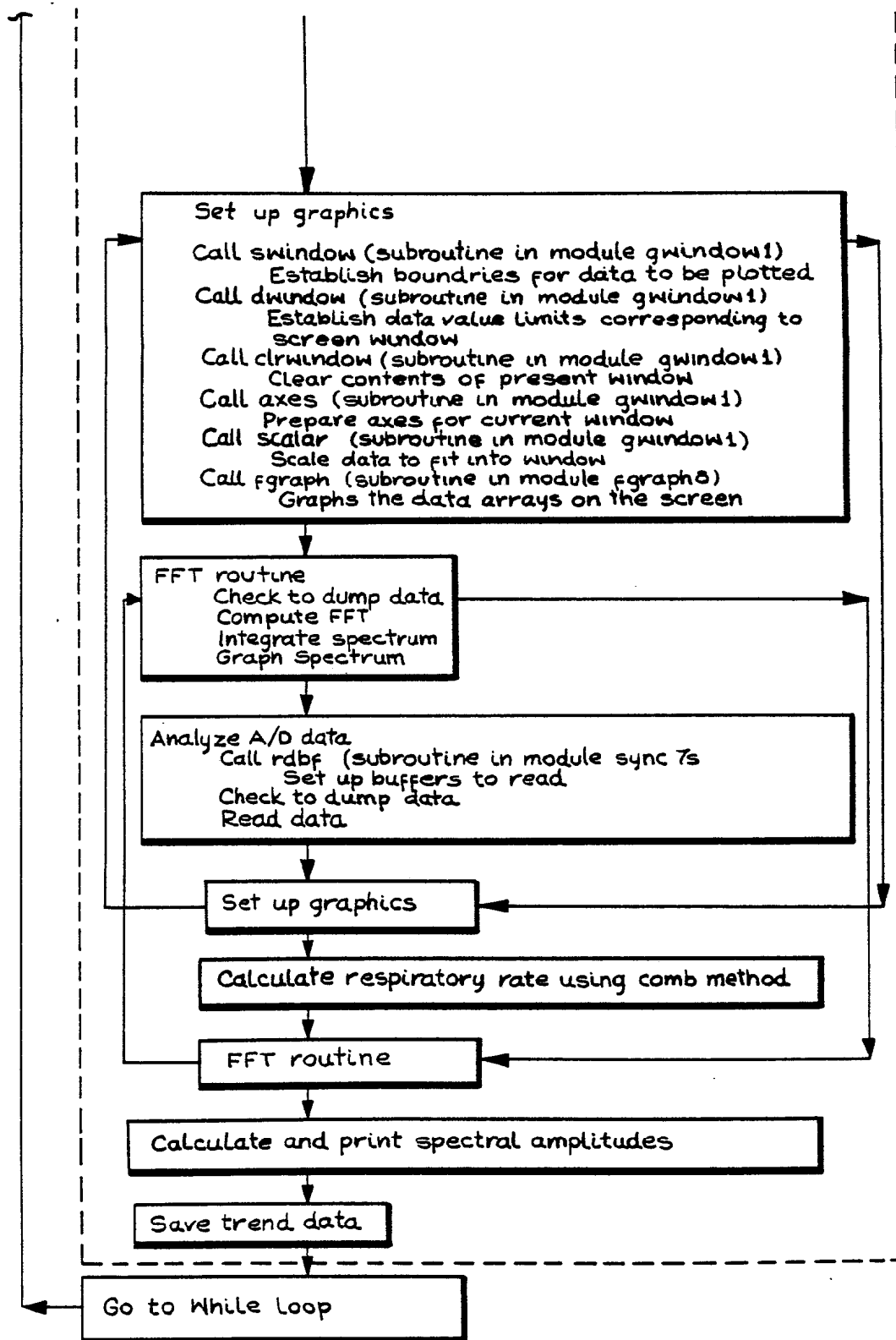


FIG. 9B

10 / 17


FIG. 10



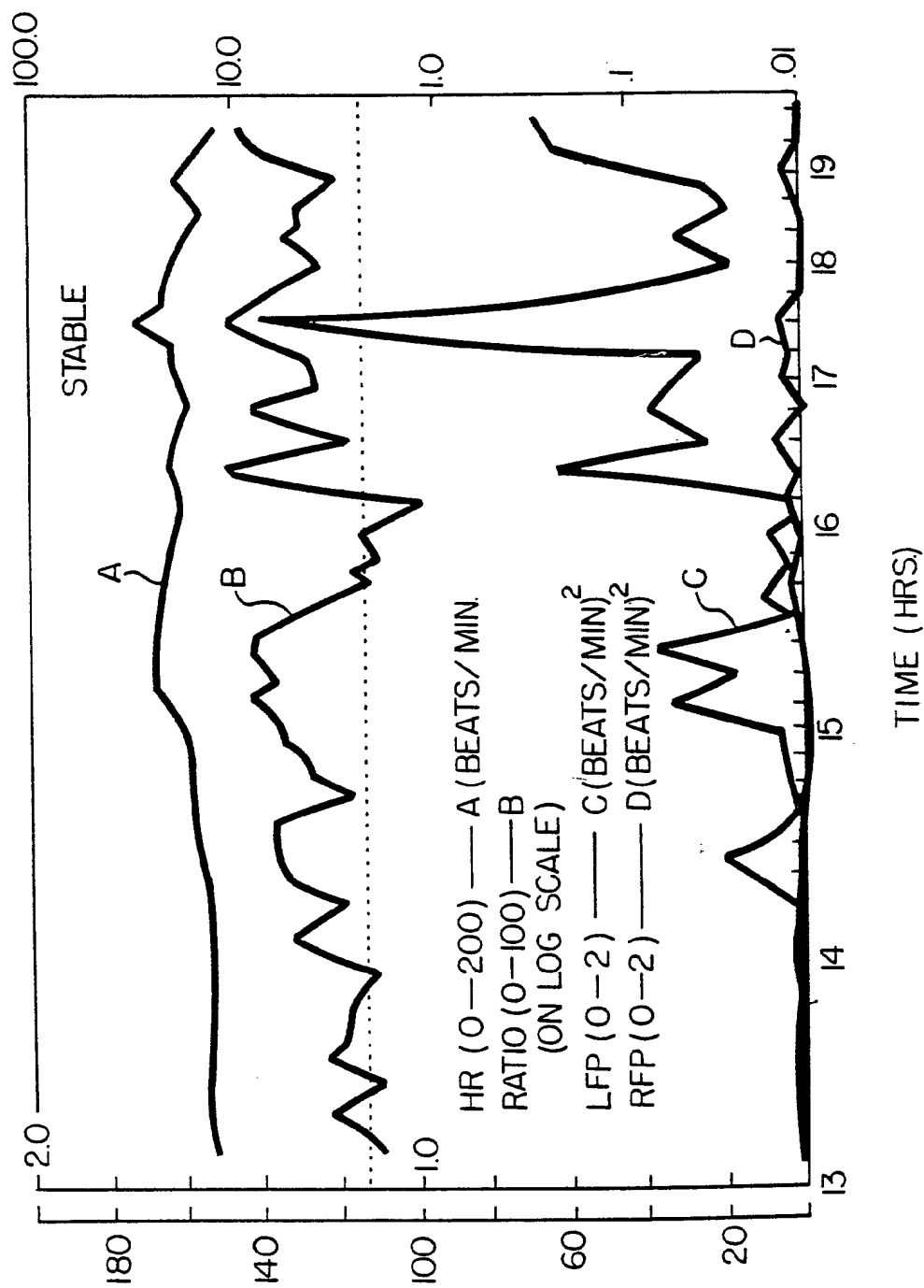

11 / 17

FIG. II

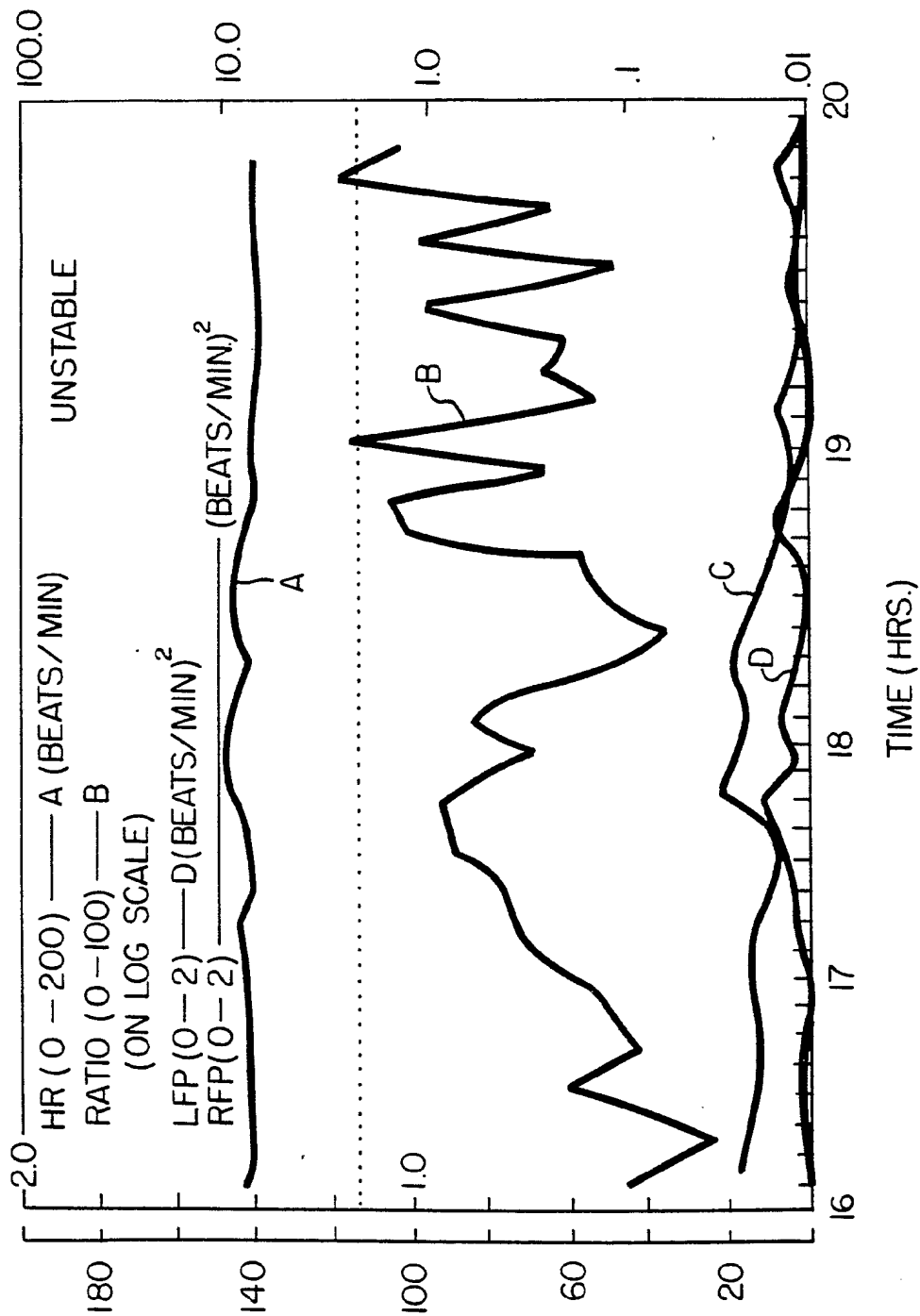

12 / 17

FIG. 12



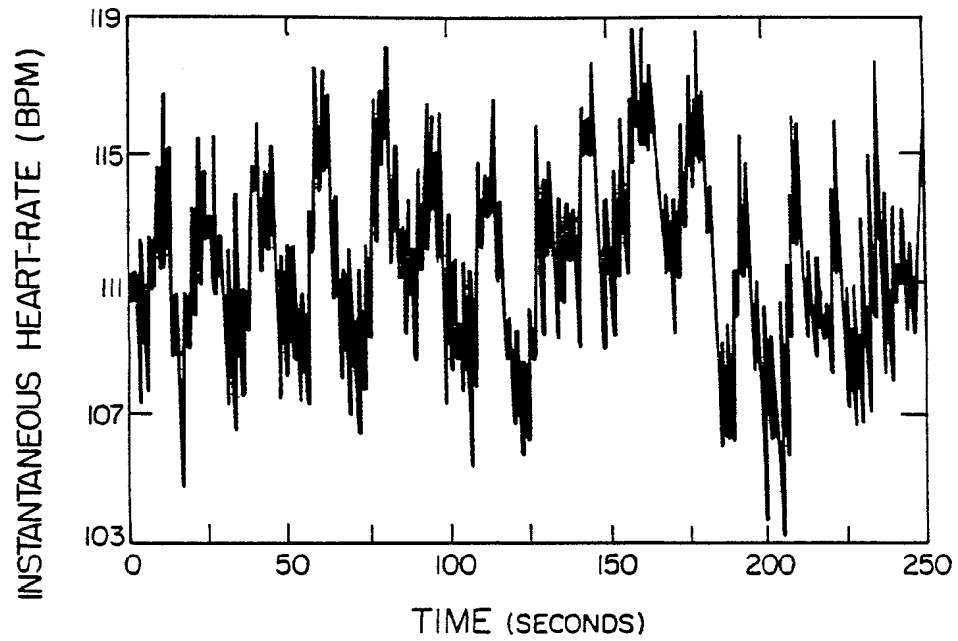


FIG. 13

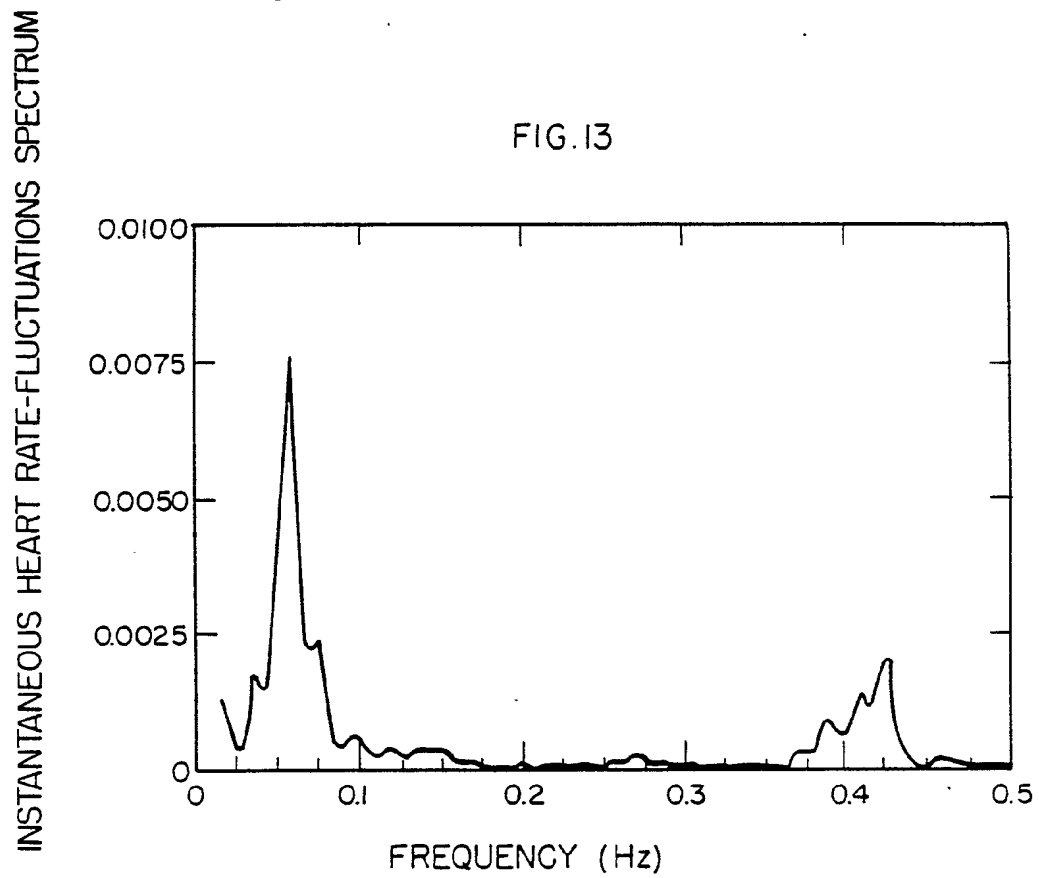


13 / 17


FIG. 14

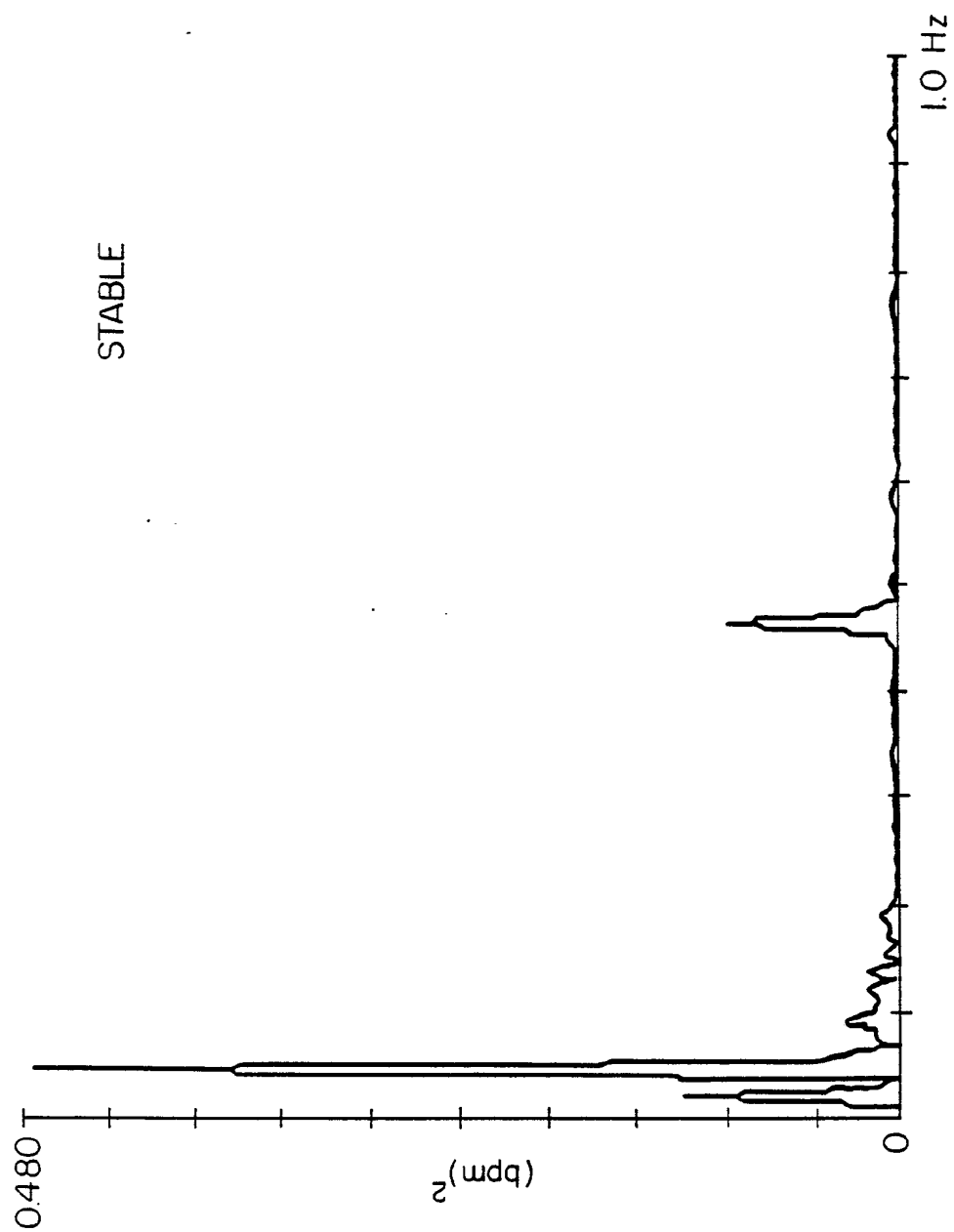

14 / 17

FIG. 15

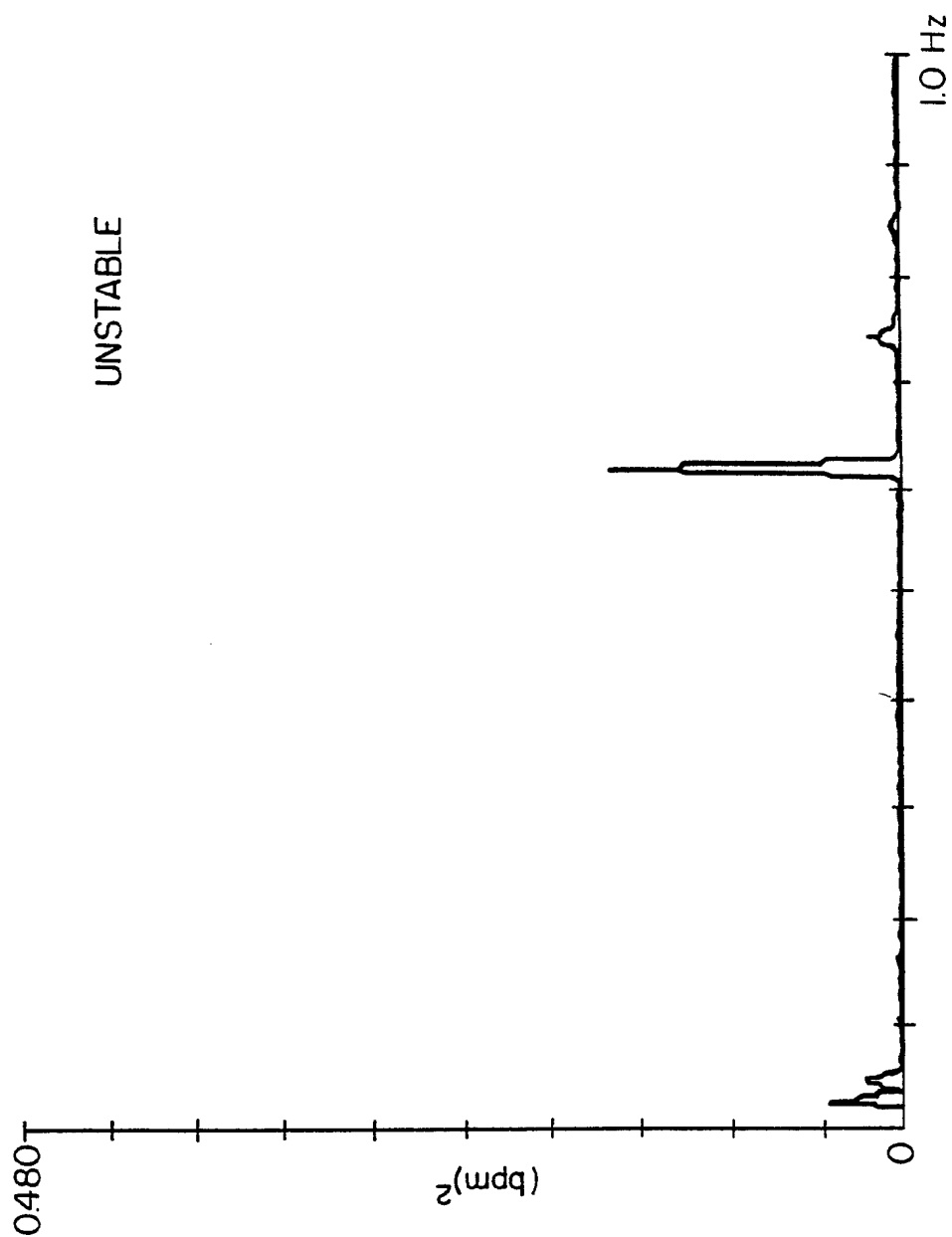

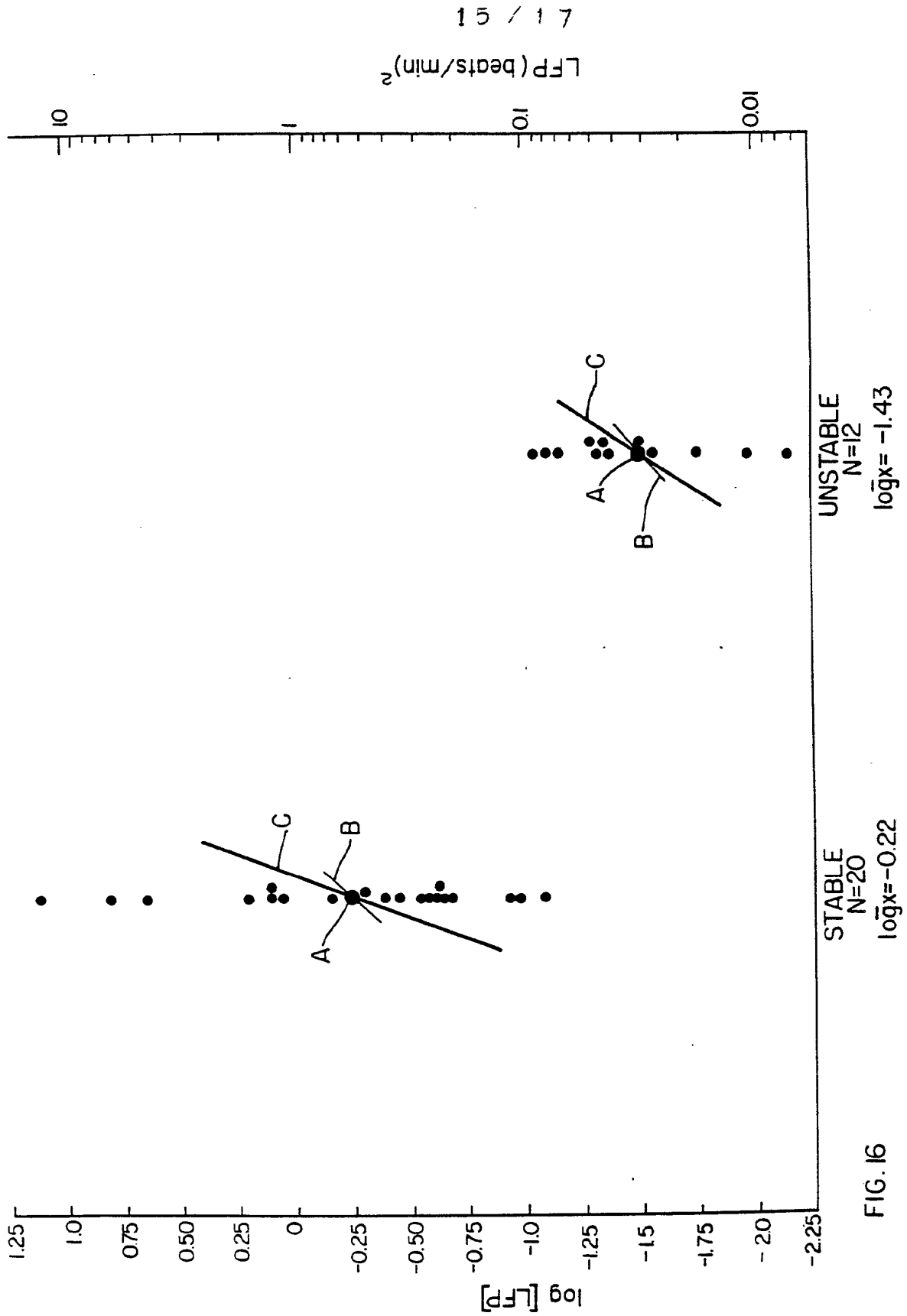

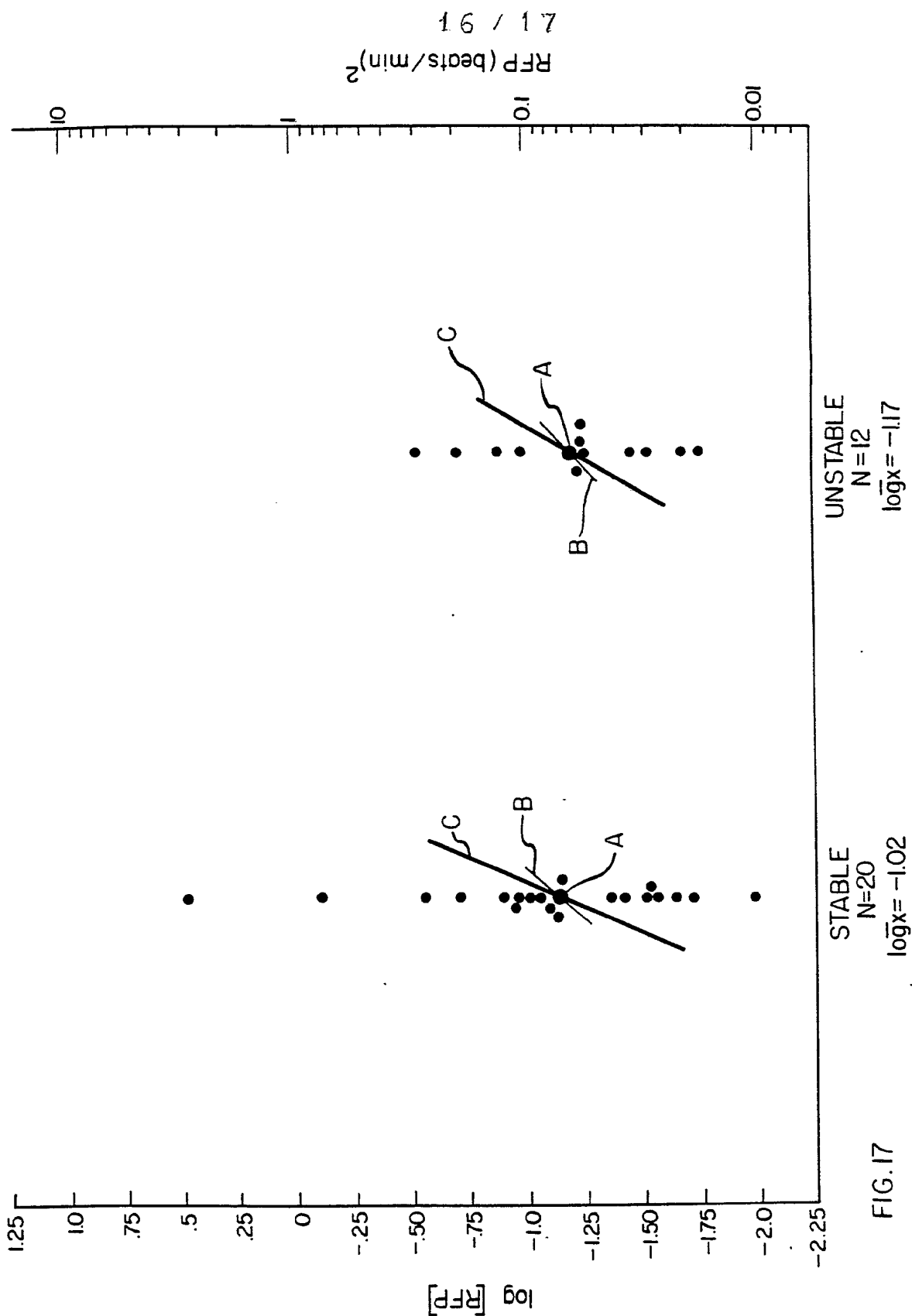


FIG. 17

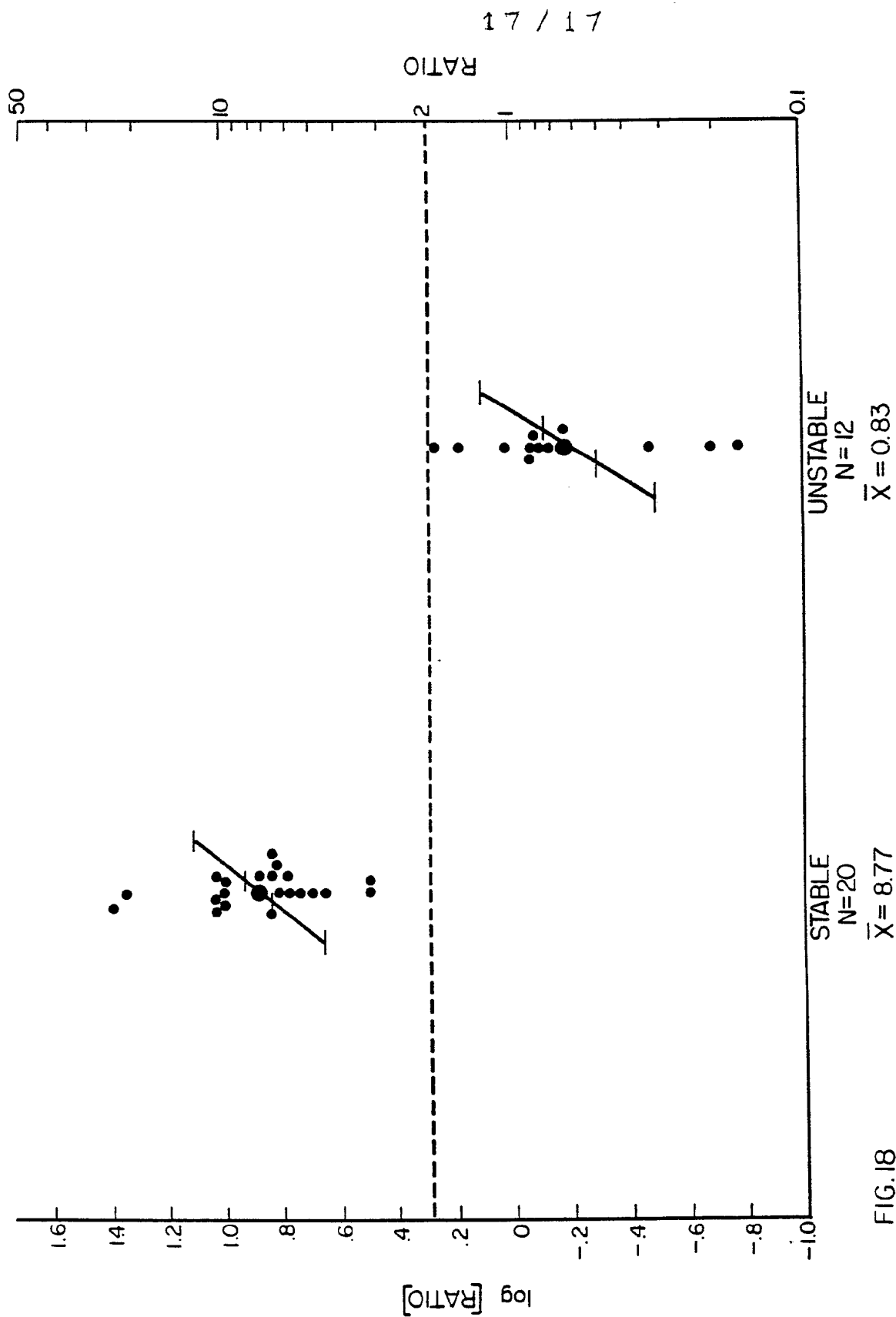


FIG. 18

# INTERNATIONAL SEARCH REPORT

International Application No PCT/US86/01193

<b>I. CLASSIFICATION OF SUBJECT MATTER</b> (if several classification symbols apply, indicate all) <sup>3</sup> According to International Patent Classification (IPC) or to both National Classification and IPC IPC (4): A61B 5/02 US Cl: 128/671		
<b>II. FIELDS SEARCHED</b>		
Minimum Documentation Searched <sup>4</sup>		
Classification System	Classification Symbols	
U.S.	128/671, 695, 696, 702, 725 364/415, 417	
Documentation Searched other than Minimum Documentation to the Extent that such Documents are Included in the Fields Searched <sup>6</sup>		
<b>III. DOCUMENTS CONSIDERED TO BE RELEVANT</b> <sup>14</sup>		
Category *	Citation of Document, <sup>16</sup> with indication, where appropriate, of the relevant passages <sup>17</sup>	Relevant to Claim No. <sup>18</sup>
A	US, A, 4,379,460 (JUDELL) 12 APRIL 1983 See entire document	
A	US, A, 4,422,458 (KRAVATH) 27 DECEMBER 1983 See entire document	
A	US, A, 4,463,764 (ANDERSON et al) 7 AUGUST 1984 See entire document	
A	US, A, 4,506,678 RUSSELL et al) 26 MARCH 1985 See entire document	
A	US, A, 4,519,395 (IRISHESKY) 23 MARCH 1985 See entire document	
A	DE, A, 2,527,475 (HOFMANN et al) 30 DECEMBER 1976 See entire document	
A	JOURNAL OF CLINICAL ENGINEERING, volume 5, number 1, issued JANUARY-MARCH 1980 (QUEST PUBLISHING CO.), J.R. PUTNAM et al, "ECG/RESPIRATION MONITOR CALIBRATOR". See entire document	
<div style="display: flex; justify-content: space-between;"> <div style="width: 45%;"> <p>* Special categories of cited documents: <sup>15</sup></p> <p>"A" document defining the general state of the art which is not considered to be of particular relevance</p> <p>"E" earlier document but published on or after the international filing date</p> <p>"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)</p> <p>"O" document referring to an oral disclosure, use, exhibition or other means</p> <p>"P" document published prior to the international filing date but later than the priority date claimed</p> </div> <div style="width: 50%;"> <p>"T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention</p> <p>"X" document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step</p> <p>"Y" document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.</p> <p>"&amp;" document member of the same patent family</p> </div> </div>		
<b>IV. CERTIFICATION</b>		
Date of the Actual Completion of the International Search <sup>1</sup>		Date of Mailing of this International Search Report <sup>2</sup>
4 AUGUST 1986		13 AUG 1986
International Searching Authority <sup>1</sup>		Signature of Authorized Officer <sup>20</sup>
ISA/US		W.E. KAMM



**FURTHER INFORMATION CONTINUED FROM THE SECOND SHEET**

A

MEDICAL AND BIOLOGICAL ENGINEERING AND  
COMPUTING, volume 20, number 3, issued  
MAY 1972, "A.J. WILSON et al, "METHODS OF  
FILTERING THE HEART-BEAT ARTIFACT FROM THE  
BREATHING WAVEFORM OF INFANTS OBTAINED BY  
IMPEDANCE PNEUMOGRAPHY".  
See entire document

**V. ☐ OBSERVATIONS WHERE CERTAIN CLAIMS WERE FOUND UNSEARCHABLE <sup>10</sup>**

This international search report has not been established in respect of certain claims under Article 17(2) (a) for the following reasons:

1. ☐ Claim numbers ..... because they relate to subject matter <sup>12</sup> not required to be searched by this Authority, namely:

2. ☐ Claim numbers ..... because they relate to parts of the international application that do not comply with the prescribed requirements to such an extent that no meaningful international search can be carried out <sup>13</sup>, specifically:

**VI. ☐ OBSERVATIONS WHERE UNITY OF INVENTION IS LACKING <sup>11</sup>**

This International Searching Authority found multiple inventions in this international application as follows:

1. ☐ As all required additional search fees were timely paid by the applicant, this international search report covers all searchable claims of the international application.

2. ☐ As only some of the required additional search fees were timely paid by the applicant, this international search report covers only those claims of the international application for which fees were paid, specifically claims:

3. ☐ No required additional search fees were timely paid by the applicant. Consequently, this international search report is restricted to the invention first mentioned in the claims; it is covered by claim numbers:

4. ☐ As all searchable claims could be searched without effort justifying an additional fee, the International Searching Authority did not invite payment of any additional fee.

**Remark on Protest**

- ☐ The additional search fees were accompanied by applicant's protest.  
☐ No protest accompanied the payment of additional search fees.